



# Smart Cupping Therapy Device


## *INSTRUCTION MANUAL*



This multifunctional smart massager is an at-home therapy device that is innovatively developed to combine traditional Chinese medicine with modern technology.

## **This device has 5 effective functions :**

1. Intelligent Cupping
2. Hot Compress
3. Red Light Pulse
4. High-Frequency Vibration
5. Magnetic Energy Field

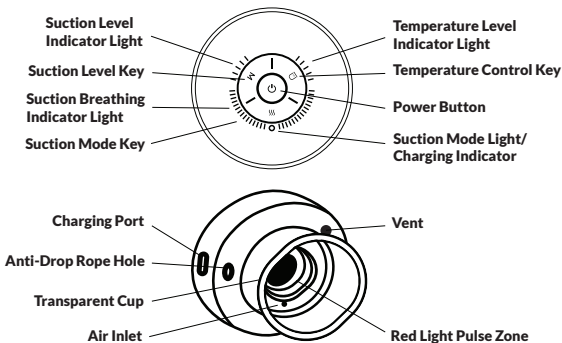

### **Power Key**

Press and hold for 1.5 seconds to power on / off

### **Suction Mode Key**

Press to switch between suction modes

### **Suction Level Key**

Press to control cupping suction level

### **Temperature Control Key**

Turn on / off red light function

## Technical Specifications

**Product Name:** BodyPro Smart Cupping Massage Therapy Device

**Product Weight:** 0.44 lbs/200g

**Product Dimension:** 3.7 in x 3.7 in x 2.83 in

**Charging Voltage:** 5V=2A

**Suction Strength:** -60~0kPa

**Charging Duration:** 2.5 hrs

**Red Light Band (Wavelength):** 650 ± 20nm

**Rated Power:** 5W

**Battery Capacity:** 3.7V/1, 800mAh

**Executive Standard:** GB4706.1-2005 GB4706.10-2008

## Preparing the Device for Use

1. To achieve optimal results, it is recommended to use wet wipes to clean and moisturize the skin before starting your therapy.
2. Place the cotton filter ring on the transparent cupping jar and re-attach the jar tightly to the device by turning it clockwise if it is unattached.
3. Attach the protection cord to the harness to prevent damage from accidental falls. You may wear the harness around your wrist during your session.
4. Make sure the device is securely attached to the transparent cupping jar and charged before turning on.



## Powering the Device On / Off

1. **Press and hold** the power button for about **1.5 seconds** to power on/power off. The device will turn off automatically after 5-12 minutes of use.
2. After powering on, its suction function will be automatically turned on. During use, **press the power button to pause**.
3. Place the device on your body. Once in contact, you will begin to feel the suction and red-light pulse therapy. The vibration massage function will automatically start as well.



## Temperature Control Button

1. Press the temperature control key while the device is powered on to adjust the heat.
2. There are 5 Levels of temperature to choose from and the default temperature is Level- 1.
3. The heating level indicator light corresponds to a different temperature level:

- ⊕ **Level 1: 38°C / 100.4° F**
- ⊕ **Level 2: 40°C / 104° F**
- ⊕ **Level 3: 42°C / 107.6° F**
- ⊕ **Level 4: 45°C / 113° F**
- ⊕ **Level 5: 48°C / 118.4° F**



## Cupping Suction Levels Button

1. Short press to adjust the strength of cupping suction (1-5 Levels). The default setting is level-1.
2. The suction intensity indicator light will correspondingly display the current strength level.
3. When the suction light indicator flashes blue and starts to move, the device is in an active working state, and the cupping suction level is displayed in real-time.
4. When you use it for the first time or if you have sensitive skin, please start with Level 1.
5. In the initial use, please pay attention to the step-by-step adjustment to find the mode and intensity that suits you.



## Using The Cupping Mode Key (M)

Press the cupping mode key (M) while the device is powered on to switch between different cupping modes.

(the default mode is "Mode 1: M1" when the machine is turned on)

### *The 4 Cupping Modes:*

- ⊕ **M1:** Vitality mode (cupping mode light is green)
- ⊕ **M2:** Soothing mode (cupping mode light is blue)
- ⊕ **M3:** Depth mode (cupping mode light is red)
- ⊕ **M4:** Professional mode (cupping mode is multi-colored)

## Additional Information

- ⊗ Short pressing the power button during therapy will stop the suction so you can remove and place it on a different part of your body. Short pressing it again will turn the suction back on.
- ⊗ Press and hold the Cupping Suction Intensity key for 1.5 seconds while the device is in use to turn the vibration on or off.
- ⊗ Press and hold the temperature level key [sss] for 1.5 seconds while the device is in use to turn the red light pulse function on or off.
- ⊗ When using this device with liquid products such as essential oils or massage rubs, make sure the liquid is applied evenly on the skin surface and is fully absorbed without any residue.
- ⊗ **Avoid looking directly at the red light**

## Charging the Device

- ⊗ Connect the Type-C charging cable to the adapter socket and into the wall, Powerbank, or computer and connect the other end to the product charging port to start charging.
- ⊗ Avoid charging the device for more than 8 hours.
- ⊗ The charging indicator light is red during charging, and green when fully charged.
- ⊗ Avoid using the device while it is charging.
- ⊗ When the battery power is too low, it will be automatically turned off after 10 warning beeps.

## Maintaining the Device

- ⊗ Wipe the body with a dry non-woven cloth. Do not use acid-alkaline cleaner, alcohol, thinner, gasoline, nail polish remover, and other solvents for cleaning.
- ⊗ Turn off the power and gently clean with a wet cloth after each use.
- ⊗ Store this product in a safe and dry place, out of reach of children.
- ⊗ Do not leave the product idle with moisture remaining on it.
- ⊗ Never use hard brushes or stainless steel wire balls to wipe the device, to avoid damaging the surface.
- ⊗ Do not submerge the device in water.
- ⊗ Do not use or store this product in places prone to vibration, exposed to direct sunlight, places with high temperature, humidity, and water, places with strong electromagnetic fields, near fireworks, etc.

# Replacing the Cotton Filter


1. Unscrew cupping jar
2. Insert cotton filter
3. Screw cupping jar tightly back in place



This machine is thoughtfully designed with a cotton filter to filter out most of the granular or liquid impurities during the use of the product to protect the internal cleaning of the fuselage as well as maintain the service life of the product. Please check and replace the filter regularly according to your frequency of use. Make sure the device is powered off before replacing the cotton filter.

## Safety Precautions

1. DO NOT use the device if you have Implanted medical electronic devices such as cardiac pacemakers.
2. DO NOT use the device if you have Life-sustaining medical electronic devices such as an artificial heart or lungs.
3. DO NOT use the device with other Medical electronic devices such as electrocardiographs.
4. All Patients with any other internal implanted medical devices should not use the device.
5. Always keep the device out of reach from children for safety and to prevent the risk of injury.
6. Avoid using the device on scars and wounds.
7. If you feel unwell during or after use, please stop using it immediately and consult a doctor.
8. Avoid using the device close to the heart to prevent the risk of heart fibrillation.
9. This product contains lithium batteries, please do not disassemble or modify it yourself, soak it in water or throw it into fire to prevent damage to the product and to avoid an explosion.
10. Avoid using the device for more than 12 minutes on one spot as it may cause a low degree burn on your skin.
11. Do not use a chipped, cracked, or damaged transparent cupping tank.
12. Consult your doctor before using the device if you suffer from any medical conditions such as an allergy to metal and red light, infectious diseases, Fever, High blood pressure, skin condition, diabetes, etc.
13. Avoid using the device if you are pregnant or menstruating.
14. Patients with malignant tumors and ulceration, traumatic fractures, varicose veins, acute suppurative inflammation should avoid using the device.
15. Avoid using the device if you suffer from severe heart disease, heart failure, respiratory failure, and severe edema.
16. Avoid using the device if you suffer from poor coagulation, a spontaneous bleeding tendency, or bleeding after injury.
17. Cupping is not suitable for people with severe skin allergies or skin infectious diseases.
18. Do not use the device on patients with high neuroticism, high fever, convulsions, manic restlessness.
19. Never use the device on facial features, private parts, or large blood vessels on the body surface.

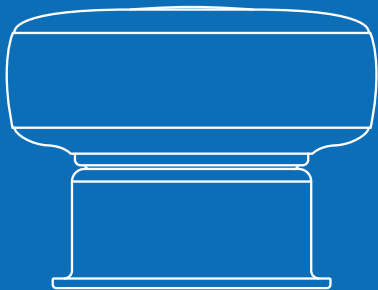

AIM FOR GREATNESS



BodyPro.us  
866.425.3209  
life@bodypro.us  
2612 Larch Lane, STE 101  
Mt. Pleasant, SC 29466