


**JC BUCKMAN**

Healthy Lifestyles

# Back Massager

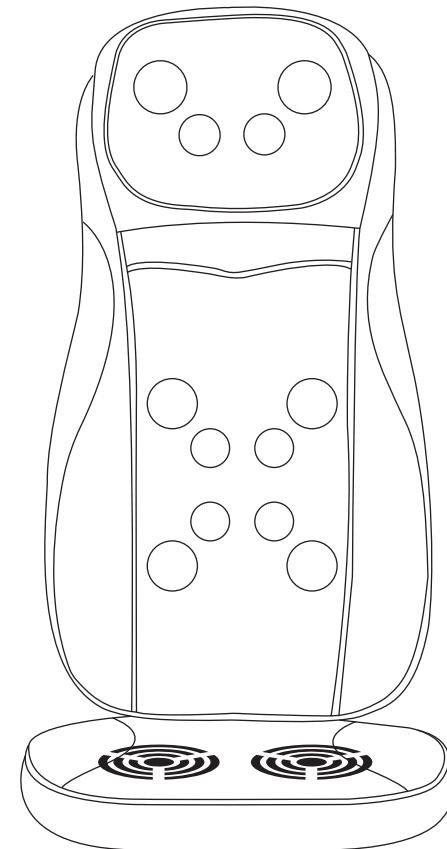

Model:OM23  
USER MANUAL

---

Dear user, firstly, thank you for using our back massager. In order to offer better services, we write this instruction manual. Please read details in this brochure carefully and pay special attention to the important Safety Instruction to ensure proper use. Please properly keep this brochure for your consultation.

All information presented is as accurate as possible but subject to change without notice, please check with the product, our company reserves the right to final explanation.

## Warning

1. Do not use this machine with unsuitable standard voltage;
2. be sure to use the plug with matching socket; .
3. Please use this product according to the user manual .
4. Please do not use this product outdoor;
5. This appliance is not intended to provide to person with body disabilities, senses and nerve defect or lack of common sense and life experience to (including children) use. Unless they get those, Who is responsible for the safety of their personnel on how to use the appliance under supervision and guidance?
6. Children must be supervised to ensure that they do not get to this equipment
7. If the power cord were damaged, in order to avoid danger, must be changed by professional personnel of the manufacturer from the maintenance department or from the similar functional departments.
8. This equipment with hot surface, the person who is not sensitive to heat must pay attention to use.

---

## Operational Environment

1. Please do not use it in a highly humid environment such as to use the products in the bathroom
2. Do not use this product in the environment of sharply changing temperature.

## Specification

Product Name: Back Massager

Model: OM23

Rated Voltage: 12V $\overline{=}$

Rated Power: 30W

## Feature

- >> Kneading and Shiatsu massage
- >> Clockwise and Anticlockwise Kneading
- >> With Airbag massage, Promote the body Relax
- >> Can be used at many places of home


## Feature



Controller



Seat Part Airbag

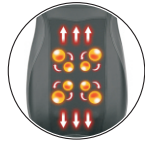

Wait Part Airbag Massage



Overheat Protection



Knead




Heating

## Care and Maintenance

1. Requests for routine maintenance at the neutral detergent, clean it use the soft cloth, do not use the corrosive cleaning agents.
2. Press the switch lightly when operating, do not overexert to avoid stress, spine item split the surface
3. Please reseave it well into the box or kept in dry places after using
4. When using, if find some symptom, please cut of the power, stop using immediately, and give to our company maintenance point, do not dismantle by yourself.
5. The massage have the temperature protetion device, it can protect your massager overheating, then please turn off the switch, stop using for 30 minutes.
6. If the fuse blew out, please stop the massage and power off.

## Back Massager

### Products Diagram

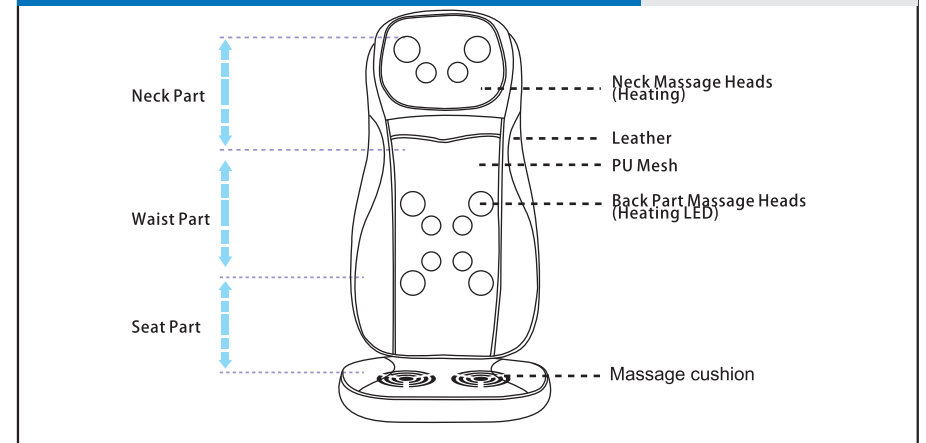

## Product simple fault analysis and processing

Phenomenon	Probable cause	Processing method
The machine doesn't work	Power supply hasn't been connected, power plug connected undesirable Whether the power switches open.	connect the power Struck the power plug well Open the switch
Suddenly shut down when using	Use for a long time overheating	Stop the machine Turn off the power
Strength is too strong or too weak	1. The machine is put angle unsuitable 2. Make too much effort, the body is too clingy massager or your temporary is not adapt to.	1. Adjust the putting angle until the dynamic is normal rate 2. Move and adjust your body until the body colose to the machine properly
Make a sound of "ZhiZhi" or "XiuXiu"	1. The parts have the mormal damage 2. The normal working voice of the massager 3. Use for a long time	No treatment

## Instructions

- On/off 

In standby mode, press , indicator goes on, chair working automatically, setting time for 20 minutes, press again, light extinguished, chair recovering to standby mode.

- Neck switch key

Manually turn on / off neck kneading

▲ In the product local working state, press this key neck massage ball local upward

▼ In the product local working state, press this key neck massage ball local downward

- Area position

▲ In the product standby or working state, press this button to turn on or off the local kneading massage function, the corresponding indicator will be lit or extinguished

▼ In the partial working state of the product, press this button to move the back part of the massage ball upward

In the partial working state of the product, press this button to move the back part of the massage ball downward

- All back massage

In the standby or working state of the product, press this button to turn on the full back massage, and then press this button to turn off the full back massage, the corresponding indicator light will be turned on or off.

- Upper back massage

In the standby or working state of the product, press this button to start the upper back massage, and then press this button to turn off the upper back massage, the corresponding indicator light will be lit or extinguished.

- Lower back massage

In the standby or working state of the product, press this button to start the lower back massage, and then press this button to turn off the lower back massage, the corresponding indicator light will be lit or extinguished.

- Heating

In the standby or working state of the product, press this key to turn on or off the heating function of the back and neck massage head, and the corresponding indicator light will be turned on or off.

## Instructions

- Seat cushion vibration

In the product standby or working state, press this key to cycle through low, frequency mid, frequency high, frequency and turn off the seat cushion vibration massage function, the corresponding indicator light will be turned on or off. (In low frequency; medium frequency and high frequency the blinking frequency of the indicator is slow, flashing, fast flashing, and no flashing. The default vibration intensity is: low, frequency)

## Controller Diagram

