MAMBA

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USER'S MANUAL

GESS®
Expert

Thank you for purchasing ultrasonic face cleansing device «Mamba»!

We reserve the right to modify the design and additions to the description of our products. Changes may be made without notice, it refers to the color and components.

Please, read the instructions for use before using the appliance. Keep the instructions in good condition for a possible appeal to it if necessary.

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Deutsch

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CONTENTS

Description and purpose	2	Preparation and carrying out of the procedure	8
Indications and contraindications.....	3	Recommendations	12
Safety precautions	5	Cleaning and care	14
Device components.....	6	Troubleshooting.....	14
Device charger	7	Disposal	15
Control	7	Specifications.....	16

DESCRIPTION AND PURPOSE

All skin types need a complex treatment, which includes cleansing, moisturizing and nourishing, lifting and toning. Ultrasonic device «Mamba» combines several types of action: ultrasonic peeling, phonophoresis, micromassage, miostimulation (EMS), roller massage with pinching function, as well as gentle cleansing using a silicone brush. The possibilities of the device make face and neckline skin care more effective.

ULTRASOUND ACTION

A course of cleansing will give your skin a fresh and well-groomed appearance. The device cleanses and narrows pores, exfoliates keratinized particles, smoothes mimic wrinkles, and tightens the facial oval.

- **Cleansing mode (Ultrasonic peeling).** An effective stimulating facial treatment for skin care at any age. The skin is being gently and slowly exfoliated and cleansed without stretching or mechanical damaging. In just a few sessions acne disappears, sebaceous glands function becomes normalized, the upper layer of the epidermis is renewed and the skin becomes noticeably softer and more elastic.
- **Nutrition mode (Nutrition/ phonophoresis).** Ultrasound makes the preparations more active. It improves the process of penetration of valuable nutrients into the deepest layers of the skin, thereby multiplying the effectiveness of cosmetic products.
- **Up mode (Lifting/Micromassage).** A deep ultrasound massage tones muscles, massages tissues at the cellular level, visibly rejuvenates the skin, and tightens the facial oval.

MIOSTIMULATION (EMS)

It stimulates muscle contractions, making the facial muscles work, it allows you to tighten skin that has lost its elasticity. The metabolic processes and blood circulation are intensified, which improves the condition of the skin. Electric muscle stimulation (EMS) has proven to be an active technique, resulting in the correction of the oval shape and strengthening of the tone of skin that has lost elasticity, reduction of «creases» and

wrinkles, strengthening of muscles, and tightening of the skin. Renewed skin becomes fresher, more elastic and tighter, wrinkles are smoothed, and the facial contour is lifted.

ROLLER MASSAGE WITH ELECTRODE HEADS.

It provides a delicate pinching massage with 3D-rollers and acupuncture effect on the skin. The rotating rollers grab the skin for effective kneading, pinching and massaging. The cooling action of the rollers is helpful in reducing swelling and puffiness. The massage tones up the muscles, thereby tightening the face and slowing down its aging.

* The massage can be done without turning the device on.

SILICONE CAP (BRUSH HEAD)

The silicone cap has soft silicone fluff, so it functions as a brush for gentle daily cleansing to keep your skin clean between ultrasound cleansings.

INDICATIONS AND CONTRAINDICATIONS

Indications for use:

- Blackheads, enlarged pores, skin inflammation;
- The presence of a keratinized layer of the skin, peeling, post-acne;
- Edema, formation of nasolabial folds;
- Disorders of the sebaceous glands: dry and dehydrated skin or excessively oily skin;
- Age-related changes: deep and mimic wrinkles, decreased skin tone, unhealthy skin color, double chin formation;
- Sagging skin;
- Scars and post-acne;

English

Deutsch

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- Loss of muscle tone;
- Face contour changing.

Contraindications:**Do not use the device:**

- on inflamed and burned skin areas;
- in case of skin allergy to cosmetics;
- in case of infectious skin lesions;
- in case of infectious diseases aggravation;
- in case of oncological diseases;
- on damaged skin areas;
- herpes and the presence of unhealed scratches;
- the presence of pathologies of the cardiovascular system;
- pregnancy and lactation period;
- diseases of the blood vessels;
- The presence of implants located under the skin or in the muscles (gold threads and other frame installations used for facelifts, metal postoperative braces, silicone). This contraindication does not apply to dental crowns, posts and braces.

Not recommended:

- To carry out the EMS procedure up to 18 years, since the procedure stimulates the collagen production, and before this age, collagen is naturally produced in a sufficient amount.
- Carry out the EMS procedure for injections (botox, dysport, fillers, gels, mesothreads).

SAFETY PRECAUTIONS

Please read the safety instructions before using the device. The safety cautions noted below are to reduce unpredicted danger and reduce loss through the safe and correct use of the product.

Attention!

- Do not use the device near flammable liquids (gasoline, thinner, spray, etc.). Doing so could result in a fire.
- Do not touch the device with wet hands. Doing so may result in electric shock and / or short circuit.
- Do not keep this apparatus near water (baths, swimming pools, etc.) Doing so may result in electric shock and / or short circuits.
- Keep the device out of the reach of children. Do not let children play with the appliance. Failure to do so may result in burns, electric shock, or other injury.
- Use only the USB cable supplied with the product.
- Remove jewelry and metal objects from your body before use.

Avoid mechanical shock

- In case of malfunction, do not try to disassemble the device yourself.
- Do not open the instrument case. In case of malfunction, contact a specialist in the service center.

It is prohibited to use this device on the following areas:


- Eyes;
- The area above your heart and thyroid gland;
- Reproductive organs;
- The area with silicone implant and gel wrinkle filling;

English

Deutsch

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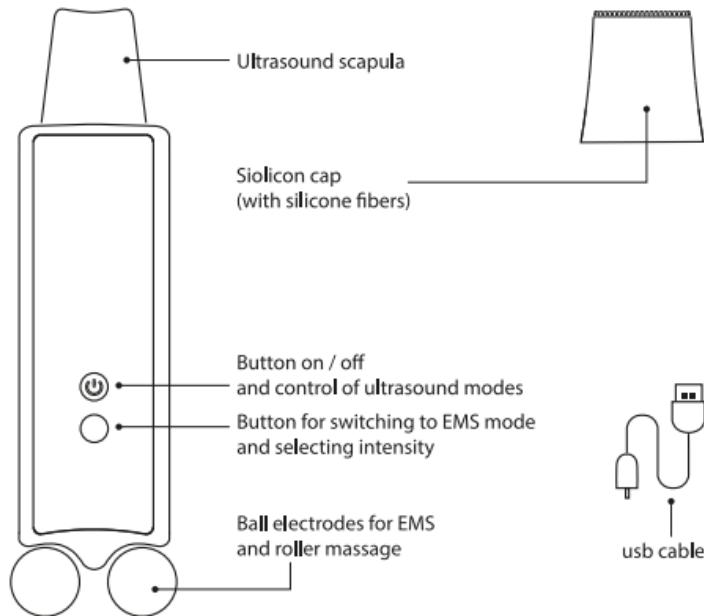DEVICE COMPONENTS

Complete set: device, silicone cap, usb cable, manual, warranty card.

DEVICE CHARGER

Make sure to charge the device fully before the first use. Full charging of the device takes no more than 2 hours (when it is completely discharged).

1. To charge the device, use the USB cable supplied with the device.
2. During charging the light indicator will be on, when the device is fully charged the light indicator will turn off.
3. The duration of charging may vary depending on operating time and conditions.
4. During long-term storage, the device will be discharged even without use. In this case, it is necessary to recharge the device at least monthly.
5. The battery charge lasts for 2.5 hours of continuous operation.

CONTROL

1. To turn the device on and off, press and hold the power button  for 2 seconds. The default setting is «Cleansing» mode.
2. To select the necessary ultrasonic mode, press the same button, the light indicator will tell you which mode is selected.
3. To switch to the EMS mode press the button  (see picture) and the miostimulation function using the electrodes balls will start working on the minimum intensity mode. There are 3 modes in total: L, M, H (low, medium, high).
4. To select the necessary intensity, press the same button, the light indicator will tell you which mode is selected.

English

Deutsch

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5. To return from EMS mode to ultrasonic modes, press the power button again 

*The device has an auto-off function after 10 minutes of operation, so there is no need to monitor the time of the procedure. If necessary, turn on the device again and continue operation.

PREPARATION AND CARRYING OUT OF THE PROCEDURE

PREPARATION

- Before starting the procedure, it is necessary to remove impurities and makeup from the skin of the face and décolleté.
- To increase the effectiveness of the ultrasound procedure, you can prepare the skin with a light scrub.
- **Ultrasound procedures** — should be performed using a water-based gel/serum/lotion designed for ultrasound procedures. The main function of using cosmetics is that ultrasound increases the activity of drugs and enhances the penetration of valuable nutrients into the deepest layers of the skin. An important parameter is that the gel / serum / lotion is an ideal conductor of ultrasonic waves.
- **EMS** — Apply a gel that acts as a conductor for current pulses on the skin or a special GEES lifting gel for all skin types.
- Perform procedures with progressive movements along the massage lines of the face, neck and décolleté.
- After the procedure, it will be helpful to apply a face mask suitable for the skin type.

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1. ULTRASONIC PEELING (mode Cleansing)

- Press the power button to select this mode.
- Clean your face in the direction of the massage lines shown in the drawing.
- Correct position of the device spatula during the ultrasonic face cleaning procedure. The spatula is positioned at an angle of 35°-45°

2. PHONOPHORESIS (mode Nutrition)

- Press the power button to select this mode.
- Turn the device with the other side of the shoulder spatula and
- perform the facial phonophoresis procedure in the direction of the
- massage lines shown in the drawing.
- Correct position of the instrument blade during the phonophoresis
- procedure. The spatula is positioned with the reverse side at an angle
- of 45°.

English

Deutsch

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3. MICROMASSAGE (mode Up)


- Press the power button to select this mode.
- Turn the device with the other side of the shoulder spatula and perform the facial micromassage procedure in the direction of the massage lines shown in the drawing.
- Correct position of the instrument blade during the procedure. The spatula is positioned with the reverse side at an angle
- of 45°.

5. EMS (Miostimulation mode)

- Press the corresponding button to select the EMS mode.
- Adjust the intensity as needed.
- Carry out miostimulation of the face and body by smoothly rolling the electrode rollers over the desired areas.

Basic schemes of massage lines:

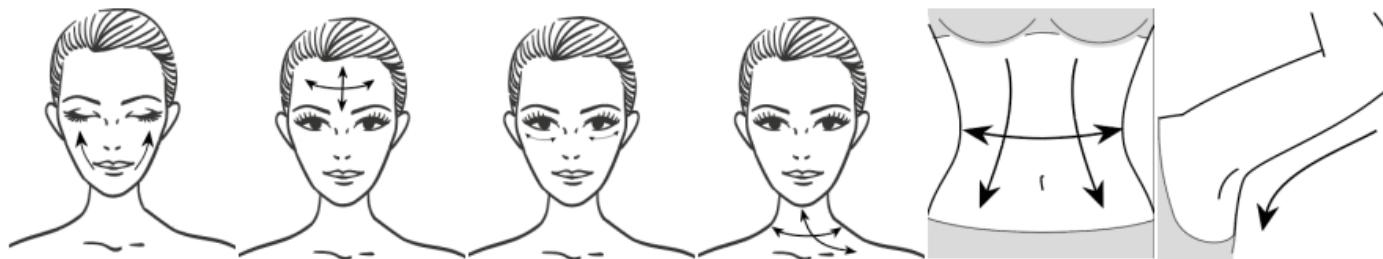Facial skin «lifting»:
from the corners of the mouth toward the outer corners of the eyes.

Forehead:
up and down,
left and right

Eyes:
around the eye area

Neck and decollete area:
up and down along the line of the neck,
avoiding the thyroid area.

Waist:
up and down
the waistcurve

Body:
vertically along
the bodyline

6. ROLLER MASSAGE

Without turning the device on, perform a roller massage of the face along the massage lines.

English

Deutsch

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RECOMMENDATIONS

Ultrasound care

- **For oily and combination skin** an intensive treatment course of 10–15 procedures, 1–2 times a week. Then an upkeep treatment — one procedure, once a month.
- **For dry and sensitive skin** an intensive treatment course of 10–12 procedures, not more than once a week. Then an upkeep treatment — one procedure, once a month.
- **For normal skin** an intensive treatment course of 10–15 procedures, 1–2 times a week. Then an upkeep treatment — one procedure, every two weeks.

EMS

Procedure time:

- Up to 30 minutes for the whole body (about 5–10 minutes for each individual site)
- The recommended frequency of use is 2–3 times a week.

Use original GEES gels they are specially designed for different skin types:

1. Cleaning gel for ultrasonic peeling. For oily and combination skin. Effectively cleans pores; gently exfoliates skin; removes oily sheen; normalizes skin acidity; moisturizes and makes skin matte and fresh. Gel has antiseptic effect and does not overdry skin.
2. Cleaning gel for ultrasonic peeling. Suitable for dry and sensitive skin. Gently cleans pores and exfoliates; starts skin refreshing; calms skin and prevents irritation or redness; moisturizes; softens; evens out skin tone.
3. For phonophoresis, microcurrents, myostimulation use original GEES lifting gel for all skin types. Smoothens deep and small wrinkles, lifts skin, prevents new wrinkles. Rejuvenates and nourishes.
4. Cleansing gel by GEES for desincrustation of all skin types. It is suitable for ultrasonic peeling and Ion cleansing.

* After the procedure, it will be helpful to apply a face mask suitable for the skin type.

English

Deutsch

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CLEANING AND CARE

Cleaning

Make sure the unit is unplugged and the OFF button is pressed before cleaning the unit. Use the soft, and dry cloth for the device cleaning.

ATTENTION

Do not use petrol, thinner or other volatile liquids for cleaning the appliance because it may damage the product serious.

Care

Make sure that the unit is turned off and the power cord is disconnected from the network. Keep it in a clean, dry place at temperatures between 5° to 45 °C and a humidity of 60%. Avoid direct sunlight on the unit.

TROUBLESHOOTING.

Problem	Cause	Elimination
The device does not turn on	1.The device is discharged2.The device is not working properly	1.Charge the device2. Call for service.
The device has shut down	1.The device is discharged2. The 10 minute mode has ended and the device automatically turned off	1.Charge the device2. Press the ON button to continue.

DISPOSAL

ATTENTION

It is not recommended to dispose of the device at the end of its service life with household waste. Disposal must be carried out through the appropriate waste collection points. If you have any questions, please contact the local municipal service responsible for waste disposal.

Do not dispose with household waste

Packaging is also a subject to recycling. Paper, wood, polyethylene and plastic must be disposed of separately in accordance with the rules of collection, accounting and disposal of packaging. The packaging is disposed of in waste paper for further recycling.

Paperboard, paper

English

Deutsch

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SPECIFICATIONS

Article	GESS-689 ems
Ultrasound frequency, kHz	25
Current strength, μ A	Pulse current: 200, 300 и 400
Pulse frequency, Hz	28
Voltage in electrodes, V	2; 2,5 and 3
Source of power	accumulator
Battery capacity, mAh	600
Rated voltage, V	3,7
Battery rated power, W	1,9
Charging current V/A	5/0.38
Charger	usb

Full charge time, h	2
Time of continuous work, h	2,5
Auto power off, min	10
Number of ultrasound modes	3
Number of EMS intensity	3
Moisture protection	IPX5
Electric shock protection class	III
Video instruction	yes, by QR
Net / gross weight, g	120/350
Packing size, mm	210 x 110 x 35

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Authorized organization by the manufacturer in Europe: GESS MARKET OÜ
Reg.code: 127106981; VAT No.: EE101758947
Room No. 219, 2nd floor, Tulika 19 (A) 10613 Tallinn, Estonia | EU
Phone: +372 51990171 info@gessmarket.com

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Product is not intended for medical use.

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Made in PRC.

GESS®

Warranty Maintenance (provided warranty provisions are complied with) includes warranty repair of the items within the warranty period or replacing them with equivalent good ones if repairing is impossible.

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gessmarket.com | gessmarket.ru | gessmarket.kz

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GESS International Trading (Shanghai) Co., LTD