

# UGI

## ELECTRIC FOOT WARMER




4 745010 327056



**EcoSapiens**



Smart textile

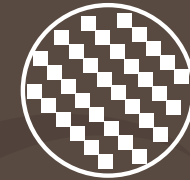

Carbon fiber



Energy Saving



Detachable remote control

# EcoSapiens

Article: S285 E  
Model: ENILTH/U/OZ/K-220v  
Size: 30x30x20 cm  
Warranty: 24 months



# UGI ELECTRIC FOOT WARMER 30x30x20 cm, EcoSapiens

## USER MANUAL

Dear customer!

Thank you for having chosen UGI infrared electric foot warmer. We hope you will enjoy it!

Please read carefully the manual before use.

The electric pad is a heating device, that is why it is important to know how to use it correctly. We give detailed information about device operation, correct usage and care.

### How Infrared Heating Pad operates (Please read this information. It is important and interesting)

The heating pad gives dry heat that promotes pain relief and muscle relaxation, positively affects joints and improves blood flow.

Infrared heating is produced by carbon fiber (not metal wire). Carbon fiber is a unique light and durable material that conducts heat. Unlike heating metal wire, carbon fiber conducts heat in the infrared range that is close to the heat temperature of the human body. Such heat is natural and safe. You will not feel dehydration and dry skin even if you spend a long time next to the heating pad. The heating pad practically does not heat up by itself. If you set the maximum heating mode the pad will stay cold. It is necessary to touch the pad in order to feel the warmth. It is important to know that Infrared heating implies warming up the body or subjects but not the air. Therefore, the device has lower power consumption in comparison with other heating pads.



### Precautions

The infrared heating pad is a household heater that requires gentle use in compliance with all recommendations in the manual. Always disconnect the device from the power-supply source if you do not plan to use it for a long time. Pay close attention when using the electric heating pad in the following conditions:

- When the pad is used by children;
- When the sensitivity of the skin is diminished; When using drugs that cause drowsiness or being under the influence of alcohol;
- In case of extreme fatigue;
- In case of diseases when the use of heat is not recommended by a consulting physician.

### PROHIBITED:

- Leave switched on the pad without supervision;
- Use the pad when it is folded or crumpled
- Stick metal or other objects in the electric heating pad;
- Disassemble insulating covers and switchers;
- Use the electric pad when the cover is torn and wire or switch are damaged
- Turn on wet electric pad;
- Rip or open the device otherwise;
- Wash and dry in a washing machine;

- Iron;
- Expose the switched-on heater to moisture;
- Put other heating devices on the pad;
- Apply to helpless people, people who are not sensitive to heat, people with disabilities and kids unsupervised.

### How to use the infrared heating pad

UGI Infrared heating pad is designed for foot warming at home. UGI is perfect for people who suffer from cold feet, have low blood circulation, or feel discomfort. Cozy electric foot pad will warm you up in cold weather or poorly heated room at home, country house and office. It is suitable for any foot size. The inner lining can be removed and washed if necessary.

Place your feet inside the heating pad. The device should be unfolded and without any wrinkles. Plug the pad into the mains.

In order to warm up your feet choose the maximum mode. The approximate warmup time is 5 minutes.

After that regulate the temperature to a comfortable one.



### Modes:

- **MIN is a minimum temperature mode.**

The heating pad operates for 60 minutes at a minimum heat level and then turns off.

- **MED is a medium temperature mode.**

The heating pad operates for 60 minutes at medium heat level and then turns off.

- **MAX is a high temperature mode.**

The heating pad operates for 60 minutes at maximum heat level and then turns off.

**When any mode is enabled, switching to another mode does not increase the total working time (60 minutes). For your safety, the device will turn off automatically after 1 hour of operation.**

If you want to continue using the device, turn it off with the button on the remote control and turn it on again. The warming up in the selected mode will repeat.

### Storage and care

- For the proper operation of the heater, store it at a temperature in the range from -15 to +25°C, in a dry room, in an individual package, without placing other items on top of it;
- Before use, inspect the product for signs of damage (if found, do not use products and contact the service center);
- To clean the electric pad, disconnect the device from the mains;
- Dry hand wash and hand wash is allowed;
- Light stains can be removed with a damp cloth or a sponge and a small amount of soap;
- Hand wash with soap is allowed only at a temperature below 30 °C and detached remote control; Further use is possible only when the blanket dries completely the natural way;
- Avoid using clothes pegs and other mechanical attachments when drying;
- It is not recommended to wash the item too often; not more than 5 times during service life.

## PRODUCT CERTIFICATE

Operating supply voltage, V \_\_\_\_\_ 220;

Europlug (for UK adapter is included)

The heating temperature in contact with human body

(in order to feel warmth, press the palm down to the device and hold it for some time):

low °C \_\_\_\_\_ 35( ± 5);

medium °C \_\_\_\_\_ 45( ± 5);

high °C \_\_\_\_\_ 55 ( ±5);

Number of temperature modes, pcs \_\_\_\_\_ 3;

Automatic shut-off timer \_\_\_\_\_ yes, 60 min;

Power consumption, Watt \_\_\_\_\_ 50;

(with rated voltage 220V)

Moisture protection level \_\_\_\_\_ IP10, protection class II;

Weight, kg \_\_\_\_\_ ±0,5;

Size, cm \_\_\_\_\_ 30\*30\*20 (a tolerance of ± 5 cm is accepted;

Wire length, m \_\_\_\_\_ 1,1 (±10 cm);

Material \_\_\_\_\_ item: 100% polyester;

Inner lining: single jersey – cotton 100%, fur - polyester 100%;

Model \_\_\_\_\_ ENILTH/U/OZ/K-220v;

Article \_\_\_\_\_ S285 E;

Color \_\_\_\_\_ brown

In a set: electric heating pad, remote control, user manual, package.

**The product is not intended for medical use.**

# EcoSapiens

Produced for EcoSapiens.

Manufacturer: GESS MARKET OÜ,

Room No. 219, 2nd floor, Tulika 19 (A) 10613 Tallinn, Estonia | EU

Place of manufacture: GME Europe Ltd., Ilzenes 2, Riga, LV-1005, Latvia

Exclusive distributor in EU is GESS MARKET OU [www.gessmarket.com](http://www.gessmarket.com)

Please send your inquiries and offers: [info@gessmarket.com](mailto:info@gessmarket.com)

The service life is 5 years.

The warranty period is 24 months.

Date of manufacture: see on the sewn-in label.

**Your opinion is important to us! You can always leave a review about the device in the store you purchased the item or send it by e-mail.**

Customer service:

E-mail: [info@gessmarket.com](mailto:info@gessmarket.com)

Tel. number: +372 51990171

Official website: [gessmarket.com](http://gessmarket.com)

