


Please keep this manual for future reference.



SHIATSU FOOT MASSAGER

with Heat & RockerTech™

User Manual

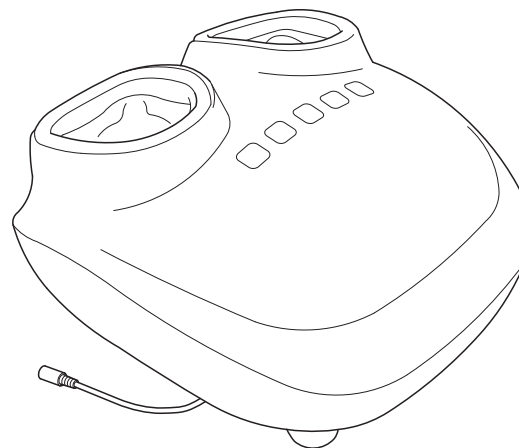

Infinity Massage Chairs
72 Stard Road
Seabrook, NH 03874
603-910-5000

www.infinitymassagechairs.com

Thank you for choosing the Infinity Shiatsu Foot Massager with Heat & RockerTech™. Please read the User Manual carefully before using the massager. Keep this booklet in a safe place for future reference.


IMPORTANT SAFETY WARNINGS:

- Please use this product in accordance with the guidelines of this manual.
- Do not allow any foreign objects to enter the clearance between the massage head and casing.
- To prevent short-circuiting, avoid getting the massager wet. Avoid contact with liquids/ water.
- If the main cable or plug is damaged, do not turn on the massager or attempt to repair the plug. Please contact Infinity for troubleshooting instructions.
- Disconnect the product from the power supply after use and before cleaning to avoid electrical injury.
- Children and disabled individuals should be supervised while using this product.
- Do not attach anything else to the massager during use.
- Do not place heavy objects on massager when not in use.
- Do not use this product within one hour after eating.
- Do not use this product immediately after there has been a sharp change in the environmental temperature.
- Do not use this product in a dusty or corrosive environment.
- Do not cover the massager with a blanket or pillow during use.
- Do not fall asleep while using this product.
- Do not use product on pets.


The following people should consult their doctor prior to using this massager:

- People with contagious skin conditions.
- People with osteoporosis or other bone and muscle conditions.
- People with heart disease and/or people with electrical medical instruments such as pacemakers.
- People with high fevers.
- Women who are pregnant.
- People with injuries or recovering from injuries.
- The massager has a heated surface. Those sensitive to heat should use this massager with caution.

MASSAGER COMPONENTS

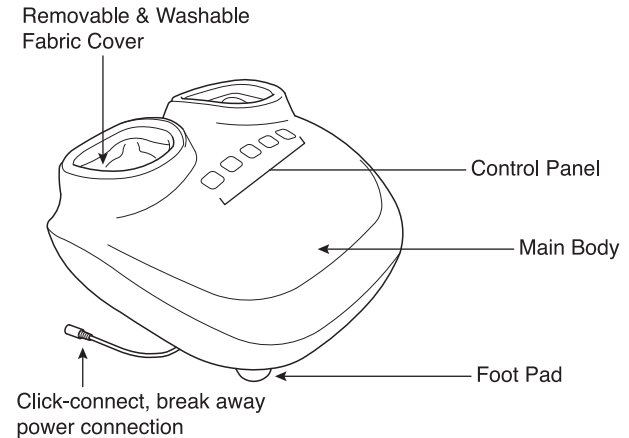

MASSAGER OPERATION

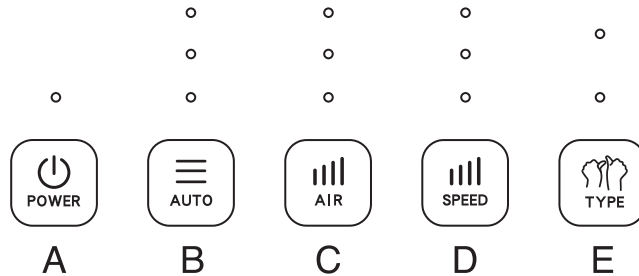

A. Power button:

- Turns massager on/off; Toggles heat function
- Press once: Massager turns on. Indicator light is red, heat is on.
- Press again: Indicator light is green, heat is off.
- Press a third time to turn massager off.

B. Auto button:

- Toggles Auto programs
- Cycle through the three available massage programs with this button.
 - 1) **Energy (Top indicator light):** Massage rollers, airbags, and RockerTech™ work quickly in concert for an energizing foot massage.
 - 2) **Renew (Middle indicator light):** Massage rollers, airbags, and RockerTech work the feet more slowly with longer holds to boost recovery.
 - 3) **Relax (Bottom indicator light):** Massage rollers are inactive while airbags and RockerTech gently soothe sore, tired feet.
- Press Auto a fourth time to enter Manual mode, in which you can adjust Air, Speed, and RockerTech settings.

C. Air button:

- Controls air pressure intensity
- Adjust the intensity level with this button. Air pressure can be adjusted during both Auto and Manual programs.
 - Mode A (Top indicator light):** High intensity
 - Mode B (Middle indicator light):** Medium intensity
 - Mode C (Bottom indicator light):** Low intensity
- Once a mode is selected, press the Air button again to change the level of intensity.

D. Speed button:

- Controls speed of RockerTech™ and massage rollers
- Adjust the speed of all functions (except the airbags) with this button. Speed can only be adjusted in Manual mode.
 - Mode A (Top indicator light):** Fast
 - Mode B (Middle indicator light):** Medium
 - Mode C (Bottom indicator light):** Slow
- Once a mode is selected, press the Speed button again to adjust.

E. Type button:

- Toggles between RockerTech only and RockerTech with massage rollers
- Select either RockerTech only or both RockerTech and massage rollers.
 - Mode A (Top indicator light):** RockerTech tapping functions alone.
 - Mode B (Middle indicator light):** RockerTech and massage rollers function together.
- Once a mode is selected, press the Type button again to toggle them.

Toggle Manual Mode Functions

If you would like to turn off the Air, Speed, or Type functions when already in a Manual massage, power the entire massager off and back on again to then select your preferred settings.

Massage Session Time & Auto Shut-off

- The massager will automatically end the massage session after 15 minutes. If another session is desired, please make a new selection.
- Multiple back-to-back uses could result in auto shut-off to prevent overheating. If this occurs, please wait at least 60 minutes to allow the motors to cool down before attempting to start a new session.

MAINTENANCE

A. Storage

- When not using the massager, store it back in the box and place it in a dry, well-ventilated, cool area.
- Do not twist wires.
- Avoid sharp objects to prevent scratching surface of the product.
- Keep out of direct sunlight and high temperatures.

B. Cleaning

- Unplug power cord before cleaning.
- Use a dry cloth only to clean the surface of the product. Avoid direct contact with water.
- Do not use corrosive cleaners such as benzene or bleach products.
- The fabric pouch can be removed for cleaning.

C. Maintenance

- Keep the product clean and away from moisture.
- Cover up the product when not in use to avoid dust.



Warning: This product has an internal temperature protection device. When the safe temperature or safe heat use duration is exceeded, the heat will be turned off automatically. Please unplug the power to allow machine to cool after use.

Infinity Technical Support

If you have any questions about maintenance or are experience technical difficulties with your massager, please contact Infinity's technical support team at 603-910-5000 for assistance.

TROUBLESHOOTING

Problem	Possible Cause	Possible Solution
Massager not responding	No power	Check power supply
	Function keys are not responding	Restart massager using On/Off button or unplug the device
Massager suddenly stops operating	No power	Check power supply
	15 minute massage session has ended	Press "Power" button for another session
	Overheated after continuous use for 45 minutes	Unplug massager and allow to cool down for 60 minutes before using it again.
	Auto shut-off due to excessive force	Unplug the massager. Lighten the pressure before turning on again.

SPECIFICATIONS

Model	SHIATSU FOOT MASSAGER with Heat & RockerTech™
Dimensions	16.14" x 15.35" x 10.28" (L x W x H)
Weight	16.6 lbs
Rated Voltage	24V \approx 2.5A
Rated Power	45W
Rated Duration	15 Minutes