


Relax Energy Table Lamp

Instruction Manual

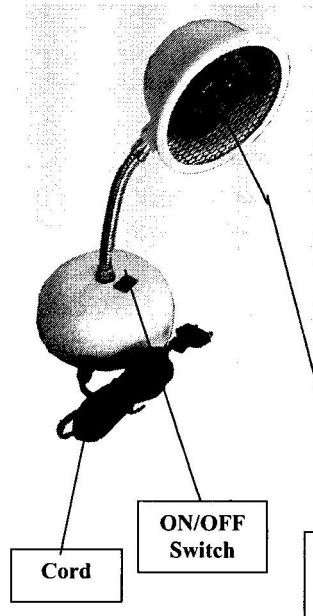

**Patented Ceramic Semiconductors
100% Safe - Unique Automatic
Temperature Control System**

Portable

Hand Free

- ⊗ Seniors, Children, and/or those who require medical care must use Relax products only under the supervision of a competent adult.
- ⊗ Seniors, Children, and/or those who are under medical care must consult with their physicians first before using the Relax Energy products.

Guidelines on usage



- The recommended use time is 15 minutes. User should adjust the use time based on individual condition.
- In the process of operation, do not touch the metal screen, as it is hot.
- Keep the eyes firmly shut when radiate on the face.
- The temperature of the metal screen shaft is between 165°F to 185°F. The distance between the radiator and the user depends on the individual reaction of acceptable heat. The recommended distance is 10 inches from your skin to the Lamp face.
- The product has a protective device for overheating. The electric current will be automatically cut-off once overheated. After cooling, the electric current will resume automatically.
- As with any other electric appliance, the same rules apply with the FIR Table Lamp.
- Do not use the FIR Table Lamp near wet areas in case of electrocution.
- There should be an accompanying adult on hand when an elderly person/chronic illness sufferer desires to use the product.
- At any point of time the user is feeling uncomfortable, stop using immediately and consult the doctor.
- User with special medical problem should consult the doctor before using.
- Turn on the switch to activate the device. After using, turn off the switch immediately so as to save energy and prolong the life span of the product.



Caution:

- Do not attempt to change or modify any part of the product.
- It is forbidden and dangerous to alter the specification or modify this appliance in any way.
- In case of any abnormal situation, please switch off the product; remove the plug and contact local dealer for RMA service.
- Follow this instruction book strictly and do not use any parts outside the specifications.

Contraindications

- **Do not** use the Relax Energy Table Lamp under the following conditions:
 - * Insensitivity to Heat
 - * While under the influence of alcohol
 - * Under doctor's instructions not to use heat lamp
 - * Suffering from a deep cut or any fresh wounds
 - * Infants are never to be use the Relax Energy Table Lamp
 - * Seniors, Children, and/or those who require medical care must use the it only under the supervision of a competent adult.
- **Those who require medical care must consult with their doctor before use.**
- **Diabetics** and those with poor circulation (such as cold feet and/or icy hands), **must never exceed 5 minutes per session.**
- Patients with cardiovascular conditions and chronic conditions should only use it under supervision of medical professionals
- The occurrence of red spots and itchiness is normal in that it is only a transitional phase in the process. When toxins are unable to be carried away by the body's blood stream, they appear as red spots on the skin. If some areas more affected than others, it is often due to an old injury. Different individuals experiences will differ and, normally in 2 weeks, the reactions will clear. **Decrease the time per session short if the red spots occur.**
- If your skin itches, avoid scratching and rest.
- Not everyone will experience dizziness. **If you feel dizzy, please stop using it immediately**

Guidelines on maintenance

- **Beware of the following safety precautions on the electric cord to avoid accident**
 - **When the appliance is in operation, do not twist, pull or blind the electric cord. Do not put any heavy object on the cord.**
 - **Use only outlet of AC 120V!**
 - **Once the cord torn or faulty, do not use the appliance and send to the company for repair.**
- **For cleaning and maintenance, use soft and damp cloth or mild detergent to wipe. Do not use volatile cleaning liquid.**
- **After a period of usage, use a hair dryer (use only cool wind) to spout out the dirt in the FIR radiators.**

Warranty

Product	FIR Table Lamp
Serial No.	
Name	
Address	
Telephone	
Fax	
E-mail	
Expiry	One year from date of purchase
Purchase Date	
Distributor stamp:	

1. We sincerely support the best service for users. Please fill in warranty card completely and retained with the appliance for future reference.
2. Should the warranty card be misplaced / lost, manufacture is under no obligation to reissue Warranty card.
3. Within the warranty period, we will not provide free repair service if under the following circumstances:
 - a) Damaged caused by improper use (not according to this operating instruction)
 - b) Damaged caused by the user' s repairs, changes or dismantle the product without consent from the management.
 - c) Damaged caused by accident, fire, flood or any natural disaster
 - d) Warranty is void if the trademark or purchased date is replaced or missing.