
tria[®]

SMOOTHBEAUTY™ LASER

Quick Start Guide

GET TO KNOW YOUR SMOOTHBEAUTY™ LASER

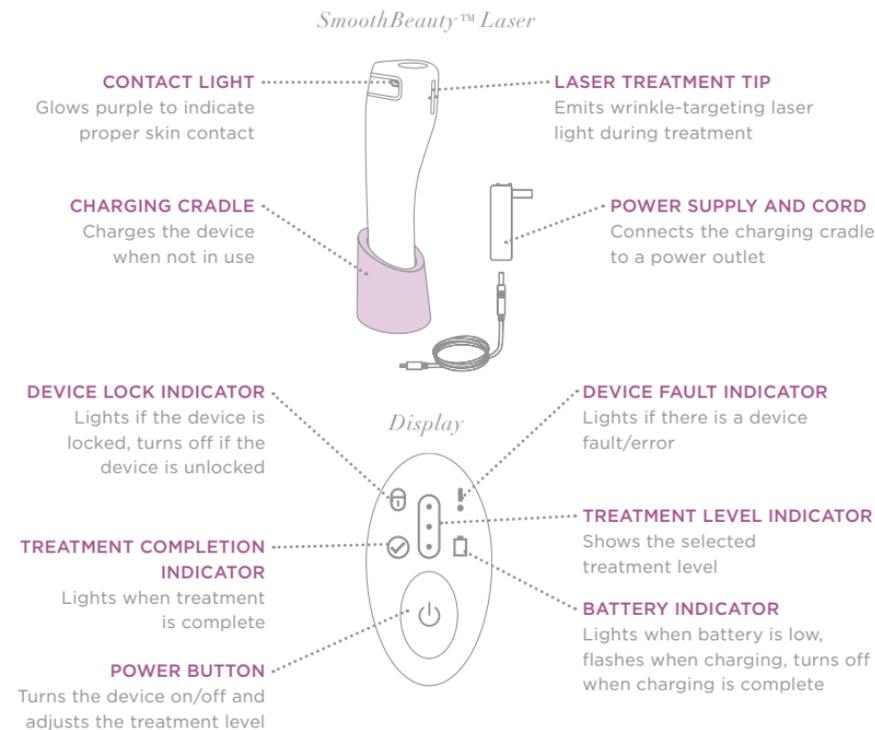TAKE OUR ADVICE

PRACTICE IN A MIRROR

Perfect your technique by using your device in front of a mirror. The purple light will glow when you are treating correctly.

INCREASE YOUR TREATMENT LEVEL FOR BEST RESULTS

Some women feel a slight stinging or warming sensation during treatment. This is normal and will diminish. To help your skin adjust, start on Training Level 1, Increase to Treatment Level 2, then Treatment Level 3 for best results.

MAKE IT A NIGHTLY BEAUTY RITUAL

For best results, treat skin in the evening to allow any redness to diminish.

WEAR SUNSCREEN

Make it a routine to apply a broad-spectrum sunscreen in the morning and reapply during the day.

 Tria Beauty, a company of Channel Investments, LLC
14796 Wicks Blvd.
San Leandro, CA 94577 USA

TRIA BEAUTY CUSTOMER CARE
877-321-8742
customer@triabeauty.com
triabeauty.com

tria[®]

Refer to the *Instructions for Use* booklet for complete information about the Tria SmoothBeauty™ Laser, including warnings and indications for use.

©2019 Tria Beauty, 14796 Wicks Blvd., San Leandro, CA 94577 USA
All rights reserved. SmoothBeauty™ is a trademark of Tria Beauty. Tria, Tria Beauty and the Tria logo are registered trademarks of Tria Beauty. Protected by patent(s). See triabeauty.com/legal.

BEFORE FIRST TREATMENT

For detailed directions and safety information, read the complete *Instructions for Use*.

CHARGE: Charge the laser until the battery indicator turns off, indicating the device is fully charged (up to 2.5 hrs). The device may arrive fully charged and not require additional recharging before first use.

HOW TO TREAT

Treat nightly by first cleansing and drying your skin. Treat with your laser, and then apply a post-treatment serum. For best results, use the Tria Priming Cleanser and Finishing Serum.


1 TURN ON DEVICE

Remove the laser from its charging cradle and press the power button. Make sure the device is unlocked.


2 SET THE TREATMENT LEVEL

Press the power button to adjust levels (Training Level 1, Treatment Level 2, or Treatment Level 3). Your SmoothBeauty™ Laser defaults to the level you last used.

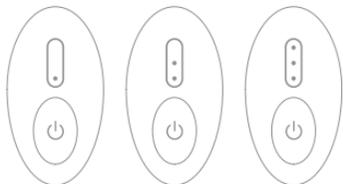3 START YOUR TREATMENT

- Divide the face into 4 zones. Start in Zone 1 by placing the laser tip flat against your skin and gliding it slowly from side-to-side.
- Cover the entire zone and maintain proper skin contact. The purple light will glow if you are treating correctly.
- Treat until you hear a beep and then progress to the next zone. The device will automatically shut off after you treat the fourth zone.

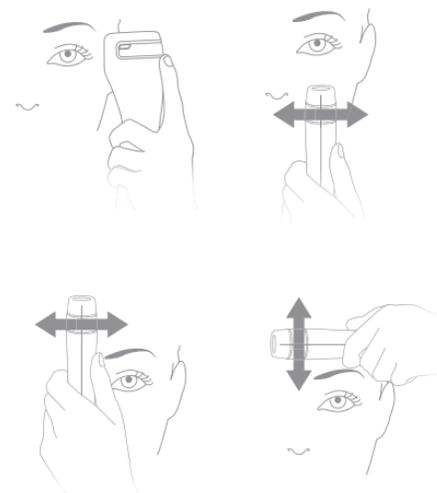
Increase levels for best results
As your skin adjusts, increase the treatment level for best results.

⚠ DO NOT use your SmoothBeauty™ Laser on your eye or eyelid or within the eye socket. Treat around the eyes only where you can feel bone underneath the skin.

TREAT THE FULL FACE BY ZONE


Treatment Zones

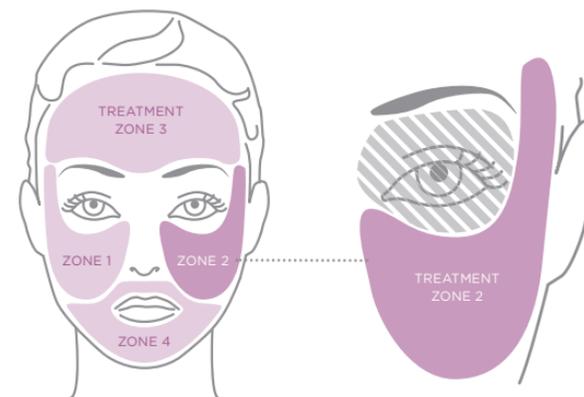Approximate Time Per Setting

Training Level 1	30 seconds per zone
Treatment Level 2	1 minute per zone
Treatment Level 3	2½ minutes per zone

Treatment Levels

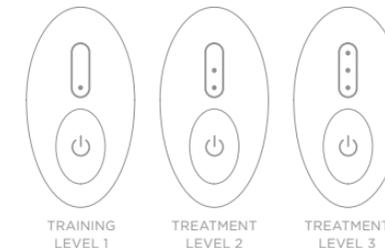GIVE YOUR WRINKLE-FIGHTING REGIMEN A BOOST

Enhance your SmoothBeauty™ Laser Experience with our Cleanser and Serum.


PRIMING CLEANSER

Prepare skin for laser treatment with this luxurious cleanser. Formulated with Willow Bark Extract to promote exfoliation and Birch Tree Sugars to help preserve moisture, our cleanser will leave your skin feeling refreshed and hydrated.

FINISHING SERUM

Our nourishing serum is formulated to be the perfect finish to your laser treatment. Sea Fennel helps to soothe and hydrate laser-treated skin, while Retinol and Biopeptides promote elasticity and firmness and help to improve the appearance of fine lines and wrinkles.

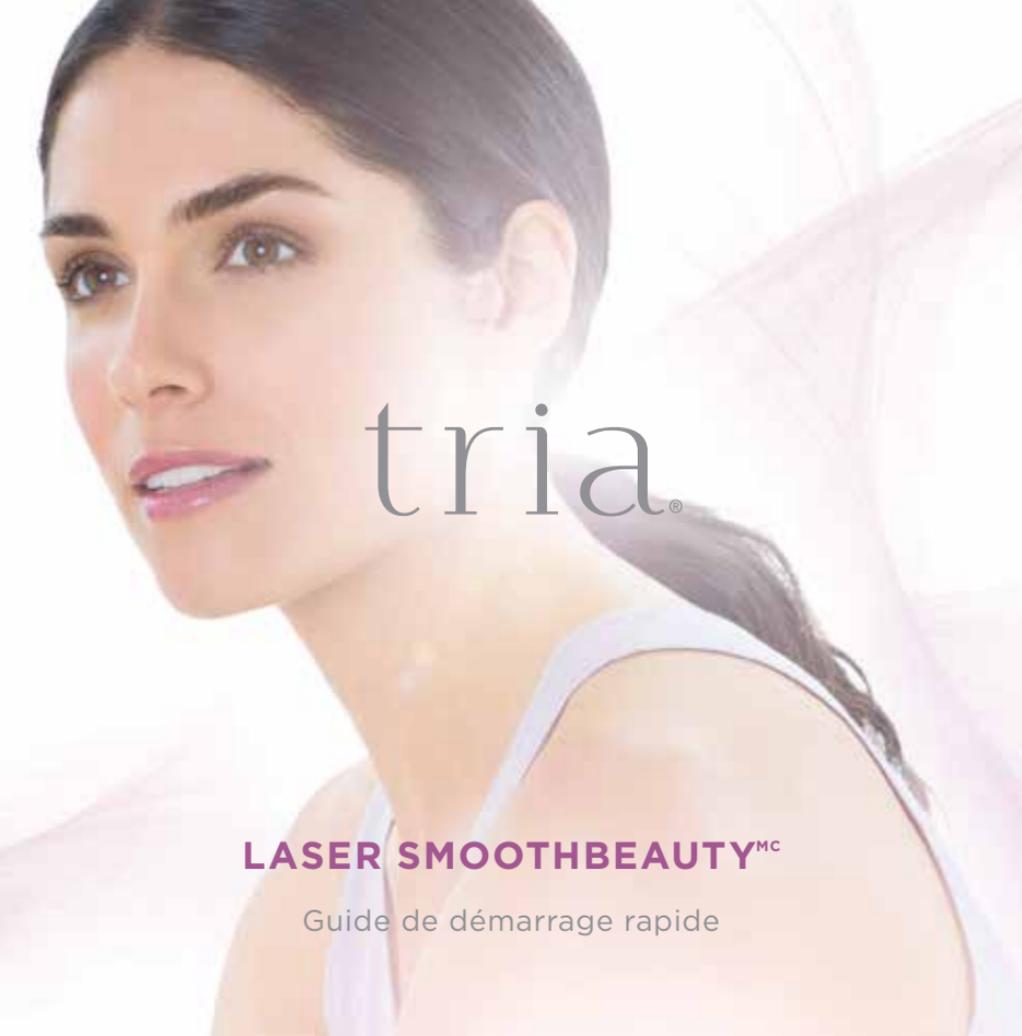For more information on Tria's wrinkle-fighting skincare products visit our website.

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