

CryoPad




1. Treatment Area
2. Show Time
3. Time Setting
4. Switch
5. Supply Hub
6. Inlet Air Area

Binging the machine on the area that you want to do the treatment, pug in power,you will hear a sound “di”, when showing ---- ---- means normal power supply

Power Button: Press once means machine on (second hand flashing and blue light on, Default time is 8 minutes), press second time, machine off(Time showing”----” and blue light off)

Time setting button: Only when machine on then you can set the time you need, press once means add 1 minute. Long press button can add time to 30

minutes fleetly but when the time reach that you want please stop pressing,otherwise the treatment time will start from 1 minute to 30 minute again and again. Treatment time is 1-30minutes.

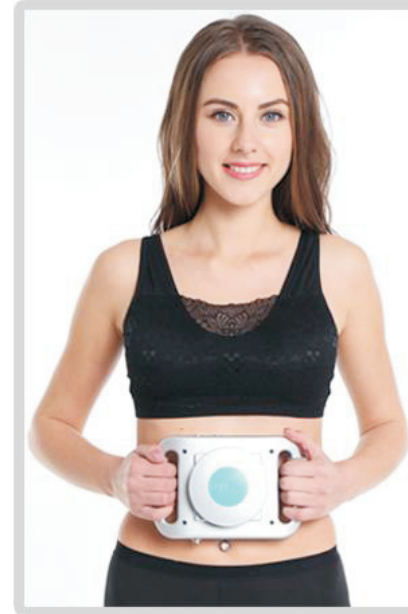

Attention:

1: Contraindication:

Raynaud's disease(Raynaud's disease is a rare disorder of the blood vessels,usually in the fingers and toes,where digits show blue discoloration secondary to cold exposure)

Pregnancy

Implanted medical device

Recent surgery in the treatment area(within 1 year)
Umbilical or inguinal hernia
Cold sensitivities
Patients with Hepatitis C

2. For Avoiding frostbite,please do not do the treatment exceed 8 minutes at first time,if patient Ok for that then lengthen treatment time suitable

3. To ensure the results, please do not block the Inlet Air Area of both side

Parameter:

Treatment area: 98x98mm

Product dimensions: 20cm(L)x13cm(W)x4.4cm(H)

Gross weight: 1.9kg

Output Voltage: 12V

Output Electric Current: 10A

Automatic working time: 20 minutes, but the treated time can be set in 20,40 and 60 minutes manually by own

Temperature: Min -2 degree(in door), the machine will as per the treatment time to control the temperature automatic

Competitive Advantage

This mini cryolipolysis machine is very light, with the good price, also can save the freight, you can do the treatment when you are listening music or enjoy you own time after the hard working. You don't need to take any more time on it. It is reasonable human design for the beauty machine.What's more it is effective and safety for the weight loss,believed it is your top option, you deserve it.

What is cryolipolysis?

Cryolipolysis is the non-invasive cooling of adipose tissue to reduce lipolysis(breaking down fat cells) without damage to other tissues. The result over two to four months is a noticeable, natural-looking reduction in fat.

How does cryolipolysis work?

Cryolipolysis is based on the principle that fat cells are more vulnerable to energy extraction (cooling) than surrounding tissues. A non-invasive applicator is attached to the targeted procedure area to extract energy from the underlying fat tissue, while protecting the skin,nerves,muscles and other tissue. The cooled fat cells undergo apoptosis(controlled cell death) and are gradually eliminated, reducing the thickness of the fat layer.

What are the results?

Clinical studies demonstrate that cryolipolysis procedure provides

noticeable,measureable fat reduction in properly selected customers over the course of two to four months after the procedure. Additional procedures may be administered two to four months after the initial procedure to achieve greater fat reduction.

How is the cryolipolysis procedure different from other procedures?

The controlled cooling of the cryolipolysis procedure targets and eliminates only fat cells. Other treatment modalities,such as lasers, radio frequency and focused ultrasound, affect fat cells and may affect other adjacent tissue in a way that is not comparable to the cryolipolysis method of cryolipolysis equipment

How long do the results last?

Customers experiencing fat layer reduction show persistent results at least six months after the cryolipolysis treatment. Longer-term studies have not yet been completed. It's known that the exposure to cooling causes the fat cells to begin a process of natural removal, which gradually reduces the thickness of the fat layer. The fat cells in the treated area are gently eliminated through the body's normal metabolism process,like expected to last about as long as fat cells removed by invasive procedures such as liposuction.

treat the area with just 1 application a month... You can treat multiple areas, but make sure that you only treat each area ONCE, then WAIT 30 days before treating it again. No need to use every day. For continuous results, use once every month!Results will be seen within 4 to 8 weeks! Simply keep using each month for best and continuous results.