Ab Tone & Shape

Tones, Tightens and Firms

Bodi-Tek®

User Guide

Warnings

Please read all instructions thoroughly before use and retain for future reference.

This appliance can be used by children aged 14 years and above and persons with reduced physical, sensory or mental capabilities or lack of experience and knowledge if they have been given supervision or instruction concerning the use of the appliance in a safe way and understand the hazards involved. Children shall not play with the appliance. Cleaning and maintenance shall not be made by children without supervision.

Do NOT use:

- If you have an electronic implant (e.g. heart pacemaker).
- In places with high humidity such as the bathroom.
- If you suffer from high/low blood pressure or from any other heart problems.
- If you are pregnant.
- Over the abdominal area during your menstrual cycle.

- On loose skinned parts of your body or any areas showing signs of: varicose veins, skin diseases (for example eczema, psoriasis, dermatitis), cuts and abrasions, irritation, warts, moles, skin tags, skin rash, hypersensitivity or sunburn, phlebitis, recent scar tissue, other skin conditions. We recommend you consult your physician before treatment if you suffer from any of the above conditions or on areas that have had cosmetic surgery.
- Over metal objects such as belt buckles or jewellery.

Get your doctor's permission before using if:


- You have a history of heart disease or cardiac condition.
- You have an orthopedic implant.
- You suffer from epilepsy, multiple sclerosis, phlebitis in its active phase (inflammation of a vein), inflamed tissue through disease or injury or varicose veins in its later stage.
- You have suffered acute trauma, fracture or are recovering from recent surgery.

If you are in any doubt about using the Ab Tone & Shape for any reason, please consult your doctor before use.

Precautions

- Keep out of reach of children and pets.
- Do NOT use if you are under the influence of alcohol.
- Do NOT use whilst driving or operating machinery.
- Do NOT touch the conductive side of the belt during operation.
- Do NOT drop, place heavy objects on the unit or immerse in water.
- Should you experience any discomfort, firstly reduce the intensity, if discomfort persists discontinue use and contact our Customer Care.
- In very rare cases, skin irritation occurs around the electrode area. If this happens discontinue use.
- Bodi-Tek will not accept responsibility if the guidelines and instructions supplied with the unit are not followed.
- This device must only be used for its intended purpose as described in this user guide.


Know your Unit

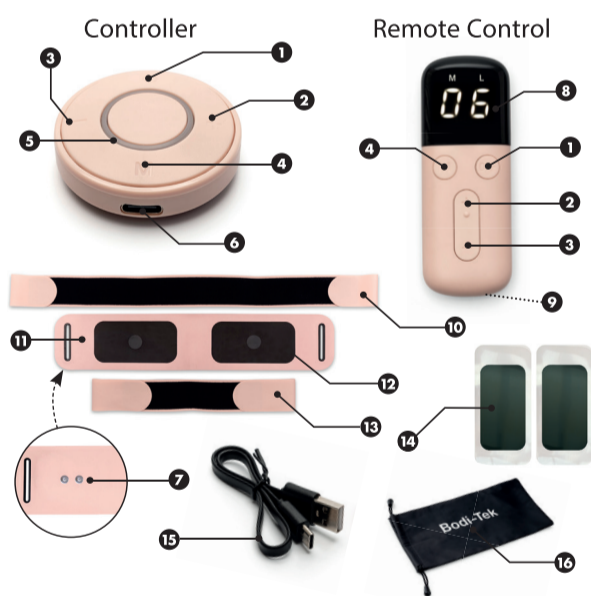- | | |
|---------------------------|-------------------------|
| 1. On/Off | 9. Battery Compartment |
| 2. (+) Intensity Increase | 10. Belt Extender Large |
| 3. (-) Intensity Decrease | 11. Belt |
| 4. 'M' Programme Button | 12. Conductive Pads |
| 5. LED Light | 13. Belt Extender Small |
| 6. USB Charging Port | 14. Gel Pads x 2 |
| 7. Magnetic Connectors | 15. USB Charging Cable |
| 8. LCD Display | 16. Protective Pouch |


How to Charge

1. Use the USB Charging Cable (included) to connect to a USB port. It takes 90 minutes to fully charge.
2. The central LED light will show a consistent red while charging.
3. When fully charged the LED light will turn green.

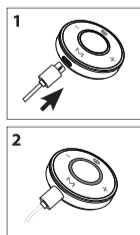Remote Control


Open the Battery Compartment on the back of the Remote Control. Insert CR2032 3V battery and close Battery Compartment.

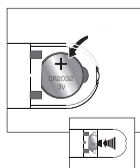To use the Remote Control, first switch the Controller on, then turn on the Remote Control. To turn off the Remote Control, press the On/Off button.

Getting Started


1. Remove any excess oil or moisture from the skin.
2. Moisten the skin. For best results use Bodi-Tek Conductive Gel (<https://www.bodi-tek.com/products/conductive-gel-85ml>). Using conductive gel gives stronger contractions while reducing the sensitivity of the skin.
3. Ensure the Controller is switched off and attach to the Belt using the Magnetic Connectors.
4. Attach one or both of the belt extenders, passing them through the end slots in the belt and secure using the hook & loop sections.
5. Remove the Gel Pads from the protective seal and attach to the Conductive Pads on the inside of the Belt. Avoid touching the Gel Pads as this will reduce tackiness.
6. Place the Belt around your waist with the Conductive Pads placed directly over the targeted muscles.
7. Press the On/Off button on the Controller. The dark blue LED indicates that programme 1 is activated. To increase or reduce the intensity press the "+" or "-" button on the Remote Control/Controller. The Controller will automatically switch off at the end of the 25 minute workout programme.
8. To change between workout programmes press the 'M' button on the Remote Control/Controller. When switching between programmes the intensity level resets to level 1.
9. Once finished, replace the plastic protective sheet and store in a cool dry place. You do not need to remove the Gel Pads from the Belt.

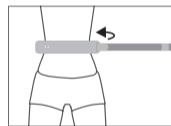Note: The Controller LED will flash for 30 seconds and then turn off if the pads are not in contact with the body.

Intensity Adjustment

Intensity levels range from 1 to 15. Press the "+" button to increase and the "-" button to decrease.

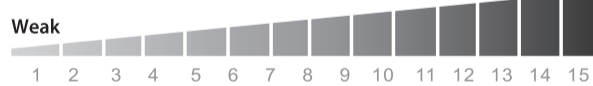
EMS Programmes

Each programme is 25 minutes. At the end of the programme the Controller will beep twice and switch off.

Programme	1	2	3	4	5	6	7
Name	Auto	Beginner	Strength 1	Strength 2	Strength 3	Tone	Relax
LED Indicator	Dark Blue	White	Green	Yellow	Pink	Light Blue	Red


* Auto Training Programme

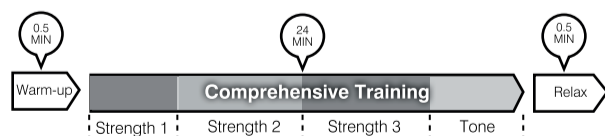Cleaning and Maintenance

Gel Pads: You may extend the life of the Gel Pads by briefly wetting the adhesive side in cold running water and allowing them to air dry (do not heat). Over-saturation with water will reduce the adhesive properties of the pads. Alternatively apply to slightly damp skin and allow 1 minute for the adhesive to anchor to the skin. Do not continue using a pad once the adhesive no longer sticks or the pad has become damaged. Important: Ensure the control unit is removed from the Belt when wetting the adhesive side.

Belt & Controller: Wipe clean with a damp cloth, allow to air dry. Never immerse the belt or controller in water or machine wash.

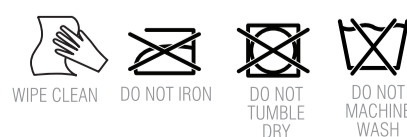General Care

- Do not place close to excessive heat.
- Do NOT disassemble or repair. Repair and service must only be undertaken by Bodi-Tek or an authorised distributor.
- There are no user serviceable parts inside the unit. If you have any problems using the product contact our Customer Care.

After Use

The skin may appear slightly red and warm after use. This should subside after 30 minutes.

Customer Care

If you have any questions, please contact us at: customerservice@bodi-tek.com

Specification

Model: BT-ABST
 Input: 5V 500mA 2.5W
 Battery: Li-Ion 3.7V 150mAh
 Dimension of main device: 44mm x 44mm
 Material: PU
 Outer: PU 40% Lycra 30% Nylon 22%
 Inner: Conductive Film PU 8%

Environment

Operating Temperature: 10–30 degrees Celsius
 Humidity: 20%–65% RH
 Atmospheric Pressure: 700–1060 hPa

Storage & Transport Temperature: 5–40 degrees Celsius
 Humidity: 10%–90% RH

Read these instructions carefully before use and then retain for future use.

This means that the unit complies with current product legislation.

This appliance is for indoor use only.

Guarantee

To register your two year guarantee online, visit: www.bodi-tek.com/warranty

Hereby, The Dezac Group declares that the radio equipment type "Non-FHSS adaptive equipment with only one antenna" is in compliance with Directive 2014/53/EU and Radio Equipment Regulations 2017.

The full text of the EU/UK declaration of conformity is available at:

<https://www.bodi-tek.com/products/ab-tone-shape>

UK Office: PO Box 17, Cheltenham Spa, England, GL53 7ET

EU Office: The Dezac Group (Europe) BV, Azura Saturnusstraat 46-62, 2132 HB Hoofddorp, Netherlands

Disposal

The Waste of Electrical and Electronic Equipment (WEEE) Directive (2012/19/EU) has been put in place to recycle products using best available recovery and recycling techniques to minimise the impact on the environment, saving energy and resources and avoiding hazardous materials from going to landfill. The product is classed as Electrical or Electronic equipment so please ensure that at the end of the life of this product it is disposed of in the correct manner in accordance with local authority requirements. It must NOT be disposed of with household waste.

Recycle Responsibly:

1. Check your local municipal guidelines
2. Separate the components and dispose of them correctly
3. Empty the packaging of its contents before recycling
4. Reduce the volume of the packaging before recycling

UK RP: The Dezac Group Ltd, PO Box 17, Cheltenham Spa, England GL53 7ET
 EU RP: The Dezac Group (Europe) B.V. Azura, Saturnusstraat 46-62, 2132 HB Hoofddorp, Netherlands
 ©The Dezac Group Ltd 2022, PO Box 17, Cheltenham Spa, England GL53 7ET.
www.bodi-tek.com

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RU RP: The Dezac Group Ltd, PO Box 17, Cheltenham Spa, Inglaterra GL53 7ET
EU RP: The Dezac Group (Europe) B.V. Azura, Saturnusstraat 46-62, 2132 HB Hoofddorp, Países Bajos
©The Dezac Group Ltd 2022, PO Box 17, Cheltenham Spa, Inglaterra GL53 7ET.