





INSTRUCTION MANUAL  
INSTRUCTION MANUAL  
INSTRUCTION MANUAL


MK9199 Plus

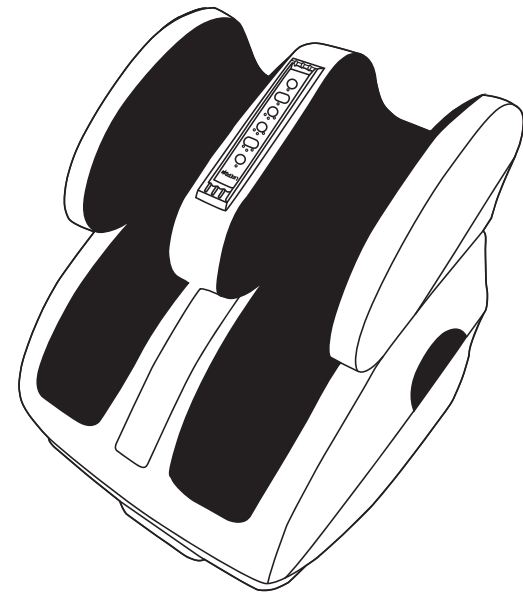

*Charming Health*

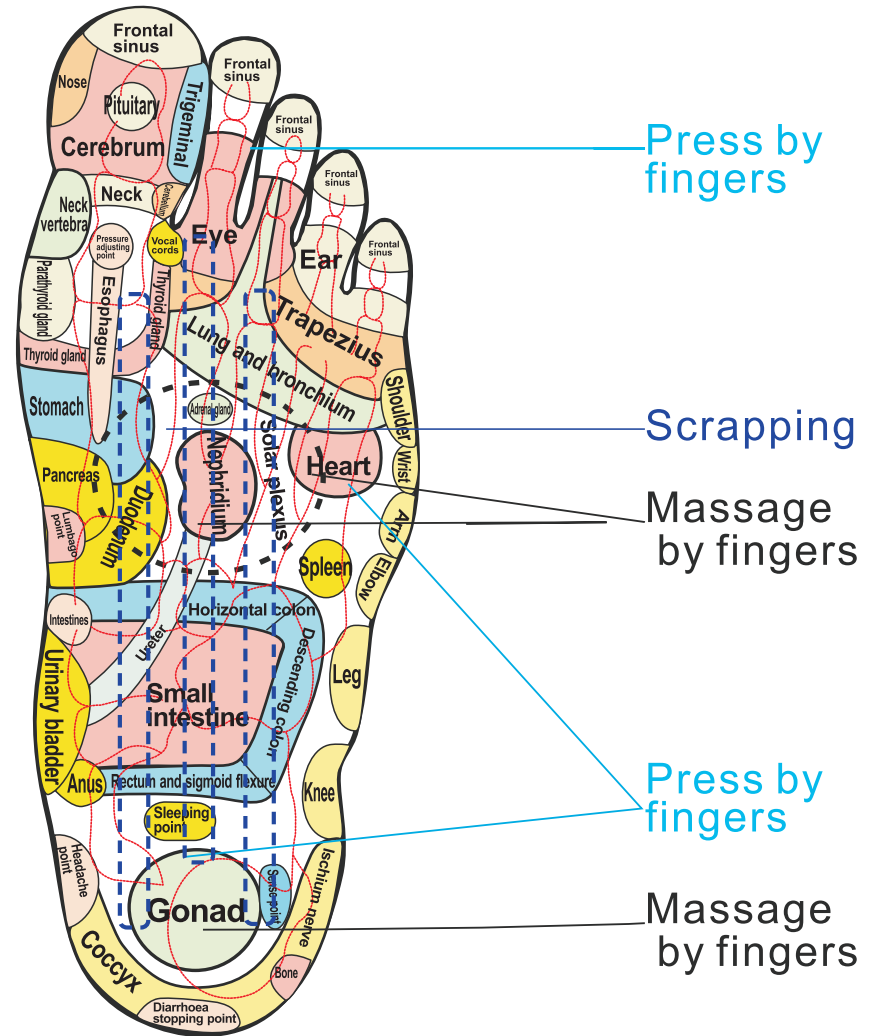

Reflecting areas on foot

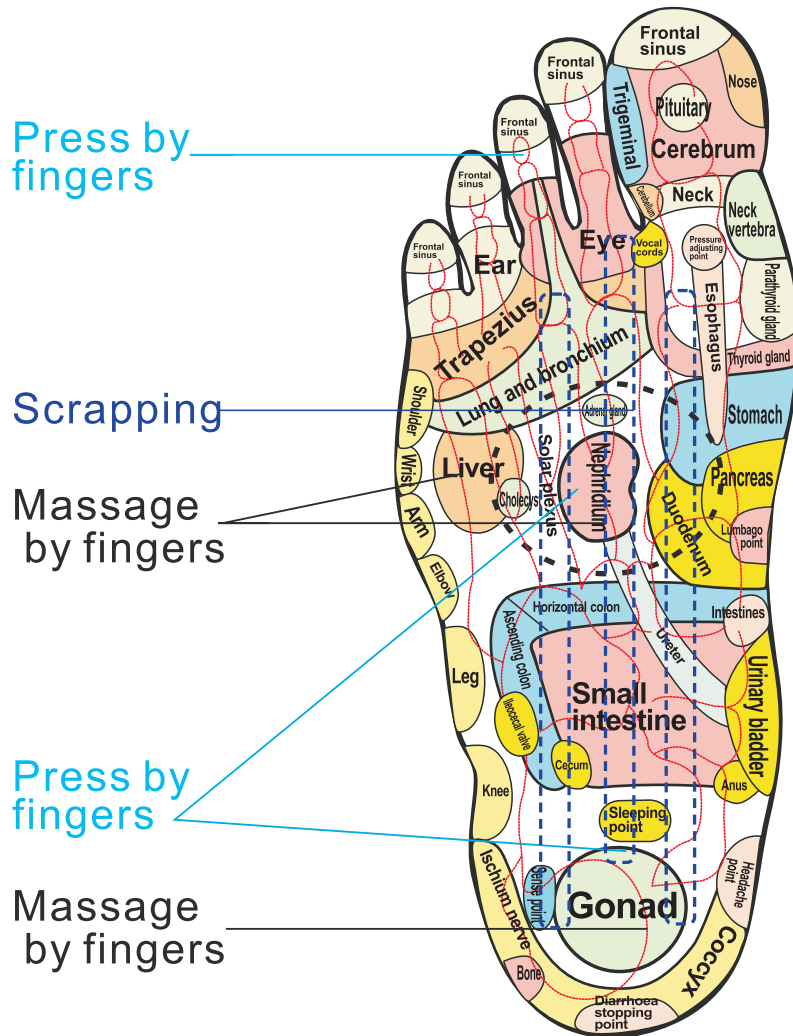

# TABLE OF CONTENTS

## Before Use

- Safety Precautions ..... 1
- Product Features ..... 2
- Preparation before Use ..... 3

## Operating Instructions

- Operating instructions ..... 4
- Product composition (Front side view) ..... 5
- Product composition (Back side view) ..... 5
- Product composition (Control Panel)
- Illustration for massage ..... 7
- Program selection ..... 8
- Changing the airbag cover ..... 9-10

## Troubleshooting/Warranty

- Maintenance/Storage/Transport ..... 11
- Storage environment ..... 12
- Troubleshooting ..... 13
- Specifications ..... 13

## Reflecting areas on foot

- Reflecting areas on foot ..... 14-15



## Safety Precautions

- Please read the “Safety Precautions” thoroughly before using massager
- The following precautions help you safely and correctly use this equipment, and avoid danger.
- The precaution items are categorized into “Warning” items or “Caution” items according to the degree of danger which may occur if such items are not abided by.
- After reading, please place this manual in a place that is easily accessible for future reading.

### Warning

Do not use power supply with unacceptable voltage	The use of unacceptable voltage may cause accident such as electric shock or fire.
When there is power outage	Please cut off the power switch and disconnect the power plug to extend the operating life of this equipment.
When there is failure with the device	Please pull out the power plug from the socket immediately and call the dealer for inquiry or service.
When not in use	Please turn off the power switch and pull out the plug from the socket.
When the power cable or plug is damaged	To avoid any danger, plug should be replaced by professional technicians from the manufacturer or its service agency or by other qualified personnel.

### Caution


The proper time of using this device is 15-30 minute. It may cause ache to muscle and nerve if you use it for a long period of time

- Do not use the device in the environment with high temperature or high humidity (places such as hot spring, bathroom, swimming pool, etc)
- Do not use the device at places with too much dust, oil or smoke.
- Please consult your doctor if you have had an operation or if your leg is injured.
- Do not use this device within 20 minutes after meals.
- In the following cases, the users should consult their doctors before using this device:  
Persons with malignant tumor, pregnant women, and those who recently had an operation done.

## Troubleshooting

Trouble	Check	Solution
Noise occur during massage	The device has been used for an extended period of time	Turn off the device and allow it to cool down for one hour
Hot surface	The device has been used for an extended period of time	Turn off the device and allow it to cool down for one hour
Control Panel has no light	Check whether the power cord is loose	Reconnect the power cord
No heating function	The heating pad connector may be loose	Unzip the Airbag Cover and inspect the connector
If the above solutions do not fix the problem, please contact iLIVING USA at 800-317-1688.		

## Specification

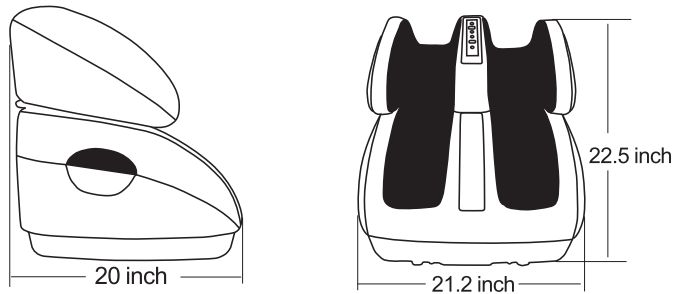

<b>Model</b>	MK9199 Plus
<b>Function</b>	Reduce the pressure of feet, promote blood circulation, eliminate tiredness, and improve sleep; 6 programs, 6 strengths and 30 combinations.
<b>Voltage</b>	AC 110V~60Hz <input type="checkbox"/> 220V~50Hz <input type="checkbox"/>
<b>Power</b>	120W
<b>Timing</b>	Each cycle is 15min. Automatic shut off occurs after 15min
<b>Security device</b>	Over current and over voltage protection
<b>Net/gross weight</b>	49 lbs / 56 lbs
<b>Dimensions</b>	20x21.2x22.5 inch

\* Subject to change without prior notice!

## Product Features

### ■ Foot Roller Function

The foot rolling function stimulates the acupressure points at the bottom of the feet to soothe stiff muscle and improve blood circulation.

### ■ More acupressure points

Flexible shiatsu, full range foot massage, and more stimulation

### ■ Full air pressure massage

When first using massager and air pressure technology, the intensity is soft, the feel is solid and the massage methods are more flexible. The intensity changes progressively and thus would be able to stimulate the acupuncture points better.

### ■ 6 Strengths, 6 Programs, and 30 Combinations

The multiple program combinations allow you to select your favorite strength and program, providing the perfect massage.

### ■ Size selection

Large: more than 40 yards; medium: 35-40 yards; small: less than 35 yards.

### ■ Far-infrared heating

The activated carbon cloth has the far-infrared heat function, which is the ideal light spectrum recognized by the whole medical profession.

The far-infrared will produce a resonance that will resonate with deeper human cells. This helps to reduce fatigue, accelerate blood circulation, activate cells, improve metabolism, and improve the immune system. It can efficiently prevent and adjuvant treats some diseases caused by microcirculatory disturbance (e.g. cold, rheumatism, rheumatoid etc.).

### ■ Removable cover

Allow for easier cleaning


### ■ Full range massage

The massage ranges from the bottom of the feet to the knees.

## Preparation before Use


### 1. The storage environment

- Put the device in dry places. Avoid exposure to direct sunlight.
- To reduce noise when using, operate the massager on carpet or rug.

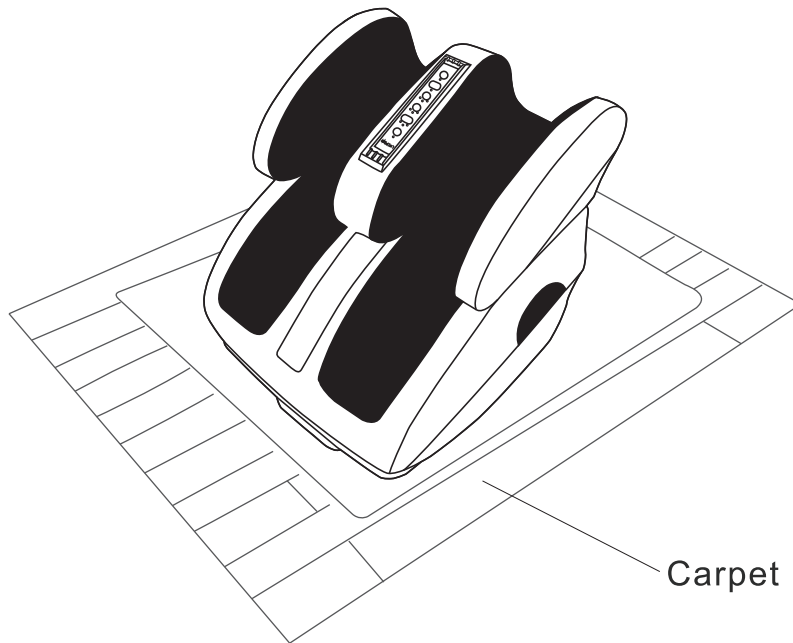

## Storage Instructions

Clean the device and put it in a place where it is dry and dust free. Do not use this device in places with too much dust and moisture. Do not put this device in places that are too close to a heating source or under direct sunlight. Do not use this device outdoor.

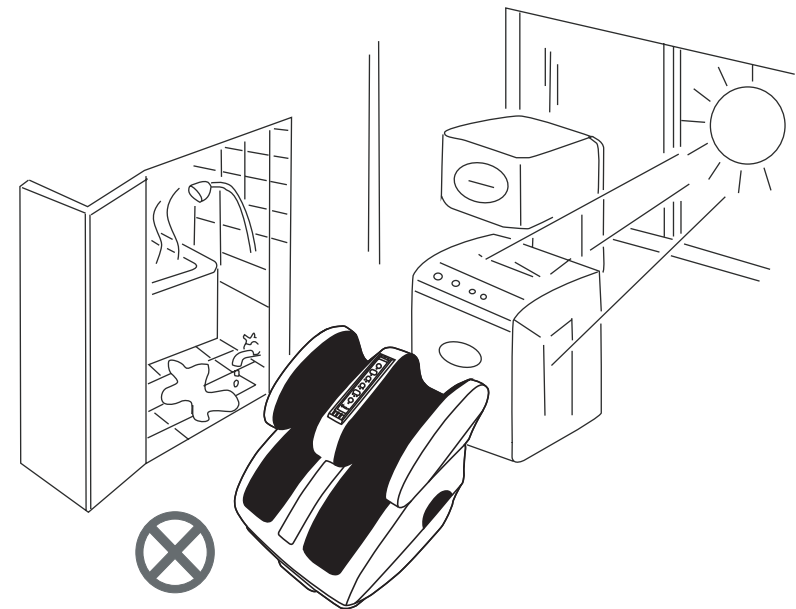

Maintenance

 Warning

Before maintenance and cleaning, please disconnect the power supply.

1. If the device is stained after using, please clean it with wet cloth soaked in neutral detergent (do not let chemicals such as detergent, alcohol or thinner dropping on the device) Then wipe with dry cloth and allow it to dry naturally. (other methods to dry it rapidly.)

Tip

\* Use for 15-30 minutes at a time for optimal effect.

Storage


The following conditions are necessary:

Temperature: 32F - 104F

Humidity: (no condensation) 30%RH – 85%RH

**Do not store the device in places with high temperature, humidity, radiated directly by sunlight, or where there are corrosion gases.**


Transport

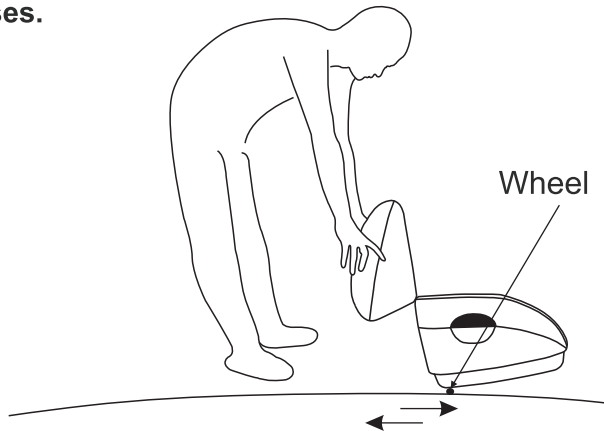

Preparation before Use


2. Connecting power supply

- (1) Insert the plug into the outlet.

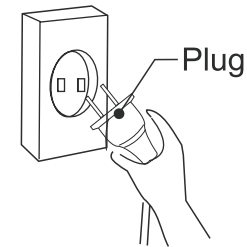

※

Make sure the cover of the power cord isn't peeled off or damaged.

※ Before inserting the plug, make sure the power switch  is in the OFF position.



Do not connect or disconnect plug with wet hands.

Operating Instructions

1. Lift the front holder before switching on.



2. Insert the plug into the outlet.



3. Switch on the power supply.


4. Press the Power  key to start the device.

5. Choose Size according to the size of your feet.

6. Choose "male/female" program to start the massage.


Product Composition

Front side view

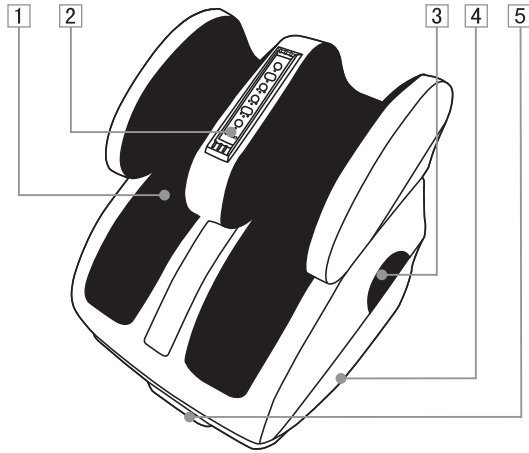

- 1 (Top/bottom) Airbag protection cover
- 2 Operating panel
- 3 Handle
- 4 Base
- 5 Front holder

Back side view

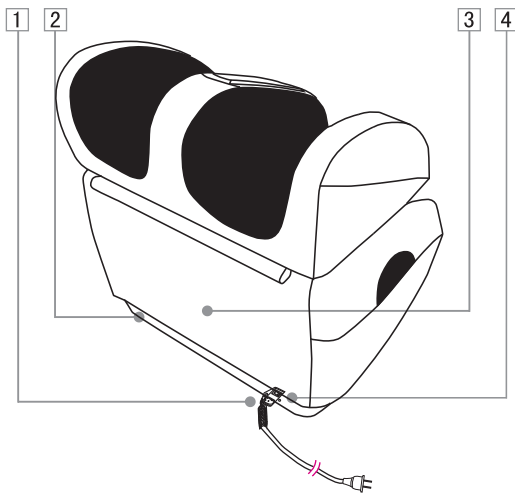


- 1 Power cord.
- 2 Wheel
- 3 Back cover
- 4 Power switch

Changing the Airbag Cover (Bottom)


To remove

1. Unzip the slide fastener.
2. Detach the carbon fiber heating cloth with VELCRO.

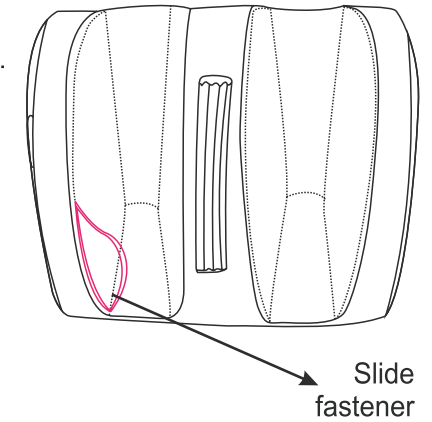

To install


1. Reconnect the male and female connector
2. Zip the slide fastener.

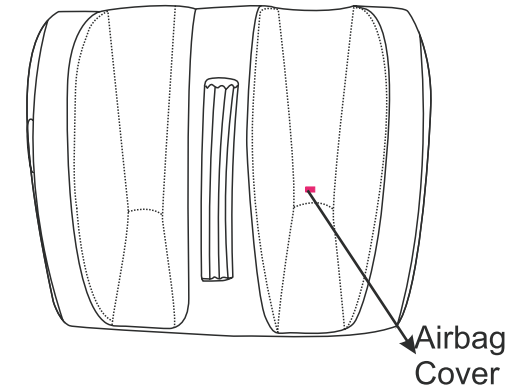




## Changing the Airbag Cover (Top)

1. Unzip the slide fastener and take out the protection cover.
2. Clean
3. Install the airbag cover (pay attention to the left and right direction)

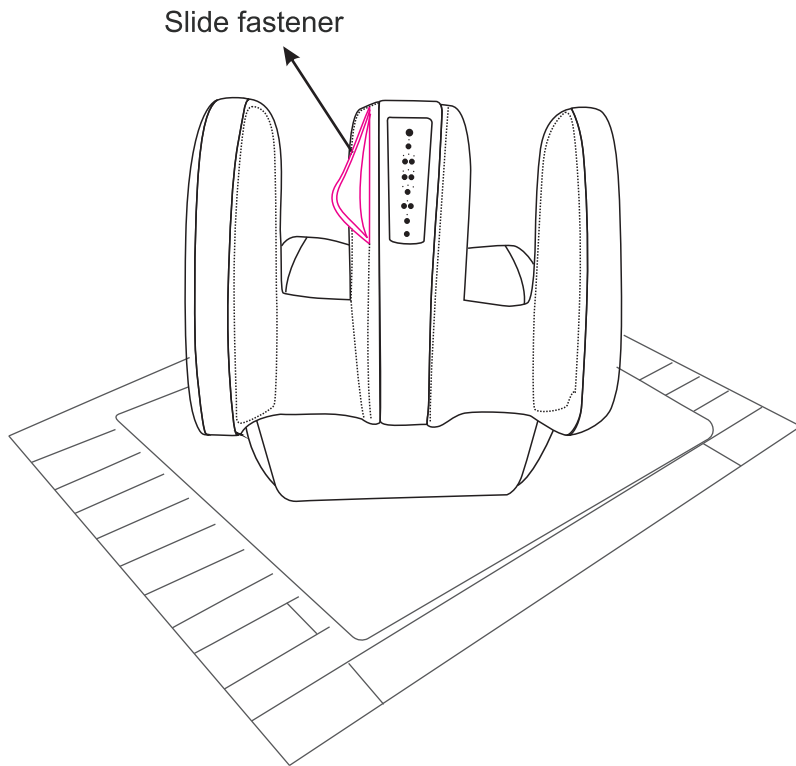


## Product Composition

### Control Panel

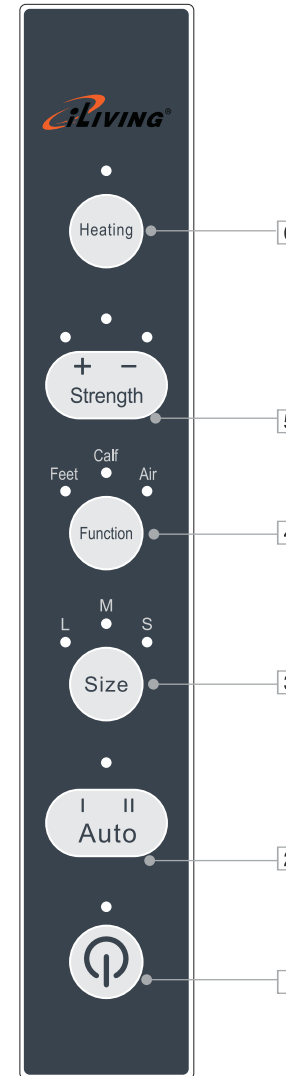


- 6 Heating**  
Press it once, it will be low-temperature hyperthermia, and the lights will be on for a long time. Press it again, it will be high-temperature hyperthermia and the lights will twinkle. Press the button for the third time, the hyperthermia function will be turned off.
- 5 Strength**  
Press once to increase one grade. There are six grades available.
- 4 Function selection**  
Press the button once to select Feet mode to start Scrapping and Shiatsu function. Press the button again to select Leg and Air Pressure massage function
- 3 Size**  
There are three different sizes.
- 2 Auto program selection**  
Press this button to select FOOT or LEG message function
- 1 On/Off switch**  
Press this button to turn on or turn off the foot massager.


Illustration for massage

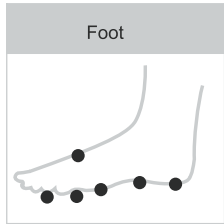


Sole reflection function

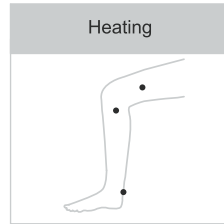

Heating function

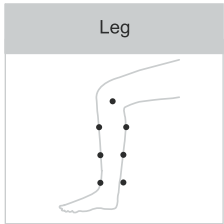


Blood circulation of legs, point relaxing function

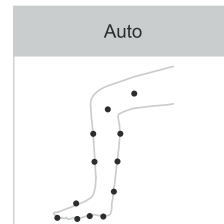

All massage points are engaged




Full foot massage using shiatsu

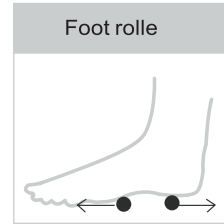

Roller moves up and down the foot to stimulate the acupunctural points

Strength Key

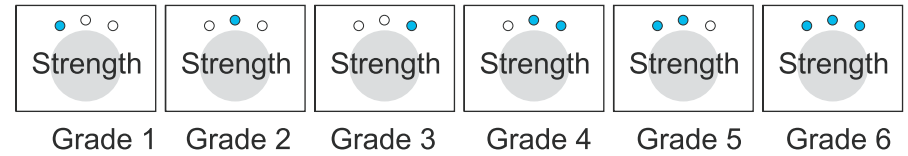

Program selection

Press SIZE button to choose size after power on the machine, then press AUTO1/AUTO2 button to choose your preferred programs

1. Press AUTO button once, the blue light stays on. AUTO 1 program (mainly for feet massage) is selected.
2. Press AUTO button again, the blue light starts flashing. AUTO 2 program is (mainly for leg massage) is selected.
3. Light indicator for FEET is on when Auto button is
4. pressed again, indicating the feet massage is starting.
5. Light indicator for CALF is on when Auto button is pressed again, indicating the calf massage is starting.
6. Light indicator for AIR is on when Auto button is pressed again, indicating air pressure is starting.