


MEMORISM PARAFFIN WAX MACHINE

AURA PWM65

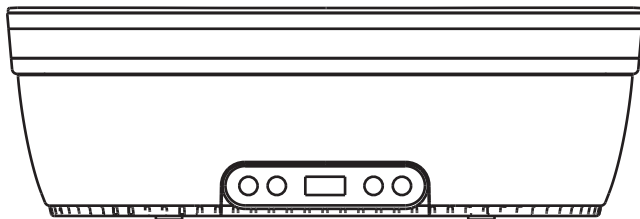

USER MANUAL

Please read this manual carefully before use and keep it for future reference.

CONTENTS

Introduction	01
Important Safety Instructions	01
Specifications	03
How to Use	04
Cleaning and Maintenance	07
Limited Warranty	08

INTRODUCTION

Congratulations on your purchase and welcome to Memorism!

To fully benefit from the support that Memorism offers please register your product at www.mymemorism.com/register.

IMPORTANT SAFE INSTRUCTIONS

When using an electrical appliance, especially when children are present, basic safety precautions should always be followed, including the following:

1. To avoid burns, electrocution, fire, or injury, never use the appliance if it has a damaged cord or plug, if the appliance has been damaged, or if the appliance doesn't work properly.
2. This product should never be used by any individual suffering from a physical ailment that would limit the user's capacity to operate the controls.
3. This product should not be used in the presence of open cuts or wounds. If you have any concerns, consult a physician or physical therapist before using it.
4. Always unplug the product when not in use.
5. Do not plug or unplug this product when your hands are wet.
6. Do not place the appliance where it can fall or be pulled into a tub or sink. Do not place in or drop into water or other liquid. Store in a cool, dry place.

7. Keep the cord away from the heated surface.
8. Do not carry this appliance by the supply cord or use the cord as a handle.
9. Never use sharp instruments to remove solidified paraffin wax from the unit.
10. Never operate the unit without paraffin wax.
11. Use this appliance only for its intended use as described in this manual.
12. This appliance can be used by children aged 8 years and above and persons with reduced physical, sensory or mental capabilities or lack of experience and knowledge if they have been given supervision or instruction concerning the use of the appliance in a safe way and understand the hazards involved.
13. Children shall not play with the appliance.
14. Cleaning and user maintenance shall not be made by children without supervision.
15. This device complies with Part 15 of the FCC Rules. Operation is subject to the following two conditions: (1) this device may not cause harmful interference, and (2) this device must accept any interference received, including interference that may cause undesired operation.

SPECIFICATIONS

Product Name: Paraffin Wax Machine

Model: PWM65

Power Supply: 110-120V

Power: 200W

Capacity: 2000ml aluminum pot + 4500ml rubber pot

Temperature Adjustable Range: Hand Care mode 86-158°F
Foot Care mode 86-194°F

Time Adjustable Range: 1-24 hour

Product Size: 15 x 8.9 x 5 inch

Product Weight (without wax): 6.6 lb

HOW TO USE

This machine has 2 separate pots to meet different demands (e.g. hands and feet). Please remove jewelry and choose the correct mode before use.

- Instructions on Foot Care Mode

1. Make sure the unit is on a level surface, and plug it in. Remove the lid and put wax in the **rubber pot**. (We recommend adding 2-3 packs of wax.)
2. Press the power icon to turn it on. The default mode is **Foot Care** when it is powered on.

- Instructions on Hand Care Mode

1. Make sure the unit is on a level surface, and plug it in. Remove the lid and put wax in the **aluminum pot**. (We recommend adding 2-3 packs of wax.)
2. Press the power icon to turn it on and long press the clock icon for 3s to switch to **Hand Care** mode.

3. Adjust the Temperature Heating, Time Heating, and Reservation Heating as follows.

① Temperature Heating

- Press +/- to adjust the heating temperature. After flashing for 3s the temperature is set and the machine starts to heat at the set temperature. The screen shows the real-time temperature.
- Upon reaching the set temperature, the machine will beep 3 times. The temperature will keep at the set temperature until powered off.

* **We recommend setting 122-131°F, which is the proper temperature to use.**

② Time Heating

- Press the clock icon to set the heating time. The screen will show "h30" as a default, which means 30 minutes. Press +/- to adjust from 1-90 minutes.
- After flashing for 3s the time is set and the machine starts to heat. The screen shows the countdown time. Under this mode, the real-time temperature will not be shown.
- After the countdown is finished, the screen displays the real-time temperature and the machine will automatically cool down. It needs around 35 minutes to cool down to the proper temperature (127°F) and beep 3 times. It will keep at this temperature until powered off.

* **We recommend setting 40 minutes to melt 2-3 packs of wax.**

③ Reservation Heating

- Press the clock icon twice to set a time reservation. The screen will show “H01” as a default, which means 1 Hour. Press +/- to adjust from 1-24 Hours. After flashing for 3s the reservation is set and the machine starts to countdown.
 - The machine will heat in the last hour of the the countdown time shown on the screen.
 - After the countdown is finished, the screen displays the real-time temperature and the machine will automatically cool down to the proper temperature (127°F) and beep 3 times. It will keep at this temperature until powered off.
4. After the wax is totally melted, **insert the heat-insulating grid when using the aluminum pot.** Spread your fingers/toes and slowly dip your hand or foot in the machine to test the temperature. Immerse your hand or foot in the wax for 2s.
 5. Gently remove your hand or foot from the machine and wait for 3s.
 6. Repeat steps 4-5 for three to six times. You can also use the brush to apply wax.
 7. Cover your hand or foot with a plastic liner. Then put your hand or foot into the thermal mitt to retain heat for 15 minutes.
 8. After 15 minutes, remove the mitt and liner and peel the hardened wax from the top of your wrist or ankle. Discard the used plastic liner and wax.
 9. When finished, power off the machine and unplug. The remaining wax can be used for next time.

CLEANING AND MAINTENANCE

Clean the pot when paraffin wax becomes discolored or dry.

- Do not scratch or damage the inside coating.
- Do not immerse the machine in water.
- Unplug the unit and cool down before cleaning.

To clean the rubber pot

1. Let the wax cool down and become solid.
2. Press upwards the bottom of the pot. One side of the wax will tilt, then you can remove the wax easily.
3. Wipe with a cloth or paper towel to remove the remaining wax.

To clean the aluminum pot

1. Let the wax cool down and become solid.
2. Plug the unit in until the wax loosens from the pot, then unplug the unit.
3. Push down firmly on one side of the wax to tilt the opposite side up. Lift out the wax and discard.
4. Wipe with a cloth or paper towel to remove the remaining wax when it is a little warm. The remaining wax is hard to clean when it is totally cool and solid.

LIMITED WARRANTY

Memorism warrants Paraffin Wax Machine PWM65 for 1 year from the date of purchase. For more information, please visit www.mymemorism.com/warranty.

What's covered?

Defects due to faulty materials and/or workmanship of the product.

What's not covered?

1. Damage caused by the owner when attempting to fix or altering the product by himself or herself.
2. Damage caused by misuse, abuse, neglect, alterations, or unauthorized repair.
3. Natural depreciation.

Warranty Period

Warranty becomes effective from the date of original purchase and will last for 1 year.

Who's covered?

This limited warranty is only effective for the original purchaser. It is non-transferable.

How to request warranty services?

To obtain the benefit of this warranty, please email support@mymemorism.com.

Our customer service team will address your issue as soon as possible!

MEMORISM



www.mymemorism.com
support@mymemorism.com

6945 Speedway Blvd Ste, H102 Las Vegas, NV 89115