150mm


Dr Dennis Gross

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Welcome to DRx SpectraLite™ BodyWare Pro, your FDA cleared, anti-aging, anti-inflammatory and anti-acne light therapy body device.

Using the high-performance power of natural light therapy to treat acne, dark spots, and aging in just 3 minutes a day - plus intermittent relief of muscle and joint pain in 14 minutes. Inspired by the same medical-grade LED technology used in Dr. Dennis Gross' dermatology offices, DRx SpectraLite™ BodyWare Pro is an at-home, wearable device that molds to fit each body zone. Anti-aging and anti-inflammatory red light helps stimulate skin's natural collagen production for firmer, smoother skin and relief of minor muscle and joint pain. Blue light kills acne-causing bacteria and prevents future body breakouts.

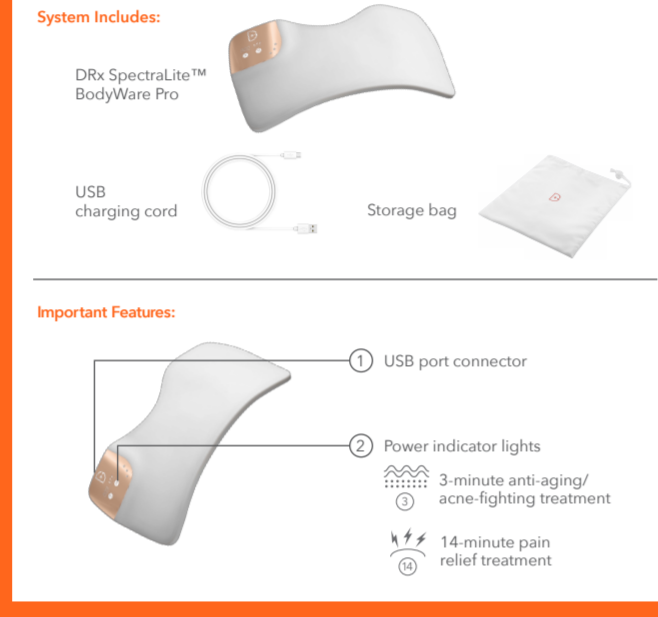
USER GUIDE AND MANUAL Includes important safety information. Read all instructions prior to use. SAVE THESE INSTRUCTIONS MANUEL D'UTILISATION Contient des informations importantes. Bien lire toutes les instructions avant utilisation. CONSERVER CES INSTRUCTIONS BHS658810.081

IMPORTANT SAFETY INFORMATION
Do not use on face or near eyes.
Do not operate device where it could come in contact with water.
Do not use on or charge device with power adapter that is damaged and/or faulty.
Do not handle device or above area of the USB-C charging cord.
Do not attempt to take device apart. Device is not serviceable by user. This will void the warranty.
Stop using device immediately if you think it is defective. Please contact Client Care at +1 888 830 SKIN.

CONTRAINDICATIONS
Please consult your physician prior to use, if:
1. You are or suspect you are pregnant.
2. You are under treatment or suspect you may have thyroid ailments.
3. You are taking sensitizing medications.
TEST YOURSELF FOR LIGHT SENSITIVITY PRIOR TO USING YOUR FULL BODY
1. Flare up on your device.
2. Place the device on your arm and conduct a full 3 minute treatment.
3. After treatment is complete, visually check the area. If the area has turned red and lasts for more than 2 hours, your skin has a sensitivity to light. Discontinue use and contact Client Care at +1 888 830 SKIN.

DEVICE OVERVIEW
DRx SpectraLite™ BodyWare Pro is a Light Emitting Diode (LED) device that harnesses the energy of 148 anti-aging red lights and 16 breakout fighting lights (red/blue) to help minimize the appearance of dark spots, enhance skin tone, and treat wrinkles, acne, and mild pain over the entire body. The device incorporates natural light energy in the Blue (415nm), Amber (630nm), Red (630nm), Deep Red (650nm), and Infrared (880nm) spectrums, that are proven to eliminate acne-causing bacteria, reduce fine lines and wrinkles, and temporary relief of minor muscle and joint aches, pains, and stiffness.
How It Works
Youthful skin is rich in collagen, which lifts and supports the skin. Collagen works deep within skin like columns of tissue to support skin's surface. Over time, skin produces less collagen, and with exposure to sunlight, healthy collagen becomes damaged. This can result in lines and wrinkles. Medical-grade red light stimulates collagen production and stimulates fine lines and wrinkles.

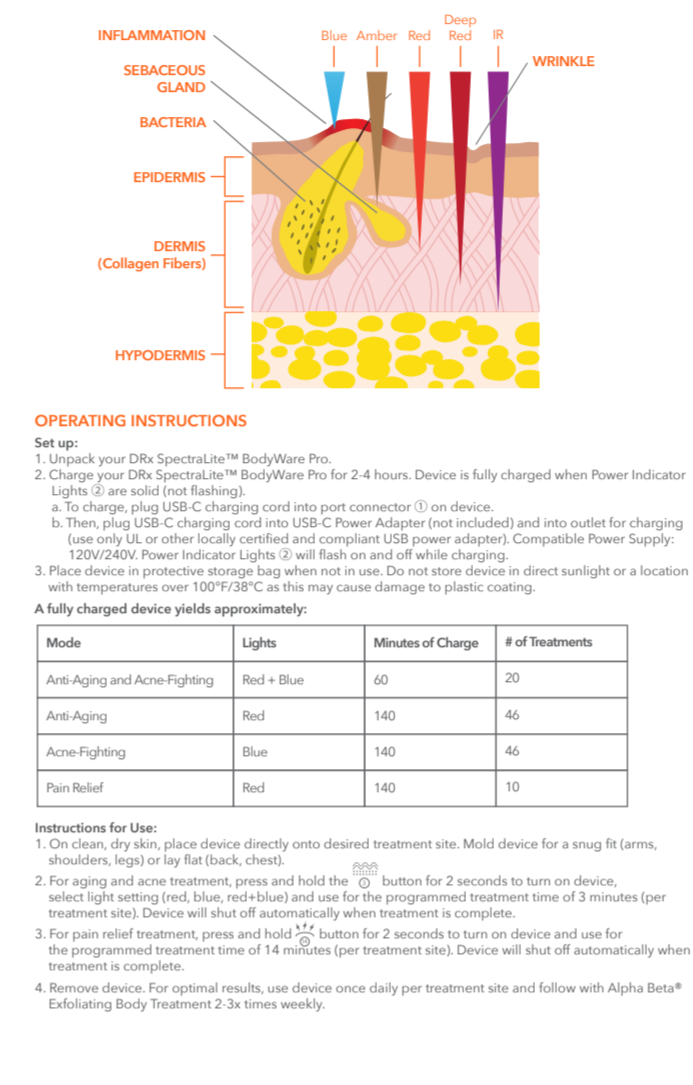OPERATING INSTRUCTIONS
Set up:
1. Unpack your DRx SpectraLite™ BodyWare Pro.
2. Charge your DRx SpectraLite™ BodyWare Pro for 2-4 hours. Device is fully charged when Power Indicator Lights (2) are solid (not flashing).
3. Plug USB-C charging cord into port connector (1) on device, and into power adapter (not included).
4. Charge device only with the supplied USB-C charging cord and a USB power adapter (not included).
5. Voltage of BodyWare Pro: 3.7V
6. Compatible Power Supply: 120V/240V
Charging:
1. Turn device off and unplug from any power source.
2. Wipe the inside of device with an alcohol swab or damp cloth to remove any residue from skin oils after each use. This is especially important if the device is being used by more than one individual.
3. Clean device by using a damp or microfiber cloth. The microfiber storage bag can be used for cleaning and polishing the body of device. Do not use alcohol or alcohol-based cleaners on exterior as health chemicals may discolor the plastic panels.
4. Allow device to dry completely prior to use.

LIMITED WARRANTY INFORMATION
We value you as a client and strive for the highest quality products to earn your satisfaction. Please use this as your reference for any customer service or warranty needs you may have when purchased within the United States. You will register your product within 30 days of purchase at www.dennisdgross.com/registration for the warranty to be effective. When purchased outside the United States, please contact your place of purchase for local warranty information. Our warranty is limited to the terms set out below.
What is covered by the limited warranty?
This 1-year limited warranty covers defects in materials and workmanship in your SpectraLite™ BodyWare Pro when purchased through an authorized partner or directly from Dr. Dennis Gross Skincare™. Your 1-year limited warranty is valid for you alone, the original purchaser, and cannot be transferred. If the product was a gift, obtain a copy of the original receipt and send to registration@dennisdgross.com, or have the original purchaser process the warranty claim on your behalf.

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