Dr Dennis Gross

DRx SpectraLite™  
**EyeCare Max Pro**

DRx SpectraLight™  
**Masque LED EyeCare Max Pro**

DRx SpectraLight™  
**Dispositivo EyeCare Max Pro**

DRx SpectraLight™  
**EyeCare Max Pro**

DRx SpectraLight™  
**EyeCare Max Pro**

**USER GUIDE AND MANUAL**

Includes important safety information. Read all instructions prior to use.

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**SPARA DESSA ANVISNINGAR**

# Welcome to DRx SpectraLite™, your state-of-the-art, anti-aging eye device.

HELLO, BRIGHT EYES! Get ready to show off the look of rested, revived, youthful eyes...in just 3 minutes a day. DRx SpectraLite™ brings the latest professional light therapy to your home in a next-generation, hands-free, wearable design. Maximizing the full spectrum of anti-aging wavelengths, DRx SpectraLite™ features 96 LED lights to reduce wrinkles, crow's feet, and elevens by stimulating natural collagen production. This is the future of your eyes, and it's looking bright.

## System Includes:

DRx SpectraLite™  
EyeCare Max Pro

Durable, silicone strap  
that easily adjusts to  
comfortably fit all head  
shapes

USB  
charging cord

Microfiber  
Storage Bag

## Important Features:

① ON/OFF button;  
Power indicator light

② USB port connector

## IMPORTANT SAFETY INFORMATION

- Do not operate device where it could come in contact with water.
- Do not use or charge device with power adapter that is damaged and/or faulty.
- Do not harshly bend or abuse the wire of the USB charging cord.
- Do not attempt to take the device apart. Device is not serviceable by user. This will void the warranty.
- Stop using the device immediately if you think it is defective. Please contact Client Care at +1-888-830-SKIN.

## Client Care

You can reach Client Care by calling +1-888-830-SKIN Monday-Friday 9:00AM-5:00PM Eastern Standard Time. Client Care can also be reached by e-mail at [customerservice@drdennisgross.com](mailto:customerservice@drdennisgross.com)

## Disclaimer

This product is not intended to diagnose, cure, or prevent any disease. Dr. Dennis Gross Skincare® makes no claims, representations, or warranties regarding the ability of this product to cure any physical, skin, or mental condition. A qualified health professional should be consulted with regards to any condition requiring medical attention.

## Product Representation

Dr. Dennis Gross Skincare® is constantly improving products with the latest technology. As a result, the actual product may differ slightly from the product shown in images on the packaging and in this user guide.

## WARNINGS

- Do not aim or look directly into the lights. Doing so may cause temporary vision discomfort. Placement of the device should ensure that no part of the device is covering the eyes.
- This device is designed for use on adults over 18 years of age and not recommended for use on infants or children.
- Please consult your doctor prior to use if you are, or suspect you are, pregnant.
- Please consult your doctor prior to use if you have, or suspect you have, cancer.
- Do not use device if it is wet or if it could fall into water while it is plugged in. If device has fallen into water, do not touch until you first unplug it.
- Do not treat if your skin is red or swollen from a previous treatment. Wait until the redness and swelling have disappeared before treating again. If you have skin redness or swelling that lasts for more than 2 hours, you may have light-sensitive skin. Stop using the device and contact a healthcare professional.
- Do not use device to treat skin that is painful to touch, raised (higher than surrounding skin), crusting, or bleeding. This may be a sign of a more serious condition and should be checked by your doctor.
- No modification of this equipment is allowed.
- This device contains a Lithium Ion Battery. Before disposal of this device, please review your local laws and requirements surrounding Lithium Ion Battery disposal. The preferred method of disposal is through recycling. You can visit <http://www.call2recycle.org/lithiumion-battery-safety/> and search for your nearest battery recycling center.

## CAUTIONS

- Discontinue use if device overheats or becomes uncomfortably hot.

- Always unplug device from outlet after charging is complete.
- Do not store device in direct sunlight or location with temperatures over 100°F (38°C) as this may cause damage to the plastic coating.
- Use only the supplied Universal USB Power Cord. Use of incorrect power cord could damage device.
- Do not charge with fast charge adapters designed to rapidly charge mobile devices and laptops. The USB C charging cord is designed to work with any standard universal USB A (5V/1 or 2A standard UL/ETL rated) power adapter.
- Do not attempt to open or repair device or power supply.
- Do not drop device on hard surfaces, as this could damage it. Using a damaged device could result in skin or eye injury.
- Do not use device on any area of the body other than around the eyes. Using the device on other areas of the body has not been tested and the risks are unknown.
- If the area of skin for which you are using the device to treat wrinkles appears to be worsening, you should cease use and check with your doctor to make sure you do not have a more serious condition.
- Contact your doctor if you have any medical questions regarding use of this device.

## CONTRAINDICATIONS

### Ask your doctor prior to using if you:

- Are taking prescription or over-the-counter drugs known to cause light sensitivity.
- Have any suspicious or cancerous lesions.

### Test yourself for light sensitivity prior to using on your face:

1. Power on your device.
2. Place the device on your arm and conduct a full 3-minute treatment.
3. After treatment is complete, visually check the area. If the area has turned red and lasts for more than 2 hours, your skin has a sensitivity to light. Discontinue use and contact Client Care at +1-888-830-SKIN.

## DEVICE OVERVIEW


DRx SpectraLite™ EyeCare Max Pro is a Light Emitting Diode (LED) device that harnesses energy to treat signs of aging around the eye area. The device incorporates natural light energy in the Amber (605nm), Red (630nm), Deep Red (660nm), and Infrared (880nm) spectrums that are proven to reduce fine lines and wrinkles.

### How it Works

Youthful skin is rich in collagen, which gives the skin lift and support. Collagen works deep within the skin like columns of tissue that support the skin's surface. Over time, skin produces less collagen, and with exposure to sunlight, healthy collagen becomes damaged. This can result in wrinkles, especially around the eyes (crow's feet and elevens). The lights from this device penetrate deep within the layers of your skin to stimulate collagen production, diminishing fine lines and wrinkles. 96 red LED lights to even tone and help firm skin below the eyes.

Easy to use, DRx SpectraLite™ is programmed with a 3-minute automatic shutoff function for professional anti-aging results with short, simple, daily treatments in the comfort of your own home.

### Light Wavelength Diagram
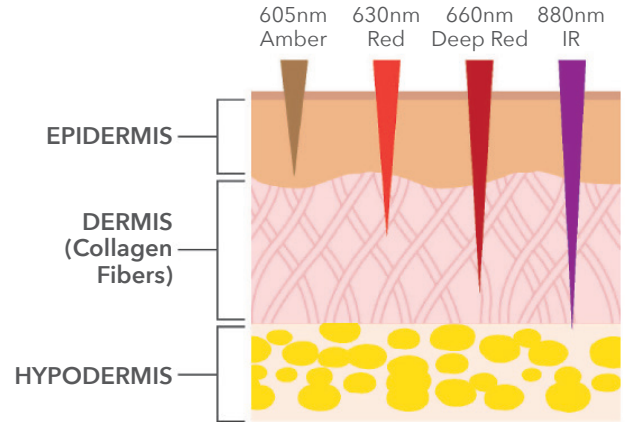
### What to Expect

Users typically see initial results within 2 weeks, with the best results occurring after a 10-week course of treatment.

## OPERATING INSTRUCTIONS

### Set up:

1. Unpack your DRx SpectraLite™.
2. Plug USB power/charging cord into port connector @ on device.
3. Plug USB power/charging cord into Universal USB Power Adapter (not included) and into outlet for charging (Use only a standard universal USB A (5V/1 or 2A standard UL/ETL rated) power adapter).
4. Charge your DRx SpectraLite™ for 2-4 hours. Device is fully charged when Power indicator light ☉ is solid (not flashing).
5. Place device in protective storage bag when not in use. Do not store device in direct sunlight or location with temperatures over 100° Fahrenheit as this may cause damage to plastic coating.
6. A fully charged device yields 30 minutes of use, approximately 10 treatments.

### Instructions for Use:

1. Wash and fully dry face prior to use.
2. Place device onto face over eyes and adjust strap to fit snugly around head.
3. Power device on. DRx SpectraLite™ features an automatic shutoff function after completion of a 3-minute treatment.
4. Remove device from face and wipe device clean (see cleaning instructions).
5. Place device into protective storage bag.
6. After treatment, apply creams or topical treatments as desired.

### Helpful Tips for Best Results:

- Use daily.
- Use device on clean, dry skin only.
- Do not use device after the application of creams that contain SPF as this could block light penetration.
- Wipe device clean after each use (see cleaning instructions).

### Keep Skin Hydrated

Immediately following treatment, apply a serum or moisturizer of your choice.

### Maintenance:

1. It is important to keep the device clean (see cleaning instructions).
2. This device has no serviceable parts. Battery is not accessible and not serviceable. User should not attempt to access or service the battery as this will damage the device and void the warranty.
3. Device includes only an ON/OFF button. All settings are set to the factory setting.
4. Charge device only with the supplied USB power/charging cord and a universal USB power adapter (not included).
5. Power indicator light will begin to flash when battery levels are low. This indicates it is appropriate to begin recharging the device.
6. Compatible Power Supply: 120V/240V

### Cleaning:

1. Turn device off and unplug from any power source.
2. Wipe the tops of the LEDs with alcohol wipes or damp cloth to remove any residue from skin oils after each use. This is especially important if the device is being used by more than one individual.
3. Clean device body by using damp or microfiber cloth. The microfiber storage bag can be used for cleaning and polishing the body of device. Do not use chemical or alcohol based cleansers on exterior as harsh chemicals may discolor the metallic panels.
4. Allow device to dry completely prior to use.

### Troubleshooting:

1. System is not turning on:
  - a. Is the device fully charged? (Power indicator light is on and solid)
  - b. Is the Power indicator light illuminated?

- c. If device still does not turn on, call Client Care at +1-888-830-SKIN.
- 2. System is not charging:
  - a. Is USB power cord firmly inserted into device? (does Power indicator light turn on?)
  - b. Try a different USB Power Adapter to ensure it is working properly.
  - c. Try a different wall outlet to determine if there is power to the outlet.
  - d. If device still does not charge, call Client Care at +1-888-830-SKIN.
- 3. System shuts off abruptly, prior to completion of 3-minute treatment:
  - a. System is designed to shut off automatically after 3 minutes.
  - b. System may shut off after multiple repeated uses if the device exceeds desired temperature. Allow device to cool for several minutes and try again.
  - c. If device continues to shut off prior to 3 minutes, please contact Client Care at +1-888-830-SKIN.
- 4. Some lights appear turned off:
  - a. The infrared (IR) light spectrum is not visible to the human eye. These lights will look as if they are not on, but they are operating properly. If the red, visible lights are on, then the device is operating as designed.

## LIMITED WARRANTY INFORMATION

Congratulations on your purchase. We value you as a client and strive for the highest quality products to earn your satisfaction. Please use this as your reference for any customer service or warranty needs you may have when purchased within the United States. You must register your product within 30 days of purchase at [www.drdenisgross.com/registration](http://www.drdenisgross.com/registration) for the warranty to be effective. When purchased outside the United States, please contact your place of purchase for local warranty information. Our warranty is limited to the terms set out below:

### ***What is covered by the limited warranty?***

This 1-year limited warranty covers defects in materials and workmanship in your DRx SpectraLite™ when purchased through an authorized partner or directly from Dr. Dennis Gross Skincare®. Your 1-year limited warranty is valid for you alone, the original purchaser, and cannot be transferred. If the product was a gift, obtain a copy of the original receipt and send to [registration@drdenisgross.com](mailto:registration@drdenisgross.com), or have the original purchaser process the warranty claim on your behalf.

### ***What is not covered by this limited warranty?***

This limited warranty does not cover: damage, defects, failure caused by abuse or misuse, improper care, impact or insertion of foreign objects, use with a voltage connector, or loss of device.

### ***How long does this limited warranty last?***

The limited warranty begins after the device has been registered and expires after one year. The warranty period is not extended if your device is replaced by Dr. Dennis Gross Skincare® during the warranty period.

### ***What will Dr. Dennis Gross Skincare® do if your device is defective?***

If you discover a defect in your DRx SpectraLite™ covered by the warranty, before the 1-year warranty period expires, contact Client Care at +1-888-830-SKIN for assistance. If we confirm a defect covered by this limited warranty, Dr. Dennis Gross Skincare® will decide whether to replace your DRx SpectraLite™ with either a new or refurbished device. If we find that your DRx SpectraLite™ is not defective, we will notify you and return the device to you. All Limited Warranty services must be authorized and performed by Dr. Dennis Gross Skincare® or a service provider authorized by Dr. Dennis Gross Skincare®.

## WARRANTY LIMITATIONS & DISCLAIMERS

This warranty gives you specific legal rights, and is in addition to your statutory rights. Your remedies for the breach of this limited warranty are limited solely to the repair and replacement of your device as set forth in this warranty statement.

To the maximum extent permitted by applicable law, Dr. Dennis Gross Skincare® disclaims all implied warranties. If applicable law does not permit the disclaimer of implied warranties, then the duration of implied warranties are hereby expressly limited to the duration of the express limited warranty set forth above.

In no event shall Dr. Dennis Gross Skincare® be liable for any special, incidental, consequential, punitive or indirect damages based upon breach of warranty, breach of contract, tort (including negligence or strict product liability) or otherwise, even if Dr. Dennis Gross Skincare® is advised or aware that such damages may arise. Such damages include, but are not limited to: loss of DRx SpectraLite™, cost of any substitute equipment, the claims of third parties or damage to property.

If the applicable law does not allow limits on warranties, or on remedies for breach in certain transactions, the limits in this paragraph and the disclaimer or exceptions from warranty may not apply. This limitation does not apply to claims for death or personal injury caused by a defect in the device or any other liability that cannot be excluded or limited by law.

These warranty terms, conditions, and limitations allocate the risks of product failure between Dr. Dennis Gross Skincare® and you, the purchaser of DRx SpectraLite™. This allocation is recognized by both parties and is reflected in the price of DRx SpectraLite™. By using DRx SpectraLite™ you acknowledge that you have read these terms, conditions, and limitations, understand them and are bound by them.

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Dr. Dennis Gross Skincare® LLC.
UK AR: Obelis UK Ltd., Oxford, OX4 6LB