Dr Dennis Gross

DRx SpectraLite™

FaceWare Pro

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Dispositivo LED Para el Rostro

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#### **USER GUIDE AND MANUAL**

Includes important safety information.  
Read all instructions prior to use.

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**SPARA DESSA ANVISNINGAR**

BH568180.044

# Welcome to SpectraLite™ FaceWare Pro, your next generation, full face, anti-aging and anti-blemish LED device.

Using the high-performance power of LED technology, help target spots, blemishes, lines and wrinkles over your entire face in as quickly as 3 minutes a day. Inspired by the same LED technology used in Dr. Gross' dermatology office, Spectralite™ FaceWare Pro is an at-home, hands-free, wearable device that uses anti-aging red light and blemish-reducing blue light to target wrinkles and spots. FaceWare Pro helps to stimulate skin's collagen production to reduce even the deepest lines, while reducing the appearance of blemishes for a smooth, beautiful, youthful-looking complexion.

## System Includes:

SpectraLite™  
FaceWare Pro

Durable, silicone,  
adjustable strap


USB  
charging cord

Storage bag

## Important Features:
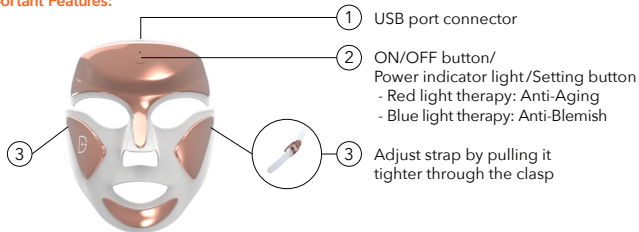
## IMPORTANT SAFETY INFORMATION

- Do not operate device where it could come in contact with water.
- Do not use or charge device with power adapter that is damaged and/or faulty.
- Do not harshly bend or abuse wire of the USB charging cord.
- Do not attempt to take device apart. Device is not serviceable by user. This will void the warranty.
- Stop using device immediately if you think it is defective. Please contact Client Care at +1-888-830-SKIN.

## Client Care

For information regarding device usage, warranty, and additional information, please contact Client Care by calling +1-888-830-SKIN Monday-Friday 9:00AM-5:00PM Eastern Time. Client Care can also be reached through e-mail at [customerservice@drdennissgross.com](mailto:customerservice@drdennissgross.com)

## Disclaimer

This product is not intended to diagnose, cure, or prevent any disease. Dr. Dennis Gross Skincare™ makes no claims, representations, or warranties regarding the ability of this product to cure any mental, physical, or skin condition. A qualified health professional should be consulted with regards to any condition requiring medical attention.

## Product Representation

Dr. Dennis Gross Skincare™ is constantly improving products with the latest technology. As a result, the actual product may differ slightly from the product shown in images on the packaging and in this user guide.

## WARNINGS

- Do not look directly into the lights. Doing so may cause temporary vision discomfort. Placement of the device should ensure that no part of the device is covering the eyes.
- This device is designed for use on adults over 18 years of age and not recommended for use on children under the age of 12.
- Please consult your doctor prior to use if you are, or suspect you are, pregnant.
- Please consult your doctor prior to use if you have, or suspect you have, cancer.
- Do not place device near water or any wet surface.

- Do not use device if it is wet. If device comes in contact with water, do not touch it until you first unplug it.
- Do not use on red or swollen skin caused by a previous session. Wait until the redness and swelling have disappeared before using the device again. If you have skin redness or swelling that lasts for more than 2 hours you may have light-sensitive skin. Stop using the device and contact a healthcare professional.
- Do not use device on open wounds.
- Do not use device to target skin that is painful to touch, raised (higher than surrounding skin), crusting, or bleeding. This may be a sign of a more serious condition and should be checked by your doctor.
- No modification of this equipment is allowed.
- This device contains a Lithium Ion Battery. Before disposal of this device, please review your local laws and requirements surrounding Lithium Ion Battery disposal. The preferred method of disposal is through recycling. For North America, you can visit <http://www.call2recycle.org/lithiumion-battery-safety/> and search for your nearest battery recycling center.

## CAUTIONS

- Discontinue use if device overheats or becomes uncomfortably hot.
- Always unplug device from outlet after charging is complete.
- Do not store device in direct sunlight or a location with temperatures over 100°F/38°C as this may cause damage to the plastic coating.
- Use only the supplied Universal USB charging cord. Use of incorrect charging cord could damage device.
- Do not attempt to open or repair device or power supply.
- Do not drop device on hard surfaces, as this could damage it. Using a damaged device could result in skin injury.
- Do not use device on any area of the body other than the face. Using the device on other areas of the body has not been tested and the risks are unknown.
- If the area of skin for which you are using the device to target wrinkles and/or blemishes appears to be worsening, you should cease use and check with your doctor to make sure you do not have a more serious condition.
- Contact your doctor if you have any medical questions regarding use of this device.

## CONTRAINDICATIONS

### Ask your doctor prior to using if you:

- Are taking prescription or over-the-counter drugs known to cause light sensitivity.
- Have any suspicious or cancerous lesions.

### Test yourself for light sensitivity prior to using on your face:

1. Power on your device.
2. Place the device on your arm and conduct a full 3-minute test.
3. After test is complete, visually check the area. If the area has turned red and lasts for more than 2 hours, your skin has a sensitivity to light. Discontinue use and contact Client Care at +1-888-830-SKIN.

## DEVICE OVERVIEW

SpectraLite™ FaceWare Pro is a Light Emitting Diode (LED) device that harnesses the energy of 100 anti-aging red lights and 60 blemish-reducing blue lights to reduce the appearance of spots, blemishes, lines and wrinkles while improving your natural skin tone. This device incorporates natural light energy in the Blue (415nm), Amber (605nm), Red (630nm), Deep Red (660nm), and Infrared (880nm) spectrums, that is proven to reduce the appearance of blemishes and helps to increase collagen production to minimize lines and wrinkles.

### How it Works


Youthful skin is rich in collagen, which lifts and supports the skin. Collagen works deep within skin like columns of tissue to support skin's surface. Over time, skin produces less collagen, and with exposure to sunlight, healthy collagen becomes damaged. This can result in lines and wrinkles. The spectrum of red lights in this device helps to increase collagen production to diminish the look of lines and wrinkles.

Balanced, healthy skin is clear of spots and blemishes. Unbalanced skin is due to excess sebum and dirt, leading to clogged pores that can become unsightly blemishes. Our blue and red lights target these spots to reduce redness, minimize the appearance of spots, and help prevent future blemishes from forming.

Easy-to-use SpectraLite™ FaceWare Pro is a state-of-the-art, anti-aging and anti-blemish device that delivers quick, daily 3-minute sessions in the comfort of your own home.

**This device contains 3 settings:**

- All Red Mode: Anti-Aging Lights
  - Helps to increase the production of collagen
  - Helps reduce fine lines and wrinkles
  - Evens skin texture
  - Helps minimize appearance of sunspots, age spots and enhances skin tone
- All Blue Mode: Anti-Blemish Lights
  - Helps reduce blemishes
  - Helps decongest skin that is blemish prone
  - Helps calm redness
- Red and Blue Mode: Both Anti-Aging Lights and Anti-Blemish Lights simultaneously

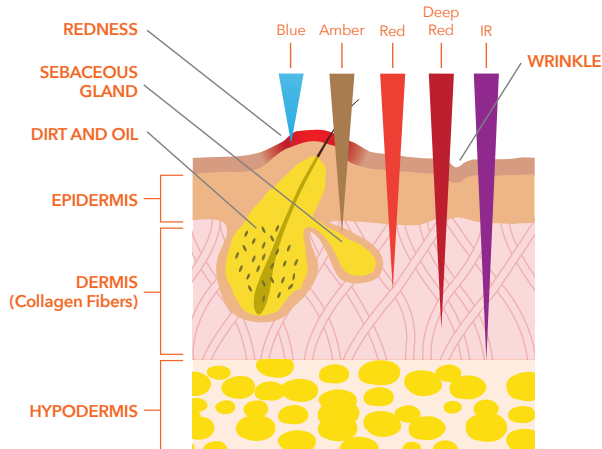**What to Expect**

Users typically see a reduction of blemishes within a few days and a clearing of spots as well as initial reduction of lines and wrinkles within 2 weeks of daily use. The best results for reduction of lines and wrinkles occur after a 10-week course of daily use.

**OPERATING INSTRUCTIONS**

**Set up:**

1. Unpack your SpectraLite™ FaceWare Pro.
2. Plug USB charging cord into port connector ① on device.
3. Plug USB charging cord into Universal USB Power Adapter (not included) and into outlet for charging (use only UL or other locally certified and compliant USB power adapter).
4. Charge your SpectraLite™ FaceWare Pro for 2-4 hours. Device is fully charged when Power Indicator Light ② is solid (not flashing).
5. Place device in protective storage bag when not in use. Do not store device in direct sunlight or a location with temperatures over 100°F/38°C as this may cause damage to plastic coating.
6. A fully charged device yields 24-30 minutes of use, approximately 8-10 sessions. When operating device multiple times, allow device to rest 5-10 minutes in between use.

**Instructions for Use:**

1. Wash and fully dry face prior to use.
2. Place device onto face and adjust straps to fit securely around head.
3. Press and hold button for 2 seconds to power on device (lights will turn on) and select device setting of your choice. SpectraLite™ FaceWare Pro will automatically shut off after completion of a 3-minute session. To turn off device prior to completion simply press and hold power button until lights turn off.
4. Remove device from face and wipe device clean (see cleaning instructions).
5. Place device into protective storage bag.
6. After completing a session, apply topical products as desired.

### Helpful Tips for Best Results:

- Use daily.
- Use device only on clean, dry skin.
- Do not use device after the application of products that contain SPF as this may block light penetration.

### Keep Skin Hydrated:

Immediately following session, apply a serum and/or moisturizer of your choice from Dr. Dennis Gross Skincare™.

### Maintenance:

1. It is important to keep device clean (see cleaning instructions).
2. This device has no serviceable parts. Battery is not accessible and not serviceable. User should not attempt to access or service the battery as this will damage the device and void the warranty.
3. Device includes only one button, which features the following:
  - a. Press and hold: Power on device in Red Light Mode
  - b. Click twice: Blue Light Mode
  - c. Click 3 times: Red Light and Blue Light Mode
  - d. Press and hold for 2 seconds to turn device off at any time.

All settings are set to the factory setting.

4. Charge device only with the supplied USB charging cord and a universal USB power adapter (not included).
5. Power indicator light will begin to flash when battery levels are low. This indicates it is appropriate to begin recharging the device.
6. Voltage of SpectraLite™ FaceWare Pro: 5V
7. Intensity of SpectraLite™ FaceWare Pro: Maximum Intensity Output:  
Blue Light: 40 mW/cm<sup>2</sup>, Red Light: 75 mW/cm<sup>2</sup>
8. Compatible Power Supply: 120V/240V

### Cleaning:

1. Turn device off and unplug from any power source.
2. Wipe the tops of the LEDs with an alcohol wipe or damp cloth to remove any residue from skin oils after each use. This is especially important if the device is being used by more than one individual.
3. Clean device body by using a damp or microfiber cloth. The microfiber storage bag can be used for cleaning and polishing the body of device. Do not use chemical or alcohol-based cleansers on exterior as harsh chemicals may discolor the metallic panels.
4. Allow device to dry completely prior to use.

### Troubleshooting:

1. System is not turning on:
  - a. Did you press and hold the on/off button for 2 seconds?
  - b. Is the device fully charged? (Power Indicator Light is illuminated and solid)
  - c. Is the Power Indicator Light illuminated?
  - d. If device still does not turn on, call Client Care at +1-888-830-SKIN or email an inquiry to customerservice@drdennissgross.com (Be prepared with device serial number on hand.)
2. System is not charging:
  - a. Is USB charging cord firmly inserted into device? (Does Power Indicator Light illuminate?)
  - b. Try a different USB Power Adapter to ensure it is working properly.
  - c. Try a different wall outlet to determine if there is power to the outlet.
  - d. If device still does not charge, call Client Care at +1-888-830-SKIN.
3. System shuts off abruptly, prior to completion of 3-minute session:
  - a. System is designed to shut off automatically after 3 minutes.
  - b. System may shut off after multiple, repeated uses if the device exceeds desired temperature. Allow device to rest for a minimum of 10 minutes and try again.
  - c. If device continues to shut off prior to 3 minutes, please contact Client Care at +1-888-830-SKIN.
4. Some lights appear turned off:
  - a. The infrared (IR) light spectrum is not visible to the human eye. These lights will look as if they are not illuminated, but they are operating properly. If the red and blue visible lights are on, then the device is operating as designed.

### LIMITED WARRANTY INFORMATION

We value you as a client and strive for the highest quality products to earn your satisfaction. Please use this as your reference for any customer service or warranty needs you may have when purchased within the United States. You must register your product within 30 days of purchase at [www.drdennissgross.com/registration](http://www.drdennissgross.com/registration) for the warranty to be effective. When purchased outside the United States, please contact your place of purchase for local warranty information. Our warranty is limited to the terms set out below:

***What is covered by the limited warranty?***

This 1-year limited warranty covers defects in materials and workmanship in your SpectraLite™ FaceWare Pro when purchased through an authorized partner or directly from Dr. Dennis Gross Skincare™. Your 1-year limited warranty is valid for you alone, the original purchaser, and cannot be transferred. If the product was a gift, obtain a copy of the original receipt and send to registration@drdennissgross.com, or have the original purchaser process the warranty claim on your behalf.

***What is not covered by this limited warranty?***

This limited warranty does not cover: damage, defects, failure caused by abuse or misuse, improper care, impact or insertion of foreign objects, use with a voltage connector, or loss of device.

***How long does this limited warranty last?***

The limited warranty begins after the device has been registered and expires after one year. The warranty period is not extended if your device is replaced by Dr. Dennis Gross Skincare™ during the warranty period.

***What will Dr. Dennis Gross Skincare™ do if your device is defective?***

If you discover a defect in your SpectraLite™ FaceWare Pro covered by the warranty before the 1-year warranty period expires, contact Client Care at +1-888-830-SKIN for assistance. If we confirm a defect covered by this limited warranty, Dr. Dennis Gross Skincare™ will decide whether to replace your SpectraLite™ FaceWare Pro with either a new or refurbished device. If we find that your SpectraLite™ FaceWare Pro is not defective, we will notify you and return the device to you. All Limited Warranty services must be authorized and performed by Dr. Dennis Gross Skincare™ or a service provider authorized by Dr. Dennis Gross Skincare™.

## **WARRANTY LIMITATIONS & DISCLAIMERS**

This warranty gives you specific legal rights and is in addition to your statutory rights. Your remedies for the breach of this limited warranty are limited solely to the repair and replacement of your device as set forth in this warranty statement.

To the maximum extent permitted by applicable law, Dr. Dennis Gross Skincare™ disclaims all implied warranties. If applicable law does not

permit the disclaimer of implied warranties, then the duration of implied warranties are hereby expressly limited to the duration of the express limited warranty set forth above.

In no event shall Dr. Dennis Gross Skincare™ be liable for any special, incidental, consequential, punitive, or indirect damages based upon breach of warranty, breach of contract, tort (including negligence or strict product liability), or otherwise, even if Dr. Dennis Gross Skincare™ is advised or aware that such damages may arise. Such damages include, but are not limited to: loss of SpectraLite™ FaceWare Pro, cost of any substitute equipment, the claims of third parties, or damage to property.

If the applicable law does not allow limits on warranties, or on remedies for breach in certain transactions, the limits in this paragraph and the disclaimer or exceptions from warranty may not apply. This limitation does not apply to claims for death or personal injury caused by a defect in the device or any other liability that cannot be excluded or limited by law.

These warranty terms, conditions, and limitations allocate the risks of product failure between Dr. Dennis Gross Skincare™ and you, the purchaser of SpectraLite™ FaceWare Pro. This allocation is recognized by both parties and is reflected in the price of SpectraLite™ FaceWare Pro. By using SpectraLite™ FaceWare Pro you acknowledge that you have read these terms, conditions, and limitations, understand them, and are bound by them.

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Dr. Dennis Gross Skincare,™ LLC.