

Product KAHUNA FOOT MASSAGER

Model FLM SOLE MATE

Distributed by AJX Distribution INC

11125 Knott Ave, Suite D

Cypress, CA 90630

Tel 1-855-646-9500

Email INFO@KAHUNACHAIR.COM

You must register warranty at www.kahunachair.com/services/

INSTRUCTION MANUAL
MODEL : FLM SOLE MATE

Please read this user manual thoroughly
prior to operating the massage chair.
Store this booklet for future reference.

KAHUNA  CHAIR®

KAHUNA FOOT MASSAGER SOLE MATE MANUAL



Content

Safety instruction	1
Product structure	4
Product feature	5
Preparation before use	6
Product usage	8
Remove cloth cover	9
Massage sketch chart	10
Product maintenance and storage	11
Troubleshooting and specification	12

Safety instruction

Thanks a lot for choosing our products. Please read and comply with the 'safety instruction' carefully: these attentions will help to use the products and to avoid risk. According to the dangerous level divide to 'warning' and 'caution'. Please keep the instruction book properly for future using.

Warning

Do not use the non-specified voltage source	Will cause the electric shock or fire
When power failure	Turn off the power switch and take off the plug to extend the products life.
When out of control	Take off the plug and contact the dealer
When not in use	Power off and disconnect the plug from socket
With breakage on power cord or attaching plug	Contact the manufacturer or the repairman to replace in case of any risk.
When in use	Do not touch the massage device cloth with hands.

Attention

It's better to use it within 30 minutes, any longer time may cause the muscular stiffness or neuropathic pain.

1. Please consult the doctor with these cases:
 - Use the cardiac pacemaker or implantable electronic medical device.
 - Suffer from heart disease.
 - Suffer from strain disease.
 - Skeletal fragility.
 - After surgery.
 - Pregnancy or shortly after childbirth.

Safety instruction


Person with osteoporosis.
Person with blood circulation problems or thrombosis.
Person with skin infection.
Person receiving treatment or assisted treatment
Person with abnormalities in the limbs.
Person with foot injury caused by accident or disease.

2. Do not use this product when there is damage on its shell or mechanic electronic failure. You should contact with local distributor or manufacturer.
3. Do not use it in wet or high temperature place or near the exothermic item such as heater or oven.
4. Do not paint medical or make-up cream on this product before using or during using it.
5. Please send the product to authorized service center for repair. Do not fix it by yourself. Or the warranty won't be effective.
6. Do not vibrate or throw it heavily.
7. Do not use it within 1 hour before or after meal.
8. Do not wind the cords on the product or bend the cords in case of cord damage.
9. Do not stand on the product or put heavy thing on it.
10. If any failure occurs, unplug the cord and contact with the manufacturer immediately.

Safety instruction

11. Electrical appliances are dangerous when close to the water. Do not use this product in damp environment, for example, do not use in the bathroom or near the faucet or swimming pool.
12. If the power cord is damaged, it must be replaced by the manufacturer or its service or personnel with relevant qualifications to avoid danger.
13. Do not damage the power cord, do not bend, fold, tear open outfit, pull, tie, bending the power cord, please do not place the power cord under the heavy objects, when move the product, do not clip the power cord in the middle of the object, using a damaged power cord may cause fire or get an electric shock, do not fall asleep when use it, do not use this product after drunk, do not use this product during a high fever.
14. The product can only be placed according to the instructions in this manual.
15. This product is not used as a medical device
16. When the product is powering up, it needs to be watched.
17. Do not use in gasoline, thinner, alcohol or bleach and other high corrosive ingredients environments, so as not to cause product discoloration, deformation, cracks, or faults.
18. This product is not intended for use by people (including children) who have a weak body, sensory or neurological deficiencies or lack of relevant experience and knowledge unless supervised and instructed by the person in charge of their safety about their use of the product.
19. Do not touch this product directly with your skin.
20. Do not let children play this product.
21. This product has a heating element, People who are insensitive to fever should pay attention
22. This product is for one person at a time
23. If you feel unwell during massage, stop using it immediately.
24. Do not insert or drop items into any openings in the product.

Product structure chart

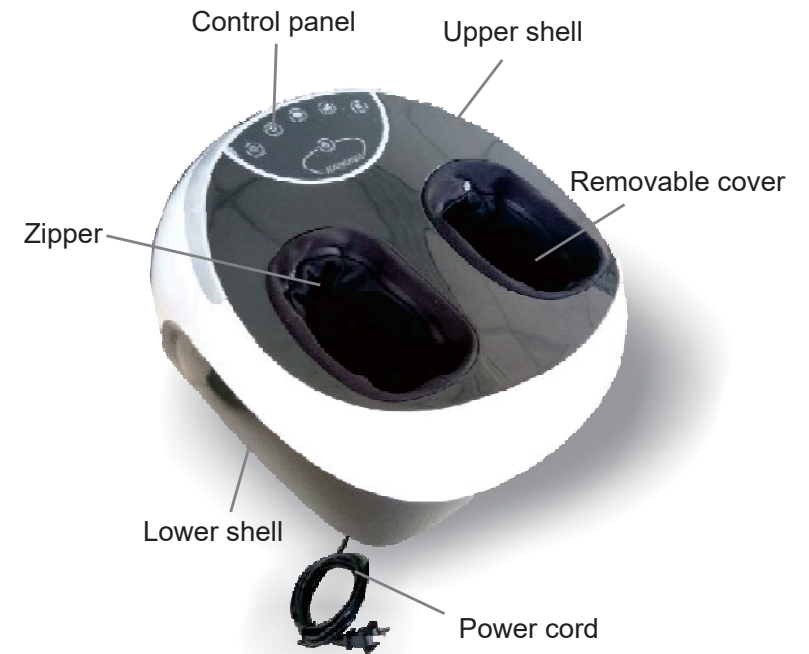

Product Feature

★ New concept, Innovation design.

New appearance design, Fluidic Sculpture design, fashion and portable.

★ Full range massage

With foot roller, scraping, vibration, air massage and heating function, It can massage the feet overall, deep acupuncture point massage, dredging meridians for whole body.

★ Heating function

Foot parts with heating function, warm massage for foot arch, wrapped carbon fiber hot compress massage for instep. It is like professional heating massage, heat therapy can eliminate fatigue, dredging meridians and promoting blood circulation.

★ Convenient for storage

Equipped with power cord collect device, it is easy for collection.

★ Removable cover

Fabric cover is removable, washable, and changeable.

★ Safety

The foot massager can deflate automatically if power off suddenly during massage.


★ Product function

Relieve foot pressure and leg pressure, improve muscle flexibility, promote blood circulation, relieve muscle pain after sport, stimulate main acupuncture point, improve body health, and relieve fatigue.

Preparation before use

1.Regarding place


Please put the machine in dry place, not exposure to sunlight long time. In order to reduce working noise, please put the machine on a carpet.

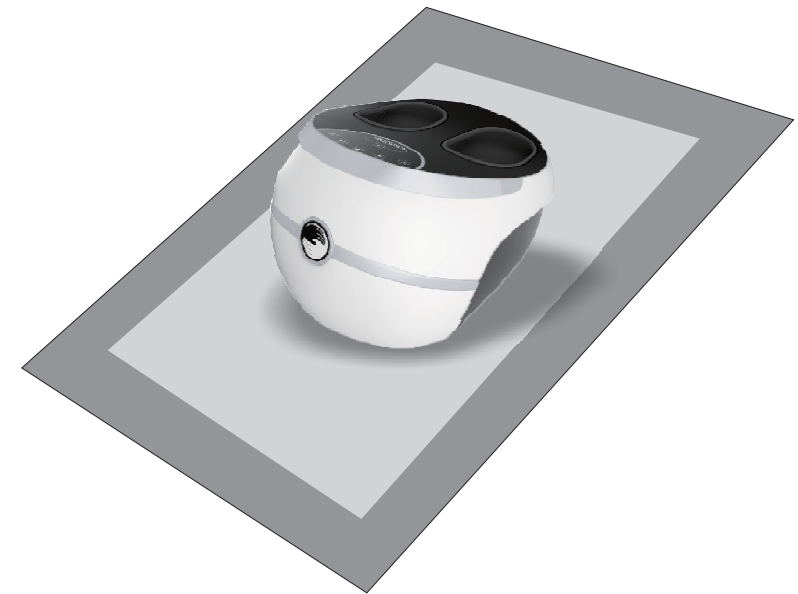


Preparation before use

2、 Method to plug in

Please insert the plug into socket

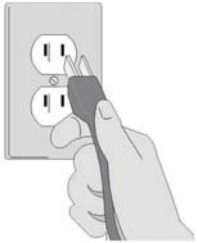

Make sure the power cord is not broken


 Forbidden to plug in or unplug with wet hand.







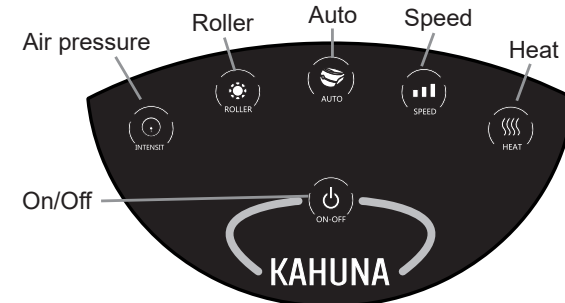

Usage

1. Please insert the plug into socket.
2. In power on situation, On/off button light will in flash status.
3. Press on/off button, start the machine.
4. Choose program for massage.

Control panel usage



1.  On/Off: Turn on the power, Press this key, the indicator on the control panel is twinkle, press it again, this massager will start, the switch indicator is always light, and goes into auto massage program 1, Air pressure massage is in the middle range, heating function is on, entering a 15-minute countdown. (When all massage functions are turned off, Countdown stop. When you start the massage again, the 15-minute countdown begins. After the 15-minute countdown, it will be turned off automatically); Press this key a third time, this machine is turned off, the switch indicator is twinkle.
 2.  Air pressure: On Auto mode, Press this key to adjust the auto massage intensity of air pressure, Low-Mid-high for your choice. On Manual mode, press this key to choose manual air pressure intensity Low-mid-High. The indicator light of relative strength (low - blue, medium - green and high - red) is bright.
 3.  Roller: In "On" Status, Press this key, Vibration function is on, the roller run clockwise, and display blue light; Press this key again, display Green light, Vibration is off, the roller run anti-clockwise; Press this key a third time, the light is off, and roller stop working.
 4.  Auto: In "On" Status, press this key, Switching between automatic mode 1 (the indicator is blue) and mode 2 (the indicator is green). During Auto mode, Roller on the sole, Vibration part, Air pressure massage is following the default program. The airbag automatic massage intensity default is mid level.
 5.  On manual mode, When Roller or Vibration function is chosen, The speed key is effective; Press this key, The speed is switching between mid (The Blue and Green indicator are bright), high level (the Blue, Green, Red indicator are bright) and Low speed (The Blue indicator is bright).
 6.  Heating: In "On" Status, Press this key, the switch indicator is light, Starting heating function. Press this key again, the indicator is off, close the heating function. Cycle by turn.
- Note: every click of the button, the window indicator will have a light indicator to indicate the success of the operation.

Remove the cloth cover

The way to remove the cloth cover

1. Open the zipper.
2. Remove the cloth cover.




The way to install the cloth cover

1. Paste the Velcro inside of the cloth cover.
2. Zip up.
3. Tidy up.



Massage sketch chart

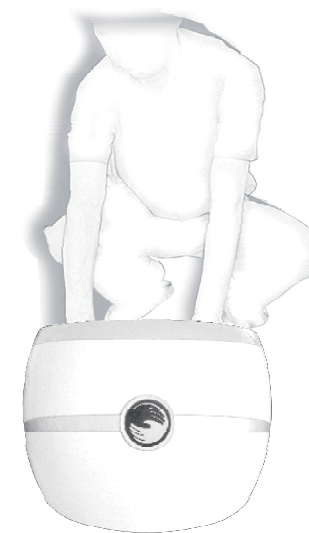

Don't put your hands in the fabric covering.

Product maintenance and storage

1. Maintenance method

Warning

When cleaning, please turn off the power switch and pull the plug out of the socket.

After using the machine, if there is stain on the machine, please wipe by the wet cloth with a neutral detergent,(Do not let chemicals such as benzoline, alcohol or diluent drop on the massager), then wipe with a dry cloth, and let it dry naturally.(Do not dry in rapid drying method)

Warm tips

The machine can achieve the ideal massage effect for about 15 minutes every time.

2. Place of storage

Please keep it in the following conditions:

Temperature : 0℃-40℃

Humidity : (No condensation) 30%RH-85%RH

Please be careful not keep it in a high temperature, humid and direct sunlight place, also don't keep it in a place with corrosive gases.

First, clean the product, then put it in a dry and clean place, do not use the product in a very dusty, greasy, humid or smoky environment. Please do not put the product too close to the heating area, do not put it in the direct sunlight area, do not use the product outdoors.

Troubleshooting and Specification

1. Troubleshooting

This product is microcomputer controlled, sensitive to voltage fluctuation. If the voltage fluctuation exceeds 10% or -15% of the limit voltage, the product may appear the following abnormal situation.

Possible trouble	Troubleshooting
The massager can't work after power is on.	1. Check the switch and wire is intact or not. 2. Check if the plug and socket connection is loose or not.
The massager can't work after power is on.	After long time work, the temperature of the product is too high, the product will turn on temperature safety status, please turn off the product for 30-50 mins at least to let the product cool down.
When working, the controller can't work and the massager stop working.	1. On timing status, if time is off, press the "on" again to restart working. 2. turn off the power and turn on again.
If the troubles are not caused by above reasons, please contact with manufacturer or local sales company.	

2. Specification

Vibrating foot massager	
Model: Sole Mate	N/W: 5. 6KG
Input Voltage: 110V	Frequency: 60Hz
Power: 60W	Time: 15mins
Packing size: 402X277X460mm	Safety class: <input type="checkbox"/>

We reserve the right for design modification. It is subject to change without notice. Please refer to the actual product.