


VIQURE

Blackhead Remover

Model: FC3001

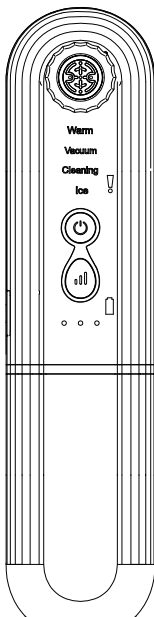

User Manual

Blackhead Remover User Manual

Thank you for choosing Blackhead Remover! Based on the principle of negative vacuum pressure, the product is designed to clean pores by absorbing deep-seated facial dirt such as blackheads and pimples through an adsorption flow. It features four modes: "Warm" mode opens pores and softens blemishes and acne; "Vacuum" mode draws dirt and oil away from the skin; "Cleaning" mode effectively cleanses the skin; and "Ice" mode shrinks pores. The separated water circulation makes it easy for you to have clean, moisturized and healthy skin by replenishing the skin with moisture and shrinking the pores while cleaning them safely and effectively. On top of that, this device also has an automatic cleaning function.

Scan the QR code below to download the user manual and detailed video tutorials!



Contents

1. Safety precautions	1
2. Overview	3
3. Get Started	6
4. Trouble shooting	12
5. Cleaning and Maintenance	15
6. Product Parameters	17
7. Configuration List	18
8. After-sales Service	19


1. Safety precautions

- This is not a medical product and is only intended for daily skin cleansing. It cannot be used for the treatment of any skin conditions. Please consult your doctor for any skin health concerns.
- Do not use this device after cosmetic surgery or rhinoplasty.
- Do not use it on areas with scabies, localized inflammation and swelling.
- Do not use this device on scars or delicate skin, such as the chin, corners of the mouth, or eyes.
- If you have sensitive skin or any skin problems, please consult your dermatologist before using this product.
- The product works with negative pressure. It is recommended to dab or gently lift and slide the suction head without staying in the same place for more than 3-5s, otherwise it will aggravate the red marks.
- It is normal to get red marks during use (see the Troubleshooting section of this manual on how to

minimize red marks), which will fade away in 30 minutes. You can speed up the process by applying ice to the area after use, or use a soothing mask if necessary.


- If serious discomfort occurs during or after use, stop using it immediately and seek medical advice.
- The device has been checked twice before delivery, a small amount of water droplets or water mist is normal, please feel free to use!
- Do not disassemble, repair or modify this device to avoid electric shock or other safety hazards.
- Do not use the device while charging.
- This product contains batteries, please keep it away from fire to avoid any explosions that could harm you or others.
- The device contains electronic components, please do not soak it in water and place it in a dry place after use.
- The health and safety of your child is of paramount importance, please ensure that the device and accessories are placed out of children's reach.

2. Overview

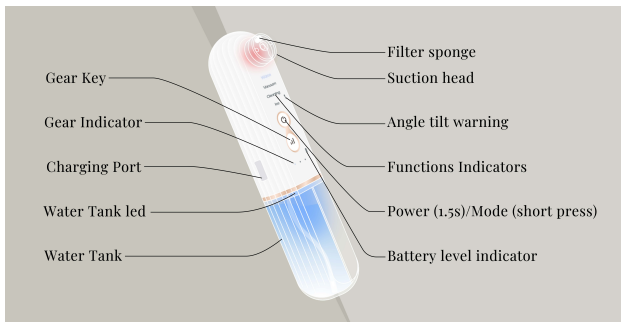

*Power/ Model: press and hold for 1.5s to turn on the device;
short press to switch modes.

*Gear Key: short press to adjust level.

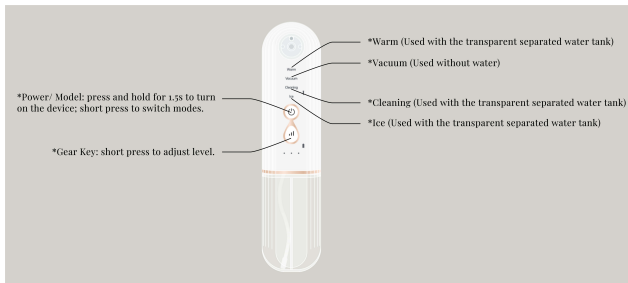

*Warm (Use with the transparent separated water tank)

This mode awakens the skin with red light and fully opens up your facial pores with a thermal treatment that softens blackheads, acne, oil and makeup residue, making it easier for later suction.

*Vacuum (without water)

This mode utilizes vacuum suction technology to thoroughly cleanse the skin by removing softened pimples, acne, oil and makeup residue.


*Cleaning (Use with the transparent separated water tank)

This mode makes use of water to take away the dirt

from the pores while moisturizing the skin.

*Ice (Use with the transparent separated water tank)


This mode delivers a blue light that calms the skin, shrinks pores, and tightens the skin.

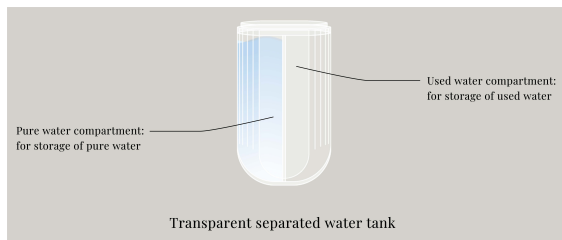

Transparent separated water tank

Pure water compartment: for storage of pure water

Used water compartment: for storage of used water

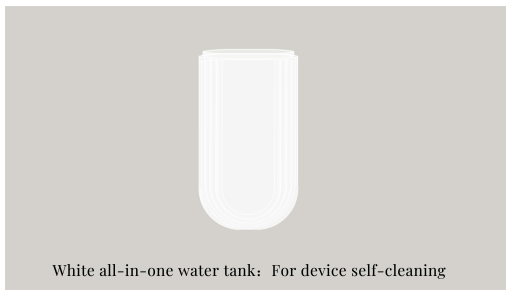


White all-in-one water tank




For device self-cleaning

3. Get Started

Step1 Installation

Rotate clockwise to remove the water tank and fill the left tank with clean water when the device is off. Insert the sponge into the water outlet of the device and the water inlet in the water tank, then rotate counterclockwise to re-install the water tank. (Installed by default)

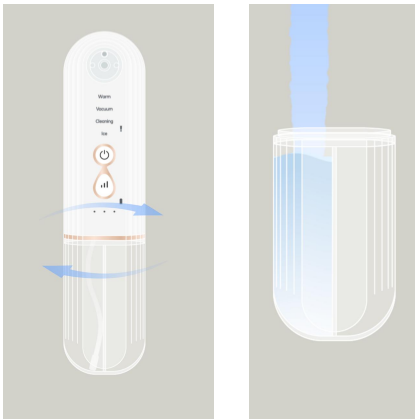


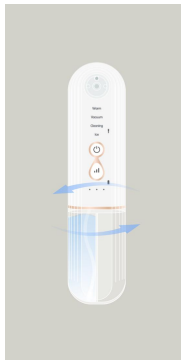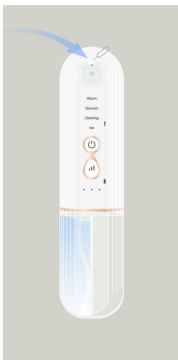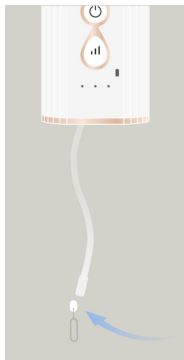


Select suction head


*Round: for cheeks and forehead

*Oval: for nose and nasal folds

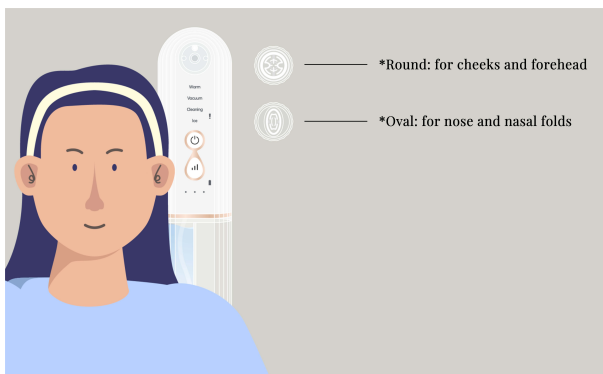


Step 2 Start cleaning

Press and hold  for 1.5s to turn on the device.

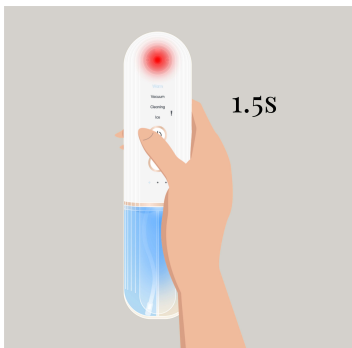

Hold the device vertically in hand during use (the device will stop working and start beeping if it tilts over 40 degrees). Press the suction head against the skin of the cleansing area (otherwise no water will be released) and clean the skin by dabbing, gently lifting and sliding, or moving in a circle without staying in the same position for more than 3-5 seconds. You can short press the "Power / Mode" button to switch modes and short press the "Gear key" to adjust the suction power.

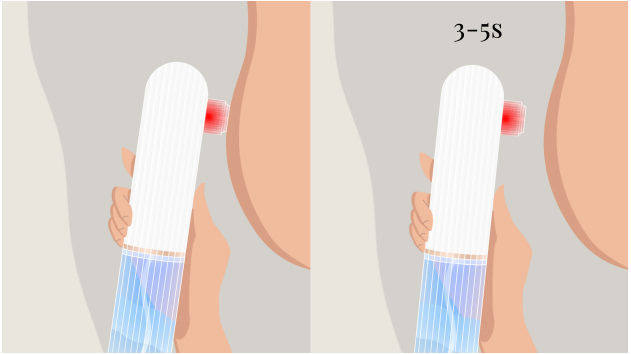

Suggested full use process

Warm (3mins) + Vacuum (3mins) + Cleaning (4mins) +
Ice (3mins) = 1 Clean Cycle


You can also customize the mode and duration of use
to your needs.

Note:


1. If you feel that the warm mode of this device can not
fully open the pores, you can use a hot towel for 5-8
minutes (especially in cold weather). This step can
effectively open the pores so that dirt can be sucked
out more easily and red marks can be reduced.

2. Please replace the fresh water in time during use, and turn off the device when doing so.

Step 3 Complete cleaning


Press and hold  to turn off the device.

Wash your skin with water. It is normal to get red marks , which will fade away in 30 minutes. You can speed up the process by applying ice to the area after use, or use a soothing mask if necessary.

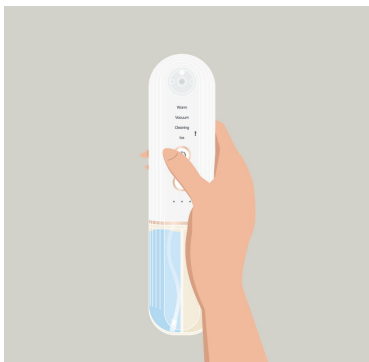

Step 4 Device self-cleaning

Install the all-in-one water tank, remove the filter sponge, and install the cleaning cap. Then switch the device to the cleaning mode and leave it through the process.

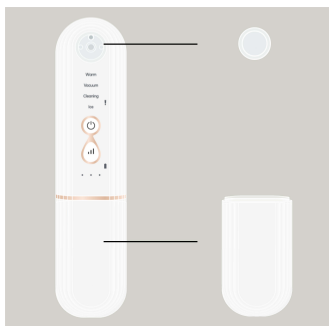

4. Trouble shooting

4.1 How to reduce red marks

1) Apply a warm compress before use to fully open pores. If you feel that the warm mode of this device can not fully open the pores, you can use a hot towel for 5-8 minutes (especially in cold weather). This step can effectively open the pores so that dirt can be sucked out more easily and red marks can be reduced.

2) First-timers are advised to try from the smallest suction power, and then increase it after the skin gets used to it to find the most suitable state.

3) Please dab or gently lift and slide suction heads without dragging hard or staying in the same place for more than 5s, otherwise it will aggravate the red marks or even hurt your skin.

4) We recommend you follow the process of Warm-Vacuum-Cleaning-Ice, choosing a lower level for Warm and Ice, and a proper level for Vacuum and Cleaning. Steps can also be reduced without

compromising the effect if you have delicate skin, for example, use a hot towel for 5-8 minutes before use to fully open the pores and then use the Warm mode for cleaning.

5) Wash your skin with water after use. You can speed up the reduction of red marks by applying ice or using a soothing mask if necessary.

4.2 Why the device does not work

1) The device can not be turned on due to low battery. Please try later after charging for 3-4 hours (The device is charging when the battery level symbol blinks). If charging fails, it may be resulted from incompatible charging devices, please use 5V-1A charging head and charging cable.

2) Press and hold the power button for over 1.5s before releasing it.

3) The tilt angel of the device is greater than 40 degrees , which triggers the protection mechanism featuring blinking light and beeping sound. Please

keep the device at less than 40 degrees.

4.3 Why the device does not release water

1) The device is designed based on negative pressure principle, which means that the water can only be pumped up when the suction head is close enough to the skin that creates negative pressure, otherwise there will be no water. (The water is from the inner circulation and only exists inside the suction head without jetting or flowing out.)


2) The sealing ring at the connection of the suction head is off, which causes leakage. Please ensure the sealing ring is installed correctly.

3) Please make sure that there is enough water in the water tank and the end of the tube is completely submerged in the water.

4) No water for Vacuum mode.


5. Cleaning and Maintenance

Suction heads

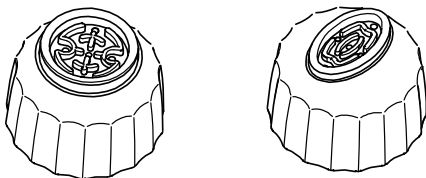

Clean the suction heads with a brush after use and store them properly. (Please change every two months)

Filter sponge

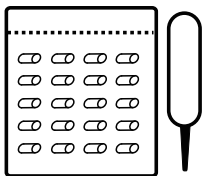

Rinse or dispose of the sponge after each use (1 week maximum use).

Water tank

Brush the water tank with a detergent (PH 6-8) after use.

Device











Having finished the self-cleaning, please wipe the device with a paper towel or wait for it to air dry. For longer service life, it is recommended not to leave the device in a humid environment as it contains delicate electronic components.

Do not expose the device to prolonged sunlight to prevent deterioration.

6. Product Parameters

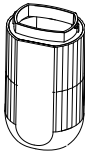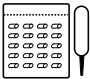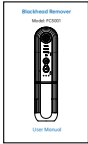Name:	Blackhead Remover	Model:	FC3001
Rated voltage:	5V-1A	Rated power:	8 W
Battery capacity:	1500 mAh	Function modes:	4
Product dimensions:	50*40*205 mm	Net Weight:	≈275 g
Execution standard:	GB4706.1-2005,GB4706.15-2008		

7. Configuration List
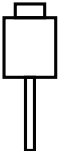
<p>Main Body * 1</p>	<p>All-in-one tank * 1</p>	<p>Separate tank * 1</p>	<p>Filter sponge Pin * 1</p>	<p>User Manual * 1</p>
				
<p>Cleaning Cap * 1</p>	<p>Round head * 1</p>	<p>Oval head * 1</p>	<p>Brush * 1</p>	<p>Charging cable * 1</p>
				

8. After-sales Service

If you have any questions, please contact us via the link or QR code below, where you can also download the user manual and detailed video tutorials.

<https://viqure.com/pages/user-faq-viqure-blackhead-remover-vacuum>



Thank you for your purchase and concern, we wish you a perfect experience!