



RF Beauty Instrument

User manual



Model: RF02

► Preface

Thank you for choosing radio frequency beauty instrument, the RF radio frequency beauty instrument converts electrical energy into heat energy , and deep into the dermis to increase the temperature, accelerate the blood circulation of dermis and subcutaneous tissue, to shrink and tighten the skin collagen , stimulate and reshape the collagen to achieve multiple cosmetic effect.

Before use it, please carefully check it and keep it properly.

In order to better use RF beauty instrument (Model: RF02) to ensure good cosmetic effect, please find below our suggestions:

To provide and ensure you more convenient after-service, when you receive RF02,you can load on the web of center of China, Input goods code to check information authenticity verification.

Contents

Contents	01-02
Product information	03-04
instrument introduction	03-06
Instrument Function Introduction	07-08
Instrument operating procedure	07-08
Mode selection	09-10
level explain	09-10
Mode declaration	09-10
Instument using step	11-13
Contraindications and safety attentions	13-14
Operation and maintenance	15-16
Contact us	17-18
Warranty	18-20

1.Product information:

- Keep away from children when using this product.
- Special caution should be taken when using this product near pacemakers, hearing AIDS, and other medical electronic devices as it may interfere with these devices.
- Non-original accessories and accessories do not have warranty qualification.
- Please do not disassemble this product by yourself. If your product has problem, please contact our designated repair department.
- When using this product, please do not throw, put, knock on the product, so as not to damage the shell and internal components.

Introductions:

Warning: Our company will not be responsible for the consequences caused by non-compliance with the above suggestions or improper use. If the content in this manual is inconsistent with your product, please refer to the product.

Our company reserves the right to amend the specification at any time without prior notice.

* Please follow the correct operation procedure to use the instrument

2.Instrument is introduced:

List: Host *1, handle probe *2, adapter *1, User manual *1



Parameters:

Product: **RF Beauty device**

Model: **RF02**

Technological type: **Radio Frequency current** rated frequency: **50/60Hz**

Rated power: **36W**


Input Voltage: **100/240V**

Instrument structure:

Front:



Display Screen:

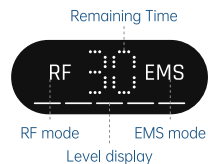

Reverse:



Handle probe:



3. Instrument Function Introduction:

- **Range of application:** Face and body
- **Applied environment:** family, beauty salon, fitness club and other occasions, but also convenient to travel to use.
- **Suitable for the crowd:** Loose skin, aging, aging fibrous tissue, lack of elasticity, dull, rough, dull skin.


Recommended usage cycle:

Operate once every 3-4 days;

For mature skin such as rough, dry and aging skin, it can be used once every 2-3 days;

For better skin substrate skin can do daily maintenance, It can be used every 1-2 weeks

4. Instrument operating procedure:

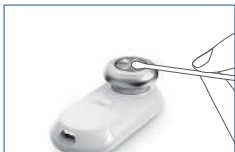

1. Please remove makeup and clean skin before use.
Clean and disinfect handle probe.



2. Connect the power cord and power it up.
The white light comes on when the electricity is switched on



3. Connect the handle probe
Note: Handle probe is not allowed to connect other Type C



4. Long press the power button to power on. When you hear a beep, it means that the RF mode is working. At this time, the display lights up and the default RF mode is 30 minutes.



5. Press the level switch button to switch level;
Choose the right level (it is recommended to try gradually from low level)
Note: if you feel overheated during nursing, please reduce the level immediately



6. Press mode switch button to switch mode;
Red light: RF mode
Blue: EMS mode
Green: RF+EMS mode

5.Mode selection:

RF Mode:

When using RF mode, the red light of the device is on, the display will show the word RF, the countdown is 30 minutes, and the default level is level 1.

EMS mode:

When using EMS mode, the blue light of the instrument is on, the display will show the word EMS, the countdown is 30 minutes, and the default level is level 1.

RF+EMS mode:

When using RF+EMS mode, the green light of the device is on, the display will show RF EMS, the countdown is 30 minutes, and the default level is level 1.

Note: after 30 minutes of countdown, the instrument will turn on white light and enter standby mode. Long press the power button when restarting

6.level explain:

- There are five level of the instrument, please choose the most comfortable one to use
- Default level is 1 after turned on
- Gear adjustment (as shown in figure)

— — — — — — — — — — — — — — —
1 level 2 level 3 level 4 level 5 level

- For the first time, please use the first gear until you get used to it, so that your skin can run in with the instrument

7.Mode declaration:

RF:

Display: RF

Features: warm feeling directly to the bottom of the skin, warm the skin deep and let the nutrients penetrate into the skin cuticle, deeply nourish the skin, make the skin moist and shiny

Function	RF	Using frequency	2-3 times every week
Time	30 minutes	Recommended cosmetic	Cool nourishing repair essence
Level	5 levels		

EMS:

Display: EMS

Myoelectric stimulation exercises facial muscles through electric wave stimulation, contracts and relaxes muscles, and achieves lifting and thinning of the face, reducing Wrinkle reduction effect

Function	EMS	Using frequency	2-3 times every week
Time	30 minutes	Recommended cosmetic	Cool nourishing repair essence
Level	5 levels		

RF+EMS:

Display: RF+EMS

Characteristics: while warming the deep skin, electromyoelectric stimulation also exercises the muscles, so as to achieve the effect of nourishing the skin, contracting and relaxing the muscles


Function	RF+EMS	Using frequency	2-3 times every week
Time	30 minutes	Recommended cosmetic	Cool nourishing repair essence
Level	5 levels		

Note: There are two electrodes on the plane of a single probe, and it will be effective if the entire area is in contact with the skin.

8.Using the step:

Facial manipulation:

During the operation, comfortable warmth is good experience feeling, the cool nourishing and repairing essence is easily absorbed by the skin. Please pay attention to supplement to keep the probe operating smoothly.

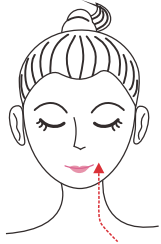

① From the upper collarbone to the lower mandible in circular motion(avoid throat)



② From the lower mandible to the preauricular in circular motion



③ From the mouth to the preauricular in circular motion



④ From the nasolabial folds to the cheekbone in circular motion(avoid temple)



⑤ From the wing of nose to the cheekbone in circular motion (avoid temple)



⑥ Sliding from the inner eye corner to the outer corner and enhances eye and eyebrow (avoid eye and upper eyelid)



⑦ From lower eyelid lifting to the corner of eye and eyebrow, hairline.




⑧ From the eyebrows (eyebrows, eyebrow end) to the mane

Tips:

1. Move at a speed of 2cm / sec.
2. If you feel too hot, please speed up a little
3. The probe should be close to the skin and move on the skin when use. Do not massage the eyelids or around the eyes with probe.

Body parts operation techniques:

The operation methods of body parts are consistent. Taking the arm as an example, the arm is divided into two parts (a) and (b). Start from the left outer arm, circle up, and then gently pull down to the origin, you can repeat the operation 3-4 times according to the time. The inner arm is operated in the same way.

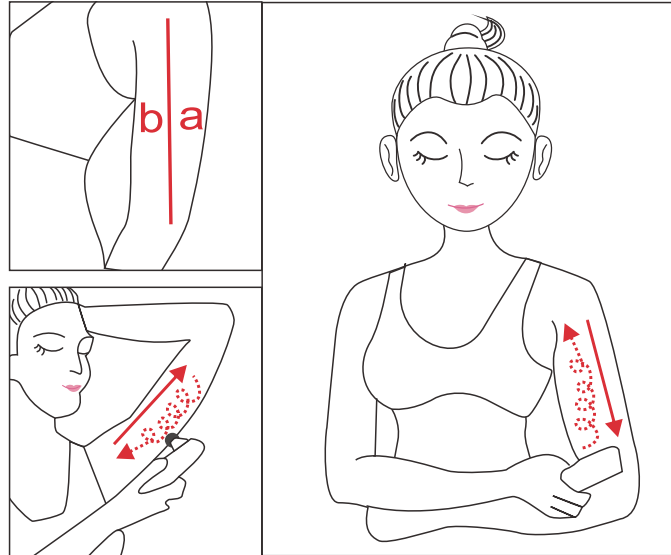

Note: all other parts can be operated by partial spiral winding.

Note: It should be used on the skin. Do not stay on the same part for a long time to avoid scalding

9. Contraindications and safety attentions

- ⊘ Banning the use for minors.
- ⊘ Pregnant women , woman's menstrual period is prohibited.
- ⊘ Heart disease, hypertension, coronary heart disease , diabetes, skin cancer, thyroid, malignancy, current sensitive group are prohibited.
- ⊘ Infections of the skin , skin burns , anti- clotting disorders , skin bleeding, or skin allergies are prohibited.
- ⊘ Don't use on the area with surgery, scars.
- ⊘ Dermatology parts , eczema, and other sensitive parts of the cornea are prohibited.
- ⊘ Skin botox and hyaluronic acid after three months and gold implantation after six months.
- ⊘ Body with hearing aids , heart pacemakers and heart machine are prohibited.

10. Operation and maintenance:

⚠ The probes of this instrument are not allowed to connect Type-C externally.

⚠ When replacing probe, please power off the instrument. can not power on the instrument until probe is well installed.

Please keep the probe clean before& after use it each time, use a soft clean cloth to wipe it, water cleaning is forbidden.

When operating, you should try to select from the lower energy ⚠ to high energy , until you feel comfortable and afford to bear.

⚠ Before use it , you should smear RF skin rejuvenation gel cream , do not use on the dry skin surface. or it will burn your skin.

⚠ During operation , the gel cream will be easily absorbed by the skin , keep wiping cream to make probe operate smoothly.

⚠ Prohibited operation on the throat , eyes , ears , eyelids , temples and other area.

⚠ Please use face probe for facial operation, please change for the body probe for body operation.

⚠ When finish operation , please avoid to putting it in the sun, and should make good protection from the sun, far away from strong and ultraviolet light.

⚠ Do not put the instrument and accessories in the liquid or wet environment.

⚠ Stored in ventilation and drying environment and away from children.

⚠ If the cases, accessories, probes are flooded accidentally, please return them to factory to check if they can be repaired for re-using, if the Instrument is flooded when connect to the power, You should cut off the power first before touching it , in order to ensure your safety.

⚠ Do not disassemble or repair the equipment by yourself.

⚠ If the equipment can not work properly , please return it with accessories to our factory, We have professional maintainer for repairing.

11.Contact us

Shenzhen Viatom Technology Co, Ltd.

●4E, Building 3, Tingwei Industrial Park, No. 6 Liufang Road, Block 67, Xin'an Street, Baoan District, Shenzhen, 518101, Guangdong, China

●E-mail: support@viqure.com

●Website: www.viqure.com

12.Warranty:

1. Warranty:

Machine warranty: from the consumer to buy the day, the machine warranty period for 1 year.

2. Warranty responsibilities:

If the instrument in the warranty period of product quality problems, in the detection of non-human factors caused by the warranty will be free; Failure caused by other abnormal use or force majeure is not covered by warranty.

3. Warranty specifications:

When users apply for free warranty service during the warranty period, they need to present valid purchase certificate (such as warranty letter).

If the product is damaged during transportation due to improper packaging when you return the instrument, no warranty or replacement will be provided free of charge.

If the following points belong to the company will not be responsible, not within the scope of warranty.

1. Failure and damage caused by fire, earthquake, third party action, customer's intentional or negligent use, misuse, use under other abnormal conditions, unpredictable events, etc.
2. Incidental damage to the action or use of the Goods (loss of business benefits, interruption of business, including loss of use of instruments, costs incurred in replacing any equipment, fixtures or services, or claims by your customer).
3. Damage caused by using methods other than the contents of this manual.
4. Circumstances beyond our control: misuse, abuse, failure to operate and maintain in accordance with our instructions, or use of accessories and accessories not sold by our company, and unauthorized disassembly and assembly of instruments, etc., resulting in problems.
5. Failure or damage caused by use in conditions other than those specified by the product, including overload.
6. Due to improper storage by users (such as rodent, liquid infiltration, etc.).

Warranty card

Product name		Type	
Custom name		Purchase date	
Shop for purchase		Contact number	
Recipient address			
Maintenance records			
Broken record			

※ Reminder: Please keep this card properly, fill in this card and send it back to the manufacturer together with the product when repairing.