
zenbed

Safety Instructions

1. Please turn off the power when you are not using this product.
2. For your safety, please turn off the power before going to bed.
3. Please avoid sticking to the skin at the same position for more than 8 minutes under high temperature and confirm the skin condition. If there is any abnormality, remove it immediately.
4. It is forbidden to be close to fire sources, corrosive liquids, and various unfavourable environments.
5. The belt cannot be washed in washing machine, clean with damp cloth.
6. Do not use when wet.
7. Do not use the product when charging.

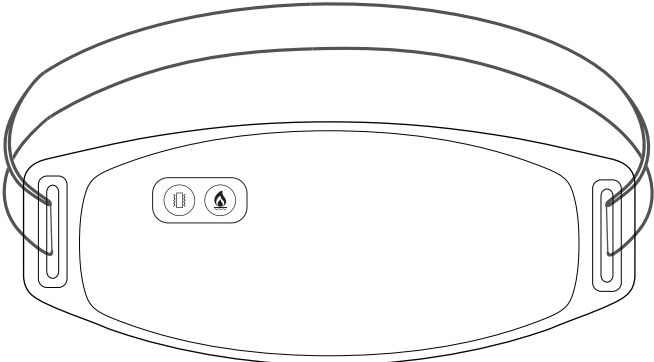
GRAPHENE HEATING BELT USER MANUAL

Please do not use if:

1. Moderate and severe cardiovascular disease, severe heart disease.
2. Skin sensory disturbance or abnormal skin.
3. Contraindicated in epilepsy patients or those who have not fully recovered from gout.
4. Pregnant women are prohibited.
5. Those with unhealed or unhealed wounds after surgery are prohibited.
6. Persons who cannot take care of themselves, infants or people who are not sensitive to heat are prohibited.

Product structure diagram

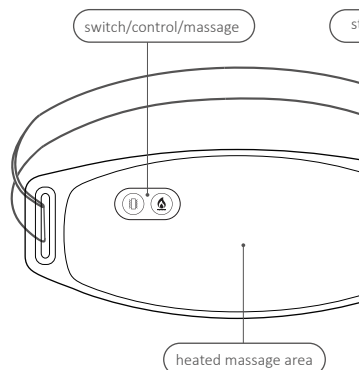Instructions for use

Bluetooth use:
The product needs to be powered on in advance (the red light is always on), then turn on the Bluetooth function of the mobile phone, scan the QR code to download the I-Health-app. Turn on the machine and enter the high temperature gear (red light), the temperature can be adjusted from high to low (65°C-45°C). Massage function 5 modes, off timer can be set between 0-120 minutes.

Not connected to Bluetooth:
Press and hold the massage button for two seconds to turn on and off the machine, the red light will be on for a long time, the massage function has 5 modes, 1 to 5 gears can be adjusted freely, short press 5 times for a cycle, 3 gears for heat (high temperature red gear 65°C, medium temperature blue gear 55°C, low temperature green gear 45°C), short press 3 times for a cycle (there is a gap before heating and massage enter the next cycle)
Power off automatically after 10 minutes without use.

Cleaning and storage

Please disconnect the product from the charging cable, turn off the heating switch and massage switch, and use a small brush to clean the surface of the product.

After brushing, put it in the room temperature area and wait until the product is completely dry before using or storing in the storage box.

When not in use, store as follows:


When storing the product, allow the product to cool down before folding.
During the storage period, do not press the item on the product to cause the product to crease.

This product is not intended for medical use in hospital.

Please keep this product out of the reach of children, away from humidity, high temperature and direct sunlight.

Maintenance

The fabric area does not support washing. When cleaning and maintenance are required, please use alcohol to remove sweat or wipe with a clean damp towel.

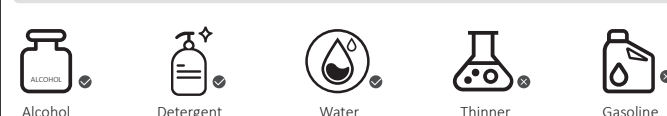Parameters

Product Name: Graphene heating belt
Item No: HEATING-BELT
Rated voltage: 5V
Rated power: 8W
Charging power input: USB-C 5V/2A

How to get rid of the product again:
This appliance has been identified in accordance with the European directive 2022/79/EG on Waste electrical and Electronic Equipment – WEEE.
The directive specifies the framework for an EU-wide value return and recycling of old appliances. Please dispose this product at the recycling site.

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