



age**LOC**<sup>TM</sup>  
*boost*<sup>TM</sup>  
USER'S MANUAL

# CONTENTS

- 1 ageLOC Boost™ Device
  - 1 Quick Start Guide
  - 1 Magnetic Charger
  - 1 Protective Treatment Surface Cover
- Treatment products designed for use with ageLOC Boost are sold separately.

## TABLE OF CONTENTS


CONTENTS	3
SAFETY INSTRUCTIONS, WARNINGS, AND PRECAUTIONS	4
BEFORE YOU BEGIN	6
CHARGING YOUR AGELOC BOOST™ DEVICE	6
USAGE INSTRUCTIONS	8
CLEANING AND CARING FOR YOUR AGELOC BOOST DEVICE	12
TROUBLESHOOTING	13
DEVICE DISPOSAL	13
MAINTENANCE	14
REPLACEMENT AND WARRANTY INFORMATION	15
PATENTS	15
DEFINITION OF SYMBOLS	15

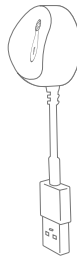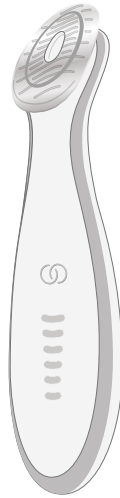



## SAFETY INSTRUCTIONS, WARNINGS, AND PRECAUTIONS


ageLOC Boost™ and Magnetic Charger are tested at pressures matching 4000 m in altitude.


This appliance contains batteries that are non-replaceable.

 **WARNING:** For the purposes of recharging the battery, only use the Magnetic Charger provided with this appliance.


 The ageLOC Boost™ device is only for use on healthy skin. Do not use if you have open sores or wounds, are prone to rosacea, have allergies to metal, or if you have problematic or overly sensitive skin. Check with your physician before using if you are pregnant, have a pacemaker or similar device, are epileptic, wear metal braces on your teeth, have a metal implant, or are ill.


 ageLOC Boost is not intended for use by persons (including children) with reduced physical, sensory, or mental capabilities or who lack experience and knowledge. Keep out of reach of children.


 Do not apply directly to lips, eyes, or eyelids. The device and treatment product are for external use only. If you experience any prolonged redness or abnormal irritation, discontinue use immediately and seek medical advice if necessary.


 To reduce the risk of shock, burn, fire, or injury, do not fully immerse your ageLOC Boost device in water. Do not place or store where it may fall or


be pulled into a bathtub, shower, sink, or toilet.


 Periodically inspect for damage; never operate the device if it has been damaged.


 ageLOC Boost contains a lithium-polymer battery. To prevent the possibility of explosion, never expose the device to heat. Do not store near a heat source, such as a radiator, fire, or heat vent. Do not leave in a hot vehicle. Check with your carrier before shipping or flying with this device.


 Rechargeable batteries, like those of ageLOC Boost, need to carry a minimal charge in order to function properly. We recommend charging your device whenever low battery is indicated.


 Note, ageLOC Boost is charged using inductive charging. Always use the provided charger when charging. Do not attempt to charge using any other inductive charger. If damaged, contact Nu Skin Support Services.

 Do not charge this device in the bathroom or near water.


 Do not use harsh chemicals or abrasives on your ageLOC Boost device.

 Do not leave ageLOC Boost in a cold environment for prolonged periods of time.

 Always place your device on a heat-resistant, stable, flat surface, especially when charging.

 Use ageLOC Boost only as directed.

## BEFORE YOU BEGIN



 Please read all safety and technical information before operating ageLOC Boost™.

## CHARGING YOUR AGELOC BOOST™ DEVICE

Fully charge ageLOC Boost before using for the first time.

Before charging, make sure both the device and the magnetic charger are completely dry.

Place the device on the magnetic charger (in the correct position, see diagram). Make sure the magnets snap into position. Plug the USB cable into a power source.

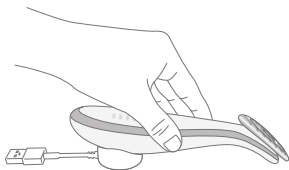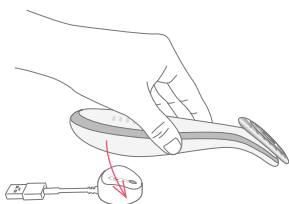

Charge ageLOC Boost completely. When used as recommended, a fully charged device will last approximately one month.

When charging begins, ageLOC Boost will vibrate once and the LED light bar will slowly cycle up—this indicates the device is charging and displays the current level of battery charge.



After two minutes, the light bars will glow on and off to indicate the current battery level. Once the device is fully charged, the lights will stop glowing on and off and the LED light bars will remain dimly lit.



ageLOC Boost™ is charged using inductive charging. Always use the provided charger to charge the device. Do not attempt to charge using any other inductive charger. If damaged, contact Nu Skin Support Services.



## USAGE INSTRUCTIONS

Use ageLOC Boost™ as part of your morning skin care routine, after cleansing and toning.


**PLEASE NOTE:** ageLOC Boost has been specifically designed to be used with ageLOC Boost Activating Serum. For your safety and comfort, we do not recommend using the device with any product other than Activating Serum.

### TREATMENT TIME


A treatment consists of two (2) one-minute, half-face treatments. After one minute, ageLOC Boost will pause briefly to allow you to apply additional ageLOC Boost Activating Serum to the treatment surface before you begin treating the second half of your face.




**STEP 1** Turn on the device by swiping your finger up over the full length of the touch surface. Lights will turn on.

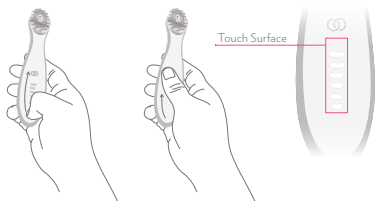


**STEP 2** Dispense three pumps of ageLOC Boost Activating Serum onto the device's treatment surface. Using the device, spread Activating Serum evenly over one half of your face.

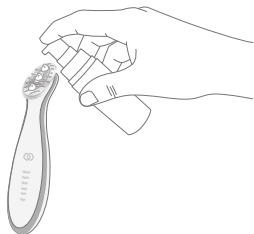

**STEP 3** Begin treatment by swiping your finger up over the full length of the touch surface. The device will gently vibrate indicating treatment is in progress. Place treatment surface on face and begin treating skin by using slow movements in an upward and outward motion for one minute.

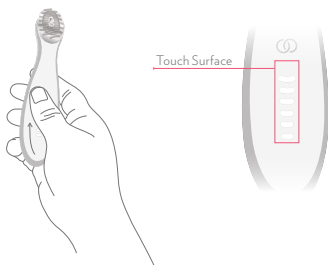


**STEP 4** After the first one-minute treatment is complete, the device will pause automatically. Dispense an additional three pumps onto the treatment surface of the device, and spread evenly over the other half of the face. Swipe up the full length of the touch surface. The device will resume gently vibrating. Continue to treat the skin using slow, upward and outward motions for the remaining minute.

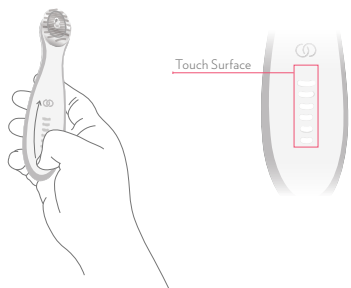


**STEP 5** The device will automatically turn off after completing a full treatment cycle. Once treatment is complete, massage excess into skin, then follow with your daily serum and moisturizer.

**TO PAUSE TREATMENT** At any time, simply swipe down over the full length of the touch surface. Swipe back up full length of touch surface to resume.

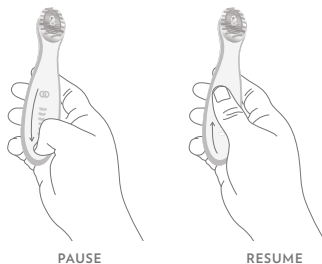


**TO STOP TREATMENT MANUALLY** during a treatment, you can stop the device at any time by swiping down the full length of the touch surface two times.

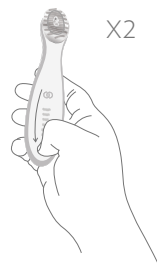

**PLEASE NOTE:** ageLOC Boost was designed to slowly glide across the skin's surface during use and should not be held in one spot. If you experience any prolonged redness or abnormal irritation, discontinue use immediately and seek medical advice if necessary.

## CLEANING AND CARING FOR YOUR AGELOC BOOST™ DEVICE

ageLOC Boost should be cleaned according to these instructions after each use to avoid the accumulation of residue. After every treatment, rinse ageLOC Boost treatment surface under running water to ensure debris and product is removed from crevices. Thoroughly wipe the device dry with a soft washcloth. Store in a dry place.

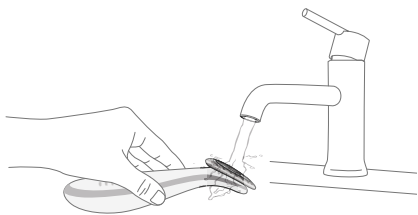

## TROUBLESHOOTING

### LOSS OF CONDUCTIVITY DURING

**TREATMENT:** If the microcurrent conductivity is temporarily lost, the device's LED light bar will flash in continuous sets of five short flashes coupled with five sharp vibrations. If this happens, apply more ageLOC Boost Activating Serum to the device's treatment surface. Swipe up along the touch surface to resume treatment.

**UNRESPONSIVE DEVICE:** Swipe up slowly from the bottom of the device along the full length of the touch surface. If the device still does not respond, charge it fully and try again. If the device still fails to respond, contact Nu Skin Support Services.

**LOW BATTERY:** The bottom light on the device will turn amber. Charge the device.

**EXTREME LOW BATTERY:** The amber light at the bottom of ageLOC Boost™ will flash repeatedly when the user attempts to start a treatment. This indicates there is not enough battery remaining to perform a treatment, and the device will turn off. Charge the device.

## DEVICE DISPOSAL

You must dispose of ageLOC Boost properly according to local laws and regulations. ageLOC Boost contains electronic components and a lithium-polymer battery, which requires disposal separate from household waste. When ageLOC Boost reaches its end of life, contact your local waste management department to learn about proper disposal and recycling options.



Dispose of Batteries Properly for Your Area

## MAINTENANCE

Do not attempt to repair the device yourself. There are no user-serviceable parts.


Do not ship damaged devices by air.

ageLOC Boost™

Maximum Output:  
5 V, 250  $\mu$ A, 0.414 Hz

 IP25

ageLOC Boost Magnetic Charger

Input: 5 V  500 mA

 IP22

In Canada, this device complies with requirements of ICES-GEN, ICES-001, and ICES-003. This charger complies with the requirements of ICES-GEN, ICES-001, and RSS-216.

**FOR THE DEVICE:** CAN ICES-003B/  
NMB-ICES003B

**FOR THE CHARGER:** CAN ICES-0  
01B/NMB-ICES001B









## REPLACEMENT AND WARRANTY INFORMATION

**LIMITED TWO-YEAR WARRANTY:** Nu Skin guarantees your device to be free from defects in materials and workmanship for a period of two years from the original date of consumer purchase. This warranty does not cover damage to the product as a result of misuse or accident, including dropping the device. If the product becomes defective within the two-year warranty period, please call your local Nu Skin Support Services.

## PATENTS

Numerous U.S. and international patents pending.

## DEFINITION OF SYMBOLS

Symbol	Definition
	Warning
	Consult Owner's Manual
	Dispose of Batteries Properly for Your Area
	Two Levels of Shock Protection
	Charge only with specified charger
	Direct Current





SOLD EXCLUSIVELY BY NU SKIN ENTERPRISES AUTHORIZED  
DISTRIBUTORS  
MFD. IN CHINA EXCLUSIVELY FOR NSE PRODUCTS, INC.

©2020 NSE PRODUCTS, INC., PROVO, UT 84601 NUSKIN.COM  
001931 94512376/1