


Dual Tapper

HANDHELD PERCUSSION MASSAGER

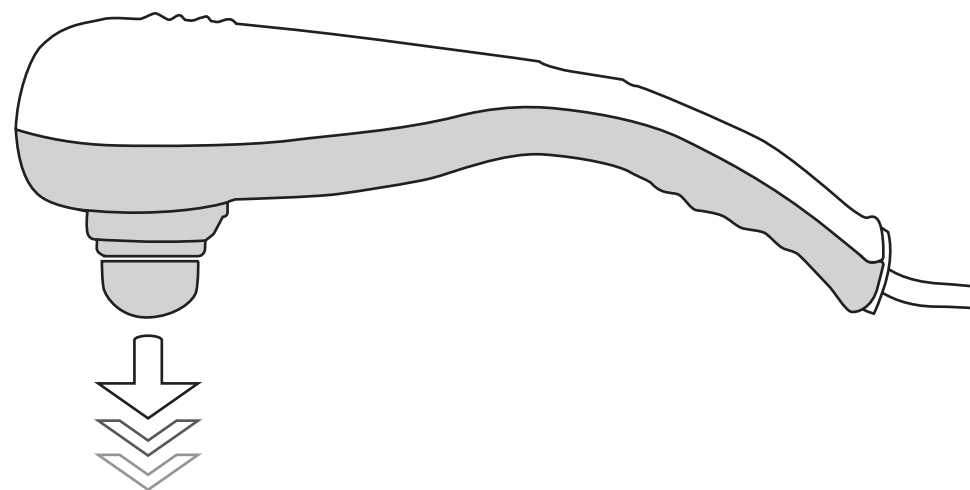

SAVE THESE INSTRUCTIONS

Thank you for purchasing the Dual Tapper.
Please read this manual carefully to ensure optimum performance and safety.

CONTENTS

Thank you for purchasing the USJ-768 Dual Tapper.
Before using the product, please read this manual carefully, paying special attention to the safety precautions. Retain this manual for future reference.

SAFETY PRECAUTIONS	2
PRODUCT FEATURES	4
PART NAMES	4
OPERATING INSTRUCTIONS	5
MAINTENANCE & STORAGE	6
SPECIFICATIONS	6
WARRANTY	7

SAFETY INSTRUCTIONS

Before using this unit for the first time, be sure to read these Safety Instructions carefully, then use the unit accordingly. These precautions are listed to prevent damage to the unit and injury to users.



To reduce the risk of injury to persons:

1. Only operate this unit under its intended usage, as described in this manual.
2. Make sure that the voltage on the massager matches that of the power outlet before connecting the plug.
3. After using the massager, always return the power switch to the (0) position, and disconnect the plug. Do the same before cleaning the unit.
4. To avoid danger, never leave the unit idle with the plug still connected to the power outlet.
5. Never sit or step on the unit.
6. Never attempt to cram material substances into any of the unit's crevices.
7. Never soak this unit in water, nor use it within humid environments, such as while bathing.
8. Do not place the unit near a radiator or any other electric appliance.
9. If the unit accidentally falls into water, shut off power immediately and disconnect its plug from the power outlet. Do not lift it out of the water too suddenly.
10. Never operate the unit while sleeping.
11. Do not use this unit right before going to sleep; the stimulation it generates may impede slumber.
12. Never use the unit under the cover of bedding, such as a blanket or pillow. Excessive heating may occur, leading to fire, electric shock, or personal injury.
13. Never operate the unit where aerosol spray products are being used, nor where oxygen is being administered.
14. Do not apply this unit around the eyes, on sensitive skin, or on hypersensitive parts of your body.
15. Do not place the hammerhead on the head, chest, elbow joint, or knee joint for extended periods of time. Doing so may damage your body.
16. Do not allow children or handicapped individuals to operate this unit.


SAFETY INSTRUCTIONS

17. Stop using the unit if you experience any discomfort. Consult a doctor.
18. If the power cord malfunctions or becomes damaged, do not attempt to repair it by yourself. Stop using the unit immediately and contact its manufacturer.
19. The recommended usage time for this unit is 15 minutes per session. Extended use may cause overheating and shorten its service life. If overheating has occurred, turn off the unit and allow it to cool before its next use. session.
20. To disconnect the machine from its power source, grasp the plug instead of pulling on the power cord.
21. Do not apply the unit within one hour after meals. Drink one cup of water before each massage; drink another two cups of water after each session.
22. Persons with the following conditions should consult a doctor before using this unit:
 - Are pregnant, have undergone an abortion, or incurred post-operative inflammation.
 - Suffer from heart disease, fractures or joint dislocations, other acute ailments, or any abnormal body conditions.
 - Currently suffer from any burns or swelling on the skin.

▶ PRODUCT FEATURES

- Convenient, stylish design
- Dual hammerheads relieve discomfort and renew vitality
- Ergonomic model draws weight to the hammer. Operate by simply holding the unit's handle to your body.
- Contoured handle adapts the unit to any body part.
- 3 interchangeable head attachments

▶ PART NAMES

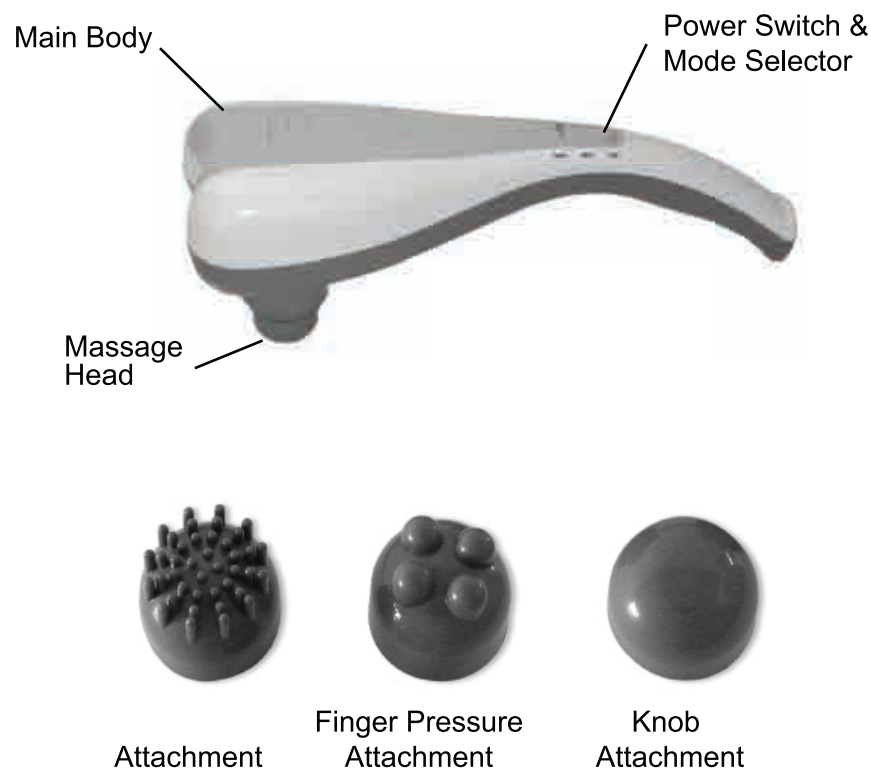

▶ OPERATION

1. Please make sure that the power switch is in the OFF position before connecting the power plug.
2. Place the dual massage heads against the body part that you wish to massage.
3. Push the power switch for LOW for low speed, or HIGH for high speed. We advise that you begin each massage session at the low speed before progressing to the high speed.
4. If desired, move the massage heads around to cover your entire target area.
5. To stop the massage, push the power switch back to the OFF position. Then pull out the plug from the power socket.
6. To prevent overheating, let the unit rest for at least 20 minutes before using it again.
7. If you feel any discomfort on your skin at any time during your massage session, stop using the unit immediately.

Massage Head Attachment Options

1. Please make sure that the power switch is in the OFF position.
2. To affix a new massage head, first ensure that no other attachments are currently placed on the two permanent massage heads.
3. To detach an existing attachment from a massage head, grasp its edges by the tips of your fingers to peel it off. If necessary, the tip of a screwdriver, the long side of a nail file, or other similarly stiff, thin, and flat object may be used in place of your fingers.
4. To attach a new massage head, place its hollow end onto one of the two permanent massage heads already on the unit. Place your palm over the attachment and press down firmly until that attachment fully covers the permanent massage head.

► MAINTENANCE & STORAGE

1. This unit requires no regular lubrication, maintenance, or repair.
2. Please keep all crevices on the unit free from dust, hair, and accumulated filth.
3. Before cleaning the unit, always disconnect the power plug from the AC outlet.
4. Gently wipe all parts of the unit with a dry cloth. Do NOT clean the unit or its power switch with a wet cloth.
5. When not using the unit for an extended period of time, coil its power cord and store the entire unit in a cool, dry place. Do not expose it to hot air or open fire.
6. Do not expose the unit to direct sunlight for an extended period of time.

► SPECIFICATIONS

Product Name: Dual Tapper
Model Number:..... USJ-768
Voltage:..... 120V AC / 60 Hz
Rated Power:..... 20W

► WARRANTY

USJ-768 DUAL TAPPER

Consumer Limited Warranty: ONE-YEAR WARRANTY

This warranty guarantees that the product purchased is to be free from defects in material and workmanship under normal usage and local voltage for a period of one year. This warranty expires in **TWO (2) years on parts, and ONE (1) year on labor** after the date of purchase.

U.S. Jaclean, Inc. (referred to as “the Warrantor”) will repair this product with new or rebuilt parts in case of defects in material or workmanship, free of charge, in the USA provided that:

- 1) The product is operated in accordance with the printed instructions supplied.
- 2) The product is used under normal conditions for which designed.
- 3) The product is not subjected to misuse, negligence, or accident.
- 4) The product receives proper care, protection, and maintenance.

This warranty is extended to the original end-user when purchased from an authorized U.S. Jaclean, Inc. distributor. A purchase receipt or other proof of date of the original purchase is required for service and parts replacement under this warranty. This warranty does not cover shipping charges for defective items if a claim is made more than 30 days after the date of purchase.

This warranty only covers failures due to defects in materials and workmanship, and does not cover normal wear or cosmetic damages. This warranty does not cover damages which occur in shipment, or failures which result from: accident, misuse (including, without limitation, failure to comply with instructions and use limitations contained in the Operating Instructions), abuse, neglect, mishandling, faulty installation, set-up adjustments, improper operation or maintenance, alteration, modification, power line surge, lightening, improper line voltage, commercial use such as hotel, office, restaurant, or other business or rental use of product, service by anyone other than from U.S. Jaclean, Inc., or damage that is attributable to “acts of God.”

Continued use of this product after discovering a defect voids this warranty. This warranty shall be governed by the state of California, United States of America.

IMPORTANT: DO NOT return defective product to the store.

Contact U.S. Jaclean, Inc. for repair services.

Please retain this Warranty Certificate, together with the original purchase receipt, and present it to our technician for inspection during service.

For assistance and more information, please contact U.S. Jaclean, Inc.



Distributed by U.S. Jaclean, Inc.
1816 West 135th Street
Gardena, CA 90249

1-800-242-2221
9AM - 5PM PST (Monday - Friday)

usjclean.com