

SHOGUN HAMMER LITE



SAVE THESE INSTRUCTIONS

Thank you for purchasing the Shogun Hammer Lite.
Please read this manual carefully to ensure optimum performance and safety.

SHOGUN HAMMER

Thank you for purchasing the USJ-894 Shogun Hammer Lite. Before using this unit for the first time, please read through all sections of this User Manual carefully, paying special attention to the Safety Instructions. Please retain this manual for future reference.

Included:

- 1x Shogun Hammer Lite
- 3x Massage heads
- 1x Type-C USB cable
- 1x Carrying case
- 1x User manual

CONTENTS

IMPORTANT SAFETY PRECAUTIONS	1
PARTS	3
OPERATION	5
CLEANING AND MAINTENANCE	6
TROUBLESHOOTING.....	6
BATTERY DISPOSAL.....	7
TECHNICAL SPECIFICATIONS.....	7
WARRANTY	8

IMPORTANT SAFETY PRECAUTIONS

Before using this unit for the first time, be sure to read these Safety Instructions carefully. To avoid burns, fire, electric shock, or other injuries, read the following instructions before using the unit.

 **DANGER** – To reduce the risk of electric shock:

1. Do not operate this unit in an excessively humid environment, such as a bathroom.
2. Keep the unit away from excessive heat, such as direct sunlight, an oven, or an open flame.

IMPORTANT SAFEGUARDS


1. Check the product carefully and read the user manual carefully before use.
2. This product will not work while charging.
3. Do not charge the product in an environment where the temperature is below 0°C/32°F or above 50°C/122°F.
4. Do not use the product in a humid environment, such as a bathroom.
5. Do not use the product on wet skin.
6. Do not use the product in an environment where the temperature is higher than 40°C/104°F.
7. Do not place any objects on the air outlet or put any objects inside the opening of the product. The blockage of the air outlet might cause the internal parts to overheat, leading to malfunction or damage.
8. If the product is moved from a low-temperature environment to a warm environment, moisture might accumulate on internal metal parts, which would influence the normal operation of the product or cause damage to it. Move the product to an environment at room temperature for 1 hour before use.
9. The default working time of the product is 10 minutes. After using it continuously for 10 minutes, press the "⏻" button to re-start.
10. Do not turn the product on and off frequently. After the product is turned off, please wait for more than 3 seconds before re-starting again.

IMPORTANT SAFETY PRECAUTIONS


11. During use, if the product stops working from overload, please wait for more than 5 seconds before using again.
12. Do not press the massage head forcefully while the product is on.
13. Do not use the product and contact our customer service in the following situations:
 - If it has a damaged adapter or charging port
 - If it is not working normally
 - If it is dropped or damaged
 - If it is dropped in water
14. Stop using the product immediately if you feel discomfort during use.
15. Do not use on or near fragile body parts such as head, neck, or heart.
16. Do not use on or near spinal vertebra or cervical vertebra without guidance from a professional doctor to avoid damage to nerves.
17. Do not use on insensitive or overly sensitive skin or areas with poor blood circulation.
18. Do not use on any hard body parts or bones. Pay attention to intensity to avoid the massage area getting bruised. Use on the soft body tissue where there won't cause pain or discomfort.
19. Do not use 1 hour before or after a meal.
20. Do not use if you have consumed alcohol.
21. People listed below should consult a doctor before using the product:
 - People who are pregnant
 - People with osteoporosis
 - People with heart diseases
 - People who have a pacemaker or other implanted medical devices
 - People with acute diseases
 - People who have trauma on the skin
 - Babies, children, or minors
 - People with physical/mental dysfunction
 - People with muscles under recovery stage

PARTS

PRODUCT DIAGRAM

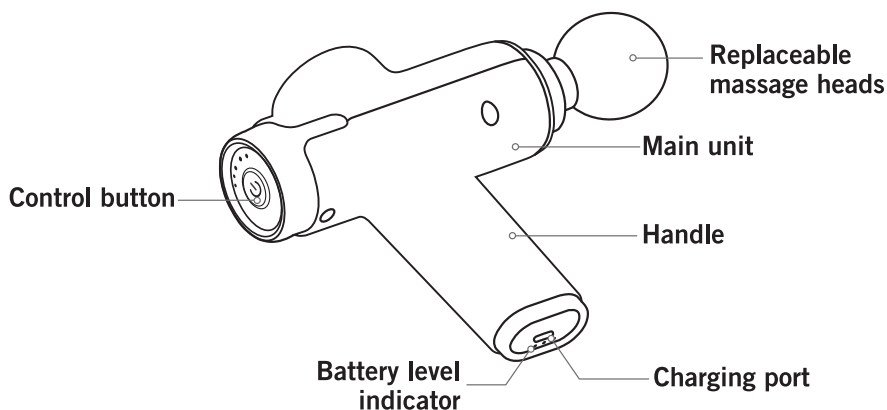

MESSAGE HEADS

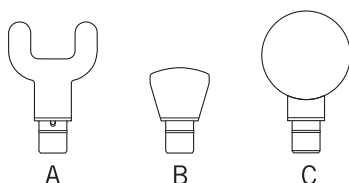

A. U-shaped massage head

- Suitable for the sides of neck, spinal column, and Achilles tendon

B. Air Cushion massage head


- Suitable for smaller muscles such as the forearms, hands, and feet

C. Ball massage head

- Suitable for large muscular parts such as the chest, back, and hips

Note: the above suggestions are for reference only. Use depending on your physical conditions or under advice from your doctor.

CONTROL PANEL

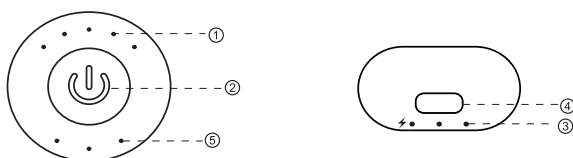

1. Intensity indicator lights

- When the product is turned on the light will illuminate
- There are 5 lights to indicate the intensity of the massager
- The more lights illuminated, the higher the intensity

2. Power button

- Press once turn on the massager
- Press again to turn up the intensity
- Press this button for 2 seconds to turn off the massager

3. Battery level indicator lights

- These 3 lights indicate the current battery level
- When the battery is fully charged, all 3 lights will remain on
- When the battery is running low, the indicator lights will flicker
- When the battery is charging the lights will also flicker

4. Charging port

- Type-c interface input charger, connect adaptor to 5v output port


5. Mode indicator lights

- These 3 indicator lights indicate the auto massage modes
- Lights will illuminate to show current auto mode selection


OPERATION

PRODUCT OPERATION

1. Fully charge the product 5 hours before first use. This product will not work while charging.
2. Select the massage head among the 3 you intend to use. Align the massage head with the massage head connector, then push firmly. Make sure that the massage head cannot be removed easily with light force.

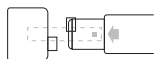

Correct installation direction

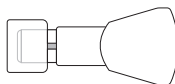

Installation finished

3. Press the “⏻” button for 2 seconds to activate the device. The battery indicator lights will display the current battery level. If the battery is fully charged, all the indicator lights will be on. If the battery is running low, the level 1 indicator light will flicker. Please charge promptly in this case.
4. When the battery indicator lights remain steadily on, press the “⏻” button to switch the speed level. Each time the button is pressed, the intensity increases by one level. There are 5 levels of intensity.
5. When the product is on the fifth speed level pressing the “⏻” button again will switch to the auto massage modes, which are recommended to use after relaxation, sports, and workouts separately. Press the “⏻” button to shift back to fixed-speed massage mode.

Automatic Massage Modes:

Option 1 switches between speed level 1, 2, and 3.

Option 2 switches between speed level 2, 3, and 4.

Option 3 switches between speed level 3 and 4.

6. Do not use the product on one part of the body for a long time, and do not leave it on the same spot for more than 1 minute.
7. Fully charge the product if it is not to be used for a while.

CLEANING AND MAINTENANCE

CLEANING AND MAINTENANCE

1. If the product is not to be used for a long time, turn off the product via ON/OFF button on the back of the handle. Then fully charge it, place it in the original package and store in a dry, cool and airy environment.
2. If the product is not to be used for a long time, it is recommended to charge it every 2 months to prevent the battery from being damaged due to lack of power.
3. Battery level should be maintained above the medium level. Regular charging and maintenance will help prolong battery life.
4. Keep the product clean. For daily maintenance, clean with a neutral detergent and a soft cloth. Do not clean using motor oil, organic solvents or; other chemicals. Do not rinse in water.
5. Avoid high pressure and avoid hard objects scratching the product surface.

TROUBLESHOOTING

Problem	Solution
The control panel does not light up	Check whether the button on the back of the handle has been turned on.
	The battery has been dead for a while. Charge the battery for at least 8 hours before re-starting.
The product is not working.	Check whether the button on the back of the handle has been turned on and if the battery level indicator lights on the control panel have lit up.
The product stops working during use.	The product has worked continuously for 10 minutes. Please re-start if you want to keep on using.
	The battery is running out. Turn off the product via the ON/OFF button on the back of the handle and then charge before use.

BATTERY DISPOSAL

BATTERY DISPOSAL

Precautions when disposing of the product:

1. This product should not be disposed of with household waste and should be disposed of in a facility for recycling and re-using.
2. The lithium battery in the product is a valuable resource that can be recycled. When the product is at the end of its service life, turn to professional to take the lithium battery out of this product and give the battery to a qualified professional organization for the recycling and regeneration of rechargeable batteries.

TECHNICAL SPECIFICATIONS

Product Name: Shogun Hammer Lite

Model Number: USJ-894

Dimensions: 7.24 x 6.46 x 2.28 in

Rated Voltage: 5V

Rated Power: 20W

Auto-off Timer: 10 minutes

LIMITED WARRANTY

USJ-894 SHOGUN HAMMER LITE

This warranty guarantees that the product purchased is to be free from defects in material and workmanship under normal usage and local voltage for a period of one year. This warranty expires in TWO (2) years on parts, and ONE (1) year on labor after the date of purchase.

U.S. Jaclean, Inc. (referred to as "the Warrantor") will repair this product with new or rebuilt parts in case of defects in material or workmanship, free of charge, in the USA provided that:

- 1) The product is operated in accordance with the printed instructions supplied.
- 2) The product is used under normal conditions for which designed.
- 3) The product is not subjected to misuse, negligence, or accident.
- 4) The product receives proper care, protection, and maintenance.

This warranty is extended to the original end user when purchased from an authorized U.S. Jaclean, Inc. distributor. A purchase receipt or other proof of date of the original purchase is required for service and parts replacement under this warranty. This warranty does not cover shipping charges for defected items 30 days after the date of purchase.

This warranty only covers failures due to defects in materials and workmanship and does not cover normal wear or cosmetic damages. The warranty does not cover damages which occur in shipment, or failures which result from accident, misuse (including, without limitation, failure to comply with instructions and use limitations contained in the Operating Instructions), abuse, neglect, mishandling, faulty installation, set-up adjustments, improper operation or maintenance, alteration, modification, power line surge, lightning, improper line voltage, commercial use such as hotel, office, restaurant, or other business or rental use of product, or service by anyone other than from U.S. Jaclean, Inc.

Continued use of this product after discovering a defect voids warranty. This warranty is governed by the state of California, United States of America, and you may also have other rights which vary from state to state.

IMPORTANT: Do NOT return merchandise to the store.
Contact U.S. Jaclean, Inc. for repair services.

Please retain this Warranty Certificate, together with the original purchase receipt, and present it to our technician for inspection during service.
For assistance and more information, please contact U.S. Jaclean, Inc.



Distributed by U.S. Jaclean, Inc.
Gardena, CA 90249

1-833-542-1918
9AM - 5PM PST (Monday - Friday)
usjaclean.com