


取扱説明書 | User Guide
音波電動歯ブラシ | Sonic Electric Toothbrush
Model: EPET005

目次

- 本体パーツ説明 01
- 仕様 01
- 使用方法 01
- 歯ブラシの使用について 02
- 歯ブラシの使用手順 03
- 安全上のご注意 04
- よくあるご質問 05
- 電寿命目安 05
- 環境保護と資源回収 05
- 保証書 06

本体パーツ説明

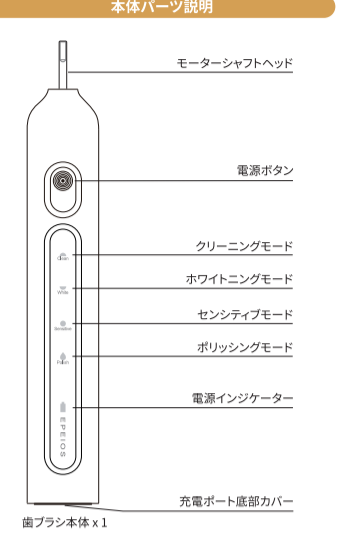

モーターシャフトヘッド
電源ボタン
クリーニングモード
ホワイトニングモード
センシティブモード
ポリッシングモード
電源インジケータ
充電ケーブル接続カバー

歯ブラシ本体 x1
専用歯ブラシヘッド
レギュラーサイズ x1
コンパクトサイズ x1
USB充電ケーブル x1
取扱説明書 x1
コンパクトサイズ x1

充電時の使用可能時間

約100分
使用時間は充電量や使用条件により異なります。

注意

- 初めて使用する際、または長期使用していない場合は、充電を開始するに充電インジケータが点灯しない場合があります。この場合は充電を続けてもう一度インジケータランプが点灯まで待ちます。
- 推奨される充電環境は室温20〜25度です。室温が上記の範囲外の場合は充電環境を調整してください。
- 充電中に本体や充電器が熱くなる場合があります。冷却してから再度充電してください。
- 充電完了後は、または他のものにさらされる場合は、充電器から取り出してください。
- 充電器は、充電完了後は電源オフにしてください。
- 充電器は充電完了後に自動的に充電を停止しますが、これは歯ブラシを初めて使用する際に充電が完了していない場合があります。この場合の充電状態は通常1〜2週間以内になります。出荷の段階では充電完了後に自動的に充電を停止する場合があります。

歯ブラシの使用手順

約30秒のブラッシング/ワイピング/タイマー機能

ブラッシングの準備

ブラッシングの準備は、歯を洗う前に歯をブラッシングすることです。歯を洗う前に歯をブラッシングすることにより、歯の汚れを落とすことができます。

モード/メモリ機能

歯ブラシモードを覚えておくには、歯ブラシを充電器から取り出し、電源をオフにしてください。電源をオフにした後、歯ブラシを充電器に挿入し、電源をオンにしてください。電源をオンにした後、歯ブラシを充電器から取り出し、電源をオフにしてください。

安全上のご注意

この製品は、子供や高齢者の手に届かない場所に保管してください。

注意


- この製品は、禁止されている手順を実行しないようにしてください。
- この製品は、禁止されている手順を実行しないようにしてください。

環境保護と資源回収

本製品はリサイクル電池を使用しています。お住まいの地域による規制がある場合は、地域の規制に従ってください。

保証書

保証期間(無料修理)は、購入日から適用されます。本製品はリサイクル電池を使用しています。お住まいの地域による規制がある場合は、地域の規制に従ってください。




取扱説明書 | User Guide
音波電動歯ブラシ | Sonic Electric Toothbrush
Model: EPET005

Contents

- Product Diagram 01
- Specifications 01
- Charging & Operating 01
- Use Instructions 02
- After Use 03
- Safety Instructions 03
- FAQs 04
- Battery Life 05
- Disposal of the Battery 05
- Warranty Card 05

Product Diagram

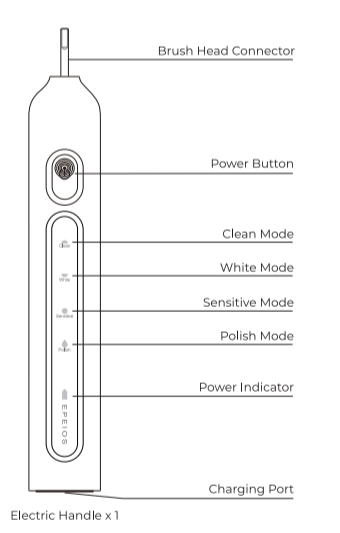

Brush Head Connector
Power Button
Clean Mode
White Mode
Sensitive Mode
Polish Mode
Power Indicator
Charging Port

Electric Handle x1
Brush Head Regular Size x1
Compact Size x1
USB Charging Cable x1
User Guide x1

Charging Indicator

The LED indicator on the handle will flash red when charging and stay solid green when fully charged.

Working time when fully charged

Approx. 100 minutes
Working time may vary due to environmental temperature or using conditions.

Note

- The charging time may vary if the product has been idle for over 6 months or during the initial use. The charging indicator may not light up immediately when charging has started. It will turn on when you continue charging the product.
- It is recommended to charge the product in temperature from 5°C to 35°C. Charging in temperatures below or above the recommended range may affect the battery life.
- It is normal for the brush handle or charging adapter to get warm during charging.

Use Instructions

This product may not function in an environment with a temperature of 5°C or below. The vibrations in cleaning mode are too intense for you to use in sensitive mode. Please keep the handle away from hair care products or hand cream.

1. Install the Brush Head

- Make sure the orientation is correct. Align the front of the handle with that of the brush head.
- Do not rotate the brush head with excessive force (to avoid damage).
- The brush head connector may move during assembly, which does not mean the product is defective.

2. Put the Brush Head into Your Mouth

- Place the brush head into your mouth first before you turn on the brush to avoid spills.

3. Press the Power Button to Turn on the Brush

- Press the Power Button gently against your teeth.
- Sonic vibration may cause itchy or tingling sensation during the initial use. These feelings will disappear after several uses.
- In cleaning mode, you should not brush too hard against your teeth, the vibration sound will be too loud and the cleaning sound may change too.

4. Set the Cleaning Mode

- Press the Power Button within 10 seconds after the brush has been turned on for over 10 seconds. Press the Power Button within 10 seconds when it has been activated to change the cleaning mode.

5. Turn off the Brush after Cleaning is Done

- Press the Power Button first before you take it out of your mouth to avoid spills.

After Use

Clean the Brush Head and the Handle
Rinse the brush head and the handle with water after use. Do not use toothpaste, powder, or tooth whitening products. (To avoid malfunctions, cracking, or discoloration.)

1. Remove the Brush Head

- Remove the brush head with care.

2. Rinse the Brush Head

- Rinse off gel from the inside and bottom of the brush head.
- Do not pull or split the bristles with excessive force to avoid shedding.
- Do not wash in water over 60°C to avoid damage.
- Do not use sharp objects like nails to clean the tongue cleaner (on EPET005 only).

3. Rinse the Handle

- Do not damage or scrape the silicone sealing ring (A) with hard objects.
- Do not immerse the handle in water for a long time to avoid warping.
- Do not wash with water over 42°C to avoid damage.

4. Dry the Brush Head and Handle

- Do not charge the product when it is wet.

Safety Instructions

To avoid electric shock, injuries, fire hazards, death or damage and property damage losses, please read and follow the below cautions strictly.

Signs

The following signs classify the levels of dangers, damage, and property loss caused by ignorance of the signs or improper use.

⚠ Danger High potential danger that may cause severe injury or death.

⚠ Warning Potential danger that may cause severe injury or death.

⚠ Attention Potential danger that may cause injury or property loss.

Safety Instructions

The below signs classify and illustrate the instructions you should follow.

⚠ Warning

Prevent accidents
Please keep this product beyond children's reach and do not let them operate this product or they may swallow the product parts to cause injury and accident.

⚠ Attention

Protect Your Teeth and Gum
Please do not use the product together with toothpaste with whitening or abrasive agents, or whitening toothpaste.
Do not use excessive force to press the brush head or the brush head's plastic part.
This product is not suitable for people without professional ability.
This product is not suitable for people with oral cavity sensitivity.
People suffering from severe periodontal disease, ulcer or gum treatments or wearing oral cavity appliances should consult a dentist before use.
Do not immerse the handle in water for a long time to avoid warping.

⚠ Notes Avoid drop or collision to avoid injury.

Specifications

Product Name	Sonic Electric Toothbrush
Model	EPET005
Input	DC 9V 50mA
Charging Type	Type-C
Battery Capacity	2200mAh
Charging Time	About 4 hours (at temperature 5°C to 35°C)
IP Rating	IPX7
Net Weight	132g±5%
Product Size	118x mm x W 28mm x H 30mm

For Household Use Only.

FAQs

Q: Why does the brush pause every 30s?
A: The brush pauses every 30s to remind you to change the brush area.

Q: Will long idle time affect the battery performance?
A: The battery performance will decrease if this product has been idle for over 3 months. Please fully charge the battery every 3 months.

Q: Will frequent charging affect the battery life?
A: The built-in battery is a lithium battery. Its lifetime will not be affected even if you charge it every time after use.

Battery Life

Approx. 3 years. If the using time has been dramatically decreased even after a full charge, the battery may be at the end of its lifetime. The battery life may vary according to use and store conditions.

Disposal of the Battery

This product contains a lithium battery. Please dispose of it as per local law and regulations.

Warranty Card

From the date of purchase, this product comes with free repair within the warranty period.

1. Please show this warranty card and your purchase receipt to claim your warranty.

a) We will not issue another warranty card if you lose your purchase receipt.
b) This warranty card only applies to the country or region of your purchase.

2. From the date of purchase and within the below period, please refer to the manufacturers under proper use.
Brush Head—Provide at a charge.
Electric Handle—Cover the Brush Head—Free for 1 year.
3. The warranty is void if one of the following items happens. But we can repair it at a charge.
A. Damages caused by failure of following the instructions, maintenance, and store conditions as stated in this user guide.
B. Damages caused by disassembly of an unauthorized repairing man.
C. Failed to provide the warranty card and valid purchase invoice.
D. Product information in this guide does not conform with that on the product or it has been altered.
E. Damages caused by force majeure (thunder, earthquake, fire hazards, flood).

4. This warranty card only applies to customers who purchased this product for personal care use and does not apply for commercial use.

Thank You for Choosing EPEIOS

Thank You for Choosing EPEIOS

