celluma SERIES

Celluma User Manual

Models PRO, LITE & ELITE

Celluma PRO

Celluma LITE

Celluma ELITE

Table of Contents

Intended Use	3
Device Description	3
Important Safety Instructions	3
Contraindications	5
Celluma Warnings & Product Safety Information	6
Preparing the Celluma for Use	7
Celluma – Protocols for Use	8
Celluma Product Assembly	9
Suggest Product Positioning	10
Celluma--Instructions for Use	12
Select and Begin a Celluma Program	12
Caring for Your Celluma	13
Celluma Accessories	14
Useful Tips and Frequently Asked Questions	14
Technical Support	15
Technical Information	15
Contact Information	17

Intended Use:

The BioPhotas Celluma® device is intended to deliver heat in the IR spectrum to provide topical heating for the purpose of elevating tissue temperature; for the temporary relief of minor muscle and joint pain, arthritis and muscle spasm; relieving stiffness; promoting the relaxation of muscle tissue; and to temporarily increase local blood circulation. The blue spectrum is intended to reduce mild to moderate acne vulgaris. The Celluma is intended to emit energy in the red and infrared region of the spectrum for use in dermatology for the treatment of periorbital wrinkles.

Device Description:

Celluma light emitting diodes (LEDs) delivers three different wavelengths (880nm, 640nm and 465nm) simultaneously to manage a variety of skin and musculoskeletal conditions. All three wavelengths work together in each of the programs to optimally address specific conditions. The benefit of using three wavelengths is that light-absorbing molecules respond to different wavelengths of light energy, resulting in a targeted increase in cellular activity at different tissue depths. The unique lightweight and flexible design of the Celluma molds closely to the body to facilitate effective energy absorption by tissue.

Important Safety Instructions:

- ❖ Read all instructions before using this device. Save these instructions for future reference.
- ❖ This manual contains important safety symbols and instructions. Please pay attention to these symbols and follow all instructions given.

Consult Accompanying Documents

The symbol in the figure above appears on the product. This symbol indicates “Follow all instructions for use”.

General Warning or Caution

The Exclamation Symbol in the figure above appears in Warning and Caution statements in this document. This symbol designates an area where personal injury or damage to the equipment is possible.

Electric Shock Symbol

The Electrical Shock Symbol in the figure above appears throughout this manual. This symbol indicates a hazard arising from dangerous voltage. Any mishandling could result in irreparable damage to the equipment, and personal injury or death.

Type B Applied Part Symbol

This symbol indicates that the portion of the product which may contact the patient is provided with type B protection as defined in IEC 60601-1.

Alternating Voltage Symbol

This international symbol implies an alternating voltage or current.

Stand-by Symbol

The Stand-by Symbol in the figure above is located on the power switch on the Celluma controller. This button switches the Celluma between “on” and “stand-by”. To turn the product “on”, please press this button. To stop treatment and place the product into a stand-by mode, press this button once again. To remove all power from the Celluma, please un-plug the power supply from the power outlet.

Keep Dry Symbol

The Keep Dry Symbol appears on the product. The Celluma is not sealed and cannot be used outdoors or where it can become wet.

Manufacturer Symbol

This symbol accompanies the name and address of the manufacturer, BioPhotas.

Serial Number Symbol

This symbol accompanies the serial number of the device.

Reference Number Symbol

This symbol accompanies the model number of the device.

WEEE Directive Symbol

This symbol on the product or on its packaging indicates that this product must not be disposed of with regular waste. Instead, it is the user responsibility to dispose of waste equipment according to the local laws. For information about where the user can drop off the waste equipment for recycling, please contact your local BioPhotas, Inc. representative.

Contraindications:

(See pages 9-10 for precautions and safety information).

NOTE: When making changes to your health and wellness routines it may be advisable to consult with your physician before proceeding.

Do not use if the patient is pregnant or breastfeeding.

Do not use on infants or children.

Do not use if the patient or administrator has epilepsy or a history of seizures.

Avoid using when the patient is taking cortisone injections or any other kind of steroid injection(s).

Avoid use if the patient is using photosensitive drugs.


Celluma® Warnings & Product Safety Information:

Read and keep these instructions.

	<p>WARNING Do not attempt to operate this equipment if there is evidence of shipping damage or you suspect the unit is damaged. Damaged equipment may present additional hazards to you. Contact technical support for advice before attempting to plug in and operate damaged equipment.</p>
	<p>WARNING Only authorized personnel may service this equipment. Contact BioPhotas for servicing or repair.</p>
 	<p>WARNING Do not lie upon or lean against the Celluma. The Celluma is not designed to withstand full body weight or extreme pressure. Do not place the Celluma under clothing, blankets, sheets or towels.</p> <p><i>Failure to follow this precaution could result in shock!</i></p>
	<p>WARNING Healthcare providers (doctors, nurses, estheticians, etc.) as well as patients should avoid looking directly into the light, and wear the iGoggles provided or use appropriate eye protection curing treatment.</p>

1. Follow instructions and heed all warnings. Install in accordance with the manufacturer's instructions.
2. **Do NOT lie on top of, or place weight on the device.**
3. **Store the Celluma flat when not in use.**
4. Make sure that the Power Supply is unplugged from wall outlet before cleaning the apparatus.
5. Do not position the device so that it is difficult to un-plug the power supply.
6. Keep Dry. Do not use this device near water.
7. Do not use this device near strong sources of electromagnetic radiation. Device may inadvertently shut off when operated in the presence of a high electromagnetic field.
8. **Do not use volatile solvents (such as toluene, alcohol, etc.) to clean the device. Such solvents may damage the housing and control surface.**
9. Clean with a soft dry cloth.
10. If the housing or control surface is seriously contaminated, use a soft cloth moistened with a **non-alcohol and non-chlorine** disinfectant to clean the device. Wring out cloth before cleaning to prevent any fluid from penetrating into the housing. Wipe the device with a clean, dry cloth after cleaning. Disinfecting wipes containing <1% alkyl dimethyl benzyl ammonium chloride are compatible with the device.
11. Do not install near any heat source such as radiators, heat registers, stoves or apparatus that produce heat.
12. Protect the power cord from being walked on or pinched particularly at plugs, convenience receptacles, or matted with extension cord or the point where they exit from the device.
13. Only use attachments/accessories specified in this manual.

14. Unplug this apparatus during lightning storms or when unused for long periods of time.
15. There are no user-serviceable parts inside this device. Refer all servicing to the manufacturer, BioPhotas, Inc.
16. Servicing is required when the apparatus has been damaged in any way, such as when the power cord is damaged, liquid has been spilled or objects have fallen into the apparatus or the apparatus has been exposed to rain or if the apparatus does not operate normally or has been dropped.
17. Operating and Storage Conditions : 3 ~ 32°C, Relative Humidity: 10~80% (non-condensing)
18. Transport Conditions: -20°C to 50°C, Relative Humidity: 20-85% (non-condensing)
19. **Do not fold, crease or crimp the Celluma beyond its design capabilities.** When shaping the Pro or Elite around a body part - such as an arm - the inward curve should not be smaller than 6" in diameter. When bending the Pro or Elite convexly (lights facing outwards) the curve should not be less than 16" in diameter. When shaping the LITE around a body part – such as an arm – the inward curve should not be smaller than 4" in diameter. When bending the LITE device convexly (lights facing outwards) the curve should not be less than 8" in diameter.

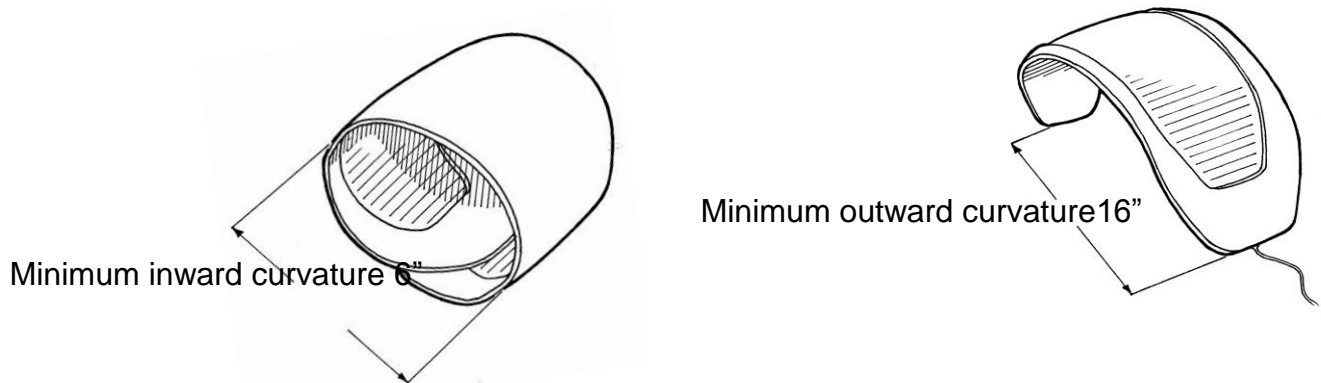***For Celluma LITE Model : Minimum inward curvature 4". Minimum outward curvature 8".**

Preparing the Celluma® for Use:

Inspect for Damage

The Celluma is carefully packaged at the factory to minimize the possibility of damage during shipping. Inspect the box for external signs of damage or mishandling. Inspect the contents for damage. If there is visible damage to the instrument upon receipt, inform the shipping company and Biophotas immediately.

BIOPHOTAS, INC.

**1000 E. Howell Ave., Unit A,
Anaheim, CA 92805 USA
(714) 838 1956**

www.biophotas.com

Celluma® – Protocols for Use:

NOTE: When making changes to your health and wellness routine, it may be advisable to consult with your physician before proceeding.

Modes of Operation

When should I use the Acne mode? Use the Celluma Acne mode at the first signs of acne breakout. The treatment consists of 12 treatment sessions spread over 4 weeks, allowing 3 treatments per week.

When should I use the Wrinkle mode? Use the Celluma Wrinkle mode for the treatment of periorbital wrinkles. For periorbital wrinkles, the treatment consists of 12 treatment sessions spread over 4 weeks, allowing 3 treatments per week.

When should I use “Aches & Pains”? Use the Aches & Pains mode for management and temporary relief of minor muscle and joint pain, stiffness, minor arthritis pain or muscle spasms and to provide a temporary increase in local blood circulation. The treatment consists of 12 treatment sessions spread over 4 weeks, allowing 3 treatments per week.

Product Assembly:

1. Connect the power supply to the Celluma console by inserting the end of the power connector into the side of the console. Use only the provided power supply, Cincon, model **TR30RAM120**.

2. Insert the power supply into an electrical outlet. The perimeter lights on the console will illuminate and circle twice.

Suggested Product Positioning for Celluma PRO:

Suggested Product Positioning for Celluma LITE:

NOTE: To position the Celluma LITE over the face, a pair of Celluma LITE Facial Rests are required. These are not included with the product and may be ordered by calling 714-838-1956 or visiting www.BioPhotas.com

Suggested Product Positioning for Celluma ELITE:

Celluma®: Instructions for Use:


User Preparation

1. Assure that the user is properly positioned and the skin area to be treated is clean and dry. Make-up and other cosmetics should be removed. Many formulas contain minerals which can deflect the light and reduce energy absorption.
2. For best results, the Celluma should be positioned in close proximity to bare skin. Place the Celluma so that the light panel is as close as possible to the skin.
3. Use the iGoggles provided with Celluma.

Select and Begin a Celluma Program:

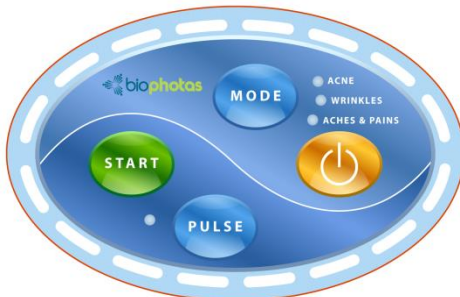1. To begin a Celluma session, press the  button. The device automatically defaults to the Acne program as indicated by the flashing blue light. Then depress the MODE button to select the desired program, if not Acne. Once a program is selected, the corresponding light and perimeter lights will turn on.

2. The Celluma may be used in either Pulse mode or without Pulse mode. Pulse mode causes pulsing lights. Each time the Celluma is turned on it defaults to the Pulse mode as shown by the amber light. To switch from Pulse mode, simply press the PULSE button once. The device will now run without the Pulse mode and the amber light will go off. To return to Pulse mode, press the PULSE button again. Selection of either mode

is based on personal preference. Scientific research has shown that cells absorb light energy better when the light is delivered in pulse mode*. If for any reason you are uncomfortable with the pulsing lights, switch Pulse mode off.

3. With the Celluma in position, press the  button to begin the selected program. The corresponding indicator light will stay lit. Illumination will increase slowly over the first 20 seconds, and remain steady.
4. All programs run for 30 minutes and may be stopped by pressing the  button. At the end of the program, the illumination will slowly decrease and the device will automatically switch off. If a program is interrupted, re-select the desired program and press the  button to begin over. For best results use the Celluma for a full 30 minutes.

* **Importance of pulsing illumination parameters in low level-light therapy.** Journal of Biomedical Optics 15_4_, 048005_July/August 2010_Daniel Barolet et al

RoseLab Skin Optics Laboratory, 3333 Graham Boulevard, Montreal, Quebec, Canada, H3R 3L5 and McGill University, Department of Medicine, Dermatology Division, Montreal, Quebec, Canada, H3A 1A1

Caring for Your Celluma®:

If your Celluma requires spot cleaning, carefully follow the instructions below.

1. Do not use volatile solvents (such as toluene, alcohol, etc.) to clean the device. Such solvents may damage the housing and control surface.
2. Wipe softly with a very slightly water moist soft towel then dry with a dry soft towel. Do not use chemicals or cleaners.
3. If the housing or control surface is seriously contaminated, or to reduce the risk of infection between users, use a soft cloth moistened with a **non-alcohol and non chlorine** disinfectant to clean the device. Wring out cloth before cleaning to prevent any fluid from penetrating into the housing. Wipe the device with a clean, dry cloth after cleaning. Disinfecting wipes containing <1% alkyl dimethyl benzyl ammonium chloride are compatible with the device.

Disposable hygiene barrier sleeves are available in Packs of 100 and may be purchased by contacting BioPhotas Customer Care at 1 (714) 838 1956 or visiting www.biophotas.com to place your order.

Celluma Accessories

The accessories listed below are available for purchase by calling **714-838-1956** or by contacting your authorized Celluma Distributor.

Item	Description	Qty	Cost
COTC-001	Celluma PRO- C/US	1 Unit	Call for pricing
CEL-001	Celluma ELITE-C/US	1 Unit	Call for pricing
CLT-001	Celluma LITE- C/US	1 Unit	Call for pricing
CHB-001	Celluma PRO Disposable Hygiene Barrier Sleeves	1 x Bag (100)	Call for pricing
CHB-002	Celluma LITE Disposable Hygiene Barrier Sleeves	1 x Bag (100)	Call for pricing
CCT-001	Celluma PRO Carrying Tote	1	Call for pricing
CCT-002	Celluma LITE Carrying Tote	1	Call for pricing
CPG-020	Celluma Protective Goggles	1 x Bag (20)	Call for pricing
CVS-001	Celluma Velcro Straps	1 pair	Call for pricing
CLS-001	Celluma LITE Facial Rests	1 pair	Call for pricing

Useful Tips and Frequently Asked Questions:

Q: What is the intended use of the Celluma?

A: The Biophotas Celluma® device is intended to deliver heat in the IR spectrum to provide topical heating for the purpose of elevating tissue temperature; for the temporary relief of minor muscle and joint pain, arthritis and muscle spasm; relieving stiffness; promoting the relaxation of muscle tissue; and to temporarily increase local blood circulation. The blue spectrum is intended to reduce inflammatory acne vulgaris. The Celluma is intended to emit energy in the red and infrared region of the spectrum for use in dermatology for the treatment of periorbital wrinkles.

Q: Does the Celluma utilize UV light?

A: No. UV light has a wavelength less than 400nm.

Q: How long does a Celluma treatment session last?

A: 30 minutes.

Q: Are there any special preparations required before using the Celluma?

The Celluma should be used as close as possible to the skin which should be free of clothing, make-up, or anything else which may deflect the light away from the skin. Do not look directly at the light, and wear the iGoggles provided or use appropriate eye protection.

Q: Why is it recommended to remove make-up?

A: Many make-up formulas contain minerals that may deflect the light. Clean, freshly washed skin is recommended for best results.

Q: How close should the Celluma be placed to the skin?

A: Position the Celluma as close to the surface of the skin as possible. The closer the Celluma is to the skin the more energy is absorbed by the cells. Do not look directly into the light, and protect your eyes with the iGoggles provided with the Celluma, or use other appropriate eye protection.

Q: Can I travel with the Celluma?

A: Yes. Celluma is lightweight and portable making it ideal for travel to or from the office or gym. For travel outside of the U.S. be sure to check for correct power source and electrical input. An international Travel Adapter Kit is provided with the power supply, and contains adapters for Europe, UK and Australia.

Q: What should I do if I cannot use the Celluma?

A: If there is a medical problem, please contact your doctor. If there is a technical problem with the Celluma, please call BioPhotas, Inc. at the number listed in Technical Support (following).

Q: Why do certain rows of diodes look like they are not emitting light?

A: Each Celluma program utilizes 3 different wavelengths of therapeutic light. For example, The Acne program utilizes blue, red and infrared wavelengths but is blue dominant. Both blue and red wavelengths are visible to the human eye. However, infrared is not visible to the human eye. The circuit on the Celluma panel which appears to emit no light is the infrared circuit. **Note:** there is one green light emitting diode in each corner of the Celluma panel.

Technical Support:

If the Celluma will not Power up	<ul style="list-style-type: none"> • Check plug and electrical socket • If the unit still has no power - call the number below and ask for support.
If Celluma does not respond to controls	<ul style="list-style-type: none"> • First turn the main power switch off • Wait five seconds • Then switch back on • If the problem persists call the number below and ask for support.
For additional issues	<ul style="list-style-type: none"> • Call the number below and ask for support.

Disposal:

WEEE Directive Symbol

This symbol on the product or on its packaging indicates that this product must not be disposed of with regular waste. Instead, it is the user responsibility to dispose of waste equipment according to the local laws. For information about where the user can drop off the waste equipment for recycling, please contact your local BioPhotas, Inc. representative.

Technical Information:**Celluma Output Specifications**

The Celluma has two output modes: Continuous and Pulsed

Red LED (640 nm): max output 71.6 microwatts/sr/m² (continuous)

Blue LED (464 nm): max output 76.4 microwatts/sr/m² (continuous)

Infrared LED (880 nm): max output 47.3 microwatts/sr/m² (continuous)

In the Pulsed mode, the LEDs cycle from 1-90% duty cycle using a proprietary algorithm. The pulse duration, repetition rate, pulse interval and pulse trains vary continuously during the treatment to achieve the therapeutic effect. In pulse mode, the maximum output is a percentage of the continuous mode as noted below:

Red Light Mode	90.75%
Acne Mode	99.73%
Aches & Pain Mode	88.09%

The Celluma is classified “Exempt” in accordance with IEC 60601-2-57:2012

Electromagnetic Emissions Guidance

CAUTION: This device meets the electromagnetic emissions requirements per IEC 60601-1-2. The Celluma was not tested for immunity to electromagnetic disturbances. If the user observes unusual device behavior, particularly if such behavior is intermittent and associated with nearby use of radio or TV transmitters, cell phones, or electro-surgical equipment, this could be an indication of electromagnetic-interference. If such behavior occurs, the user should try to move the interfering equipment further from this device.

WARNING: Use of accessories other than those specified or provided by with the Celluma could result in increased electromagnetic emissions and result in improper operation.

WARNING: The Celluma should not be used adjacent to or stacked with other electrical equipment. If adjacent or stacked use is necessary, observe the Celluma and the other electrical equipment to make sure they operate as intended. Make sure the Celluma operates correctly when it is used near other electronic devices. Portable and mobile radio frequency (RF) communications equipment can affect electrical equipment. Medical equipment needs special precautions in regard to electromagnetic compatibility (EMC) and needs to be installed and put into service according to the EMC information supplied in the table below.

Guidance and Manufacturer's Declaration—Electromagnetic Emissions

The Celluma is intended for use in the electromagnetic environment specified below. The customer or the user of the Celluma should make sure it is used in such an environment

Emissions Test	Compliance	Electromagnetic Environment—Guidance
RF emissions CISPR 11	Group 1	The Celluma uses RF energy only for its internal functions. Therefore, its RF emissions are low and are not likely to cause any interference in nearby electronic equipment The Celluma is suitable for use in all establishments, including domestic establishments and those directly connected to the public low-voltage power supply network that supplies buildings used for domestic purposes.
RF Emissions CISPR 11	Class B	
Harmonic Emissions IEC 61000-3-2	Not Applicable	
Voltage Fluctuations/ Flicker Emissions IEC 61000-3-3	Complies	

Guidance and Manufacturer's Declaration—Electromagnetic Immunity

The Celluma is intended for use in the electromagnetic environment specified below. The customer or the user of the Celluma should make sure it is used in such an environment

Immunity Test	IEC 60601 test level	Compliance level	Electromagnetic Environment—Guidance
Electrostatic Discharge (ESD) IEC 61000-4-2	±6 kV Contact ±8 kV Air	±6 kV Contact ±8 kV Air	Floors should be wood, concrete, or ceramic tile. If floors are covered with synthetic material, the relative humidity should be at least 30%.
Electrical Fast Transient/Burst IEC 61000-4-4	±2 kV for Power Supply Lines ±1 kV for Input/Output Lines	±2 kV for Power Supply Lines ±1 kV for Input/Output Lines	Mains power quality should be that of a typical commercial or hospital environment
Surge IEC 61000-4-5	± 1 kV Line(s) to Line(s) ±2 KV Line(s) to Earth	± 1 kV Line(s) to Line(s) ±2 KV Line(s) to Earth	Mains power quality should be that of a typical commercial or hospital environment
Voltage Dips, Short Interruptions, and Variations on Power Supply Lines IEC61000-4-11	<5 % 230V (>95 % dip in 230V) for 250 cycles 40 % 230V	<5 % 230V (>95 % dip in 230V) for 250 cycles 40 % 230V	Mains power quality should be that of a typical commercial or hospital environment. If the user of the Celluma requires

	(60% dip in 230V) For 5 cycles 70 % 230V (30 % dip in 230V) for 25 cycles <5 % 230V (>95 % dip in 230V for 5 sec	(60% dip in 230V) For 5 cycles 70 % 230V (30 % dip in 230V) for 25 cycles <5 % 230V (>95 % dip in 230V for 5 sec	continued operation during power mains interruption, it is recommended that the Celluma be powered from an uninterruptable power supply or a battery.
Power Frequency (50/60 hz) Magnetic Fields IEC 61000-4-8	3A/m	3A/m	Power frequency magnetic fields should be at levels characteristic of a typical location in a typical commercial or hospital environment.
Note 1: U_T is the AC mains voltage prior to application of the test level.			

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