

GLO Science  
**PRO POWER**  
At Home Teeth Whitening Kit



USER MANUAL  
Bluetooth-Enabled Smileware

# User Manual Contents

Product Components.....p 03

Illuminating Heat Technology.....p 04

GLO Vials .....p 05

Device Operation .....p 06

Quick Start Guide.....p 08

App Set Up .....p 10

FAQs .....p 12

Cleaning and Care.....p 14

Ingredients.....p 15

Warranty .....p 16

Electrical Specifications.....p 17


GLO Good.....p 18

## General Warnings and Cautions

Please read this entire manual before using this product. Not doing so may result in an unsafe condition and/or an ineffective experience.

# This is Whitening Made Modern

## PRO POWER At Home Teeth Whitening Kit Contents

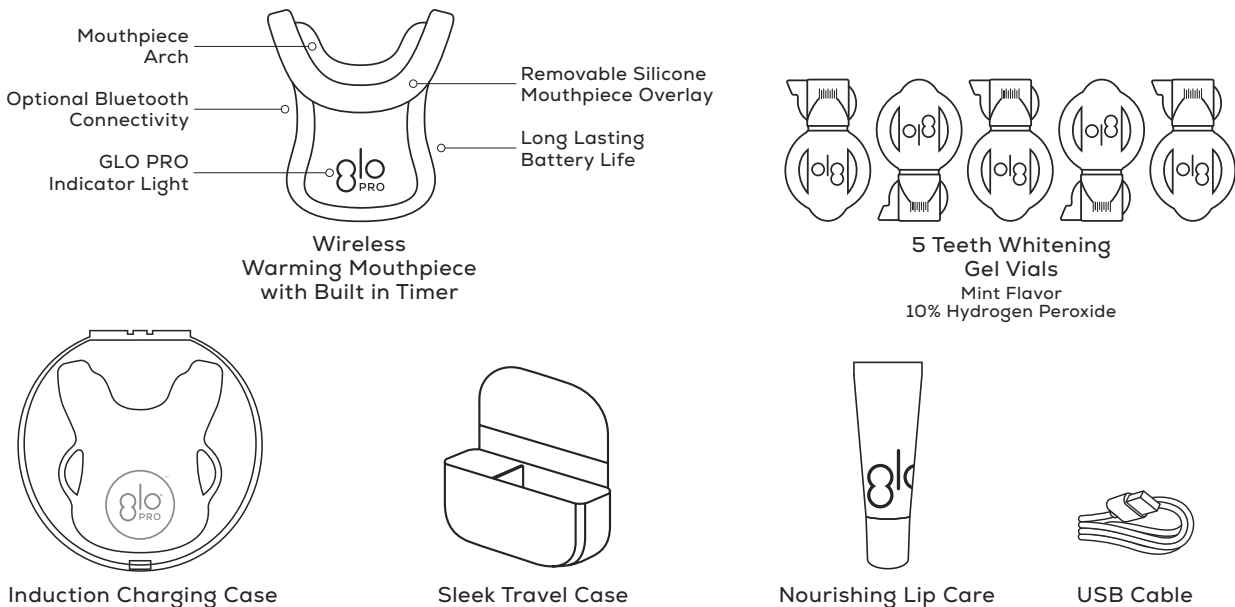

## We're Heating Things Up


Congratulations! You made a smart decision when it comes to safely and effectively whitening your teeth at home.

We all know that LED lights look fun, but they simply aren't functional when it comes to whitening your teeth. You need **warmth** to get the whitest smile in the shortest amount of time.

Our exclusive **Illuminating Heat Technology™** is key to delivering the best and safest whitening results at home - without sensitivity. The GLO PRO POWER Teeth Whitening Device produces the perfect amount of warmth to amplify our clean, specially formulated pro-strength whitening gel.

Often imitated. Never duplicated. It's the dental professional recommended **best** and **safest** way to whiten your teeth and maintain your results at home.

Available only from your dentist.

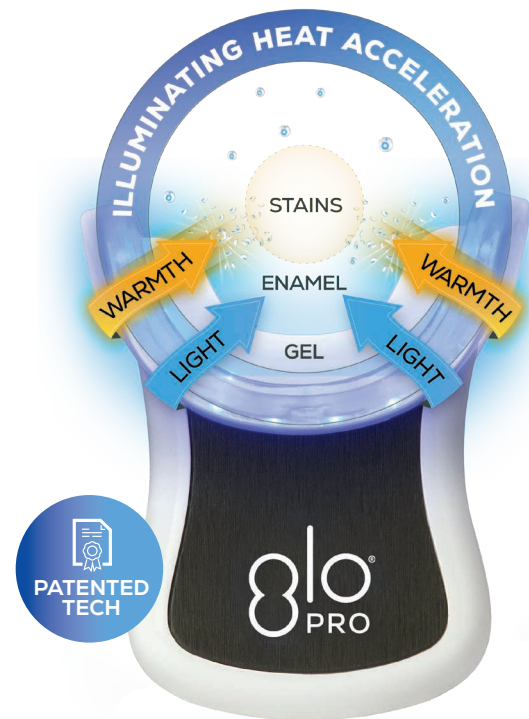

## And Keeping It Clean

Whole body health begins with what we put in our mouths so it's important to choose carefully. That's why we take special care formulating our whitening gel using only clean ingredients that are proven safe and effective.

In fact, we designed our recyclable GLO Vials to ensure optimal freshness, efficacy and ease of application for our clinically proven, hydrogen peroxide based whitening formulations.

### DID YOU KNOW?

*Decades of clinical experience and research confirm hydrogen peroxide as the gold standard active ingredient for safe and effective teeth whitening. Still, many other companies choose to use inferior active ingredients like carbamide peroxide because they are less expensive and don't require specially designed packaging to maintain their efficacy. Unlike hydrogen peroxide which dissolves cleanly into oxygen and water, lesser ingredients can leave you with potentially harmful breakdown byproducts like urea and ammonia in your mouth.*

Pro-strength for home use available only from your dentist.




SAFE FOR CROWNS, BRIDGES, AND VENEERS

## Device Operation




The induction charging case provides the perfect place to store your device while keeping it charged so it is always ready for use. The device should **ONLY** charge in the charging case provided in your kit.

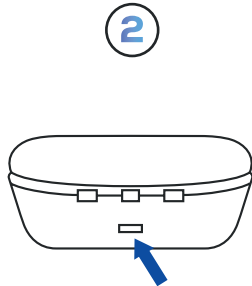


Attach small end of USB cable to back of charging case. Attach other end of cable into live active USB outlet.




Place device into charging case. GLO PRO indicator light will begin to blink. Device is fully charged when GLO PRO indicator light turns off (~5 hours).

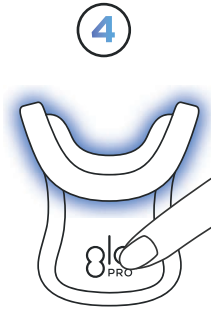

Hold finger on the GLO PRO indicator light until unit powers ON (indicator light will be white and will vibrate). Mouthpiece will not illuminate but device will be powered ON.

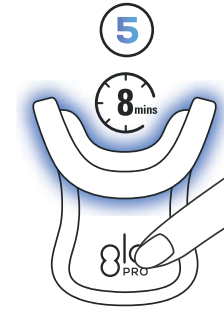

Tap GLO PRO indicator light once again to begin 8 minute session (mouthpiece will illuminate to indicate active session).

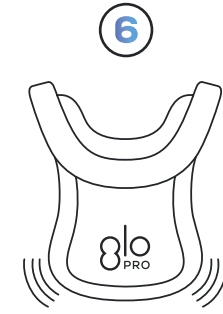

When the full 8 minute session is almost complete, device will vibrate and indicator light will begin blinking before mouthpiece lights shut off, indicating session is complete. The device will remain powered ON until GLO PRO indicator is tapped again or until automatic shut off after 90 seconds of inactivity (or after 15 minutes of inactivity when paired with App).

## Powering On & Off

### Powering ON

**TAP & HOLD ONCE** to turn device **ON**.

**TAP ONCE AGAIN** to illuminate mouthpiece and begin timed 8 minute session.

### Powering OFF


**TAP ONCE** to turn mouthpiece lights **OFF**. Device will vibrate and mouthpiece lights will turn off automatically when 8 minute whitening session is complete.

**TAP ONCE AGAIN** to turn device **OFF**. Device will shut off automatically after 90 seconds of inactivity (or after 15 minutes of inactivity when paired with App).


## GLO PRO Indicator Light Key

- **BLINKING** indicator light indicates **CHARGING**
- **SOLID** indicator light indicates device is **ON**
- **BLUE** indicator light indicates device is **PAIRED** with the APP
- **RED BLINKING** indicator light indicates battery needs to **RECHARGE**


## Quick Start Guide To Your 16 Minute Smile Makeover

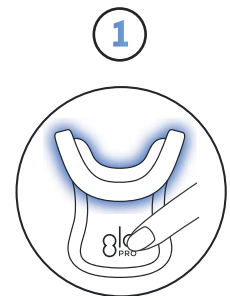


**Power ON**  
charged device.

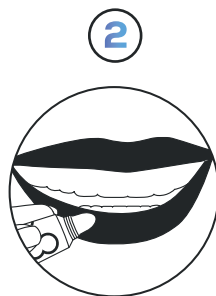


**Squeeze**  
small bead of gel  
(~1/2 pea size) onto  
brush tip applicator  
after removing  
GLO Vial cap.

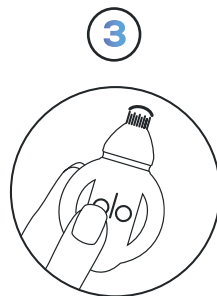

**Remove**  
GLO Vial cap and  
squeeze small  
bead of gel  
(~1/2 pea size)  
onto brush tip  
applicator.

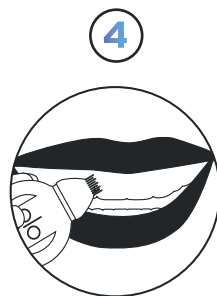

**Brush**  
thin layer of  
whitening gel on  
top and bottom  
teeth, squeezing up  
fresh gel for each  
arch and ensuring  
gel contact with all  
visible tooth  
surfaces.

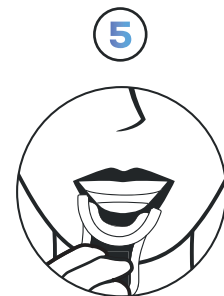

**Insert**  
mouthpiece into  
mouth and tap GLO  
PRO indicator light to  
illuminate and begin  
8 minute session -  
maintain a relaxed jaw  
and swallow as usual.

**REPEAT 2X/day (16 minutes total) for 4 days**

## Tips For Best Results




Brush and floss your  
teeth before you GLO.  
Whitening works best on  
a clean surface.



Wait 20 minutes after  
final application before  
eating and drinking.



Brush twice daily, topping your  
favorite toothpaste with ECO  
Balance to keep your smile  
healthy, happy, and bright.

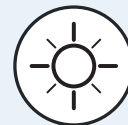

For full whitening treatment,  
complete two 8 minute applications  
in a row (16 minutes total) each  
day for four days. Your teeth will  
get whiter with each application.  
Depending on your desired results,  
you may continue whitening for  
up to four additional days.



For maintenance, one or  
two 8 minute applications  
(8-16 minutes total) each  
month will keep your  
smile stain-free and  
shining bright.




Get a complete in office teeth  
whitening treatment at  
your dentist once a year and  
a touch up after every  
hygiene visit.




Keep a GLO Vial handy to use  
after consuming foods or drinks  
that can stain your teeth. Simply  
brush a layer of gel on your  
teeth immediately after eating  
or drinking to instantly remove  
stains before they set and  
freshen your breath.

## App Set Up

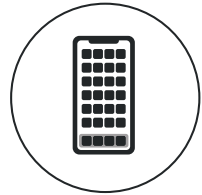


Download the GLO Whitening App from Google Play or the App Store.

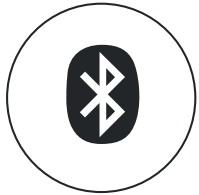

Turn your phone's Bluetooth ON.  
(Settings > Bluetooth)



Hold your finger on the GLO PRO indicator light until logo lights up WHITE indicating unit is powered ON.

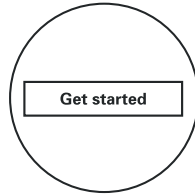


Open the App and select 'Get started' to begin device set up.

Select your device from the dropdown menu.

Select 'Continue.'



Pairing is complete when the GLO PRO indicator light turns BLUE.

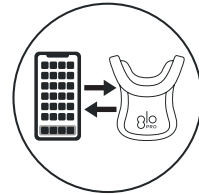

Your phone must remain physically close enough to your teeth whitening device to maintain Bluetooth connectivity.

## Your Personal Whitening Assistant

Optional Bluetooth connectivity to the GLO Whitening App is designed to enhance your whitening experience and includes these additional features:



Schedule treatments & get reminders



Preview results & track progress



Battery power indicator



Visual countdown timer



Pause session option



Restart session option

## FAQs

### RESULTS & EXPECTATIONS

#### How white can I get my smile and how long will it take?

Your teeth will get whiter with each application, so the more you GLO, the whiter your teeth will get. Keep in mind that everyone's smile is unique and there are many factors involved including age, genetics, and lifestyle, so some people's teeth will whiten faster and some will take more time.

#### Will my teeth look too white or unnatural?

No – you are in complete control of how white your teeth get so you can achieve your perfect share of white.

#### What if I want even faster results?

The fastest GLO whitening is an in office whitening treatment at your dental practice, which gets you your whitest white in only 16 minutes. Ask your dentist today!

### SPECIAL CIRCUMSTANCES

#### Can I use the GLO PRO POWER device with braces?

We don't recommend using the device while you have braces but you can use GLO Vials on their own to whiten during orthodontic treatment. The device is perfect to use as soon as you get your braces off for a straight, white smile!

#### Will GLO work on crowns, bridges, veneers and other dental work?

Yes, GLO will safely lift stains and return any unnatural tooth surfaces to their original shade.

### WHITENING GEL

#### Can I swallow when the gel is on my teeth?

Yes! Swallow as you normally would when the gel is on your teeth and/or the mouthpiece is in your mouth. The whitening gel is designed to stay put on your teeth until it evaporates and dissolves.

#### Do I need to apply the GLO Lip Care before each session?

We recommend applying the lip care to keep your lips nourished and hydrated while whitening.

#### What if my teeth have been sensitive to other whitening products?

GLO whitening gel is specially formulated to be sensitivity free. It is free of alcohol and carbamide peroxide and is designed to stay put on teeth and off gums and soft tissue.

### DEVICE

#### Do I bite down on the mouthpiece when it is in my mouth?

Yes, bite down gently to keep the mouthpiece in place and relax your lips and jaw.

#### How long will the GLO PRO POWER Device last when fully charged?

When fully charged, your device will have sufficient power for at least 16 applications (two complete 4 day treatments of 2 applications per day).

#### What do I do if my mouthpiece will not light up?

Make sure your device has been fully charged in an active charging case that was plugged into an active power source. Hold your finger on the GLO PRO indicator light for a second or two until the light illuminates and then tap again to illuminate the arch. If you require additional assistance, please contact us at [info@gloscience.com](mailto:info@gloscience.com)

#### Will the light from the mouthpiece affect my skin?

No, there is no UV light component. The light is 470 nm wavelength (blue) light generated by LEDs inside the mouthpiece which indicate that your 8 minute session is active.

#### Is GLO registered with the FDA?

Yes, the GLO PRO POWER At Home Teeth Whitening Device is a Class 1 Medical Device and has been registered with the FDA.

**We're here to help!**  
Please contact us at [info@gloscience.com](mailto:info@gloscience.com)

## Cleaning & Care



### Silicone Mouthpiece Overlay

- After each application, wipe any gel residue off silicone mouthpiece overlay with clean, dry cloth or tissue.
- After you have completed your consecutive applications for the day, remove silicone mouthpiece overlay from device mouthpiece arch.
- Rinse overlay and wipe or drip dry.
- Do not replace mouthpiece overlay onto device until both are completely dry, inside and outside.



### Device

- Device can be wiped clean with a damp cloth.
- Dry before returning to charging case.
- Avoid wetting the circuit on the mouthpiece arch.
- **DO NOT SUBMERGE THE DEVICE IN ANY LIQUID.**



### Charging Case

- Charging case can be wiped clean with a damp cloth.
- Avoid wetting the port on the back of the case.
- **DO NOT SUBMERGE THE CASE IN ANY LIQUID.**

### Maintenance & Calibration

Other than cleaning and charging, the PRO POWER device and charging case do not require any maintenance or calibration.

## Ingredients



### GLO Teeth Whitening Gel Vials

Our whitening gel is specially formulated to be safe, effective, and sensitivity-free! It is also vegan, gluten-free, alcohol-free and does not contain carbamide peroxide. Plus, it features a refreshing mint flavor you will love!

**Ingredients:** Distilled Water, Polyethylene Glycol, Hydrogen Peroxide 10%, Potassium Nitrate, Glycerin, Edetate Disodium, Sodium Saccharin, Sodium Fluoride, Sodium Hydroxide, Flavor.



### GLO Lip Care

You are in for a treat! GLO Lip Care features a superior blend of vitamins and antioxidants to moisturize, nourish, repair, and protect your lips. Featuring African shea butter, aloe, COQ10, and vitamin E to leave your lips soft and smooth. Perfect for everyday use and before any whitening treatment.

**Ingredients:** Ricinus Communis (Castor) Seed Oil, Beeswax, Mentha Piperita (Peppermint) Leaf, Stearalkonium Hectorite, Butyrospermum Parkii (Shea) Butter\*, Vanilla Planifolia Fruit Extract, Menthol, Alcohol, Aloe Barbadensis Leaf Juice\*, Copernicia Cerifera (Carnauba) Wax, Tocopherol, Ammonium Glycyrrhizate, Ubiquinone. \*Organic

**Warnings:** If irritation (such as redness, swelling, soreness) of the gums or the mouth occurs, discontinue use and consult a dentist. Products containing peroxides are not recommended for use by children under 12 years of age or by pregnant or breastfeeding women. Use for periods of longer than 14 days is to be only under the supervision of a dentist. Avoid swallowing the cosmetic or part thereof. Avoid contact of the product with the eye. Avoid direct contact of the active surface of the tooth whitening product with the gums and/or salivary flow.

**Cautions:** Consult with a dental professional before use if you have untreated dental cavities or any form of gum disease. If gum irritation occurs, stop using the product and contact your dental professional. Avoid getting the gel on fabric or clothing as it may lead to discoloration.

GLO Teeth Whitening Gel Vials are Made in the USA.

GLO Lip Care is Made in the USA.

GLO PRO POWER Device is Designed in USA and Made In China.

GLO Science, Inc.  
10 W. 37th Street  
New York, NY 10018  
855-GLO-5976 (855-456-5976)  
www.gloscience.com  
info@gloscience.com



# One Year Limited Warranty

GLO Science™ Inc. commits to the purchaser that we will replace your GLO Science Device and/or components within one year from the date of purchase, for any manufacturing defect in the PRO POWER Device.

To report defective product and request an exchange, please call 855 7GLOPRO (855 745 6776). If we determine that your device should be returned to GLO for inspection, we will provide you with a prepaid return shipping label.

If the product is a model no longer being made and is unavailable for exchange, we will happily replace the product with the most recent model.

**The 1 Year Limited Warranty does not cover any of the following:**

- Defects, damage or other conditions caused, in whole or in part, by misuse, abuse, neglect, negligence, alteration, accident, freight damage, improper repair, improper maintenance, improper, abnormal or excessive use and/or failure to obtain a replacement in a timely manner, as prescribed by GLO Science, Inc.
- Defects or damages that are a result of an act of nature (such as a flood or earthquake) or acts of war or terrorism.
- Damage resulting from improper use of parts or accessories that are not compatible with device and/or components, including but not limited to the breakage or fracture of the device and/or components by any reason other than by manufacturing defect.

- Ordinary wear and tear on the device and/or components.
- Defects or damages as a result of any modification of device and/or components.
- Cosmetic damage, including scratches or dents.
- Accessories or parts that are not those of GLO Science.
- Products purchased, new or used, from unauthorized distributors or through unauthorized distribution channels.
- Damage resulting from use of the device with a voltage source other than the voltage source supplied by GLO Science.
- GLO Science's only obligation under this warranty is the replacement of defective material that comprises the PRO POWER Device. GLO Science shall not be liable and hereby disclaims any direct, special, indirect , incidental, exemplary or consequential damages or delays or losses under any other legal theory, including but not limited to, loss of use, loss of revenue, or loss of profits.

**Warranty Disclaimer**

This warranty is GLO Science's only warranty and is in lieu of all other warranties, express or implied. GLO Science makes no implied warranties of any kind, including any implied warranties or merchantability or fitness for a particular purpose.


Register your device at [gloscience.com/register](https://gloscience.com/register)

**Electrical Specifications:**

NOTE: This equipment has been tested and found to comply with the limits for a Class B digital device, pursuant to Part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in residential installation. This equipment generates, uses and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation. If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures: Reorient or relocate the receiving antenna. Increase the separation between the equipment and receiver. Connect the equipment into an outlet on a circuit different from that to which the receiver is connected. Consult the dealer or an experienced radio/TV technician for help.

Input: 100-240 VAC, 50-60 Hz, 0.2A  
Output: 5 VDC, 1A  
Light output: 10/0.085W 470 nm Blue LEDs

Recommended storage conditions: 5° - 25° C




 **Battery Disposal Instructions:** Batteries should be recycled or disposed of as per state and local guidelines.

GLO® (stylized) and GLO Science™ are Registered Trademarks of GLO Science™, Inc. PRO POWER, PRO POWER+, Smileware, The Professional Whitening Experts, Illuminating Heat Technology, We All Speak Smile, GLO Vial, G Vial, GLO-ing, G.L.O and Guided Light Optics are Trademarks/Service marks of GLO Science, Inc. All rights reserved.

Patents			
US D636,074	US D765,255	US 9,827,174	US D863,590
US 8,029,278	US 9,572,645	US 8,267,609	US D888,961
US 8,398,324	US D795,501	US 10,201,408	US 10,918,470

## GLO For Good

**Our mission is simple** – to empower more people with healthy, confident smiles. To be able to smile freely and without pain or shame makes it possible to laugh, to sing, and to connect with others. Everyone deserves to share their smiles with the world.

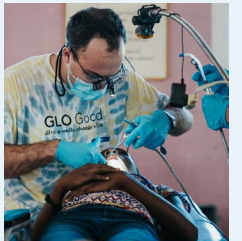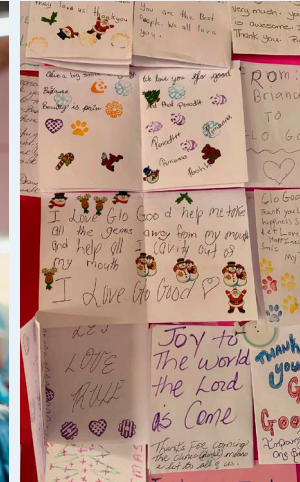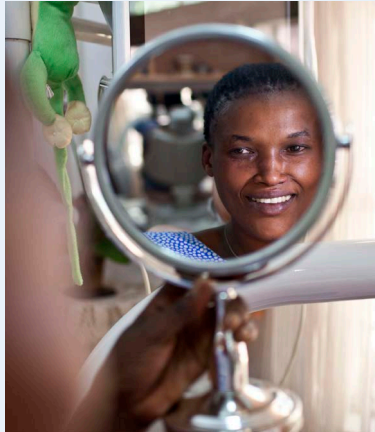

Your purchase helps provide free dental care to those without access through the GLO Good Foundation.



**GLO SCIENCE®**  
professional whitening

**TIME TO SMILE**

ITEM NUMBER  
658113

**[gloscience.com/glopro](https://gloscience.com/glopro)**  
[info@gloscience.com](mailto:info@gloscience.com)

06102022