



COVIDIEN

User's Guide

ForceTriad™

Energy Platform

User's Guide

ForceTriad™

Energy Platform

For use with software version 3.6x

Part Number: 1063547

Preface

This guide and the equipment it describes are for use only by qualified medical professionals trained in the particular technique and surgical procedure to be performed. It is intended as a guide for using the Covidien ForceTriad Energy Platform only. Additional technical information is available in the *ForceTriad Energy Platform Service Manual*.

Equipment covered in this manual

ForceTriad Energy Platform with software version 3.6x

Conventions Used in this Guide

Warning

Indicates a potentially hazardous situation which, if not avoided, could result in death or serious injury.

Caution

Indicates a hazardous situation which, if not avoided, may result in minor or moderate injury.

Notice

Indicates a hazard which may result in product damage.

Important

Indicates an operating tip or maintenance suggestion.

Limited Warranty

Covidien warrants each covered product listed below to be free from defects in material and workmanship for normal use and service for the period(s) set forth below. Covidien's obligation under this warranty is limited to the repair or replacement, at its sole option, of any product, or part thereof, which has been returned to it (or its authorized distributor) within the applicable time period shown below after delivery of the product to the original purchaser, and which examination discloses, to Covidien's satisfaction, that the product is defective. This limited warranty does not apply to any product, or part thereof, which has been repaired or altered in a way so as, in Covidien's judgment, to affect its stability or reliability, or which has been subjected to misuse, neglect, or accident.

The warranty periods for Covidien products are as follows:

ForceTriad™ Energy Platform	One year from date of shipment
Electrosurgical Generators	One year from date of shipment
Cool-tip™ RFA Generator	One year from date of shipment
Evident™ MWA Generator	One year from date of shipment
LigaSure™ Vessel Sealing System	One year from date of shipment
LigaSure™ Reusable Instruments	One year from date of shipment
Mounting Fixtures (all models)	One year from date of shipment
Footswitches (all models)	One year from date of shipment
Valleylab™ Argon Gas Delivery Unit II	One year from date of shipment
RapidVac™ Smoke Evacuator	One year from date of shipment
LigaSure™ Sterile Single Use Items	Sterility only as stated on packaging
Cool-tip™ Sterile Single Use Items	Sterility only as stated on packaging
Sterile Single Use Items	Sterility only as stated on packaging
Patient Return Electrodes	Shelf life only as stated on packaging
All purchased or supplemental software programs or updates	180 days from delivery

Notwithstanding any other provision herein or in any other document or communication, Covidien's liability with respect to this limited warranty and the products sold hereunder shall be limited to the aggregate purchase price for the products sold to the customer. This limited warranty is non-transferable and runs only to the original purchaser of the covered

Limited Warranty

product(s). There are no warranties which extend beyond the terms hereof. Covidien disclaims any liability hereunder or elsewhere in connection with the sale of products and for any form of indirect, tort, or consequential damages.

This limited warranty and the rights and obligations hereunder shall be construed under and governed by the laws of the State of Colorado, USA. The sole forum for resolving disputes arising under or relating in any way to this limited warranty is the District Court of the County of Boulder, State of Colorado, USA.

Covidien reserves the right to make changes in covered products built or sold by it at any time without incurring any obligation to make the same or similar changes to equipment previously built or sold by it.

THE OBLIGATION TO REPAIR OR REPLACE A DEFECTIVE OR NONPERFORMING PRODUCT IS THE SOLE REMEDY OF THE CUSTOMER UNDER THIS LIMITED WARRANTY. EXCEPT AS EXPRESSLY PROVIDED HEREIN, COVIDIEN DISCLAIMS ALL OTHER WARRANTIES, WHETHER EXPRESS OR IMPLIED, ORAL OR WRITTEN, WITH RESPECT TO PRODUCTS, INCLUDING WITHOUT LIMITATION ALL IMPLIED WARRANTIES, WARRANTIES OF MERCHANTABILITY OR FITNESS FOR A PARTICULAR PURPOSE.

Software License

COVIDIEN AG, Tyco Healthcare Group LP, and its affiliate EbD (collectively called "COVIDIEN" herein) own the entire right, title, and interest in and to all of the computer programs and all portions thereof, and associated documentation (collectively, the "Software") provided to Customer as may be installed in the Products and equipment addressed herein or provided separately, and it has the sole right to grant licenses hereunder.

The evaluation allowance herein and any ultimate price paid by Customer for the products incorporating the Software include as a portion of that evaluation allowance, or price, a license fee granting Customer only the rights set forth in this Software License. Customer further acknowledges and agrees that the Software is owned exclusively by COVIDIEN. The Software is licensed to be used on only one computing device or Product, and a valid license must be purchased for each computing device on which the Software is installed.

Single User License Grant: COVIDIEN grants to Customer a limited, nonexclusive, non-sublicensable, nontransferable and revocable license to use the Software, exclusively at Customer's location as identified by Customer as the ship-to location of the Product, solely in machine-readable object code form only on a single central processing unit owned or leased by Customer or otherwise embedded in equipment provided by COVIDIEN, and for the sole purpose of Customer's internal business purpose in the operation of the Product or equipment purchased from, other otherwise provided by, COVIDIEN or its affiliates.

Except to the extent expressly authorized in this Software License or by law, Customer shall not and shall not cause any third party to: (i) decompile, disassemble, or reverse engineer the Software; (ii) modify or create any derivative works (including, without limitation, translations, transformations, adaptations or other recast or altered versions) based on the Software, or alter the Software in any way; (iii) merge the Software with any other software or product not supplied by Supplier; (iv) use, copy, sell, sublicense, lease, rent, loan, assign, convey or otherwise transfer the Software except as expressly authorized by the Agreement; (v) distribute, disclose or allow use of the Software, in any format, through any timesharing service, service bureau, network or by any other means, to or by any third parties; (vi) remove or modify any copyright, confidential and/or proprietary markings, legends or restriction which are in the Software originally supplied to Customer; or (vii) violate any obligations with regard to COVIDIEN's Confidential Information. To the extent that Customer is expressly permitted by applicable mandatory law to undertake any of the activities listed in the preceding sentence, Customer will not exercise those rights until Customer has given COVIDIEN thirty (30) days written notice of Customer's intent to exercise any such rights unless an order of a government agency of competent jurisdiction will not so allow.

Except for the limited license rights expressly granted in this Software License, COVIDIEN reserves all rights in and to the Software and any modifications thereto and derivations thereof, including, but not limited to, all title, ownership, intellectual property rights and all other rights and interests. Customer will own only the hardware or physical media on which the Software is stored or processed, if any.

Customer agrees that the Software, including the specific design and structure of individual programs, constitute confidential information and trade secrets of COVIDIEN, whether or not the programs may be copyrighted or copyrightable, and/or patented or patentable. Customer agrees not to disclose, provide, or otherwise make available such confidential information, trade secrets or copyrighted material in any form to any third party. Customer agrees that it will make the Software available only to employees, contractors, or consultants with a need to know, who are obligated to comply with all license restrictions contained in this Software License Agreement and to maintain the secrecy of the Software and all other Confidential Information. Customer is responsible for the compliance of all users with these obligations.

Customer may, from time to time, request that COVIDIEN incorporate certain features, enhancements or modifications into the Software. COVIDIEN may, in its sole discretion, undertake to incorporate such changes and distribute the Software so modified to all or any of COVIDIEN's customers. All such error corrections, bug fixes, patches, updates or other modifications provided to COVIDIEN shall be the sole property of COVIDIEN.

This Software License is effective until terminated. Customer may terminate this License at any time by destroying all copies of Software including any documentation. This License will terminate immediately upon notice from COVIDIEN if Customer fails to comply with any provision of this License or any supplier agreement. COVIDIEN may terminate the Software licenses granted herein and exercise all available rights by giving written notice, effective immediately, if within ten (10) business days of Customer's receipt of a reasonably detailed written request to cure, Customer has not cured all breaches of this License's limitations or restrictions. Upon such termination, Customer will immediately pay all undisputed fees outstanding, cease use of all Software, return or delete, at COVIDIEN's request, all copies of the Software in Customer's possession, and certify compliance with all of the obligations herein to COVIDIEN in writing.

Limited Warranty: COVIDIEN represents and warrants to Customer that the Software will perform substantially as described in COVIDIEN's then current documentation for such Software for the longer of (a) the remaining warranty applicable to the product with which such Software was delivered (not to exceed one year) or (b) ninety (90) days from the date such Software was shipped or first made available to Customer for electronic download from COVIDIEN's service site. If you notify COVIDIEN of defects during the warranty period, COVIDIEN will replace the Software or, at its option, refund the purchase price. Your remedy for breach of this limited warranty shall be limited to replacement or refund and shall not encompass any other damages. No dealer, distributor, agent or employee of COVIDIEN is authorized to make any modification or addition to the warranty and remedies stated above.

Notwithstanding these warranty provisions, all of COVIDIEN's obligations with respect to such warranties shall be contingent on Customer's use of the Software in accordance with this Agreement and in accordance with COVIDIEN's instructions as provided by COVIDIEN in the documentation, as such instructions may be amended, supplemented, or modified by COVIDIEN from time to time. COVIDIEN shall have no warranty obligations with respect to any failures of the Software which are the result of accident, abuse, misapplication, extreme power surge or extreme electromagnetic field.

This warranty does not apply to any damages, malfunctions, or non-conformities caused to or by: (i) Customer's use of Software in violation of the license granted under the Agreement or in a manner inconsistent with any provided documentation; (ii) use of non-

COVIDIEN furnished equipment, software, or facilities with its equipment or Products; (iii) Customer's failure to follow COVIDIEN's installation, operation, repair or maintenance instructions; (iv) Customer's failure to permit COVIDIEN timely access, remote or otherwise, to Products; (v) failure to implement all new Updates to Software provided under the Agreement; (vi) Products or equipment with their original manufacturer's serial numbers altered, defaced or deleted; (vii) Products or equipment that have been altered, serviced or modified by a party other than COVIDIEN; or (viii) Software that has been subjected to abnormal physical or electrical stress, misuse, negligence or accident by Customer or a third party.

DISCLAIMER: EXCEPT AS SPECIFIED IN THIS WARRANTY, ALL EXPRESS OR IMPLIED CONDITIONS, REPRESENTATIONS, AND WARRANTIES INCLUDING, WITHOUT LIMITATION, ANY IMPLIED WARRANTY OF MERCHANTABILITY, FITNESS FOR A PARTICULAR PURPOSE, OR ARISING FROM A COURSE OF DEALING, USAGE, OR TRADE PRACTICE, ARE HEREBY EXCLUDED TO THE EXTENT ALLOWED BY APPLICABLE LAW.

IN NO EVENT WILL EITHER PARTY BE LIABLE FOR ANY LOST REVENUE, PROFIT, OR DATA, OR FOR SPECIAL, INDIRECT, CONSEQUENTIAL, INCIDENTAL, OR PUNITIVE DAMAGES HOWEVER CAUSED AND REGARDLESS OF THE THEORY OF LIABILITY ARISING OUT OF THIS SOFTWARE LICENSE EVEN IF SUCH PARTY HAS BEEN ADVISED OF THE POSSIBILITY OF SUCH DAMAGES. IN NO EVENT SHALL ONE PARTY'S LIABILITY TO THE OTHER PARTY, WHETHER IN CONTRACT, TORT (INCLUDING NEGLIGENCE), OR OTHERWISE, EXCEED THE PRICE PAID OR TO HAVE BEEN PAID BY CUSTOMER. THE FOREGOING LIMITATIONS SHALL APPLY EVEN IF THE ABOVE-STATED WARRANTY FAILS OF ITS ESSENTIAL PURPOSE. SOME STATES DO NOT ALLOW LIMITATION OR EXCLUSION OF LIABILITY FOR CONSEQUENTIAL OR INCIDENTAL DAMAGES.

U.S. Government Rights. The Software is a "commercial item" developed exclusively at private expense, consisting of "commercial computer software" and "commercial computer software documentation" as such terms are defined or used in the applicable U.S. acquisition regulations. The Software is licensed hereunder (i) only as a commercial item and (ii) with only those rights as are granted to all other customers pursuant to the terms and conditions of this License. Customer shall not use, duplicate, or disclose the Software in any way not specifically permitted by this License. Nothing in this License requires COVIDIEN to produce or furnish technical data for or to Customer.

If any provision of this Agreement shall be held by a court of competent jurisdiction to be illegal, invalid or unenforceable, the remaining provisions shall remain in full force and effect.

This License Agreement contains the entire understanding and agreement between the parties respecting the Software. This Agreement may not be supplemented, modified, amended, released or discharged except by an instrument in writing signed by each party's duly authorized representative. All captions and headings in this Agreement are for purposes of convenience only and shall not affect the construction or interpretation of any of its provisions. Any waiver by either party of any default or breach hereunder shall not constitute a waiver of any provision of this Agreement or of any subsequent default or breach of the same or a different kind.

The construction and performance of this Agreement will be governed by the laws of the State of Colorado without reference to its choice of law principles. The parties hereby submit to the jurisdiction of the courts of the State of Colorado.

Contents

Preface	ii
Conventions Used in this Guide	ii
Limited Warranty	iii
Software License	v

Chapter 1. Overview and General Features

ForceTriad Energy Platform Front Panel.	1-2
Introduction	1-3
Indications for Use.	1-3
System Conventions.	1-4
Touchscreens	1-4
Common Symbols	1-4
Power Modes	1-6
Monopolar Modes.	1-7
Bipolar Modes	1-7
LigaSure Mode.	1-8

Chapter 2. Patient and Operating Room Safety

General	2-2
Setting Up the System.	2-2
Fire/Explosion Hazard	2-4
Energy Platform.	2-5
Active Instruments.	2-5
Implanted Electronic Devices (IEDs)	2-6
After Surgery	2-6
Monopolar	2-7
Patient Return Electrodes	2-7
Inadvertent Radio Frequency (RF) Burns	2-8
Laparoscopic Procedures	2-9
Bipolar.	2-10
LigaSure	2-10
LigaSure in Laparoscopic Procedures	2-11
Servicing	2-11
Shunt Cords	2-11
Conductive Fluid In the Surgical Site	2-12

Chapter 3. System Setup

Setup.	3-2
Before Startup.	3-2

Powering Up the ForceTriad Energy Platform	3-2
System Functions	3-2
Adjusting Display Brightness	3-2
Activation Log	3-2
Service Display	3-3
Restore	3-3
Setup	3-3
Demo Mode	3-9

Chapter 4. Monopolar Function

Front Panel Monopolar Features	4-2
Rear Panel Monopolar Features	4-3
Monopolar Quick Setup Instructions	4-4
Monopolar Function Overview	4-4
Monopolar Power Output Modes	4-4
Monopolar Footswitch.	4-5
Standard Monopolar Footswitching	4-5
Software-Enabled Monopolar 1 Footswitching.	4-6
Return Electrodes – REM Contact Quality Monitoring System	4-7
Patient Return Electrode Considerations.	4-7
How the REM System Works.	4-8
Patient Return Electrode Setup	4-9
Monopolar Electrodes	4-10
Standard Monopolar Mode Functionality	4-11
Valleylab Mode Functionality.	4-12
Using a Valleylab Mode Instrument.	4-13
Manual Power Control Functionality	4-14
Accessory Port Functionality	4-15

Chapter 5. Bipolar Surgery

Front Panel Bipolar Features	5-2
Rear Panel Bipolar Features.	5-3
Bipolar Quick Setup Instructions	5-3
Bipolar Function Overview	5-4
Bipolar Power Output Modes.	5-4
Footswitch	5-4
Bipolar Electrode Function	5-5
Autobipolar Function.	5-6
Activation Impedance	5-6

Chapter 6. LigaSure Tissue Fusion

Front Panel LigaSure Features.	6-2
--	-----

Rear Panel LigaSure Features	6-3
LigaSure Quick Setup Instructions.	6-3
LigaSure Function Overview	6-4
LigaSure Receptacles	6-5
Footswitch	6-6
Start System	6-6
LigaSure Instruments	6-7
Reusable-Instrument Assembly.	6-7
Connecting LigaSure Instruments to the Energy Platform	6-7
LigaSure Settings	6-9
Changing the Energy-Delivery Setting.	6-9
Activating the LigaSure Instrument	6-9
Alert Situations	6-10
After Surgery	6-11
Disconnect the instruments	6-11
Reprocessing Instruments	6-12
Chapter 7. Troubleshooting	
General Troubleshooting Guidelines	7-2
REM Alarms	7-2
Correcting a REM-Alarm Condition	7-2
Correcting Malfunctions	7-3
System Alarms	7-8
Chapter 8. Maintenance and Repair	
Responsibility of the Manufacturer	8-2
Routine Maintenance and Periodic Safety Checks.	8-2
Cleaning	8-3
Product Service.	8-3
Returning the Energy Platform for Service	8-3
Adjustment to Factory Specification (Calibration).	8-4
Software Upgrades	8-4
Covidien Technical Service.	8-4
Chapter 9. Technical Specifications	
Performance Characteristics	9-2
General.	9-2
Dimensions and Weight	9-2
Operating Parameters.	9-3
Transport and Storage	9-3
Internal Memory	9-3

Activation Tone	9-4
Alarm Tone	9-4
REM Contact Quality Monitor	9-5
Autobipolar	9-6
Duty Cycle	9-8
Low Frequency (50/60 Hz) Leakage Current	9-8
High Frequency (RF) Leakage Current	9-8
Input Power	9-9
Power Cord Specification	9-9
Input Frequency.	9-9
Input Current.	9-10
Backup Power	9-10
Equipotential Ground Connection.	9-10
ECG Blanking	9-10
Standards and IEC Classifications	9-11
Symbols.	9-12
Class I Equipment (IEC 60601-1)	9-13
Type CF Equipment (IEC 60601-1)/Defibrillator Proof.	9-13
Liquid Spillage (IEC 60601-2-2:2006 Clause 44.3 and IEC 60601-2-2:2009 Clause 201.11.6.3)	9-13
Voltage Transients (Emergency Energy Platform Mains Transfer).	9-13
Electromagnetic Compatibility (IEC 60601-1-2 and IEC 60601-2-2).	9-14
Output Characteristics	9-15
Maximum Output for Bipolar, Monopolar, and LigaSure Modes.	9-15
Available Power Settings in Watts.	9-16
Output Waveforms	9-21
Output Power vs. Resistance Graphs.	9-22
Monopolar Graphs	9-22
Bipolar Graphs.	9-29

Chapter 1

Overview and General Features

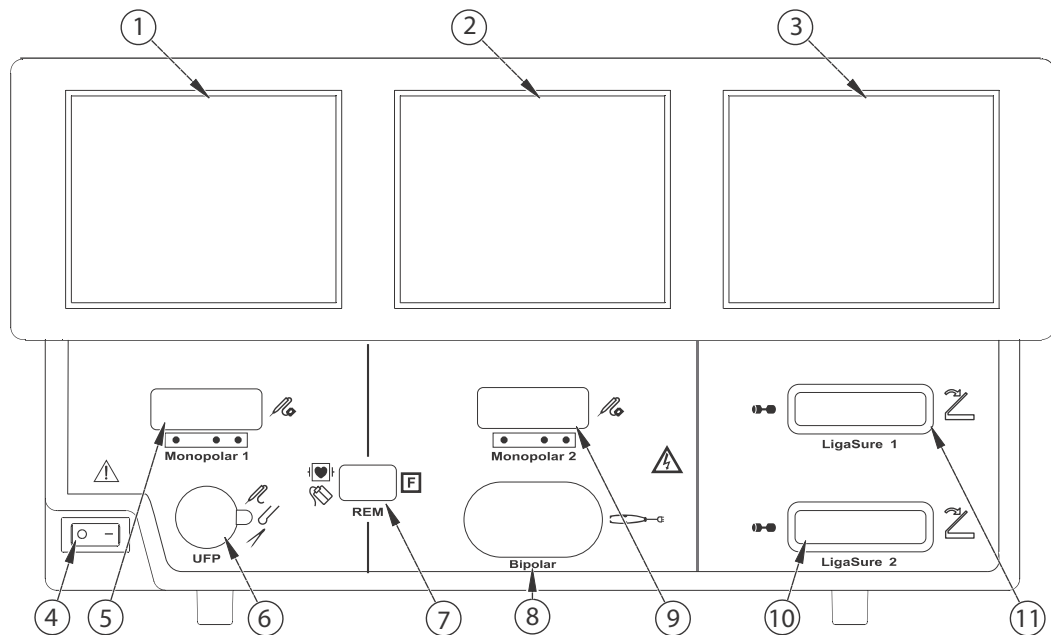
This chapter provides an overview of the features and functions of the ForceTriad Energy Platform.

Caution

Read all warnings, cautions, and instructions provided with this system before use.

Read the instructions, warnings, and cautions provided with electro-surgical instruments before use. Specific instructions for electro-surgical instruments are not included in this manual.

ForceTriad Energy Platform Front Panel



- ① Monopolar 1 and Accessory Touchscreen
- ② Monopolar 2 and Bipolar Touchscreen
- ③ LigaSure and System-Tray Touchscreen
- ④ Power Switch
- ⑤ Monopolar 1 Instrument Receptacle
- ⑥ Universal-Footswitching-Accessory Receptacle
- ⑦ REM-Patient-Return-Electrode Receptacle
- ⑧ Bipolar-Instrument Receptacle
- ⑨ Monopolar 2 Instrument Receptacle
- ⑩ LigaSure 2 Receptacle
- ⑪ LigaSure 1 Receptacle

Introduction

The ForceTriad Energy Platform is designed to provide RF energy for monopolar and bipolar surgical applications and tissue-fusion applications. It features three touchscreen user interfaces, and has the ability to automatically detect handsets and configure the energy platform accordingly. Safety and diagnostic functionality include automatic fail-safe functions.

The Covidien ForceTriad Energy Platform, patient return electrodes, and active instruments are designed to work as a system. Covidien offers a selection of patient return electrodes and active instruments that are fully compatible with this energy platform. When considering other manufacturers' patient return electrodes and/or active instruments, customers should seek detailed user instructions and warning information from the manufacturer.

Indications for Use

The ForceTriad is a full-featured electrosurgical generator intended for open and laparoscopic surgical procedures where the surgeon requires electrosurgical cutting (resecting, dividing, or separating), coagulation (hemostasis, coagulating, or desiccating), or vessel sealing (sealing or fusing). The generator is intended for use in general, laparoscopic, and gynecologic surgical procedures where ligation of vessels, pulmonary vasculature, or lymph vessels, is desired. The system creates a vessel ligation (seal) by the application of bipolar electrosurgical RF energy (coagulation) to vessels interposed between the jaws of the device. The generator can also be used with standard bipolar devices where bipolar cutting or coagulation is desired or with compatible resectoscopes for endoscopically controlled removal (resection) or coagulation of tissue using 0.9% NaCl solution (saline) as the irrigation medium.

The indications for use include general (including urologic, thoracic, plastic, and reconstructive), laparoscopic, and gynecological procedures where electrosurgical cutting and coagulation of tissue, and sealing (fusion) of vessels, including pulmonary vessels, and tissue bundles is performed, including such procedures as bowel resections, hysterectomies (both vaginal and abdominal), laparoscopic cholecystectomies, laparoscopically assisted vaginal hysterectomies, gall bladder procedures, Nissen fundoplication, adhesiolysis, oophorectomy, etc. The devices can be used on vessels (arteries, veins, pulmonary arteries, pulmonary veins, lymph) up to 7mm and bundles as large as will fit in the jaws of the instruments.

The LigaSure tissue fusion function has not been shown to be effective for tubal sterilization or tubal coagulation for sterilization procedures. Do not use this function for these procedures.



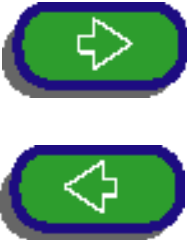

For warnings and specific contraindications for instruments used with the ForceTriad Energy Platform, refer to each instrument's instructions for use (IFU).









System Conventions






Touchscreens

The ForceTriad Energy Platform features a user-friendly interface with three touchscreens that allow the user to control system functions. The active touchscreen or touchscreens will illuminate, and the unavailable touchscreens will dim.

Common Symbols

Symbol	Name	Description
	Page Up/Page Down	Scroll through blocks of options that cannot be displayed on a single screen.
	Up/Down	Pressing once increases/decreases the associated value or moves highlighted selection up/down one line. Pressing and holding scrolls up/down.
	Next/Back	Progresses/regresses to the next screen.
	Back Space	Regresses one character.

Symbol	Name	Description
	Bipolar Mute On/Off	Turn on/off the audio tones produced by the system that indicate the increase or decrease of current during a bipolar procedure.
		
	Monopolar Footswitching On/Off	Turn on/off monopolar footswitching capability for Monopolar 1.
		
	LigaSure Hand-Activation On/Off	Turn on/off hand activation on the LigaSure 1 or 2 control panel. See <i>Hand-Activation On/Off Button</i> on page 6-8.
		
	Cancel	Cancel current screen and returns to the previous screen.
	Enter	Accept and initiate current selections.

Symbol	Name	Description
	System Tray	The system tray contains controls that allow you to access and adjust system settings including screen brightness and main menu options as well as a connection indicator.
	Brightness	Adjusts screen brightness to the next of the two available brightness settings. When maximum brightness is reached, the next selection resets to the least bright setting.
	Wrench	Displays the main menu, which provides user-selected options for language, appearance, and operation.
	Connection Indicator	Indicates active communication with another system such as Valleylab Exchange Remote Software System or a third party system.
	Errors Disabled	This icon on a yellow background overlays the screen when error warnings have been disabled using the service menu. The energy platform will not sound an alarm or give error conditions when this symbol is activated. Touching the screen removes the icon for five seconds.

Note: Additional information on symbols may be found in Chapter 9, *Technical Specifications*.

Power Modes

As a safety feature, simultaneous activation of multiple instruments is not possible on the ForceTriad Energy Platform.

Monopolar Modes

The system produces five modes of power output.

Notice

To provide expected functionality from a hand piece, proper insertion is required. Refer to the orientation drawing near the receptacles for proper insertion orientation.

Cut Modes

Pure cut provides a clean, precise cut in any tissue with little or no hemostasis.

Blend cut is a conventional blended waveform that provides slower cutting with simultaneous hemostasis.

Valleylab Mode

Valleylab mode is a unique combination of hemostasis and dissection that allows the user to slow down for more hemostasis and speed up for faster dissection. Thermal spread is equal to or less than cut or blend modes.

Coag Modes

Fulgurate coagulates tissue by sparking from the active electrode, through air, to the patient tissue.

Spray delivers wider fulguration; penetration is shallower and the affected tissue area is larger than with the Fulgurate mode.

Bipolar Modes

Three bipolar modes are available:

Low delivers precision and fine control over the amount of desiccation.

Standard is a conventional bipolar output at low voltage.

Macro (macrobipolar) may be used for bipolar cutting or rapid coagulation. Power remains constant over a wide range of tissue types.

Autobipolar

The autobipolar feature senses tissue impedance between the two bipolar electrodes, then uses the impedance information to automatically start or stop bipolar RF energy delivery. Optionally, the user may choose between footswitch start and auto start, or program a delay between auto start and RF activation.

Note: When using autobipolar mode, the tissue in the grasp of the bipolar device must have an impedance less than 1000 Ω . The activation impedance safety feature will not deliver RF power to the tissue if it is not within the specified range. This is a factory-set value that cannot be reset by the user.

LigaSure Mode

The LigaSure tissue-fusion mode can be used on arteries, veins, pulmonary vasculature, and lymphatics up to and including 7 mm in diameter and tissue bundles. This system provides precise energy delivery and electrode pressure to vessels for a controlled time period to achieve a complete and permanent fusion of the vessel lumen. The system has been designed to produce minimal sticking, charring, or thermal spread to adjacent tissue.

Warning

Do not attempt to fuse lung tissue with LigaSure mode or instruments.

LigaSure Instruments

The LigaSure instruments that complete the ForceTriad tissue-fusion system include multiple reusable and single-use instruments for open and laparoscopic procedures. Each reusable instrument requires a corresponding single-use electrode. The LigaSure function is only available when using Covidien LigaSure instruments.

Chapter 2

Patient and Operating Room Safety

The safe and effective use of electrosurgery depends to a large degree upon factors solely under the control of the operator. There is no substitute for a properly trained and vigilant surgical team. It is important that the operating instructions supplied with this or any electrosurgical equipment be read, understood, and followed.

Electrosurgery has been used safely in millions of procedures. Before starting any surgical procedure, the surgeon should be trained in the particular technique and surgical procedure to be performed, should be familiar with the medical literature related to the procedure and potential complications, and should be familiar with the risks versus the benefits of utilizing electrosurgery in the procedure.

General

Setting Up the System

Warning

Electric Shock Hazard Connect the system power cord to a properly grounded power receptacle. Do not use power plug adapters.

Fire Hazard Do not use extension cords.

Patient Safety Use the energy platform only if the power-up self-test has been completed as described in this manual, otherwise inaccurate power outputs may result.

Caution

When using a smoke evacuator in conjunction with the ForceTriad Energy Platform, set the system volume control at a level that ensures that the activation tones can be heard.

Connect only Covidien-approved footswitches. Using footswitches from other manufacturers may cause equipment malfunction.

Warning

Hazardous Electrical Output This equipment is for use only by trained, licensed physicians.

Do not use electrosurgical equipment unless properly trained to use it in the specific procedure being undertaken. Use of this equipment without such training can result in serious, unintended patient injury, including bowel perforation and unintended, irreversible tissue necrosis.

Always use the lowest power setting that achieves the desired surgical effect. The active electrode should be utilized only for the minimum time necessary in order to lessen the possibility of unintended burn injury. Accidental and unintended burn injury has occurred during procedures in small surgical fields and on small appendages. Pediatric applications and/or procedures performed on small anatomic structures may require reduced power settings. The higher the current flow and the longer the current is applied, the greater the possibility of unintended thermal damage to tissue, especially during use on small structures.

Do not wrap the instrument cords or patient return electrode cords around metal objects. This may induce currents that could lead to shocks, fires, or injury to the patient or surgical team.

Electric Shock Hazard Do not connect wet instruments to the energy platform. Ensure that all instruments and adapters are correctly connected and that no metal is exposed at any connection points.

Confirm proper power settings before proceeding with surgery. If the proper power settings are not known, set the power to a low setting and cautiously increase the power until the desired effect is achieved. If increased power settings are requested, check the patient return electrode and all instrument connections before major power-setting adjustments.

Warning

Contact between the active electrode and any metal will greatly increase current flow and can result in unintended surgical effect.

While using electrosurgery, the patient should not be allowed to come into direct contact with grounded metal objects (e.g., surgical-table frame, instrument table, etc.). If this is not possible during certain procedures (e.g., those in which noninsulated head frames are used), use extreme caution to maximize patient safety:

- Use the lowest power setting that achieves the desired effect.
- Place the patient return electrode as close to the surgical site as possible.
- Place dry gauze between the patient and the grounded object if possible.
- Continually monitor the contact point(s).
- Do not use metal needle monitoring electrodes.

Caution

Read all warnings, cautions, and instructions provided with this energy platform before using.

Read the instructions, warnings, and cautions provided with electrosurgical instruments before using. Specific instructions for electrosurgical instruments are not included in this manual.

For surgical procedures where the current could flow through delicate parts of the body, the use of bipolar techniques may be desirable in order to avoid unwanted coagulation.

Examine all instruments and connections to the system before using. Ensure that the instruments function as intended. Improper connection may result in arcs, sparks, instrument malfunction, or unintended surgical effects.

Do not turn the activation tone down to an inaudible level. The activation tone alerts the surgical team when the energy platform is delivering RF energy.

A non-functioning ForceTriad Energy Platform may cause interruption of surgery. A backup system should be available for use.

Studies have shown that smoke generated during electrosurgical procedures can be potentially harmful to patients and the surgical team. These studies recommend adequately ventilating the smoke by using a surgical-smoke evacuator or other means.¹

Inadvertent activation may occur while installing, removing, or bending electrodes. Ensure that the instrument cord is not connected to the ForceTriad Energy Platform or that the system is OFF.

1. U.S. Department of Health and Human Services. National Institute for Occupational Safety and Health (NIOSH). Control of Smoke from Laser/Electric Surgical Procedures. HAZARD CONTROLS, Publication No. 96-128, September, 1996

Notice

Connect the power cord to a properly grounded power receptacle having the correct voltage. Otherwise, product damage may result.

Important

If required by local codes, connect the energy platform to the hospital equalization connector with an equipotential cable.

Fire/Explosion Hazard

Warning

Danger: Explosion Hazard Do not use electrosurgery in the presence of flammable anesthetics.

Fire Hazard Do not place active instruments near or in contact with flammable materials (such as gauze or surgical drapes). Electrosurgical instruments that are activated or hot from use can cause a fire. When not in use, place electrosurgical instruments in a safety holster or safely away from patients, the surgical team, and flammable materials.

Fire Hazard Sparking and heating associated with electrosurgery can be an ignition source. Keep gauze and sponges wet. Keep electrosurgical electrodes away from flammable materials and oxygen (O₂) enriched environments.

Use of electrosurgery in O₂ rich environments increases the risk of fire. Therefore, take measures to reduce the O₂ concentration at the surgical site.

Avoid enriched O₂ and nitrous oxide (N₂O) atmospheres near the surgical site. Both O₂ and N₂O support combustion and may result in fires and burns to patients or surgical personnel.

If possible, stop supplemental oxygen at least one minute before and during use of electrosurgery.

Do not activate the energy platform until flammable vapors from skin-preparation solutions and tinctures have dissipated.

Avoid the accumulation of naturally occurring flammable gases that may accumulate in body cavities such as the bowel.

Prevent pooling of flammable fluids and the accumulation of flammable or oxidizing gases or vapors under surgical drapes or near the surgical site.

Tissue buildup (eschar) on the tip of an active electrode may create embers that pose a fire hazard, especially in oxygen-enriched environments. Keep the electrode clean and free of all debris.

Facial and other body hair is flammable. Water soluble surgical lubricating jelly may be used to cover hair close to the surgical site to decrease flammability.

Verify that all anesthesia circuit connections are leak free before and during use of electrosurgery.

Warning**Fire Hazard During Oropharyngeal Surgery**

Verify endotracheal tubes are leak free and that the cuff seals properly to prevent oxygen leaks.

If an uncuffed tube is in use, pack the throat with wet sponges around the uncuffed tube, and be sure to keep sponges wet throughout the procedure.

Question the need for 100% O₂ during oropharyngeal or head and neck surgery.

If necessary, scavenge excess O₂ with separate suction.

Energy Platform**Warning**

Each instrument receptacle on this energy platform is designed to accept only one instrument at a time. Follow the instructions provided with electro-surgical instruments for proper connection and use.

Caution

Do not stack equipment on top of the energy platform or place the energy platform on top of electrical equipment. This is an unstable configuration and does not allow for adequate cooling.

Provide as much distance as possible between the energy platform and other electronic equipment (such as monitors). Do not cross or bundle electronic-device cords. This energy platform may cause interference with other electronic equipment.

Active Instruments**Caution**

Read the instructions, warnings, and cautions provided with electro-surgical instruments before using. Specific instructions for electro-surgical instruments are not included in this manual.

Inspect instruments and cords for breaks, cracks, nicks, and other damage before every use. If damaged, do not use. Damaged instruments or cords may result in injury or electrical shock to the patient or surgical team.

Use only instruments that can withstand the maximum output (peak) voltage for each output mode as listed in Chapter 9, *Technical Specifications*. Using an instrument with a voltage rating that is lower than the maximum output voltage may result in injury to the patient or the operator, or damage to the instrument.

All Covidien instruments have voltage ratings that are greater than the maximum output voltages in the ForceTriad Energy Platform and are thus fully compatible.

Information on voltage ratings for non-Covidien instruments should be obtained from the instrument's manufacturer.

Implanted Electronic Devices (IEDs)

IEDs include, but are not limited to, pacemakers, neurostimulators, implantable cardioverter defibrillators (ICDs), ventricular assist devices (VAD), spinal cord stimulators, cochlear implants, infusion pumps and bone growth stimulators.

Warning

Use electrosurgery and tissue fusion with caution in the presence of internal or external pacemakers. Interference produced by the use of electrosurgical devices can cause a pacemaker to enter an asynchronous mode or can block the pacemaker effect entirely. Consult the pacemaker manufacturer or hospital cardiology department for further information when use of electrosurgery or tissue-fusion appliances is planned in patients with cardiac pacemakers.

If the patient has an implanted electronic device (IED), contact the IED manufacturer for instructions before performing an electrosurgical or tissue-fusion procedure. Electrosurgery or tissue fusion may cause multiple activations of ICDs, or interfere with the intended function of other IEDs.

After Surgery

Warning

Electric Shock Hazard Always turn off and unplug the energy platform before cleaning.

Caution

Do not reprocess, reuse or resterilize instruments labeled "disposable" or "single use only."

Notice

Do not clean the energy platform with abrasive cleaning or disinfectant compounds, solvents, or other materials that could scratch the panels or damage the energy platform.

Monopolar

Warning

Simultaneously activating suction/irrigation and electrosurgical current may result in increased arcing at the electrode tip, burns to unintended tissues, or shocks and burns to the surgical team.

Some surgeons may elect to “buzz the hemostat” during surgical procedures. It is not recommended, and the hazards of such a practice probably cannot be eliminated. Burns to the surgeon’s hands are possible. To minimize the risk take these precautions:

- Do not buzz the hemostat with a needle electrode.
- Do not lean on the patient, the table, or the retractors while buzzing the hemostat.
- Activate cut rather than coag. Cut has a lower voltage than coag.
- Firmly grasp as much of the hemostat as possible before activating the energy platform. This disperses the current over a larger area and minimizes the current concentration at the finger tips.
- Buzz the hemostat below hand level (as close as possible to the patient) to reduce the opportunity for current to follow alternate paths through the surgeon’s hands.
- Use the lowest power setting possible for the minimum time necessary to achieve hemostasis.
- Activate the energy platform after the instrument makes contact with the hemostat. Do not arc to the hemostat.
- When using a coated- or nonstick-blade electrode, place the edge of the electrode against the hemostat or other metal instrument.

Patient Return Electrodes

Warning

Do not attempt to use patient return electrodes that disable the REM system. The ForceTriad Energy Platform REM system will function correctly only with contact quality monitoring (CQM) split-style patient return electrodes. Any other patient return electrode products may cause patient injury or product damage.

The safe use of monopolar electrosurgery requires proper placement of the patient return electrode. To avoid electrosurgical burns beneath the patient return electrode, follow all directions provided with the product.

Do not cut a patient return electrode to reduce its size. Patient burns due to high current density may result.

A patient return electrode is not necessary in bipolar or LigaSure procedures.

To avoid patient burns, ensure that the patient return electrode firmly and completely contacts the skin. Always check the patient return electrode periodically and after the patient is repositioned and during procedures involving long periods of activation.

Warning

Use of duty cycles greater than 25% (10 seconds active followed by 30 seconds inactive) will increase the risk that heat build-up under a return electrode may be high enough to injure the patient. Do not continuously activate for longer than one minute.

Notice

Capacitive pads and other non-CQM patient return electrodes may not work with the ForceTriad Energy Platform.

Important

A statement of compatibility from the CQM patient return electrode manufacturer should be obtained prior to the use of a non-Covidien CQM patient return electrode.

Inadvertent Radio Frequency (RF) Burns

Warning

Electrodes and probes used with monitoring, stimulation, and imaging devices (or similar equipment) can provide a path for high frequency current even if the electrodes or probes are isolated at 50 Hz-60 Hz, insulated, and/or battery operated.

Do not use needles as monitoring electrodes during electrosurgical procedures. Inadvertent electrosurgical burns may result.

To reduce the risk of an inadvertent electrosurgical burn at the electrode or probe site, place the electrode and/or probe as far away as possible from the electrosurgical site and/or patient return electrode. Protective impedances (resistors or RF inductors) installed in the monitoring leads may reduce the risk of such burns. Consult the hospital biomedical engineer for further information.

In some circumstances, the potential exists for alternate site burns at points of skin contact (e.g., between the arm and the side of the body). This occurs when electrosurgical current seeks a path to the patient return electrode that includes the skin-to-skin contact point. Current passing through small skin-to-skin contact points is concentrated and may cause a burn. This is true for ground referenced and isolated output electrosurgical energy systems.

To reduce the potential for alternate site burns, do one or more of the following:

- Avoid skin-to-skin contact points, such as fingers touching leg or knee touching knee when positioning the patient.
- Place insulation, such as dry gauze or towel, between contact points to ensure that contact does not occur.
- Position the patient return electrode to provide a direct current route between the surgical site and the return electrode which avoids skin-to-skin contact areas.
- In addition, place patient return electrodes according to the manufacturer's instructions.

Laparoscopic Procedures

Warning

For laparoscopic procedures, be alert to these potential hazards:

- Laparoscopic surgery may result in gas embolism due to insufflation of gas in the abdomen.
- The electrode tip may remain hot enough to cause burns after the electrosurgical current is deactivated.
- Inadvertent activation or movement of the activated electrode outside of the field of vision may result in injury to the patient.
- Localized burns to the patient or physician may result from electrical currents carried through conductive objects (such as cannulas or scopes). Electrical current may be generated in conductive objects through direct contact with the active electrode, or by the active instrument (electrode or cable) being in close proximity to the conductive object.
- Do not use hybrid trocars that have a non-conductive locking anchor placed over a conductive sleeve. For the operative channel, use all-metal or all-plastic systems. At no time should electrical energy pass through hybrid systems. Capacitive coupling of RF current may cause unintended burns.
- When using laparoscopic instrumentation with metal cannulas, the potential exists for abdominal-wall burns to occur due to direct electrode contact or capacitive coupling of RF current. This is most likely to occur in instances where the energy platform is activated for extended periods at high power levels inducing high-current levels in the cannula.
- Ensure that the insulation of single-use and reusable laparoscopic instrumentation is intact and uncompromised. Compromised insulation may lead to inadvertent metal-to-metal sparking and neuromuscular stimulation and/or inadvertent sparking to adjacent tissue.
- Do not activate electrodes while in contact with other instruments as unintended tissue injury may occur.

Do not activate the energy platform in an open-circuit condition. To reduce the chances of unintended burns, activate the energy platform only when the active electrode is near or touching the target tissue.

- Use the lowest power setting that achieves the desired surgical effect and use a low-voltage waveform (Pure Cut, Blend, or Valleylab mode) to lessen the potential for the creation of capacitive currents.
- Carefully insert and withdraw active electrodes from cannulas to avoid possible injury to the patient or damage to the devices.

Covidien recommends against the use of laparoscopic surgery on pregnant patients.

Bipolar

Caution

Bipolar instruments must be connected to the bipolar instrument receptacle only. Improper connection may result in inadvertent system activation.

LigaSure

Warning

LigaSure instruments are intended for use ONLY with the ForceTriad Energy Platform and the LigaSure vessel sealing system. Use of these instruments with other Covidien generators or with generators produced by other manufacturers may not result in electrical output for which these instruments were designed and thus may not result in the desired clinical effect.

If the seal-complete tone has not sounded, an optimal seal may not have been achieved. Reactivate the RF energy until a seal-complete tone is heard.

The LigaSure tissue-fusion function has not been shown to be effective for tubal sterilization or tubal coagulation for sterilization procedures. Do not use this function for these procedures.

Use caution during surgical cases in which patients exhibit certain types of vascular pathology (atherosclerosis, aneurysmal vessels, etc.). For best results, apply the seal to unaffected vasculature.

Do not activate the energy platform in the LigaSure mode until the tissue-fusion instrument has been applied with the proper pressure. Activating the energy platform before this is done will result in an improper seal and may increase thermal spread to tissue outside the surgical site.

Tissue fusion requires the application of RF energy and pressure from the instrument. Tissue to be sealed must be firmly grasped between the instrument jaw electrodes. Tissue in the jaw hinge or outside the instrument jaw will not be sealed even if thermal blanching occurs.

Do not use LigaSure instruments on vessels in excess of 7 mm in diameter.

LigaSure instruments that require single-use electrodes must be used with the correct electrode type. Use of these instruments with any other electrodes could result in injury to the patient or surgical team, or cause damage to the instrument.

Conductive fluids (e.g, blood or saline) in direct contact with LigaSure instruments or in close proximity may carry electrical current or heat, which may cause unintended surgical effects or burns.

Caution

Energy based devices, such as electrosurgical pencils or ultrasonic scalpels, that are associated with thermal spread should not be used to transect seals.

Avoid placing fingers in the handle ratchet mechanism or between the ring handles or jaws as applicable depending on the type of instrument. Injury to the user may result.

LigaSure in Laparoscopic Procedures

Warning

For laparoscopic procedures, be alert to these potential hazards:

- The external surfaces of the LigaSure instrument jaws may remain hot enough to cause burns after the RF current is deactivated.
- Inadvertent activation or movement of the activated LigaSure instrument outside of the field of vision may result in injury to the patient.
- Do not activate the instrument while the instrument jaws are in contact with, or in close proximity to, other instruments including metal cannulas, as localized burns to the patient or physician may occur.
- Do not activate the LigaSure function in an open-circuit condition. Activate the energy platform only when the instrument is near or in direct contact with the target tissue to reduce the possibility of unintended burns.
- Carefully insert and withdraw LigaSure instruments from cannulas to avoid possible damage to the devices and/or injury to the patient.

Servicing

Warning

Electric Shock Hazard Do not remove the energy platform cover. Contact qualified personnel for service.

Notice

Refer to this system's service manual for maintenance recommendations, and function and output-power verification procedures.

Shunt Cords

Warning

Some surgical instruments (e.g., colonoscopes) may allow substantial leakage current that could burn the surgeon. If the instrument manufacturer recommends the use of a shunt cord (s-cord) to direct the current back to the energy platform, an E0507-B adaptor must also be used. To avoid a REM alarm, REM Polyhesive patient return electrode with the E0507-B adaptor must be used.

Conductive Fluid In the Surgical Site

Warning

When this energy platform is used in procedures where conductive fluid (saline or lactated Ringers) is introduced into the surgical site for distention or to conduct RF current, higher than normal currents (greater than one amp) may be produced. In this situation, use one or more adult-size return electrode. Do not use return electrodes labeled for children, infants, babies, neonatal use, or pediatric use.

Use of duty cycles greater than 25% (10 seconds active followed by 30 seconds inactive) will increase the risk that heat build-up under a return electrode may be high enough to injure the patient. Do not continuously activate for longer than one minute.

Chapter 3

System Setup

This chapter describes how to set up the energy platform, turn it on, and configure system settings.

Caution

Read all warnings, cautions, and instructions provided with this system before use.

Read the instructions, warnings, and cautions provided with electro-surgical instruments before use. Specific instructions for electro-surgical instruments are not included in this manual.

Setup

Before Startup

1. Verify the system is off by pressing the power switch off (O).
2. Place the energy platform on a flat, stable surface such as a table, platform, boom system, or ForceTriad cart. Refer to the procedures for your local institution or your local codes.
3. Plug the system power cord into the rear panel receptacle.
4. Plug the system power cord into a grounded power receptacle.

Note: Do not plug into a power strip or extension cord.

Powering Up the ForceTriad Energy Platform

1. Turn on the system by pressing the power switch on (I). Observe the following during the power-up self-test:
 - The Covidien logo appears on all three screens.
 - A blue status bar indicates activity.
 - An hourglass icon indicates activity after the status bar disappears.
 - The system revision code appears on the center screen.
 - A tone sounds upon completion of self-test.
2. If the system does not pass the power-up self-test and an error code is displayed, refer to Chapter 7, *Troubleshooting*.

System Functions


Adjusting Display Brightness



The ForceTriad screens have two levels of brightness. Touch the brightness icon on the right side of the right touchscreen to adjust the display brightness.

Activation Log

The Activation Log allows the user to view the last 1000 activations and REM alerts.

1. Touch the wrench icon  on the right side of the right touchscreen. The main-menu display will appear in the left touchscreen.

2. Touch **Activation Log** in the main menu. The activation log will appear on the center touchscreen.
3. Touch the single up or down arrow to the right of the activation log to scroll through the log one line at a time.
4. Touch the green arrow button on the bottom-right corner of the main menu screen to return the ForceTriad Energy Platform to the previous setup configuration. The last settings will be displayed.

Service Display

Refer to the *ForceTriad Energy Platform Service Manual* for complete service instructions.


Restore

Select the **Restore** button in the main menu to restore the ForceTriad Energy Platform to the previous setup configuration. The touchscreens will display the last settings entered prior to shutting the system off.

Setup

The **Setup** menu allows the user to change the language that the system touchscreens display, set the time and date, and access the features menu.

Language Setup

1. Touch the wrench icon  on the right side of the right touchscreen. The main menu display will appear in the left touchscreen.
2. Touch **Setup** in the main menu. The setup display will appear in the left touchscreen.
3. Touch **Language** in the setup menu. A list of languages will appear in the left touchscreen.
4. Touch the single up or down arrow to the right of the list to scroll through the list one line at a time.

or

Touch the double up or down arrow to scroll through the list one page at a time.


5. Touch the desired language. A confirmation box will appear and request the user to confirm that a language change is desired.
6. To proceed with the language change, touch the green check-mark button. The language will be activated and the confirmation box will close.

or

To reject the language change, touch the red 'X' button. The language setting will return to the previously selected language.

7. Touch the green arrow button to return to the setup menu.
8. Touch the green arrow button below the Setup menu to return to the main menu.

Time and Date Setup

1. Touch the wrench icon  on the right side of the right touchscreen. The main-menu display will appear in the left touchscreen.
2. Touch **Setup** in the main menu. The setup display will appear in the left touchscreen.
3. Touch the **Time and Date** button in the setup menu. The time-and-date display will appear in the left touchscreen.
4. Touch the desired numeric field (minutes, seconds, month, day, or year) to select that field.
5. Touch the up or down arrow next to the time or date row to adjust the selected numeric field.
6. Touch and hold the arrow to increase the number once per second. After four seconds, the numbers will increase once per 100 milliseconds.

Touch the green check-mark button to store the date and time information and return to the setup menu.

or

Touch the red 'X' button to return the time and date to the previous settings and return to the setup menu.
7. Touch the green arrow button below the setup menu to return to the main menu.

Features Menu

The features menu displays software features and applications that can be enabled or disabled at the system level.

- **Autobipolar**—Configures automatic activation and cessation of energy from the Auto tab.
- **Mono 1 Footswitching** —Enables footswitching control of both handswitching and footswitching devices attached to the Monopolar 1 receptacle.
- **Other features or applications**—Other features and applications may be on this menu based on special configurations or purchased applications


The default software settings for autobipolar mode and monopolar 1 footswitching are enabled or disabled at the system level from the Features menu. Once enabled, the feature is available on the affected screen and can be turned on and off.

Autobipolar and monopolar settings selected on the Features menu determine the options available from the Bipolar and Monopolar touch screens. The settings can be turned on and off locally from buttons available on the affected screens.

Autobipolar

Autobipolar mode on the Features menu enables and disables autobipolar function at the system level, controlling whether autobipolar function is available on the Bipolar tab of the center screen.


Enable/Disable Autobipolar at the system level:

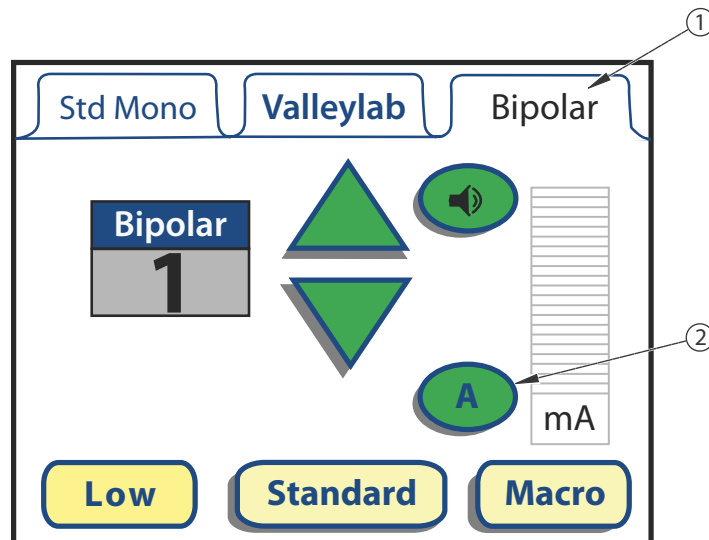
1. Touch the wrench icon  on the right side of the right touchscreen. The main-menu display will appear in the left touchscreen.
2. Touch **Setup** in the main menu. The setup display appears in the left touchscreen.
3. Touch **Features** in the setup menu. Available options appear in the left touchscreen. The factory default for all features is disabled.
4. To enable, touch **AutoBipolar**. A check appears in the accompanying box.



To disable, touch **AutoBipolar**. The check is cleared from the box.

5. Touch the green arrow buttons below the Features menu and Setup menu to return to the main menu.

When the AutoBipolar option is enabled, a green  button is added to the Bipolar tab in the center screen.



① Bipolar Tab (This label changes to Auto if autobipolar mode is turned on)

② Autobipolar Mode Button (shows autobipolar mode is enabled but not turned on)

Turning Autobipolar Mode On and Off

Autobipolar mode can be turned on and off on the Bipolar tab if it is enabled from the Features menu.

Note: If enabled on the Features menu, the default setting for autobipolar mode is off.



To turn autobipolar mode on, touch the  button on the **Bipolar** screen. The **Bipolar** tab changes to **Auto**, presenting the options and settings for autobipolar mode. See *Autobipolar Function* on page 5-6 for instructions for using autobipolar mode.

To turn autobipolar mode off, touch the green **Bipolar** button on the **Auto** tab. The **Auto** tab changes to **Bipolar**.

Mono 1 Footswitching

The Monopolar 1 footswitching option on the Features menu enables and disables footswitching on the Monopolar 1 Std Mono and Valleylab tabs. The Accessories Port tab (**Acc. Port**) is not affected. If enabled on the Features menu, footswitching can then be turned on and off on the Std Mono and Valleylab tabs as well as on the *This feature is not enabled touchscreen* as described in the following note.

Enabling and Disabling Monopolar 1 Footswitching:

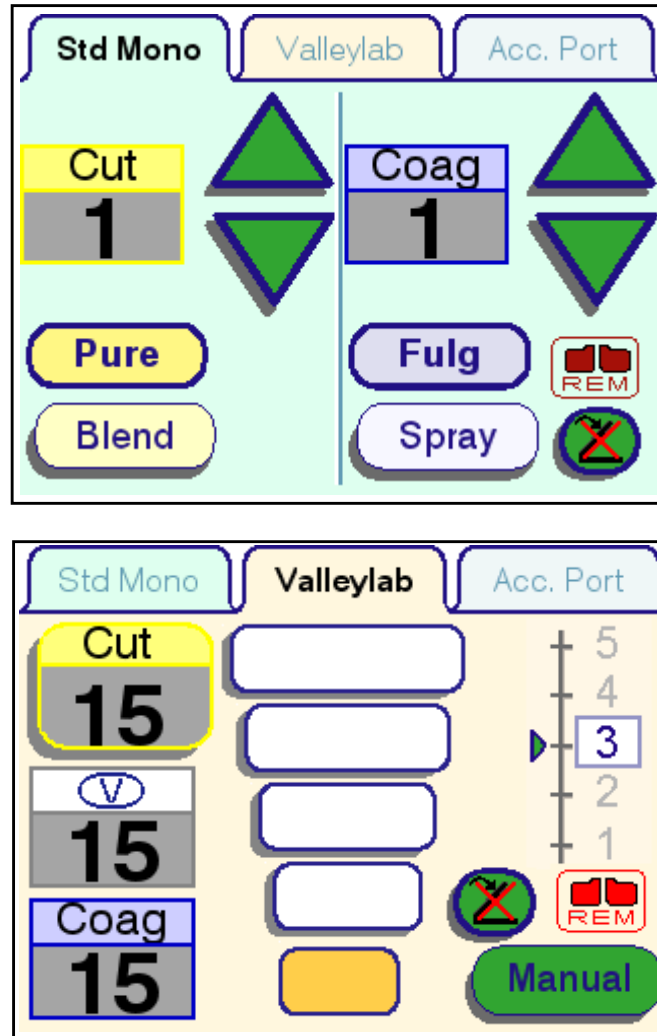
1. Touch the wrench icon  on the right side of the right touchscreen. The main-menu display appears in the left touchscreen.
2. Touch **Setup** in the main menu. The setup display appears in the left touchscreen.
3. Touch **Features** in the setup menu. Available options appear in the left touchscreen. The default setting for all features is disabled.
4. To enable the feature, touch **Mono 1** . A check appears in the accompanying box.



To disable, touch **Mono 1** . The check is cleared from the box.

5. Touch the green arrow buttons below the features menu and setup menu to return to the main menu.


When Monopolar 1 footswitching is enabled in the Features menu, the footswitching button with an "X" appears on the **Std Mono** and **Valleylab** tabs in the left screen.




Turning Monopolar 1 Footswitching On and Off:

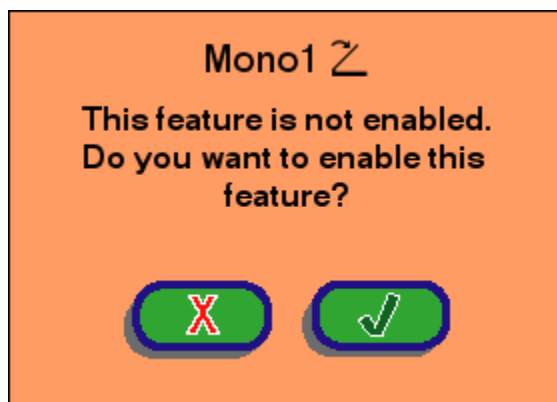
Monopolar 1 footswitching can be turned on and off on the **Std Mono** and **Valleylab** tabs as needed if it is enabled on the Features menu.

Note: If enabled on the **Features** menu, the default setting for Monopolar 1 footswitching is off, as indicated by the red "X" over the button.

To turn Monopolar 1 footswitching on, touch the footswitch button with the "X"  on the **Std Mono** or **Valleylab** tab on the left touchscreen. The red "X" is removed from the button when turned on.

To turn Monopolar 1 footswitching off, touch the footswitch button  on the **Std Mono** or **Valleylab** tabs on the left touchscreen. The red "X" appears on the button when turned off.


Note: If Monopolar 1 footswitching has not been enabled, the following touchscreen appears. Enable Monopolar 1 footswitching by selecting the button with the green check mark.



Demo Mode

Warning

Demo mode is intended for demonstration purposes only. Demo mode is not intended for clinical use.

Touch the wrench icon  on the right side of the right touchscreen. The main menu display appears in the left touchscreen.

Enable Demo Mode

1. In the main menu, the Demo-mode button will display 'ENTER DEMO' if the system is not in Demo mode. Touch the Enter Demo-mode button to begin Demo mode. The system operating displays will appear in all the touchscreens with the words 'DEMO MODE: NOT FOR CLINICAL USE' on all three screens.


DEMO MODE:
Not For Clinical Use

- Note:** Touching the Demo-mode screen will remove it briefly from all touchscreens.
2. Proceed with any practice or demonstration scenarios. While in Demo mode, the REM alarm and the dual-instrument error alarm are deactivated but RF power will still be delivered.

Note: In Demo mode the energy platform will not sense instrument type, so the appropriate tab must be selected manually for the connected instrument.

3. To exit Demo mode, either turn the system off and restart it, or follow the steps in the *Exit Demo Mode* section that follows.

Exit Demo Mode

1. Touch the wrench icon  on the right side of the right touchscreen. The main menu display will appear in the left touchscreen.
2. In the main menu, the Demo-mode button will display 'EXIT DEMO' if the system is in Demo mode. Touch the Exit Demo button in the main menu to exit the Demo mode. The system touchscreens will display the last settings entered during the Demo mode.

Chapter 4

Monopolar Function

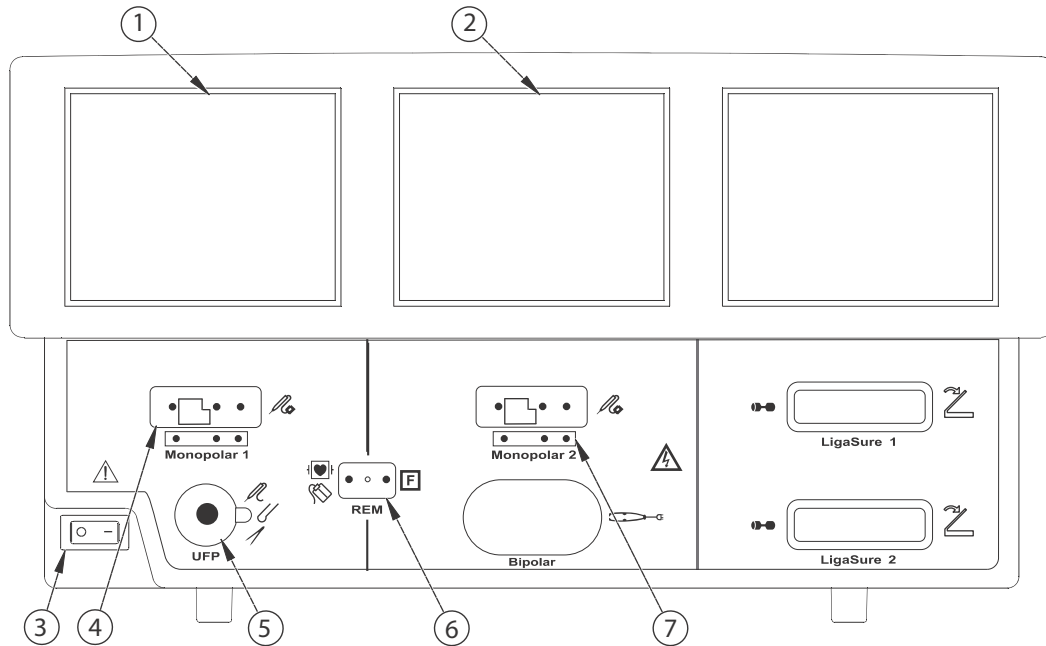
This chapter describes the monopolar surgery features of the ForceTriad Energy Platform.

Caution

Read all warnings, cautions, and instructions provided with this system before use.


Read the instructions, warnings, and cautions provided with electro-surgical instruments before use. Specific instructions for electro-surgical instruments are not included in this manual.

Front Panel Monopolar Features



- ① Monopolar 1 and Accessory Touchscreen
- ② Monopolar 2 Touchscreen
- ③ Power Switch
- ④ Monopolar 1 Instrument Receptacle
- ⑤ Universal Footswitching Accessory Receptacle (may require an E05021 or E0017 active adapter)
- ⑥ REM Patient Return Electrode Receptacle
- ⑦ Monopolar 2 Instrument Receptacle

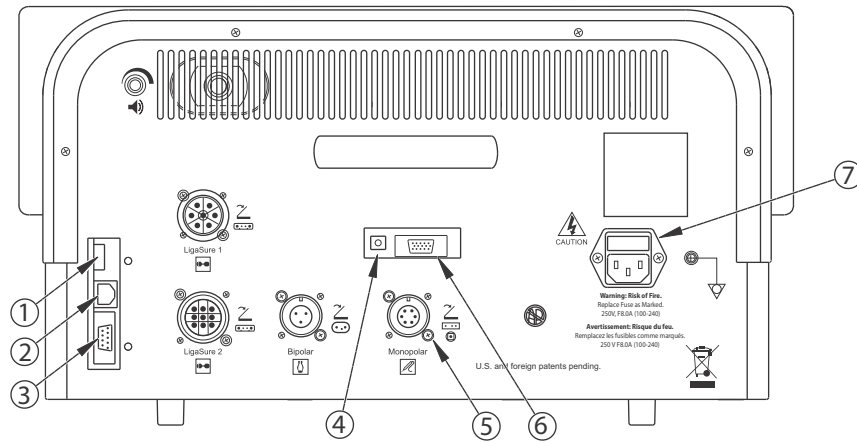Rear Panel Monopolar Features



- ① USB Port
- ② Ethernet Port
- ③ RS232 Port
- ④ Interlink Cable Port
- ⑤ Monopolar Footswitch Receptacle (requires adapter to connect standard four-pin monopolar footswitch)
- ⑥ Expansion Port
- ⑦ Fuse Port

Monopolar Quick Setup Instructions

If familiar with the ForceTriad Energy Platform, follow this abbreviated procedure to set up the system for monopolar surgery.

If not familiar with the ForceTriad Energy Platform, refer to the following sections in this chapter for detailed instructions.

1. Plug the system power cord into the rear panel receptacle.
2. Plug the system power cord into a grounded power receptacle.
3. Turn on the energy platform and verify that the self-test is successfully completed.
4. If using a footswitch, connect it to the monopolar footswitch receptacle on the rear panel. This may require an adapter to connect a standard four-pin monopolar footswitch.
5. Apply the patient return electrode to the patient and connect it to the patient return electrode receptacle on the front panel.
6. Connect the instrument to the appropriate instrument receptacle on the front panel. Align the connector pins with the dots below each Monopolar receptacle to ensure full functionality of the instrument.
7. Verify or change the mode and power settings.

Monopolar Function Overview

Monopolar Power Output Modes

The ForceTriad Energy Platform produces two cut modes, Pure and Blend; one Valleylab mode; and two coag modes, Fulgurate and Spray.

Warning

Electric Shock Hazard

- Do not connect wet instruments to the system.
- Ensure that all instruments and adapters are correctly connected and that no metal is exposed at any connection point.

Connect instruments to the proper receptacle. Improper connection may result in inadvertent instrument activation or other potentially hazardous conditions. Follow the instructions provided with electro-surgical instruments for proper connection and use.

The instrument receptacles on this energy platform are designed to accept only one instrument at a time.

Caution

Read the instructions, warnings, and cautions provided with electrosurgical instruments before use. Specific instructions are not included in this manual.

Inspect instruments and cords (especially reusable instruments and cords) for breaks, cracks, nicks, and other damage before every use. If damaged, do not use. Failure to observe this precaution may result in injury or electrical shock to the patient or surgical team.

Notice

When attaching a monopolar instrument, align the instrument's connector pins with the three dots printed on the ForceTriad face plate below the monopolar receptacles. Reversing the connector limits the instrument's functionality and will damage the connector.

Monopolar Footswitch

Monopolar footswitches control instruments connected to the accessory receptacle or the Monopolar 1 receptacle. The two receptacles are controlled by the left panel.

The ForceTriad Energy Platform provides hardware and software options for monopolar footswitching.

Standard footswitching devices attached to the universal footswitching accessory receptacle on the front panel can be controlled with a Covidien footswitch attached to the Monopolar footswitch receptacle on the back panel.

Software footswitching control can be enabled for both footswitching and handswitching instruments when Monopolar 1 footswitching is enabled from the features menu. (See page 3-7.) Monopolar 1 footswitching also accommodates a Valleylab-mode enabled three-pedal footswitch to provide Valleylab-mode functionality to all instruments.

Caution

Connect only Covidien-approved footswitches. Using footswitches from other manufacturers may cause equipment malfunction.

Standard Monopolar Footswitching

Standard monopolar footswitching can be set up using compatible accessories:

- Footswitching instrument
- Two-pedal footswitch
- Footswitch adapter

Monopolar 1 Footswitching can be set up using compatible accessories:

- FT0510 Monopolar Cord
- FT5003 ForceTriad Three-Pedal Footswitch

Setup:

1. If needed, attach the footswitch adapter to the Monopolar footswitching receptacle on the back panel.
2. Plug the footswitch connector into the adapter, or if an adapter is not needed, directly into the Monopolar footswitching receptacle.
3. If needed, attach the E05021 or E0017 adapter to the universal footswitching accessory receptacle.
4. Plug the footswitching instrument cord into the universal footswitching accessory receptacle (or the adapter, if one was attached in the previous step).

Software-Enabled Monopolar 1 Footswitching

All devices attached to the Monopolar 1 receptacle—whether handswitching or footswitching accessories—can be controlled with a footswitch when the software is enabled.

Valleylab-mode functionality can also be provided to footswitching devices when using a FT0510 Monopolar Cord and FT6003 ForceTriad Three-Pedal Footswitch, both of which are Valleylab-mode enabled.

1. Enable Monopolar 1 footswitching from the features menu as directed beginning on page 3-7.
2. Attach a FT6003 Three-Pedal Footswitch or standard footswitch to the Monopolar footswitch receptacle on the back panel:
 - If using the software-enabled Monopolar 1 footswitching, attach the FT0510 cord to the instrument and the Monopolar 1 accessory receptacle.
 - Or attach any monopolar instrument to the Monopolar 1 accessory receptacle on the front panel.
3. Touch the **Std Mono** tab, or touch the **Valleylab** tab if an FT6003 Three-Pedal Footswitch is attached.

Important

Cut and coag footswitch activation is also available at this point from the Valleylab tab for ForceTraverse™ instruments controlled by a two-pedal footswitch.

4. Enable Monopolar 1 footswitching by touching the footswitch button displaying a red "X."

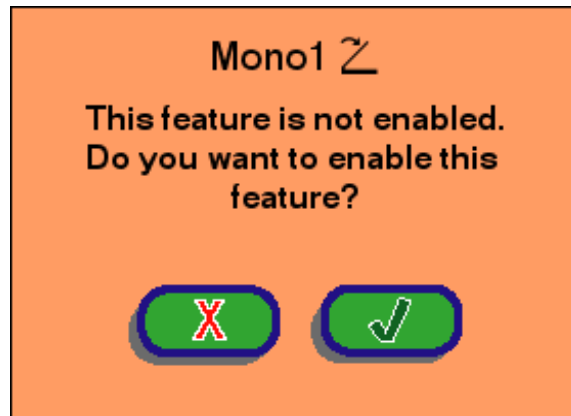


When enabled at the local level, the red "X" is removed.



5. Select settings for the procedure.
6. Step on the appropriate footswitch to deliver monopolar power.

If a Monopolar 1 footswitching has not been enabled at the system level from the features menu or at the local level from the mode tab, the following notice appears.



7. Enable Monopolar 1 footswitching by touching the button with the green check mark.

Return Electrodes – REM Contact Quality Monitoring System

Notice

Only contact quality monitoring system patient return electrodes can be used with the ForceTriad Energy Platform.

Patient Return Electrode Considerations

Warning

It is not possible to foresee what combination of current and duty cycle may be safely used in every situation, such as when higher currents and/or longer duty cycles are used on procedures such as tissue lesioning, tissue ablation, tissue vaporization, and procedures where conductive fluid is introduced into the surgical site. Under these conditions a greater risk may exist that the heating under a fully applied return electrode may be high enough to injure the patient.

When using a Covidien energy platform or a patient return electrode during these types of surgical procedures, the user should seek written guidance in the form of detailed user instructions from the manufacturer of the active accessory regarding the currents and duty cycles that can be expected. In some instances, the application of additional patient return electrodes may help mitigate the increased risk.

During monopolar electrosurgery, a patient return electrode is always required to safely recover the current that flows through the patient's body and return it to the energy

platform. A reduction in surface area contact or poor conductivity between the patient and the return electrode can cause the current to become concentrated, potentially resulting in burns at the return electrode site.

During a surgical procedure, the amount of current delivered during a given time period determines the amount of heating that occurs under the return electrode. REM Polyhesive patient return electrodes are designed for use during conventional electrosurgical procedures and duty cycles (on time compared to off time). Users should consult Chapter 9, *Technical Specifications* for the recommended maximum duty cycle specifications.

How the REM System Works

The ForceTriad Energy Platform uses the REM contact quality monitoring system to monitor the quality of electrical contact between the patient return electrode and the patient. The REM system is designed to reduce the risk of burns at the return electrode site. A non-REM return electrode is not to be used with the ForceTriad Energy Platform.

Compatible REM patient return electrodes, such as the E7507 REM Polyhesive Adult Patient Return Electrode, are available for use with the ForceTriad energy platform.

The REM system continuously measures the resistance at the return electrode site and compares it to a standard range of safe resistance (between 5 Ω and 135 Ω), thus eliminating intermittent false alarms that could result from small changes in resistance.

The REM system also adapts to individual patients by measuring the initial contact resistance between the patient and the patient return electrode and lowering the baseline resistance if the contact resistance drops.

A REM alarm sounds and the system stops producing output power when **either** of the following occurs:

- The measured resistance is below 5 Ω or above 135 Ω , the limits of the standard range of safe resistance.
- An increase in contact resistance is greater than 40% from the baseline measurement.

Patient Return Electrode Setup



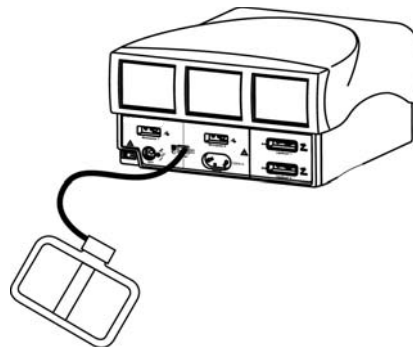
The REM indicator icon appears on the Standard Monopolar, Valleylab, and Accessory Port displays.

Warning

The safe use of monopolar electrosurgery requires proper placement of the patient return electrode. To avoid electrosurgical burns beneath the patient return electrode, follow all directions on the product package for proper return electrode placement and use.

Do not cut a patient return electrode to reduce its size. Patient burns due to high-current density may result.

1. Place the patient return electrode on the patient. Refer to the patient return electrode instructions for proper return electrode placement.
2. Connect the REM Polyhesive patient return electrode plug to the patient return electrode receptacle on the energy platform.



The REM indicator icon on the touchscreen illuminates red to indicate that the REM Polyhesive patient return electrode is disconnected from the energy platform or improperly applied to the patient.

The REM indicator icon on the touchscreen illuminates green when the system senses that the REM Polyhesive patient return electrode is properly connected to the energy platform and patient.

REM Alarm (Visual and Audible)

If the REM system senses an alarm condition, the REM indicator flashes red and yellow, emits two beeps, and discontinues RF energy delivery. When the alarm condition has been corrected, the indicator illuminates green. For detailed instructions to correct REM alarms, refer to Chapter 7, *Troubleshooting* and the troubleshooting flow chart in the REM Polyhesive patient return electrode *Instructions for Use*.



① REM Alarm Indicator

When a REM alarm occurs, the large red and yellow REM icon is displayed for a few seconds and then disappears. The small green REM icon turns red.

Monopolar Electrodes

Connect a monopolar instrument to the Monopolar 1 or Monopolar 2 instrument receptacle on the front of the energy platform. See the Covidien instrument's *Instructions for Use* (IFU) for which receptacle to use.

Notice

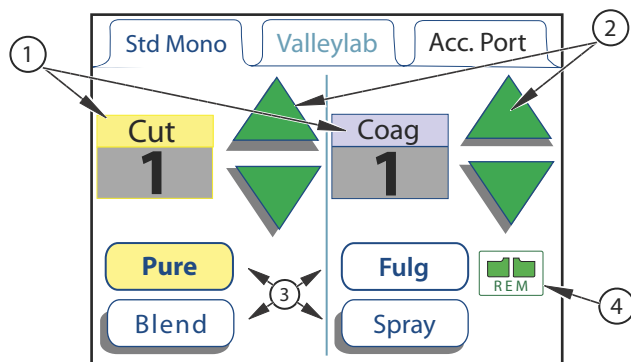
When attaching a monopolar instrument, align the instrument's connector pins with the three dots printed on the ForceTriad face plate below the monopolar receptacles. Reversing the connector limits the instrument's functionality and will damage the connector.

If a two-button instrument is connected to the energy platform, proceed to the *Standard Monopolar Mode Functionality* section that follows.

If a Valleylab-mode enabled three-button instrument is connected to the energy platform, proceed to *Valleylab Mode Functionality* on page 4-12.

Note: The Monopolar 1 and Monopolar 2 touchscreens can each control only one instrument at one time. If more than one instrument is attached under one touchscreen, an error message appears on the touchscreen and both devices are inactive until one of the devices is removed. See *Instrument Combination Not Allowed* on page 7-9.

Standard Monopolar Mode Functionality



- ① Power Output Mode Display (Watts)
- ② Green Power Control Arrows
- ③ Power Output Mode Waveforms
- ④ REM Indicator

When a two-button electrocautery instrument is attached to either the Monopolar 1 or Monopolar 2 receptacle, the ForceTriad Energy Platform detects the instrument type and displays the Standard Monopolar tab on the touchscreen. The Standard Monopolar tab allows the user to control the power mode and power output level within the energy platform interface.

1. Select the power output mode waveform by pressing the associated button at the bottom of the tab. The waveforms available in the cut mode are Pure and Blend. In the coag mode, the waveforms available are Fulgurate or Spray.
2. Set the power to the desired output level by pressing either the green up or down arrows. Power output is displayed in watts.
3. **Activate the cut mode** by pressing the yellow button on the electrocautery instrument. The cut display will illuminate yellow and a tone will sound for the duration of the activation.

Activate the coag mode by pressing the blue button on the electrocautery instrument. The coag display will illuminate blue and a tone will sound for the duration of the activation.

Valleylab Mode Functionality

Covidien instruments featuring the Valleylab mode are specialty devices some of which enable the surgeon to control ForceTriad Energy Platform output from the sterile field using a slider control on the instrument.

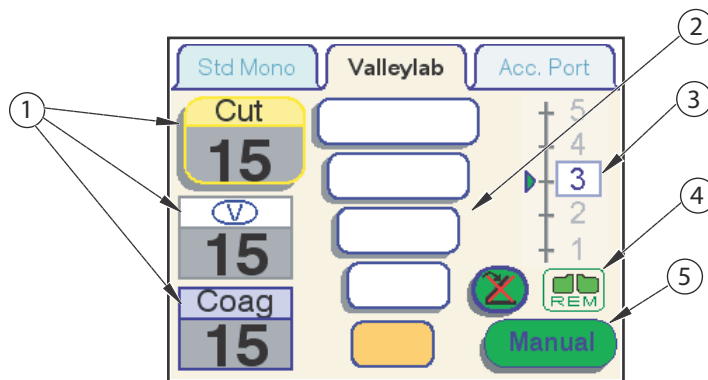
These output modes, if available on the instrument, are selected at the handset with the following buttons:

- The yellow cut button enables a cutting function.
- The clear Valleylab button enables a hemostasis function while providing dissection.
- The blue coag button enables a coagulation function.

If available on the instrument, a dual slider control switch adjusts power output in all available modes.

When a Valleylab-mode instrument's Smart connector is attached to either the Monopolar 1 or Monopolar 2 receptacle, the ForceTriad Energy Platform detects the instrument type and displays one of the following:

- Valleylab tab on the corresponding touchscreen:
 - Including a slider if the instrument has a slider
 - In manual mode if there is no slider (See *Manual Power Control Functionality* on page 4-14.)
- If the instrument can only use the Monopolar 2 receptacle but is plugged into Monopolar 1, an error message appears. See the *Invalid Instrument Alert for Monopolar 2* on page 7-10.



- ① Power Output Modes (Watts)
- ② Power Bar Indicators
- ③ Slider Position Indicator
- ④ REM Indicator
- ⑤ Manual Selector Button

Power Bars: The five gold bars in the center of the Valleylab tab represent the five power bars available for the particular instrument attached to the energy platform receptacle. The system automatically selects the default power bar setting for the particular instrument. Power bars can only be changed at the touchscreen interface on the energy platform.

Note: Refer to the individual instrument instructions for power-bar output in watts.

Slider Position (if available on the instrument): The energy platform detects the current position of the instrument's slider switch, and the slider position indicator on the right side of the Valleylab tab displays this slider position. Slider position can only be changed by the instrument user in the sterile field.

Using a Valleylab Mode Instrument

1. Select the desired power zone by touching the corresponding bar on the Valleylab-tab touchscreen. The bar that was touched, along with all the bars below it, will illuminate gold and a brief double tone will sound. Power output is displayed in watts. The power zone cannot be changed during instrument activation.

Warning

The slider, if available on the instrument, increases and decreases power output. Verify slider position prior to activation.

2. Activate power output by pressing the desired button on the instrument.
 - Activate the **cut mode** by pressing the yellow button on the electro-surgical instrument. The cut display will illuminate yellow and a tone will sound for the duration of the activation.
 - Activate the **Valleylab mode** by pressing the clear button on the electro-surgical instrument. The Valleylab display will illuminate white and a tone will sound for the duration of the activation.
 - Activate the **coag mode** by pressing the blue button on the electro-surgical instrument. The coag display will illuminate blue and a tone will sound for the duration of the activation.
3. If the instrument has a slider, the power output can be changed while in the sterile field by adjusting the slider position on the electro-surgical instrument. A double tone will sound when slider position is changed. Slider position cannot be changed while RF energy is being delivered.

If the instrument has no slider, use manual controls on the touchscreen. See *Manual Power Control Functionality* on page 4-14.

Cut Mode Disable


As a safety feature, power output can be disabled in the cut mode.

1. Disable the cut mode by pressing the yellow cut output mode display box. A '--' will replace the numeric digits in the cut box.
2. Re-enable the cut mode by pressing the yellow cut output-mode display box. The cut box will display the power setting of the current instrument slider position. The cut mode is also re-enabled when the energy platform is restarted.

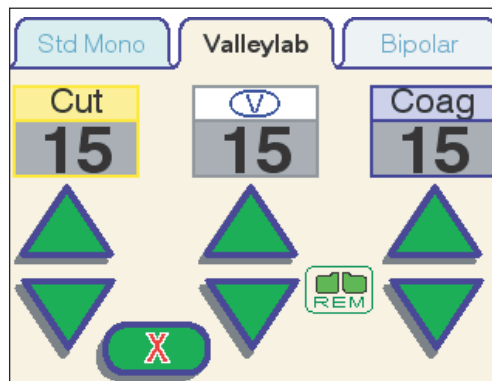
Manual Power Control Functionality

Manual mode allows the user to operate the energy platform outside the preset power ranges available on the Valleylab tab. In manual mode, if there is a slider on the electro-surgical instrument, it is disabled, and power can only be set at the touchscreen interface.


Note: Manual mode is the only available mode for instruments with no slider.

1. To set the system to manual mode for an instrument with a slider, press the green Manual button  on the Valleylab tab touchscreen. One of two manual-control displays appear on the touchscreen depending on which Valleylab-mode instrument is used.

An instrument with a slider displays this screen when placed in manual mode:



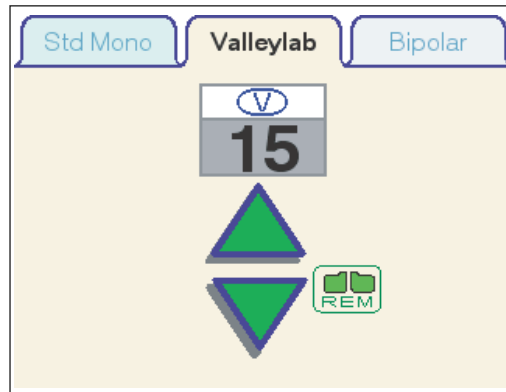An instrument with Valleylab mode but no slider displays this screen:




2. Adjust the power output for cut, Valleylab mode, and coag, if available, by pressing the associated green up or down buttons on the energy platform touchscreen.
3. For instruments with a slider, return to the sterile-field control mode by touching the "X" button. The Valleylab mode control screen re-appears on the touchscreen, and the power zone and slider position resets based on the current instrument configuration.

Accessory Port Functionality

The **Acc. Port** tab show accessory port options. Instruments with 8 mm pins connect directly to the accessory port on the ForceTriad Energy Platform. Instruments with pin diameters less than 8 mm may require the use of an E05021 or E0017 adapter. Remove the adapter when not in use.

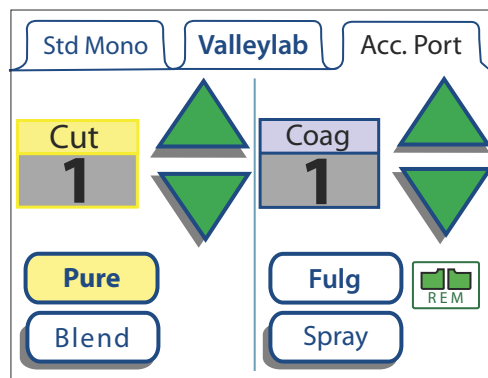When a single pin electrosurgical instrument is attached to the universal footswitching accessory port, the ForceTriad Energy Platform detects the instrument and displays the Accessory Port tab on the touchscreen.



The Accessory Port tab allows the user to control the power mode and power output level at the system interface for any footswitching instrument connected.

1. Select the power output mode waveform by pressing the associated button. The waveforms available in the cut mode are Pure and Blend. In the coag mode, the waveforms available are Fulgurate or Spray.
2. Set the power to the desired output level by pressing the green up and down arrow. Power output is displayed in watts.
3. Activate the cut mode by stepping on the cut or yellow pedal on the footswitch. The cut display will illuminate yellow and a tone will sound for the duration of the activation.

Activate the coag mode by stepping on the coag or blue pedal on the footswitch. The coag display will illuminate blue and a tone will sound for the duration of the activation.

Chapter 5

Bipolar Surgery

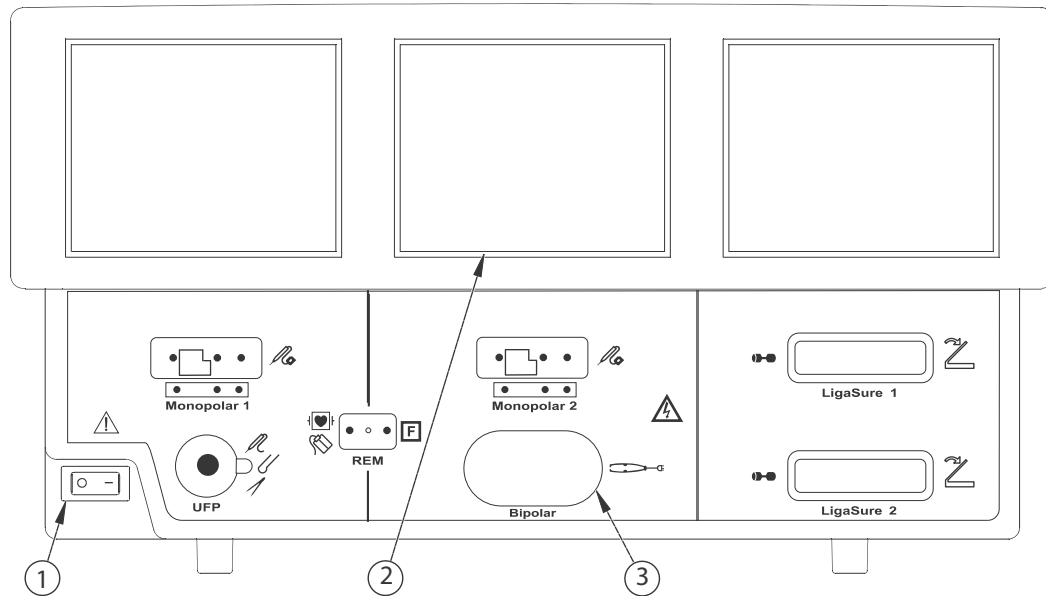
This chapter describes the bipolar-surgery features of the ForceTriad Energy Platform.

Caution

Read all warnings, cautions, and instructions provided with this system before use.


Read the instructions, warnings, and cautions provided with electro-surgical instruments before use. Specific instructions for electro-surgical instruments are not included in this manual.

Front Panel Bipolar Features



- ① Power Switch
- ② Bipolar Touchscreen
- ③ Bipolar Instrument Receptacle

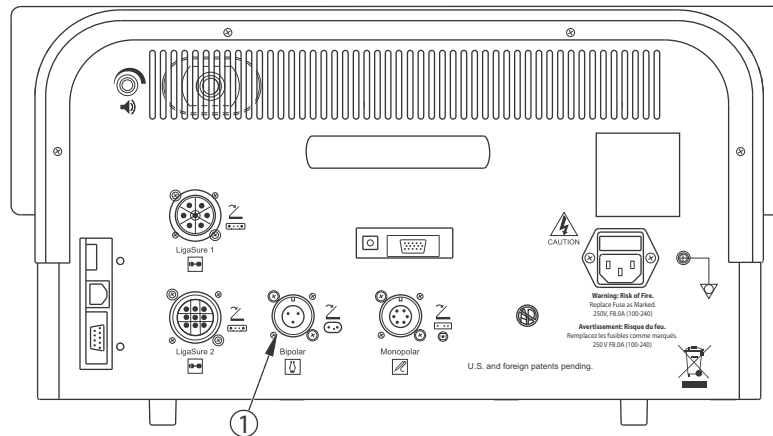Rear Panel Bipolar Features



- ① Bipolar Footswitch Receptacle

Bipolar Quick Setup Instructions

If you are familiar with the ForceTriad Energy Platform, follow this abbreviated procedure to set up the system for bipolar surgery.

If you are not familiar with the ForceTriad Energy Platform, refer to the following sections in this chapter for detailed instructions.

1. Plug the system power cord into the rear panel receptacle.
2. Plug the system power cord into a grounded wall receptacle.

Important

Do not plug into a power strip or extension cord.

3. Turn on the system and verify that the self-test is successfully completed.
4. If using a footswitch, connect it to the bipolar footswitch receptacle on the rear panel.

Warning

Use only Covidien footswitches such as the E6009 or E6019. Use of other manufacturer's footswitches is not recommended.

5. Connect the instrument to the bipolar instrument receptacle on the front panel.
6. Verify or change the mode and power settings.

Bipolar Function Overview

Delicate tissue requires less energy to desiccate. The ForceTriad Energy Platform provides low-voltage, continuous current for faster desiccation without sparking.

The possibility of sparking increases as desiccated tissue dries and becomes more resistant to energy flow. The system protects against sparking by limiting the bipolar voltage at relatively high levels of tissue resistance.

Bipolar Power Output Modes

The ForceTriad Energy Platform produces three bipolar modes: Low, Standard, and Macro.

Warning

Electric Shock Hazard

- Do not connect wet instruments to the energy platform.
- Ensure that all instruments and adapters are correctly connected and that no metal is exposed at any connection point.

Connect instruments to the proper receptacle. Improper connection may result in inadvertent instrument activation or other potentially hazardous conditions. Follow the instructions provided with electro-surgical instruments for proper connection and use.

The instrument receptacles on this system are designed to accept only one instrument at a time.

Caution

Read the instructions, warnings, and cautions provided with electro-surgical instruments before use. Specific instructions are not included in this manual.

Inspect instruments and cords (especially reusable instruments and cords) for breaks, cracks, nicks, and other damage before every use. If damaged, do not use. Failure to observe this precaution may result in injury or electrical shock to the patient or surgical team.

Footswitch

The ForceTriad Energy Platform in bipolar mode can accommodate a three-pin, single-pedal bipolar footswitch.

If footswitching is to be used with a bipolar instrument, attach the bipolar footswitch connector plug to the bipolar footswitch receptacle on the rear panel.

Warning

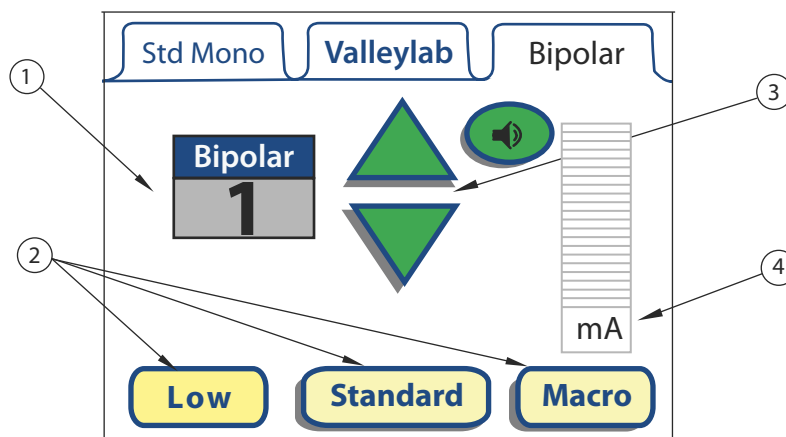
Use only Covidien footswitches such as the E6009 or E6019. Use of other manufacturer's footswitches is not recommended.

Bipolar Electrode Function

1. Connect a bipolar instrument to the bipolar instrument receptacle on the front panel.

Note: The Bipolar or Monopolar touchscreens can control only one instrument at one time. If more than one instrument is attached below this touchscreen, an error message appears on the touchscreen. Both devices are inactive until one of the devices is removed.

When a bipolar instrument is connected, the Bipolar tab appears on top on the center touchscreen.



- ① Power Output Modes (Watts)
 - ② Power Output Mode Waveforms
 - ③ Green Power Control Arrows
 - ④ Virtual Ammeter
2. Select the power output mode waveform by pressing the associated button at the bottom of the tab. The waveforms available in the bipolar mode are Low, Standard, and Macro.
3. Set the power to the desired output level by pressing the green up and down arrows. Power output is displayed in watts.

Virtual Ammeter

The virtual ammeter on the bipolar tab displays the current being delivered during bipolar-instrument activation. The ammeter registers current between 1 and 1,050 milliamps. An audio tone sounds to indicate increases and decreases in current delivery. A mute button allows the user to silence the ammeter tone, but not the activation tone.

4. Activate bipolar mode:


- With handswitching forceps, close the forceps firmly.
- With footswitching forceps, step on the single-pedal footswitch.

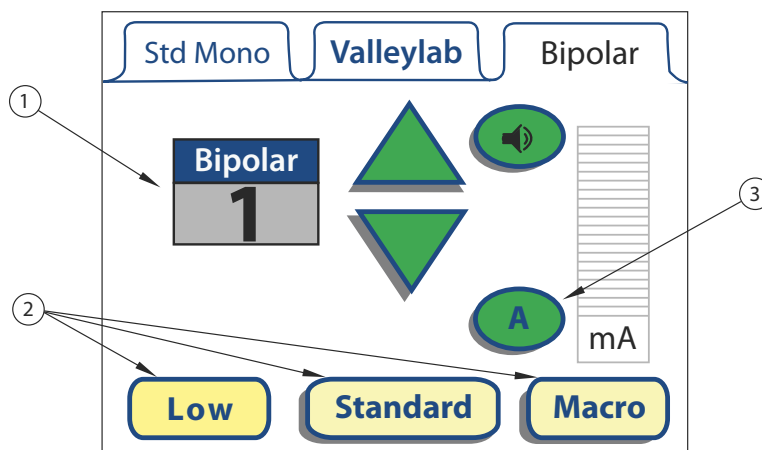
An activation tone sounds, and the delivered current is displayed on the virtual ammeter.


Autobipolar Function

The ForceTriad Energy Platform is equipped with an autobipolar feature that allows the user to configure the system for automatic activation and cessation of bipolar energy.

Autobipolar mode must be enabled in the feature menu before the autobipolar function can be used. Refer to *Autobipolar* on page 3-6 for instructions on how to enable autobipolar mode.

When the autobipolar mode is enabled, a green button “A” button  appears on the Bipolar tab.



- ① Power Output Modes (Watts)
- ② Power Output Mode Waveforms
- ③ “A” Button  (AutoBipolar) Mode

Activation Impedance


The impedance of the tissue grasped by the bipolar device must be within a factory-set range of less than 1000 Ω to activate. The activation impedance range is a safety feature that prevents power delivery if the grasped tissue is not within the anticipated range.

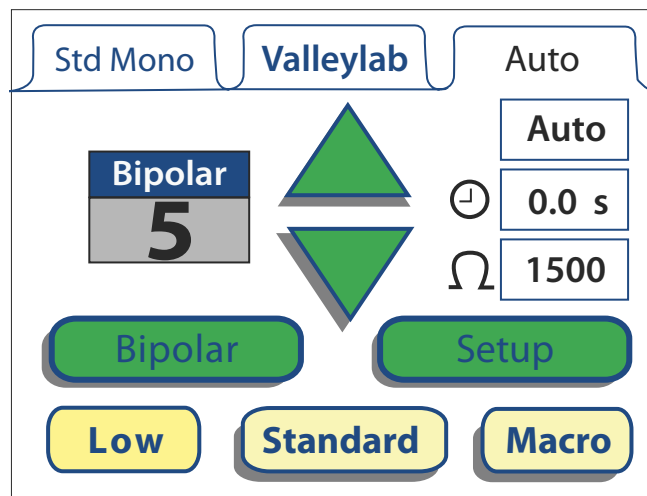
Tissues such as scar tissue, moles, and others with low blood flow are likely to have impedance greater than 1000 Ω . While grasping a high-impedance tissue, the ForceTriad will not activate regardless of customer-set preferences.

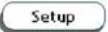
Note: The autobipolar function requires the use of these Covidien reusable footswitching bipolar cords: E0020V, E0021S, E0022W, E360150, or E360150L.


Warning

Use of different Covidien cord models or cords from other manufacturers may not achieve proper electrical output for this device, thereby failing to produce the desired clinical effect. For example, Autobipolar activation/deactivation settings may not work properly using cords other than those specified by Covidien.

1. Press the "A" button . The tab title will change from Bipolar to Auto, and the screen below will appear.



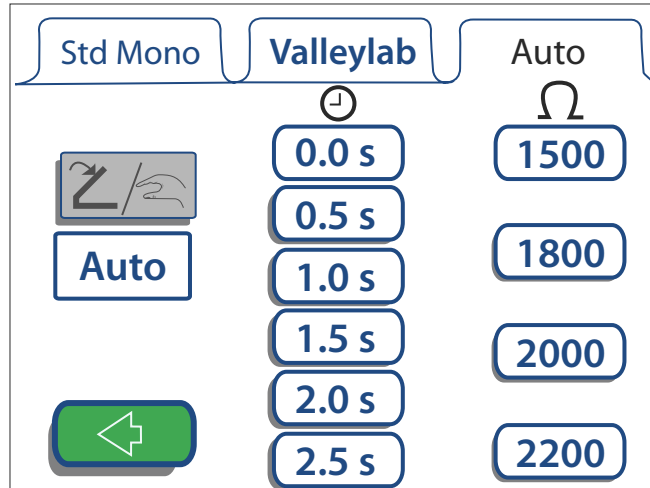2. Select the power-output mode waveform by touching the associated button at the bottom of the tab. The waveforms available in the autobipolar mode are Low, Standard, and Macro.
3. Set the power to the desired output level by touching the green up and down arrows. Power output is displayed in watts. The minimum setting for autobipolar is 5 W.
4. To change the autobipolar activation parameters, touch the green Setup button . The setup display will appear in the Auto tab.



5. To restrict RF activation to the footswitch, touch the footswitch/handswitch button .
- or**
- To enable RF activation without needing to depress the footswitch, touch the Auto button .
6. Set the desired RF output delay by touching one of the six delay times available under the timer symbol .
7. Set the desired impedance level at which RF energy will be discontinued by touching one of the four impedance values under the impedance symbol .
8. Touch the green back arrow button to return to the autobipolar activation display. The settings selected in the setup display will appear in the autobipolar display.
9. Activate autobipolar by closing the instrument forceps tines firmly or, if using a footswitch, by stepping on the footswitch pedal. An activation tone will sound.

Important

Autobipolar activation will occur only if the tissue in the forceps is in the activation impedance less than 1000 Ω. The device will not activate if the tissue has an impedance outside this range.

10. To return the energy platform to bipolar functionality, touch the green Bipolar button. The bipolar tab will replace the autobipolar tab.

or

Turn the system off. The next time the system is turned on, it will default to the bipolar function, and the bipolar tab will be visible.

Chapter 6

LigaSure Tissue Fusion

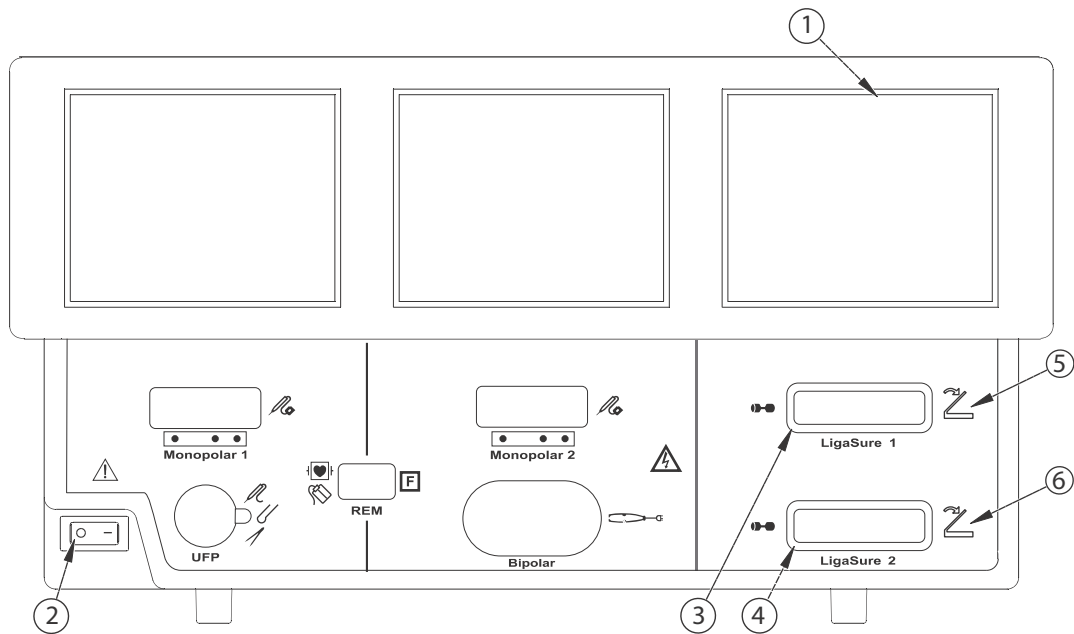
This chapter describes how to set up and operate the LigaSure tissue-fusion function of the ForceTriad Energy Platform.

Caution

Read all warnings, cautions, and instructions provided with this system before use.


Read the instructions, warnings, and cautions provided with electro-surgical instruments before use. Specific instructions for electro-surgical instruments are not included in this manual.

Front Panel LigaSure Features



- ① LigaSure Touchscreen
- ② Power Switch
- ③ LigaSure 1 Receptacle (purple)
- ④ LigaSure 2 Receptacle (orange)
- ⑤ Footswitch Symbol (purple)
- ⑥ Footswitch Symbol (orange)

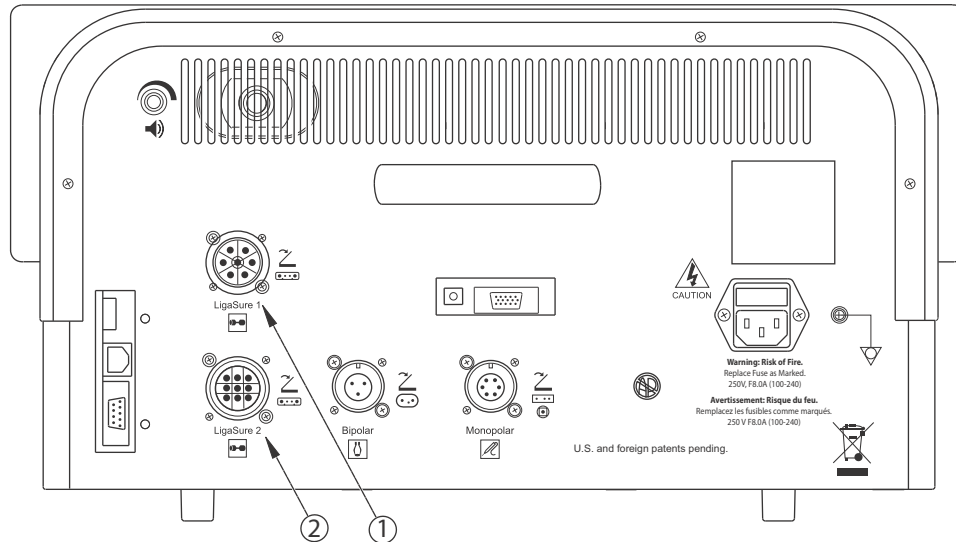Rear Panel LigaSure Features



- ① LigaSure 1 Footswitch Receptacle (purple)
- ② LigaSure 2 Footswitch Receptacle (orange)

LigaSure Quick Setup Instructions

If you are familiar with the ForceTriad Energy Platform, follow this abbreviated procedure to setup the system for LigaSure tissue fusion.

If you are not familiar with the ForceTriad Energy Platform, refer to the following sections in this chapter for detailed instructions.

1. Plug the system power cord into the rear panel receptacle.
2. Plug the system power cord into a grounded wall receptacle.

Important

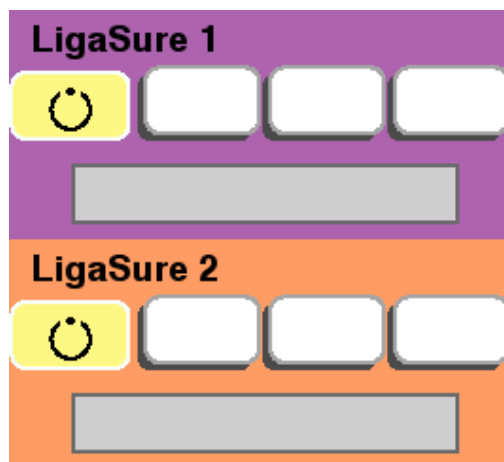
Do not plug into a power strip or extension cord.

3. Turn on the system and verify that the self-test has successfully completed.
4. If using a footswitch, connect it to the appropriate LigaSure footswitch receptacle on the rear panel.
5. Connect the instrument or instruments to the LigaSure instrument receptacles on the front panel.
6. Verify the bar setting.

LigaSure Function Overview

The LigaSure tissue-fusion mode can be used on arteries, veins, and lymphatics up to and including 7 mm in diameter and tissue bundles. This system provides precise energy delivery and electrode pressure to tissues for a controlled time period to achieve a complete and permanent fusion of tissues and vessel lumens. The system has been designed to produce minimal sticking, charring or thermal spread to adjacent tissue.

The LigaSure touchscreen is divided into two functional parts: the LigaSure 1 control panel, which controls instruments connected to the LigaSure 1 receptacle; and the LigaSure 2 control panel, which controls instruments connected to the LigaSure 2 receptacle.




Two LigaSure instruments can be set up at one time through this touchscreen, but only one instrument can be activated at one time.

During instrument activation, the status bar in the corresponding LigaSure control panel illuminates blue, and a tone sounds for the duration of energy delivery.

Important

If using an instrument that causes the hand-activation button to appear on the LigaSure control panel, see *Hand-Activation On/Off Button* on page 6-8.

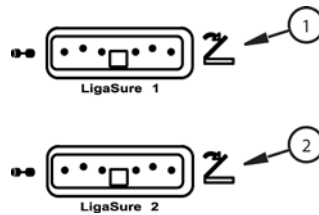LigaSure Receptacles



- ① Purple-outlined receptacle and purple footswitch symbol
- ② Purple-outlined receptacle and orange footswitch symbol

LigaSure 1 Receptacle

The LigaSure 1 receptacle is located directly below the LigaSure touchscreen and is outlined in purple with a purple footswitch icon to the right of it. This receptacle accepts all LigaSure instruments and can read either the dot patterns or bar codes on the Smart connector. Instruments attached to the LigaSure 1 receptacle are controlled from the upper, purple section of the LigaSure touchscreen.

LigaSure 2 Receptacle

The LigaSure 2 receptacle is located directly below the LigaSure 1 receptacle under the LigaSure touchscreen and is outlined in purple with an orange footswitch icon to the right of it. This receptacle accepts all LigaSure instruments and can read either the dot patterns or bar codes on the Smart connector. Instruments attached to the LigaSure 2 receptacle are controlled from the lower, orange section of the LigaSure touchscreen.

Warning

Electric Shock Hazard

- Do not connect wet instruments to the energy platform.
- Ensure that all instruments are correctly connected and that no metal is exposed at any connection point.

Connect instruments to the proper receptacle. Improper connection may result in inadvertent instrument activation or other potentially hazardous conditions. Follow the instructions provided with LigaSure instruments for proper connection and use.

The instrument receptacles on this system are designed to accept only one instrument at a time.

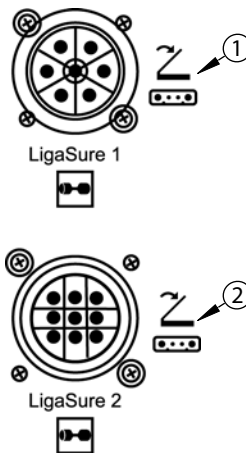
Caution

Read the instructions, warnings, and cautions provided with LigaSure instruments before use. Specific instructions are not included in this manual.

Caution

Inspect instruments and cords (especially reusable instruments and cords) for breaks, cracks, nicks, and other damage before every use. If damaged, do not use. Failure to observe this precaution may result in injury or electrical shock to the patient or surgical team.

Footswitch



① Purple

② Orange

The ForceTriad Energy Platform in LigaSure mode can accommodate two single-pedal LigaSure footswitches.

If the LigaSure instrument connected to the LigaSure 1 instrument receptacle is to be activated with a footswitch, attach the purple, seven-pin, LigaSure footswitch connector plug to the purple LigaSure 1 footswitch receptacle on the rear panel.

If the LigaSure instrument connected to the LigaSure 2 instrument receptacle is to be activated with a footswitch, attach the orange, nine-pin, LigaSure footswitch connector plug to the orange LigaSure 2 footswitch receptacle on the rear panel.

To use a footswitch with an instrument that causes the hand-activation button to appear on the LigaSure control panel, turn off hand activation (handswitching). See *Hand-Activation On/Off Button* on page 6-8.

Start System

1. Plug the system power cord into the rear panel receptacle.
2. Plug the system power cord into a grounded wall receptacle.

Important

Do not plug into a power strip or extension cord.

3. Turn on the system by pressing the power switch on (I).

LigaSure Instruments

Reusable-Instrument Assembly

Warning

LigaSure instruments that require single-use electrodes must be used with the correct electrode type. Use of these instruments with any other electrodes could result in injury to the patient or surgical team, or cause damage to the instrument.

To prepare the reusable LigaSure instruments to be used for the procedure, refer to the following steps for general preparation. To find detailed directions for each instrument, please refer to the individual instrument instructions that accompany each instrument electrode.

1. Slip the base of the disposable electrode's white shaft onto the retaining post on the instrument ring handle.
2. Snap the body of the white electrode shaft onto the instrument handle. The white shaft of the electrode must be completely flush on the reusable-instrument shaft.
3. Snap each electrode into the appropriate instrument jaw, matching electrode curvature to jaw curvature. Insert the proximal pin first. Verify that there is no gap between the electrode and the instrument jaw.

Important

Bent or broken electrode pins will not function properly and may result in an alert situation. In this case, the electrode must be discarded.

4. Gently ratchet the instrument closed on a folded 4x4 to ensure the electrodes are properly seated in the instrument jaws.

Connecting LigaSure Instruments to the Energy Platform

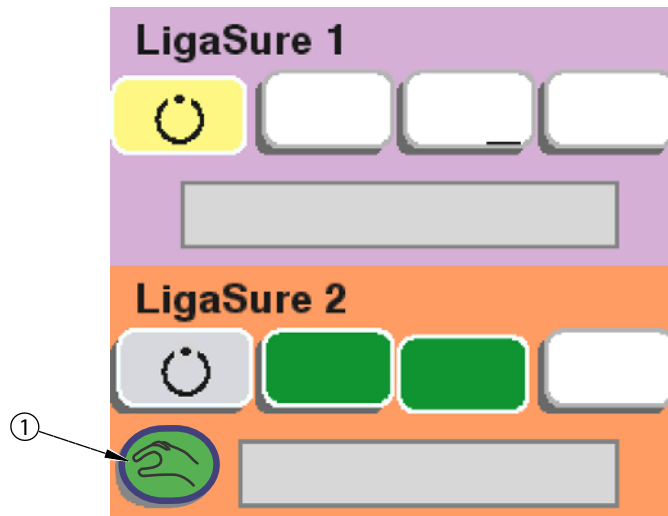
Connect the Smart connector to the LigaSure 1 or LigaSure 2 receptacle on the front panel of the ForceTriad Energy Platform. The system detects and sets the appropriate bar setting in the display. If settings have been entered in the LigaSure touchscreen prior to connecting a LigaSure instrument, these settings will be reset. If applicable, refer to the following sections:

- *Hand-Activation On/Off Button*
- *Invalid Instrument*

Hand-Activation On/Off Button

Some LigaSure instruments cause the hand-activation button to appear in the lower left corner of the LigaSure control panel. The hand-activation button toggles between the on and off setting.

For example, the LF1212 Curved, Small Jaw, Open Sealer/Divider, is a LigaSure instrument that can be used with hand-activation on or off. When hand-activation is off, a footswitch is used to initiate energy delivery.



① Hand-Activation-On Button

In this example, the hand-activation-on button appears in the lower left corner of the LigaSure 2 control panel. Touching the hand-activation-on button toggles to the hand-activation-off button:



Invalid Instrument

If the ForceTriad Energy Platform does not recognize the attached instrument, the status bar will display the words "INVALID INSTRUMENT." Refer to the following steps to resolve the issue.

1. Confirm that a LigaSure instrument is in use.
2. Reconnect the instrument using firm pressure to insert the instrument into the LigaSure 1 or LigaSure 2 instrument receptacle.
3. If "INVALID INSTRUMENT" continues to appear in the status bar, use a new LigaSure instrument or electrode.

LigaSure Settings

Changing the Energy-Delivery Setting

Warning

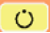
Confirm proper power or intensity settings before proceeding with surgery.

The green bars on the LigaSure vessel-fusion display panel represent different levels of desiccation. Two green bars are the default setting for all LigaSure tissue-fusion instruments. This setting allows the energy to precisely affect the target tissue, resulting in permanent tissue fusion with minimal thermal spread to the surrounding tissues.

Occasionally the surgeon may encounter tissue or vessels that are fused more effectively with a one- or three-bar fusion cycle. By selecting one bar, the surgeon can expect a more gentle and typically longer fusion cycle, which can be potentially more effective in thinner tissue bundles and smaller, isolated vessels. By selecting three bars, the surgeon can expect a longer fusion cycle, which can be potentially more effective in thicker tissue bundles.

Notice

The one- and three-bar seal cycles may result in tissue sticking due to longer desiccation periods.

1. The instrument setting can be adjusted by touching one of the three setting buttons on the respective LigaSure 1 or LigaSure 2 control panel. The button touched along with the buttons to the left of it become green, and the standby button turns grey.
2. As a safety feature, the LigaSure mode can be set on standby pressing the standby button  until the instrument is needed.
When on standby, no energy will be delivered through the LigaSure instrument. Attempting to activate the instrument will sound a single, short tone.
3. To bring the energy platform out of standby, press a desired bar setting. Previous bar settings will not be stored in standby.

Activating the LigaSure Instrument

1. Activate the LigaSure instrument either by pressing and holding the activation button on the instrument or by stepping on and holding the footswitch pedal. During instrument activation, the status bar in the corresponding LigaSure control panel illuminates blue, and an activation tone sounds for the duration of energy delivery.
2. When the end tone is heard, release the activation button or footswitch pedal. In the case of an alert condition, refer to the following section.

Alert Situations

A pulsed tone will sound when an alert condition occurs, and the LigaSure touchscreen will display an alert message that instructs the user on the corrective actions to take. When an alert condition occurs, energy delivery will be interrupted, but will be available immediately after the alert condition has been corrected.

The two alert conditions are:

- Check Instrument
- Reactivate

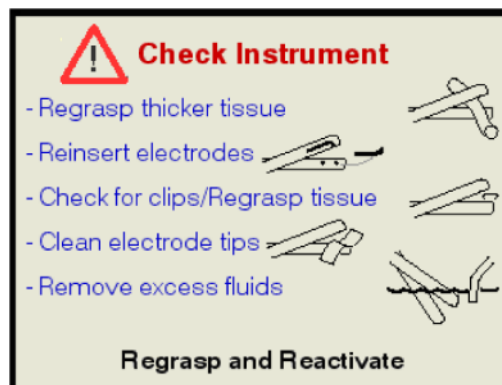
Check Instrument

A **six-pulsed** tone will sound when the Check Instrument screen is displayed.

If this message appears, the user should:

1. Release the footswitch pedal or activation button.
2. Open the instrument jaws and inspect for a successful seal.
3. Follow the suggested corrective action on the Check Instrument screen.

If possible, reposition the instrument and regrasp tissue in another location, then reactivate the seal cycle.



Regrasp thicker tissue – Thin tissue; open the jaws and confirm that a sufficient amount of tissue is inside the jaws. If necessary, increase the amount of tissue and repeat the procedure.

Reinsert electrodes – Electrodes may have become dislodged from the instrument.

Check for clips / Regrasp tissue – Avoid grasping objects, such as staples, clips or encapsulated sutures, in the jaws of the instrument.

Clean electrode tips – Use a wet gauze pad to clean surfaces and edges of instrument jaws.

Remove excess fluids – Pooled fluids around the instrument tip; minimize or remove excess fluids.

Reactivate

A **four-pulsed** tone will sound when the Reactivate screen is displayed.

If this message appears the user should:

1. Release the footswitch pedal or hand switching button.
2. Reactivate the seal cycle without repositioning the instrument.



Endpoint not reached – Additional time and energy are needed to complete the fusion cycle.

Seal cycle interrupted – The seal cycle was interrupted before completion. The handswitch or footswitch was released before the end tone activated.

After Surgery

Disconnect the instruments

1. Turn off the energy platform.
2. Disconnect all instruments from the front panel.
 - If the instrument is single use only (disposable), dispose of it according to the procedures for your institution.
 - If the instrument is reusable, clean and sterilize it according to the manufacturer's instructions for use.
3. Disconnect and store any footswitches used.

Reprocessing Instruments

Clean the Reusable LigaSure Instrument

1. Remove and dispose of single-use electrodes.
2. Wipe all surfaces with a cleaning agent and a damp cloth.
3. Follow the procedures approved by your health-care facility.
4. Soak in an enzymatic cleaning agent, such as Prolystica[®] or Enzol[®], according to the manufacturer's instructions.
5. Scrub all surfaces with a soft brush. It is important that the jaw surfaces and instrument electrode holes are cleaned of blood and tissue to ensure proper electrode assembly.
6. Rinse with water and dry with a soft cloth.

Sterilization Parameters

The hinges on reusable LigaSure instruments are extremely tight and require longer sterilization times to ensure steam penetration into the hinge.

Steam Sterilization - Wrapped

Temperature	Type	Sterilize Time	Dry Time
270° - 280° F (132° - 138° C)	Prevac	10 min.	20 min.
270° - 280° F (132° - 138° C)	Gravity	15 min.	30 min.
250° - 268° F (121° - 131° C)	Gravity	30 min.	30 min.

Steam Sterilization - Unwrapped

Temperature	Type	Sterilize Time	Dry Time
270° - 280° F (132° - 138° C)	Prevac	10 min.	1 min.
270° - 280° F (132° - 138° C)	Gravity	15 min.	1 min.

The instructions provided above have been validated by Covidien as being capable of preparing LigaSure instruments for reuse that are *not* labeled for single use only. It remains the responsibility of the processor to ensure that sterilization is performed using equipment, materials, and personnel that will achieve the desired results.

This requires validation and routine monitoring of the process. Any deviation by the processor from the instructions provided should be properly evaluated for effectiveness and potential adverse consequences.

Chapter 7

Troubleshooting

Caution

Read all warnings, cautions, and instructions provided with this system before use.

Read the instructions, warnings, and cautions provided with electro-surgical instruments before use. Specific instructions for electro-surgical instruments are not included in this manual.

General Troubleshooting Guidelines

If the ForceTriad Energy Platform malfunctions, check for obvious conditions that may have caused the problem:

- Check the system for visible signs of physical damage.
- Make sure the fuse drawer is tightly closed.
- Verify that all cords are connected and attached properly.
- If an error code is displayed on the touchscreens, note the code along with all information on the error screen and use the information in this chapter to proceed.
- If the problem is still unclear, turn the ForceTriad Energy Platform off, wait a minute or two, then turn it back on.

If the malfunction persists, the system may require service. Contact your institution's biomedical engineering department. Covidien Technical Service may be contacted in the ways indicated on page 8-4.

REM Alarms

If the ForceTriad Energy Platform does not sense the correct impedance for the connected REM Polyhesive patient return electrode, monopolar energy will be disabled, the REM symbol will illuminate red and enlarge on both the center and left touchscreen displays, and an alarm tone will sound twice. The REM symbol will return to its smaller size but will remain red, and RF energy will remain disabled until the REM alarm is corrected.

When you correct a REM alarm condition, the system is enabled and the REM-alarm indicator illuminates green.

Covidien recommends the use of REM Polyhesive patient return electrodes. Return electrodes from other manufacturers may not provide proper impedance to work correctly with the ForceTriad Energy Platform.

Correcting a REM-Alarm Condition

To correct a REM-alarm condition, follow these steps:

1. Inspect the return electrode plug and cord. If you find evidence of cracks, breaks, or other visible damage, replace the return electrode and/or the cord.
2. Verify that the patient return electrode cord is correctly connected to the energy platform.
3. Verify that the return electrode is in good contact with the patient. Follow the instructions for use provided with the REM Polyhesive patient return electrode.
4. If the REM alarm persists it may be necessary to use more than one patient return electrode. Refer to the troubleshooting flow chart in the REM Polyhesive patient return electrode instructions for use.

Note: For more detailed suggestions, refer to the troubleshooting instructions provided in the instructions for use that came with the REM Polyhesive patient return electrodes.

Correcting Malfunctions

If a solution is not readily apparent, use the table below to help identify and correct specific malfunctions. After you correct the malfunction, verify that the system completes the self-test as described in Chapter 3, *System Setup*.

Situation	Possible Cause	Solution
Abnormal neuromuscular stimulation (<i>stop surgery immediately</i>)	Metal-to-metal sparking	Check all connections to the energy platform, patient return electrode, and active electrodes.
	Can occur during coag	Use a lower power setting for the Fulgurate and the Spray modes.
	Abnormal 50 Hz-60 Hz leakage currents	Contact your biomedical engineering department or a Covidien technical service representative for assistance.
Energy platform does not respond when turned on	Disconnected power cord or faulty wall outlet	Check power cord connections (energy platform and wall outlet). Connect the power cord to a functional outlet.
	Faulty power cord	Replace the power cord.
	Fuse drawer is open or fuses are blown.	Replace the blown fuse(s). Close the fuse drawer. Refer to the <i>ForceTriad Energy Platform Service Manual</i> .
	Internal component malfunction	Use a backup energy platform. Contact your biomedical engineering department or a Covidien technical service representative for assistance.

Situation	Possible Cause	Solution
System is on, but did not complete the self-test	Software malfunction	Turn off, then turn on the system.
	Internal component malfunction	<p>Note the code along with all information on the error screen. Note the number and refer to <i>System Alarms</i> on page 7-8.</p> <p>Use a backup energy platform. Contact your biomedical engineering department or a Covidien technical service representative for assistance.</p>
Energy platform is on and instrument is activated, but system does not deliver output	Malfunctioning footswitch or handswitching instrument	<p>Turn off the energy platform. Check and correct all instrument connections.</p> <p>Turn on the energy platform. Replace the instrument if it continues to malfunction.</p>
	Power is set too low	Increase the power setting.
	An alarm condition exists	<p>Note the code along with all information on the error screen. Note the number and refer to <i>System Alarms</i> on page 7-8.</p> <p>In case of a REM alarm, refer to <i>Correcting a REM-Alarm Condition</i> on page 7-2.</p>
	Internal component malfunction	Contact your biomedical engineering department or a Covidien technical service representative for assistance.
	System does not detect tissue-fusion instrument	Firmly insert the LigaSmart connector into the appropriate receptacle on the energy platform front panel. Ensure the vessel fusion touchscreen indicates that it has detected the instrument.

Situation	Possible Cause	Solution
	System does not detect monopolar instrument	Firmly insert the Smart connector into the appropriate receptacle on the energy platform front panel. Ensure the monopolar touchscreen indicates that it has detected the instrument.
	System does not detect bipolar instrument	Firmly insert the connector into the appropriate receptacle on the energy platform front panel. Ensure the bipolar touchscreen indicates that it has detected the instrument.
CHECK INSTRUMENT screen appears, a six-pulsed tone sounds, and RF output is disabled	Excessive tissue/eschar on electrode tips or jaws	Clean electrode tips and jaws with a wet gauze pad.
	Electrodes have come loose from the instrument jaws Electrode pins may have been compromised or bent during assembly to the instrument and may need to be replaced	Re-insert the electrode into the instrument jaws making sure that all the electrode pins are firmly seated.
	Metal or other foreign object is grasped within jaws	Avoid grasping objects, such as staples, clips, or encapsulated sutures in the jaws of the instrument.
	Tissue grasped within jaws is too thin	Open the jaws and confirm that a sufficient amount of tissue is inside the jaws. If necessary, increase the amount of tissue and repeat the procedure.
	Pooled fluids around instrument tip	Minimize or remove excess fluids.
REACTIVATE screen appears, a four-pulsed tone sounds, and RF output is disabled	The seal cycle was interrupted before completion. The handswitch or footswitch was released before the end tone activated. Additional time and energy are needed to complete the fusion cycle	Reactivate the seal cycle without removing or repositioning the instrument.

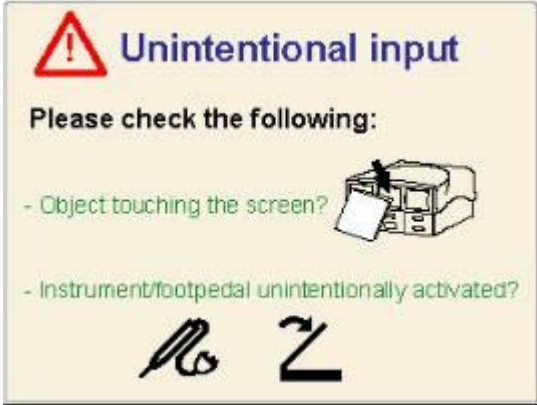

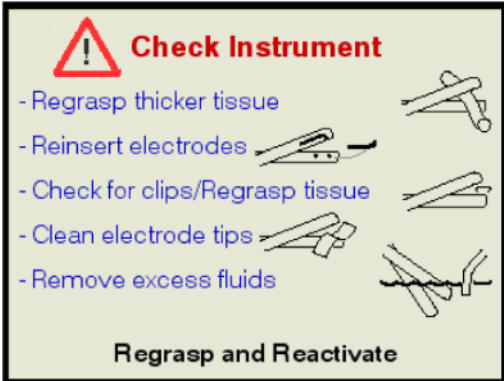
Situation	Possible Cause	Solution
Continuous monitor interference	Malfunctioning monitor	Replace the monitor.
	Faulty chassis-to-ground connections	<p>Check and correct the chassis ground connections for the monitor and for the energy platform.</p> <p>Check other electrical equipment in the room for defective grounds.</p>
	<p>Electrical equipment is grounded to different objects rather than a common ground. The energy platform may respond to the resulting voltage differences between grounded objects.</p>	<p>Plug all electrical equipment into line power at the same location. Contact your biomedical engineering department or a Covidien technical service representative for assistance.</p>
Interference with other devices only when the energy platform is activated	Metal-to-metal sparking	<p>Check all connections to the energy platform, patient return electrode, and instruments.</p>
	High settings used for fulguration	Use lower power settings for fulguration.
	Electrically inconsistent ground wires in the operating room	Verify that all ground wires are as short as possible and go to the same grounded metal.
	If interference continues when the energy platform is activated, the monitor is responding to radiated frequencies.	<p>Ask your biomedical engineering department to check with the manufacturer of the monitor.</p> <p>Some manufacturers offer RF choke filters for use in monitor leads. The filters reduce interference when the energy platform is activated and minimize the potential for an electrosurgical burn at the site of the monitor electrode.</p>




Situation	Possible Cause	Solution
Pacemaker interference	Intermittent connections or metal-to-metal sparking	<p>Check the active and patient return electrode cord connections.</p> <p>It may be necessary to reprogram the pacemaker.</p>
	<p>Current traveling from active to return electrode during monopolar electrosurgery is passing too close to pacemaker</p>	<p>Consult the pacemaker manufacturer or hospital cardiology department for further information when use of electrosurgical appliances is planned in patients with cardiac pacemakers.</p> <p>Use bipolar instruments, if possible.</p> <p>If you must use a monopolar instrument, place the patient return electrode as close as possible to the surgical site. Make sure the current path from the surgical site to the patient return electrode does not pass through the vicinity of the heart or the site where the pacemaker is implanted.</p> <p>Always monitor patients with pacemakers during surgery and keep a defibrillator available.</p>
Internal Cardiac Defibrillator (ICD) activation	ICD is activated by energy platform	Stop the procedure and contact the ICD manufacturer for instructions.




System Alarms





Most system alarms require some action on your part to correct the condition; however, some are corrected automatically. Use the following list to determine how to correct an alarm condition.

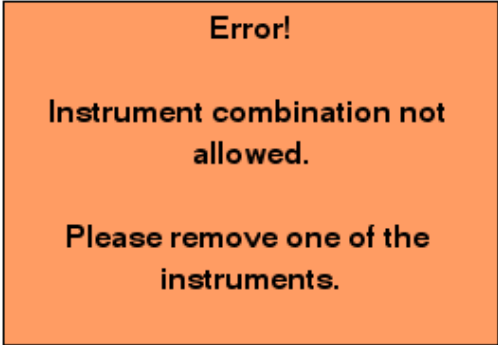After correcting the alarm condition, verify that the system completes the self-test as described in Chapter 3, *System Setup*.

Description or Screen	Solution
 <p>The screen displays a red warning triangle with an exclamation mark, followed by the text "Unintentional input". Below this, it says "Please check the following:" and lists two items: "- Object touching the screen?" and "- Instrument/footpedal unintentionally activated?". There is an illustration of a hand touching a screen and two footpedal symbols.</p>	<p>Check for stuck instrument buttons and footswitches.</p> <p>Check that nothing is bumping against the touchscreens.</p>
<p>REM Alarm</p>  <p>The icon consists of a yellow square with rounded corners, containing two red rectangular shapes representing the REM (Residual Energy Monitor) sensors, with the letters "REM" in red below them.</p>	<p>Refer to the <i>REM Alarms</i> on page 7-2 in this chapter.</p>
<p>LigaSure Check Instrument Alert</p>  <p>The screen displays a red warning triangle with an exclamation mark, followed by the text "Check Instrument". Below this, it lists five items: "- Regrasp thicker tissue", "- Reinsert electrodes", "- Check for clips/Regrasp tissue", "- Clean electrode tips", and "- Remove excess fluids". There are illustrations of various surgical instruments and electrode tips. At the bottom, it says "Regrasp and Reactivate".</p>	<p>Refer to <i>Alert Situations</i> on page 6-10.</p>

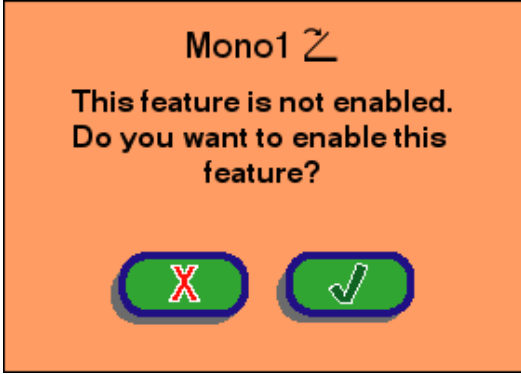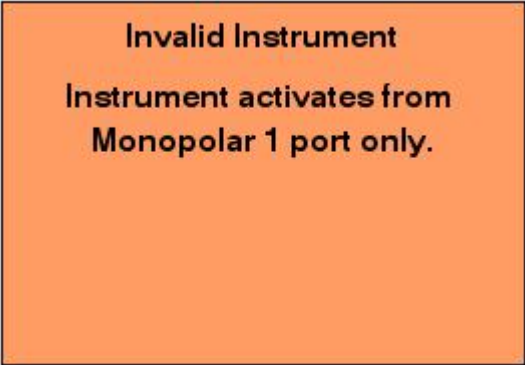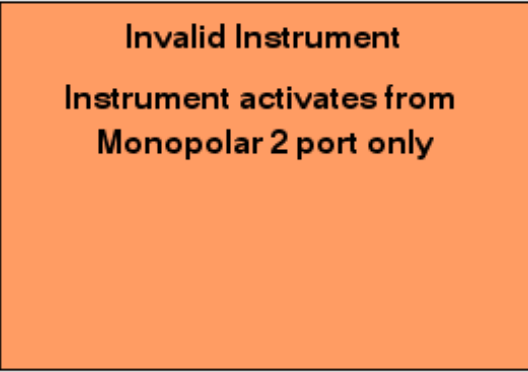Description or Screen	Solution
<p>LigaSure Endpoint Not Reached Alert</p> 	<p>Refer to <i>Alert Situations</i> on page 6-10.</p>
<p>LigaSure Advance™ Error Alert</p> 	<p>Check the two plugs on the instrument to ensure they are securely plugged into the correct receptacles. For LigaSure Advance instruments, the LigaSmart™ connector should be attached to either LigaSure receptacle, and the Smart connector should be attached to the Monopolar 2 receptacle under the center touchscreen. See <i>Monopolar Function</i> on page 4-1.</p>
<p>Instrument Combination Not Allowed</p> 	<p>The Monopolar 1 and Monopolar 2 touchscreens can each control only one instrument at one time. If more than one instrument is attached under one touchscreen, the following error message appears on the touchscreen and both devices are inactive until one of the devices is removed. See <i>Monopolar Electrodes</i> on page 4-10.</p>

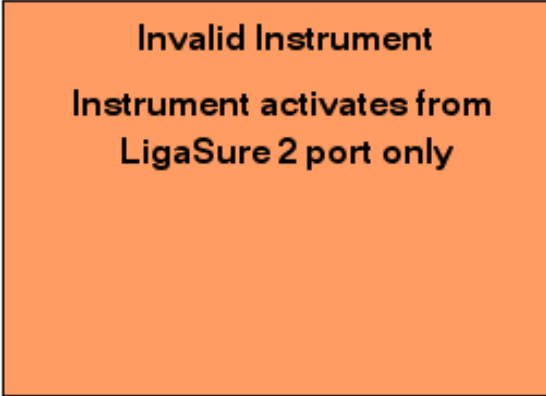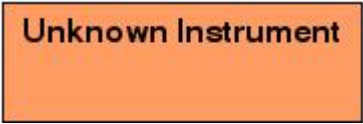Description or Screen	Solution
<p>Monopolar 1 Footswitch Not Enabled Alert</p> 	<p>Monopolar 1 Footswitching version of feature-not-enabled dialog box. Enable Monopolar 1 footswitch by touching the green-check button before attempting to use the footswitch again.</p>
<p>Invalid Instrument Alert for Monopolar 1</p> 	<p>Disconnect the instrument and attach it to Monopolar 1 port.</p>
<p>Invalid Instrument Alert for Monopolar 2</p> 	<p>Disconnect the instrument and attach it to the Monopolar 2 port. See <i>Valleylab Mode Functionality</i> on page 4-12.</p>

Description or Screen	Solution
<p data-bbox="321 401 711 428">Invalid Instrument Alert for LigaSure 2</p> 	<p data-bbox="906 401 1287 457">Disconnect the instrument connector and attach it to the LigaSure 2 port.</p>
<p data-bbox="321 886 597 913">Unknown Instrument Alert</p> 	<p data-bbox="906 886 1287 1031">Disconnect the unknown instrument and attach an instrument with a barcode recognized by the software version installed on the energy platform.</p>
<p data-bbox="321 1117 565 1144">Invalid Instrument Alert</p> 	<p data-bbox="906 1117 1287 1232">The instrument inserted into the port is not a valid instrument for that port. Refer to the instrument's instructions for use for compatibility instructions.</p>
<p data-bbox="321 1299 521 1327">Calibration Needed</p> 	<p data-bbox="906 1299 1287 1444">The ForceTriad Energy Platform needs to be calibrated. Consult the <i>ForceTriad Energy Platform Service Manual</i> for instructions regarding recalibration.</p>

Chapter 8

Maintenance and Repair

This chapter presents the following information:

- The manufacturer's responsibility
- Routine maintenance
- Returning the energy platform for service
- Service centers

Caution

Read all warnings, cautions, and instructions provided with this system before use.

Read the instructions, warnings, and cautions provided with electro-surgical instruments before use. Specific instructions for electro-surgical instruments are not included in this manual.

Responsibility of the Manufacturer

Covidien is responsible for the safety, reliability, and performance of the energy platform only if all of the following conditions have been met:

- Installation and setup procedures in this manual are followed.
- Assembly, operation, readjustments, modifications, or repairs are carried out by persons authorized by Covidien.
- The electrical installation of the relevant room complies with local codes and regulatory requirements, such as IEC and BSI.
- The equipment is used in accordance with the Covidien instructions for use.

For warranty information, refer to the *Limited Warranty* section in this manual.

Routine Maintenance and Periodic Safety Checks

Notice

Refer to the energy platform service manual for maintenance recommendations and function and output power verification procedures.

When should the energy platform be checked or serviced?

Covidien recommends that the energy platform be inspected by qualified service personnel at least once a year. This inspection should include adjusting the system to factory specifications.

When should the power cord be checked or replaced?

Check the power cord before each use of the system or at the intervals recommended by your institution. Check the power cord for exposed wires, cracks, frayed edges, or a damaged connector. Replace damaged cords.

When should the fuses be replaced?

An internal component malfunction can damage the fuses. The system fuses may need to be replaced if the system fails the self-test or if the system stops functioning, even though it is receiving power from a wall outlet. Refer to the *ForceTriad Energy Platform Service Manual* for instructions.

Cleaning

Warning

Electric Shock Hazard Always turn off and unplug the energy platform before cleaning.

Notice

Do not clean the energy platform with abrasive cleaning or disinfectant compounds, solvents, or other materials that could scratch the panels or damage the energy platform.

1. Turn off the system and unplug the power cord from the wall outlet.
2. Thoroughly wipe all surfaces of the energy platform and power cord with a damp cloth and mild cleaning solution or disinfectant. The energy platform will withstand the effects of cleaning over time without degrading the enclosure or display quality.

Product Service

Covidien recommends that all ForceTriad systems be returned to the manufacturer for all service requirements. If any service is required without returning the system to the manufacturer, Covidien recommends that only qualified personnel service the ForceTriad system.

Covidien defines qualified personnel as a person with electrosurgical equipment repair experience, such as biomedical personnel, and/or individuals who have taken official Covidien training courses.

Returning the Energy Platform for Service

Before returning the energy platform, call a Covidien sales representative for assistance. If you are instructed to send the energy platform to Covidien, do the following:

1. Obtain a return authorization number.

Call the Covidien Technical Service Center (see page 8-4) to obtain a Return Authorization Number. Have the following information ready before the call:

- Hospital/clinic name/customer number
- Telephone number
- Department/address, city, state, and zip code
- Model number
- Serial number
- Description of the problem
- Type of repair to be done

2. Clean the energy platform.
See the previous section, *Cleaning*.
3. Ship the energy platform.
 - a. Attach a tag to the energy platform that includes the return authorization number and the information (hospital, phone number, etc.) listed in step 1.
 - b. Be sure the energy platform is completely dry before packing it for shipment. Package it in its original shipping container, if available.
 - c. Ship the energy platform, prepaid, to the Covidien Service Center.

Adjustment to Factory Specification (Calibration)

Covidien recommends that only qualified personnel calibrate the ForceTriad Energy Platform. The energy platform incorporates automatic calibration where possible to reduce the required equipment and manual steps.

Software Upgrades

ForceTriad software upgrades are available over the internet directly from Covidien by using Valleylab Exchange. For more information, go to <http://www.valleylab.com> and click the link to Valleylab Exchange Software Update System.

Covidien Technical Service

For service, contact Covidien Technical Service or your ForceTriad sales representative. Contact a Covidien technical service representative by telephone, email, or through the Internet:

- USA and Canada: 1-800-255-8522 Option 2
- International: 1-303-476-7996
- Email: valleylab.technicalservice@covidien.com
- Internet: www.BioMedConnect.com

Chapter 9

Technical Specifications

All specifications are nominal and subject to change without notice. A specification referred to as "Typical" is within $\pm 20\%$ of a stated value at room temperature (77° F / 25° C) and a nominal line input voltage.

Caution

Read all warnings, cautions, and instructions provided with this system before use.

Read the instructions, warnings, and cautions provided with electro-surgical instruments before use. Specific instructions for electro-surgical instruments are not included in this manual.

Performance Characteristics

General

Output configuration	Isolated output
Cooling	Natural and forced convection, and fan
Display	Three touchscreens
Connector ports	LED illuminated Smart connector readers
Mounting	<ul style="list-style-type: none">• ForceTriad Energy Platform cart (FT900), Universal Mounting cart (UC8009), and/or the UC8010 Overshelf• Operating-room boom systems• Any stable, flat surface such as a table or cart top

Dimensions and Weight

Width	18 in. (45.7 cm)
Depth	20 in. (50.8 cm)
Height	10 in. (25.4 cm)
Weight	30 lb. (13.6 kg)

Operating Parameters

Ambient temperature range	50° F to 104° F (+10° C to +40° C)
Relative humidity	30% to 75% non-condensing
Atmospheric pressure	700 millibars to 1060 millibars
Warm-up time	If transported or stored at temperatures outside the operating temperature range, allow one hour for the energy platform to reach room temperature before use.

Transport and Storage

Ambient temperature range	-22° F to 149° F (-30° C to +65° C)
Relative humidity	25% to 85% (non-condensing)
Atmospheric pressure	500 millibars to 1060 millibars
Duration of storage	The ForceTriad Energy Platform may be stored indefinitely. If the energy platform is stored for over one year, the memory battery must be replaced and the system must be recalibrated as explained in the <i>ForceTriad Energy Platform Service Manual</i> .

Internal Memory

Nonvolatile, battery-backed RAM	Battery type – Lithium Battery life – 120 mAh
Storage capacity	256 KB

Activation Tone

The audio levels stated below are for activation tones (cut, Valleylab, coag, bipolar, and LigaSure modes) and alarm tones (REM and system alarms) at a distance of one meter.

Volume (adjustable)	45 dBA to 65 dBA
Frequency	Cut – 660 Hz Valleylab – 800 Hz Coag – 940 Hz Bipolar – 940 Hz LigaSure – 440 Hz
Duration	Continuous while the system is activated

Alarm Tone

Volume (not adjustable)	> 65 dBA
Frequency	REM – 660 Hz Reactivate/Regrasp, Check Instrument – Two tones High = 985 Hz, Low = 780 Hz Seal Complete – 985 Hz Error/System Alert – Beep tone = 1421 Hz
Duration	REM – Two 1/2 second tones separated by 1/2 second for each REM event Reactivate/Regrasp – Four 175 ms tones High, low, high, low Check Instrument – Six 175 ms tones High, low, high, low, high, low Seal Complete – Two 175 ms tones separated by 175 ms for each seal-complete event Error/System Alert – Three 250 ms tones separated by 250 ms for each error/system-alert event

REM Contact Quality Monitor

Interrogation frequency	80 kHz \pm 10 kHz
Interrogation current	< 100 μ A
Interrogation voltage	< 12 V RMS

Acceptable Resistance Range

REM resistance measurements are \pm 10% during RF activation and \pm 5% when RF output is not activated.

REM Polyhesive patient return electrode: 5 Ω to 135 Ω or up to a 40% increase in the initial measured contact resistance (whichever is less).

If the measured resistance is outside the acceptable range(s) noted above, a REM fault condition occurs.

REM Alarm Activation

REM Polyhesive patient return electrode: When the measured resistance exceeds the standard range of safe resistance (below 5 Ω or above 135 Ω) or when the initial measured contact resistance increases by 40% (whichever is less), the REM alarm indicator enlarges and flashes red and yellow, a tone sounds twice, and RF output is disabled. The indicator remains illuminated red and yellow until the user correct the condition causing the alarm. Then, the indicator illuminates green and RF output is enabled.

Autobipolar

The ForceTriad Energy Platform is equipped with an autobipolar feature that allows for automatic activation of bipolar energy.

Note: The autobipolar electrode function requires the Reusable Footswitching Bipolar Cord: E0020V, E0021S, E0022W, E360150, or E360150L.

Warning

Use of different Covidien cord models or cords from other manufacturers may not achieve proper electrical output for this device, thereby failing to produce the desired clinical effect. For example, Autobipolar activation/deactivation settings may not work properly using cords other than those specified by Covidien.

The autobipolar specifications are:

Interrogation frequency	80 kHz \pm 10 kHz
Interrogation current	< 100 μ A
Interrogation voltage	< 12 V RMS
Minimum power	5 W
Activation impedance	0 Ω to 1000 Ω
Deactivation impedance	User selectable: 1500 Ω , 1800 Ω , 2000 Ω or 2200 Ω
Keying delay	User selectable in 500 ms increments from 0 sec. to 2.5 sec.

Measurement Accuracy

Inactive		
± 5% of Full Scale activation impedance while keying inactive		
Active		
Mode – BP Low		
Load/Power	< 30 W	≥ 30 W
1—500 Ω	± 20% or ±25 Ω (Whichever is greater)	± 20% or ±25 Ω (Whichever is greater)
501—1000 Ω	± 40%	± 20%
1001—2500 Ω	+100%/-50%	± 20%
> 2500 Ω	Reads > 2200 Ω	Reads > 2200 Ω
Mode – BP Standard		
Load/Power	< 50 W	≥ 50 W
1—500 Ω	± 20% or ±25 Ω (Whichever is greater)	± 20% or ±25 Ω (Whichever is greater)
501—1000 Ω	± 40%	± 20%
1001—2500 Ω	+100%/-50%	± 20%
> 2500 Ω	Reads > 2200 Ω	Reads > 2200 Ω
Mode – BP Macro		
Load/Power	All power levels	
1—2500 Ω	± 20% or ±25 Ω (Whichever is greater)	
> 2500 Ω	Reads > 2200 Ω	

Duty Cycle

Under maximum power settings and rated load conditions, the ForceTriad Energy Platform is capable of operating a duty cycle of 25%, defined as 10 seconds active and 30 seconds inactive, in any mode for a period of 4 hours.

Caution

Use of duty cycles greater than 25% (10 seconds active followed by 30 seconds inactive) will increase the risk that heat build-up under a return electrode may be high enough to injure the patient. Do not continuously activate for longer than one minute.

Low Frequency (50/60 Hz) Leakage Current

Enclosure source current, ground open	< 300 μ A
Source current, patient leads, all outputs	Normal polarity, intact ground: < 10 μ A Normal polarity, ground open: < 50 μ A Reverse polarity, ground open: < 50 μ A Mains voltage on applied part: < 50 μ A
Sink current at high line, all inputs	< 50 μ A

High Frequency (RF) Leakage Current

	Measured with leads recommended by Covidien	Measured directly at the system terminals
Bipolar RF leakage current	< 69.2 mA RMS	< 69.2 mA RMS
Monopolar RF leakage current	< 150 mA RMS	< 100 mA RMS
LigaSure leakage	<132 mA RMS	< 100 mA RMS

Input Power

100–120 Volt	220–240 Volt
Maximum VA – 1056 VA	Maximum VA – 2080 VA
Input mains voltage, full regulation range 90–132 VAC	Input mains voltage, full regulation range 208–264 VAC
Input mains voltage, operating range 85–132 VAC	Input mains voltage, operating range 170–264 VAC
Mains line frequency range (nominal) 50 Hz to 60 Hz	Mains line frequency range (nominal) 50 Hz to 60 Hz
Fuses (2) – 5 mm x 20 mm 8 A, 250 V fast blow	Fuses (2) – 5 mm x 20 mm 8 A, 250 V fast blow

Power Cord Specification

This system is factory equipped with a 110 VAC hospital grade NEMA 5-15 power cord. Should the AC power cord need to be replaced to match another plug configuration, the replacement plug/cable/receptacle configuration must meet or exceed the following specifications:

- 100-120 VAC**
 - Cable - SJT16/3, IEC color code, maximum length 15 ft. (5 m)
 - Plug - minimum 10 A - 125 VAC
 - Unit receptacle - IEC female, minimum 10 A - 125 VAC
- 220-240 VAC**
 - Cable - H05VV3G1.0 VDE, maximum length 15 ft. (5 m)
 - Plug - minimum 6 A - 250 VAC
 - Unit receptacle - IEC female, minimum 6 A - 250 VAC

Important

Contact your local Covidien representative for alternative internationally approved power-cord options.

Input Frequency

The ForceTriad Energy Platform operates within specification at all line-input frequencies between 48 Hz and 62 Hz. The user does not need to reconfigure the ForceTriad Energy Platform for different line frequencies.

Input Current

The ForceTriad Energy Platform draws no more than 10 A at input voltages between 100 V and 240 V.

Backup Power

The ForceTriad Energy Platform retains all user programmed features, calibration, and statistical data when switched off and unplugged. The ForceTriad Energy Platform operates within specification when switched over to a supplied-line power by hospital backup systems.

Equipotential Ground Connection

An equipotential ground connection is provided to allow connection of the ForceTriad Energy Platform to ground.

ECG Blanking

An ECG blanking port is provided to signal other devices that the ForceTriad Energy Platform is active. The receptacle is a 2.5 mm mono jack. It is electrically isolated from the internal ground referenced electronics with the shell electrically connected to the chassis for ESD protection.

Standards and IEC Classifications

The ForceTriad Energy Platform meets all pertinent clauses of the IEC 60601-1 second edition and IEC 60601-2-2 third edition.



ATTENTION

Consult accompanying documents



The energy platform output is floating (isolated) with respect to ground.



DANGER

Explosion risk if used with flammable anesthetics



To reduce the risk of electric shock, do not remove the cover. Refer servicing to qualified service personnel.



Unit produces non-ionizing radiation

ETL CLASSIFIED



Intertek

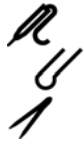
3164410
Conforms to: UL STD 60601-1
Certified: CSA STD C22.2 NO. 601.1

Classified with respect to electrical shock, fire, and mechanical hazards only in accordance with UL standard 60601-1; certified to CSA standard C22.2 No. 601.1.

Symbols



Monopolar instrument receptacle



Monopolar footswitching receptacle



Bipolar instrument receptacle



LigaSure related receptacle or footswitch



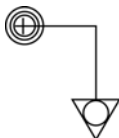
Footswitch



REM patient return electrode receptacle



Volume adjustment for activation tones



Equipotential grounding point



Equipment should not be disposed in trash

Class I Equipment (IEC 60601-1)

Accessible conductive parts cannot become live in the event of a basic insulation failure because of the way in which they are connected to the protective earth conductor.

Type CF Equipment (IEC 60601-1)/Defibrillator Proof



This ForceTriad Energy Platform provides a high degree of protection against electric shock, particularly regarding allowable leakage currents. It is type CF isolated (floating) output and may be used for procedures involving the heart.

This energy platform complies with the IEC 60601-1:1988 +A1:1991 + A2:1995 and IEC 60601-1:2005 specifications for “defibrillator proof” designation and IEC 60601-2-2:2006 and IEC 60601-2-2:2009.

Liquid Spillage (IEC 60601-2-2:2006 Clause 44.3 and IEC 60601-2-2:2009 Clause 201.11.6.3)

The ForceTriad Energy Platform is constructed so that liquid spillage in normal use does not wet electrical insulation or other components which when wetted are likely to adversely affect the safety of the equipment.

Voltage Transients (Emergency Energy Platform Mains Transfer)

The ForceTriad Energy Platform continues to operate normally with no errors or system failures when transfer is made between line AC and an emergency system-voltage source. (IEC 60601-1:1988 + A1:1991 + A2:1995 clause 49, IEC 60601-1:2005 clause 11.8, IEC 60601-2-2:2006 clause 51.101, and IEC 60601-2-2:2009 clause 201.11.8)

Electromagnetic Compatibility (IEC 60601-1-2 and IEC 60601-2-2)

The ForceTriad Energy Platform complies with the appropriate IEC 60601-1-2 and 60601-2-2 specifications regarding electromagnetic compatibility.

Notice

The ForceTriad Energy Platform requires special precautions regarding EMC and needs to be installed and put into service according to the EMC information provided in the ForceTriad Energy Platform service manual.

Portable and mobile RF communications equipment can affect the ForceTriad Energy Platform. Refer to the EMC information provided in the ForceTriad Energy Platform service manual.

The system should not be used adjacent to or stacked with equipment other than specified in the *ForceTriad Energy Platform User Guide and Service Manual*. If adjacent or stacked use is necessary, the system should be observed to verify normal operation in the configuration in which it will be used.

The system intentionally applies RF energy for diagnosis or treatment during activation. Observe other electronic medical equipment in the vicinity during the system activation for any possible adverse electromagnetic effects. Ensure adequate separation of electronic medical equipment based on observed reactions.

The use of accessories, other than specified in the *ForceTriad Energy Platform User Guide and Service Manual*, may result in increased emissions or decreased immunity of the system.

The ForceTriad Energy Platform meets the following requirements:

ESD Immunity (IEC 60601-1-2 Sub-Clause 36.202 and IEC 61000-4-2)

Radiated Immunity (IEC 60601-1-2 sub-clause 36.202.2 and IEC 61000-4-3)

Electrical Fast Transient/Burst (IEC 60601-1-2 sub-clause 36.202.3.1 and IEC 61000-4-4)

Surge Immunity (IEC 60601-1-2 sub-clause 36.202.3.2 and IEC 61000-4-5)

Emissions (IEC 60601-1-2 sub-clause 36.201.1, IEC 60601-2-2 sub-clause 36 and CISPR 11 Class A)

Harmonic distortion (IEC 60601-1-2 sub-clause 36.201.3.1 and IEC 61000-3-2)

Conducted disturbances (IEC 60601-1-2 sub-clause 36.202.6 and IEC 61000-4-6)

Power frequency magnetic fields (IEC 60601-1-2 sub-clause 36.202.8.1 and IEC 61000-4-8)

Voltage dips, short interruptions and variations (IEC 60601-1-2 sub-clause 36.202.7 and IEC 61000-4-11)

Output Characteristics

Maximum Output for Bipolar, Monopolar, and LigaSure Modes

Power readouts agree with actual power into rated load to within 15% or 5 W, whichever is greater.

Caution

To avoid injury to the patient or surgical team, use only instruments rated for use at, or greater than, the maximum peak voltages listed below. For example, bipolar instruments must have voltage ratings of 250 V peak or greater, as shown in the "Open Circuit Peak Voltage (max)" column.

Mode	Open Circuit Peak Voltage (max)	Open Circuit P-P Voltage (max)	Rated Load (max)	Power (max)	Crest Factor*	Duty Cycle
Bipolar						
Low	250 V	500 V	100 Ω	95 W	1.42	N/A
Standard	175 V	350 V	100 Ω	95 W	1.42	N/A
Macro	250 V	500 V	100 Ω	95 W	1.42	N/A
Monopolar Cut						
Cut	1050 V	2100 V	300 Ω	300 W	1.42	N/A
Blend	1485 V	2970 V	300 Ω	200 W	2.7	50%
Valleylab (HWD)	2365 V	4730 V	300 Ω	200 W	4.3	25%
Monopolar Coag						
Fulgurate	3050 V	6100 V	500 Ω	120 W	5.55	6.5%
Spray	3625 V	7250 V	500 Ω	120 W	6.6	4.6%
LigaSure	287.5 V	575 V	20 Ω	350 W	1.42	N/A
LigaSure Test	147.5 V	295 V	20 Ω	190 W	1.42	N/A

* An indication of a waveform's ability to coagulate bleeders without a cutting effect.

Available Power Settings in Watts

Autobipolar (All Modes)

5 W to 40 W available in 1 W increments

				5	6	7	8	9	10
11	12	13	14	15	16	17	18	19	20
21	22	23	24	25	26	27	28	29	30
31	32	33	34	35	36	37	38	39	40

45 W to 95 W available in 5 W increments

45	50	55	60	65	70	75	80	85	90
95									

Bipolar (All Modes)

1 W to 40 W available in 1 W increments

1	2	3	4	5	6	7	8	9	10
11	12	13	14	15	16	17	18	19	20
21	22	23	24	25	26	27	28	29	30
31	32	33	34	35	36	37	38	39	40

45 W to 95 W available in 5 W increments

45	50	55	60	65	70	75	80	85	90
95									

Monopolar Cut

1 W to 40 W available in 1 W increments

1	2	3	4	5	6	7	8	9	10
11	12	13	14	15	16	17	18	19	20
21	22	23	24	25	26	27	28	29	30
31	32	33	34	35	36	37	38	39	40

45 W to 95 W available in 5 W increments

45	50	55	60	65	70	75	80	85	90
95									

100 W to 300 W available in 10 W increments

	100	110	120	130	140	150	160	170	180
190	200	210	220	230	240	250	260	270	280
290	300								

Monopolar Blend

1 W to 40 W available in 1 W increments

1	2	3	4	5	6	7	8	9	10
11	12	13	14	15	16	17	18	19	20
21	22	23	24	25	26	27	28	29	30
31	32	33	34	35	36	37	38	39	40

45 W to 95 W available in 5 W increments

45	50	55	60	65	70	75	80	85	90
95									

100 W to 200 W available in 10 W increments

	100	110	120	130	140	150	160	170	180
190	200								

Valleylab

1 W to 40 W available in 1 W increments

1	2	3	4	5	6	7	8	9	10
11	12	13	14	15	16	17	18	19	20
21	22	23	24	25	26	27	28	29	30
31	32	33	34	35	36	37	38	39	40

45 W to 95 W available in 5 W increments

45	50	55	60	65	70	75	80	85	90
95									

100 W to 200 W available in 10 W increments

	100	110	120	130	140	150	160	170	180
190	200								

Monopolar Coag

1 W to 40 W available in 1 W increments

1	2	3	4	5	6	7	8	9	10
11	12	13	14	15	16	17	18	19	20
21	22	23	24	25	26	27	28	29	30
31	32	33	34	35	36	37	38	39	40

45 W to 90 W available in 5 W increments

45	50	55	60	65	70	75	80	85	90
95									

100 W to 120 W available in 10 W increments

100	110	120
-----	-----	-----

Output Waveforms

Tissue Sensing Technology, an automatic adjustment, controls all modes. As tissue resistance increases from zero, the energy platform outputs constant current followed by constant power followed by constant voltage. The maximum output voltage is controlled to reduce capacitive coupling and video interference and to minimize sparking.

Bipolar

Low	472 kHz sinusoid continuous
Standard	472 kHz sinusoid continuous
Macro	472 kHz sinusoid continuous

Monopolar Cut

Cut	472 kHz sinusoid continuous
Blend	472 kHz bursts of sinusoid, recurring at 26.21 kHz intervals. 50% duty cycle.

Valleylab

Valleylab	472 kHz bursts of sinusoid, recurring at 28.3 kHz intervals. 25% duty cycle.
------------------	--

Monopolar Coag


Fulgurate	472 kHz damped sinusoidal bursts with a repetition frequency of 30.66 kHz. 6.5% duty cycle.
Spray	472 kHz damped sinusoidal bursts with a randomized repetition centered at 21.7 kHz. 4.6% duty cycle.

Output Power vs. Resistance Graphs

Monopolar Graphs


Pure Cut

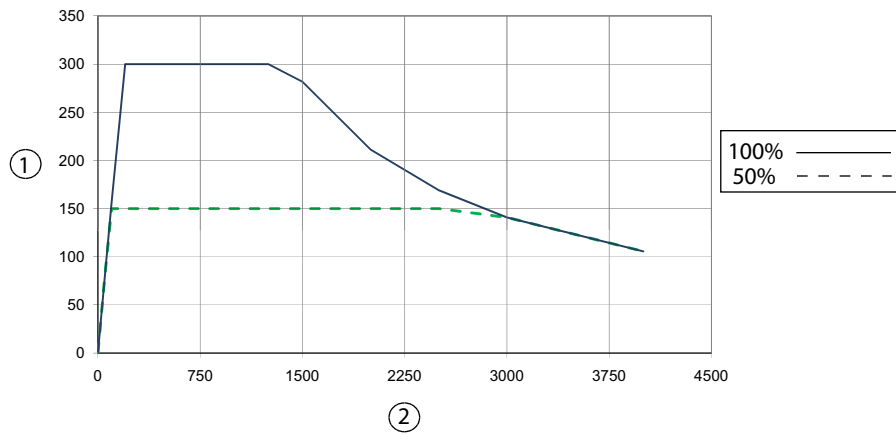Output power versus impedance for Pure cut power




- ① Output power (watts)
- ② Load impedance (ohms)

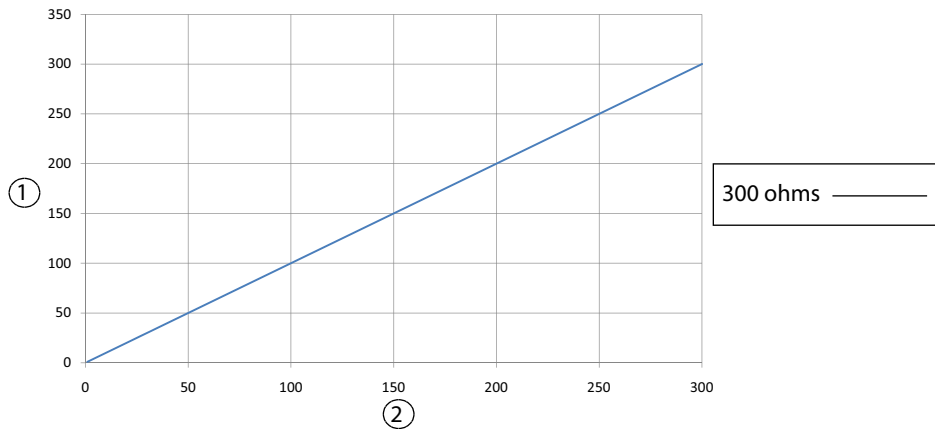Output power versus power setting for Pure cut power



- ① Output power (watts)
- ② Power setting

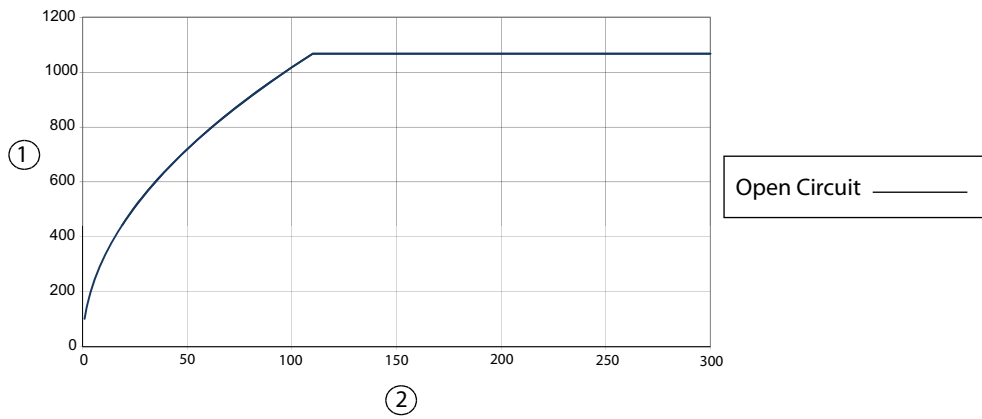
Peak voltage versus power setting for Pure cut power



- ① Peak voltage
- ② Power setting


Blend

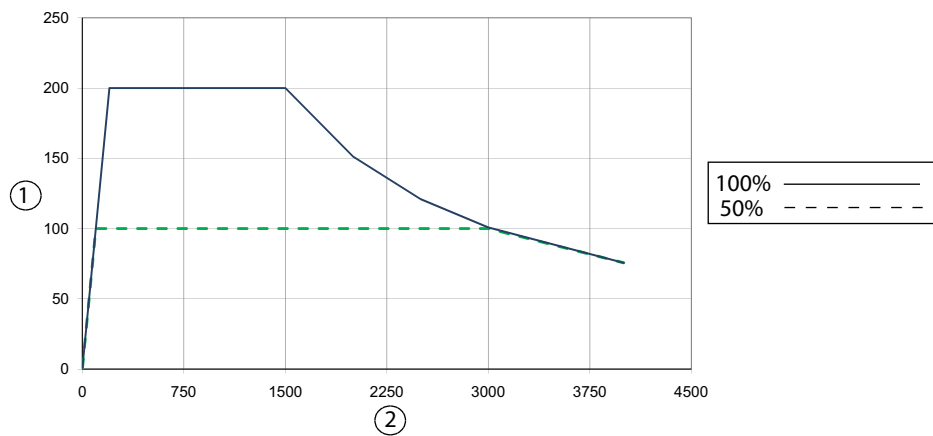Output power versus impedance for Blend power



- ① Output power (watts)
- ② Load impedance (ohms)

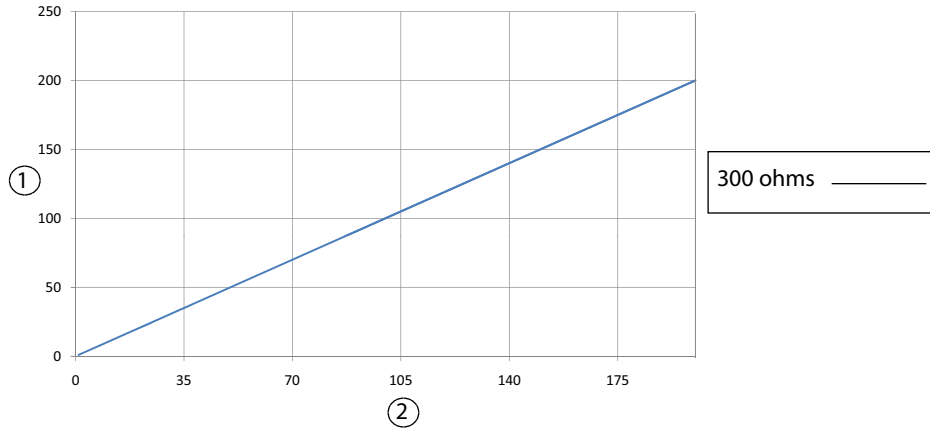
Output power versus power setting for Blend power



① Output power (watts)

② Power setting

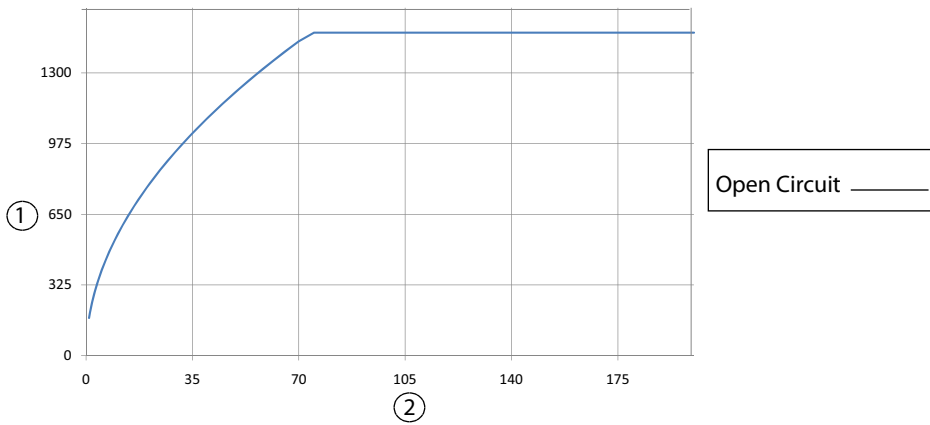Peak voltage versus power setting for Blend power




① Peak voltage

② Power setting


Fulgurate

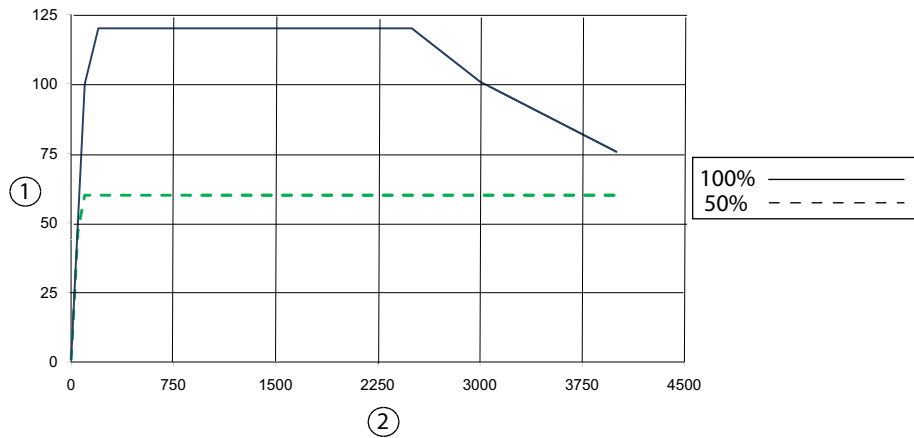Output power versus impedance for Fulgurate power




- ① Output power (watts)
- ② Load impedance (ohms)

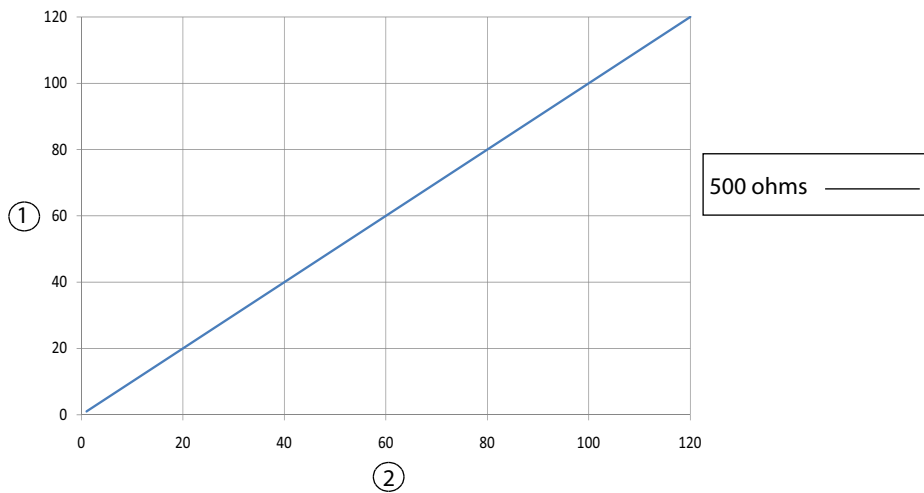Output power versus power setting for Fulgurate power



- ① Output power (watts)
- ② Power setting

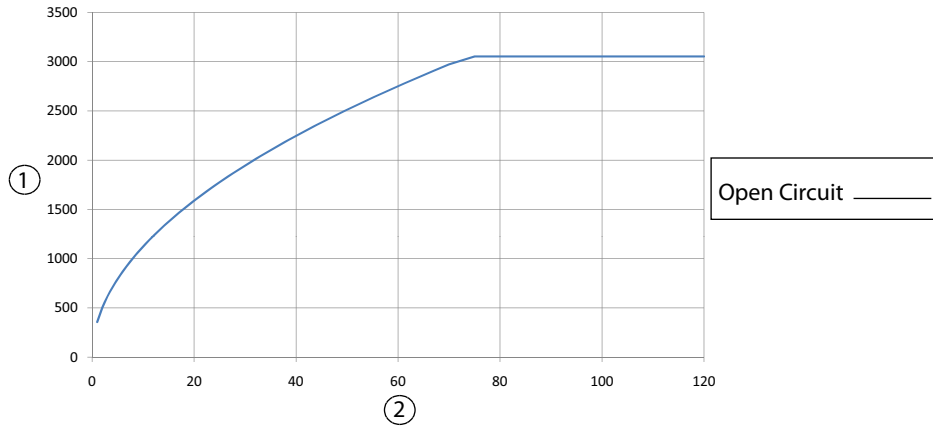
Peak voltage versus power setting for Fulgurate power



- ① Peak voltage
- ② Power setting


Spray

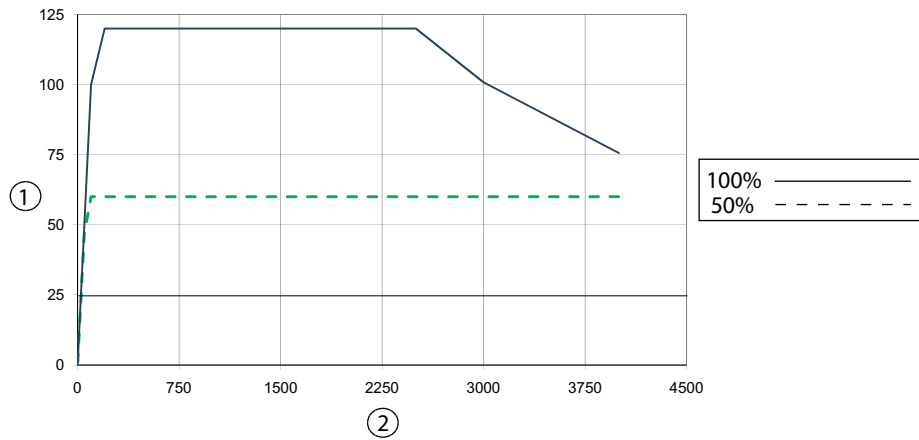Output power versus impedance for Spray power




- ① Output power (watts)
- ② Load impedance (ohms)

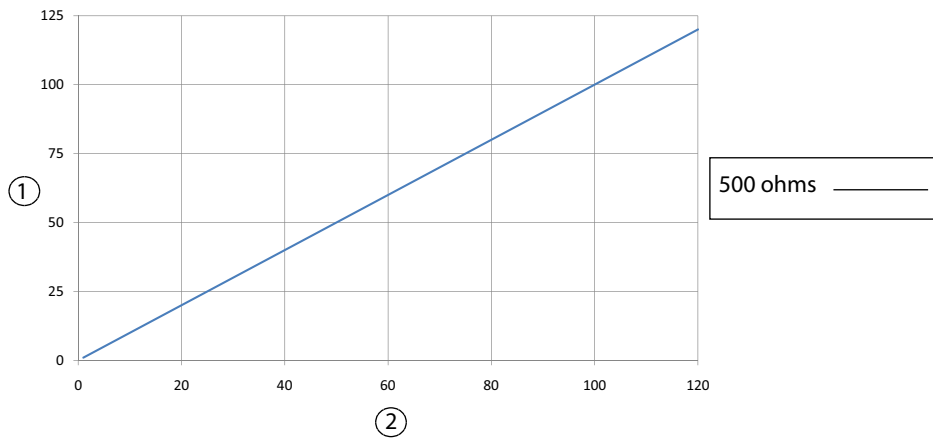Output power versus power setting for Spray power



- ① Output power (watts)
- ② Power setting

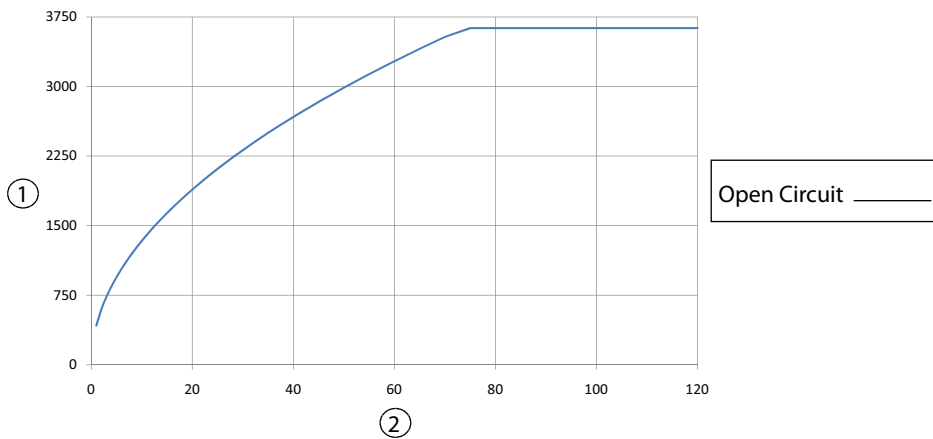
Peak voltage versus power setting for Spray power



- ① Peak voltage
- ② Power setting


Valleylab

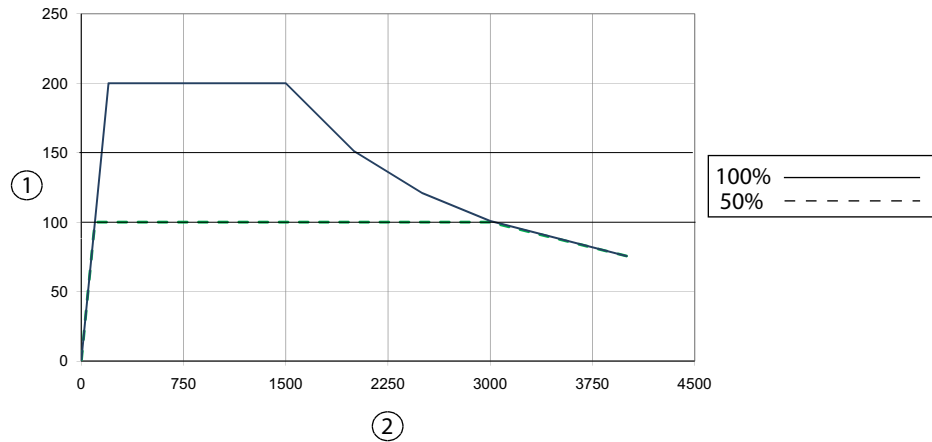Output power versus impedance for Valleylab power




- ① Output power (watts)
- ② Load impedance (ohms)

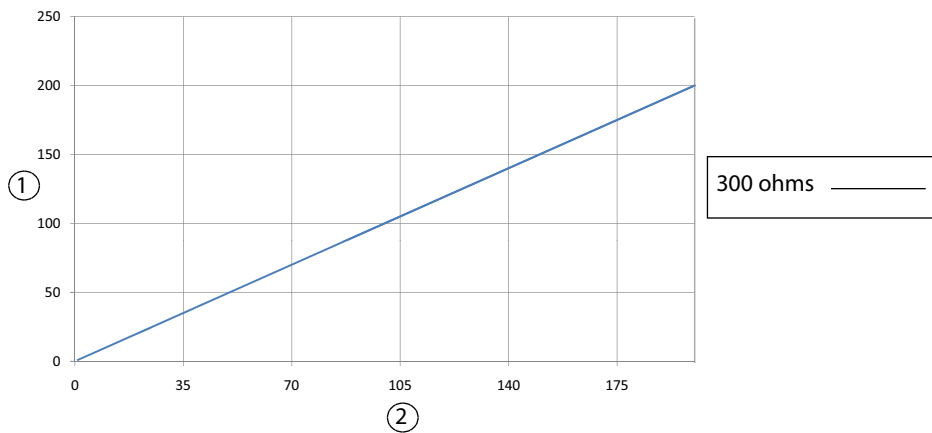Output power versus power setting for Valleylab power



- ① Output power (watts)
- ② Power setting

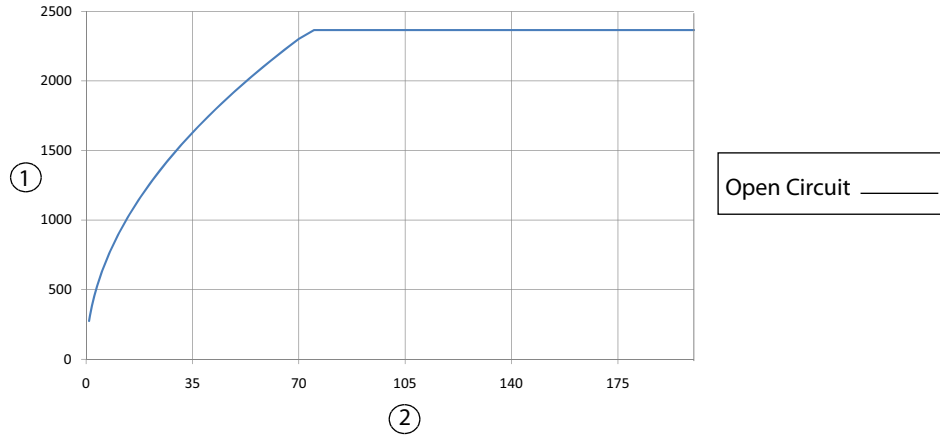Peak voltage versus power setting for Valleylab power




- ① Peak voltage
- ② Power setting

Bipolar Graphs


Bipolar Low

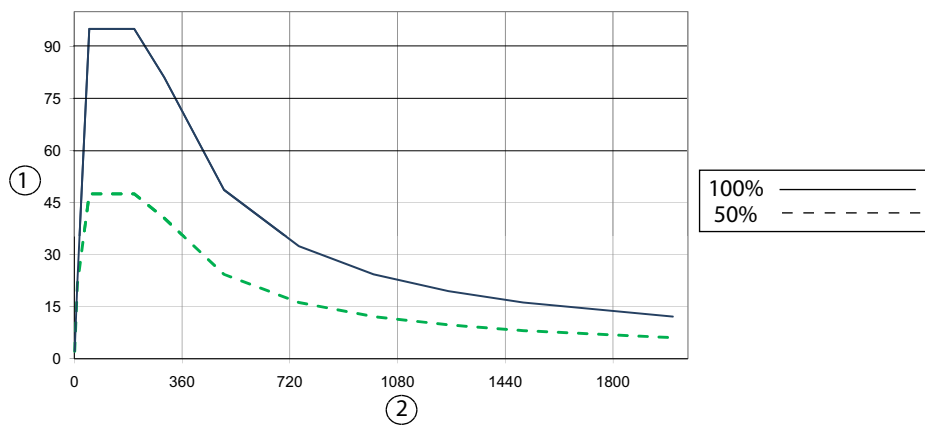Output power versus impedance for Bipolar Low power



- ① Output power (watts)
- ② Load impedance (ohms)

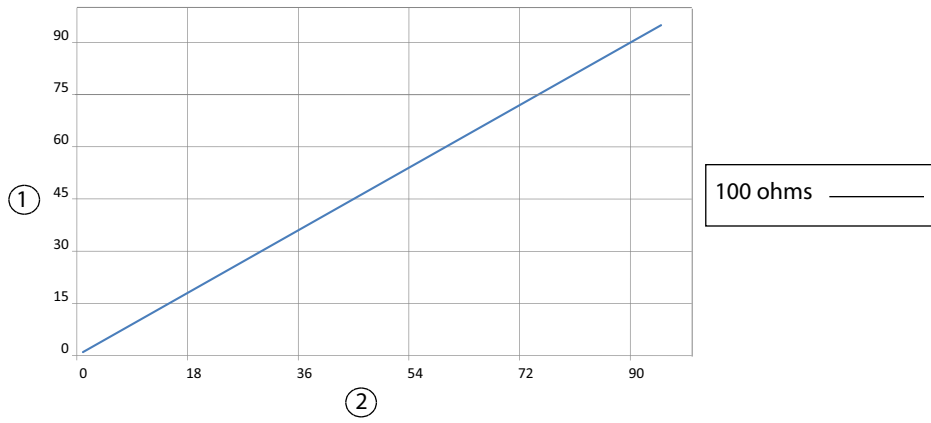
Output power versus power setting for Bipolar Low power



① Output power (watts)

② Power setting

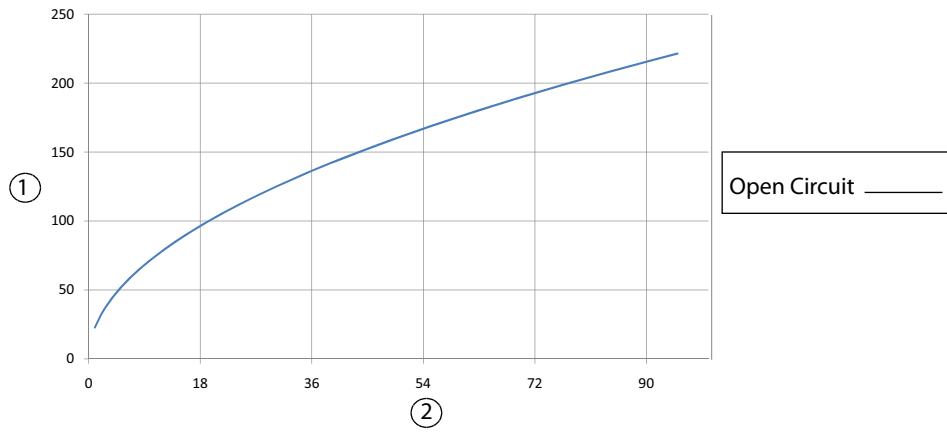Peak voltage versus power setting for Bipolar Low power




① Peak voltage

② Power setting


Bipolar Standard

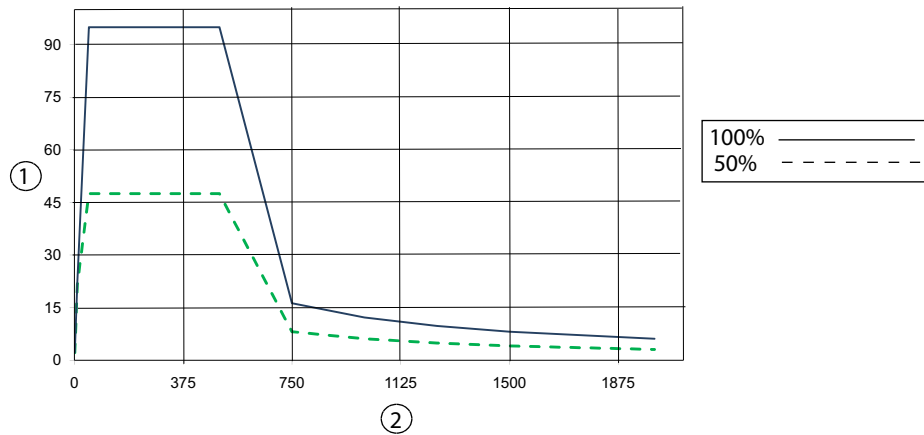Output power versus impedance for Bipolar Standard power




- ① Output power (watts)
- ② Load impedance (ohms)

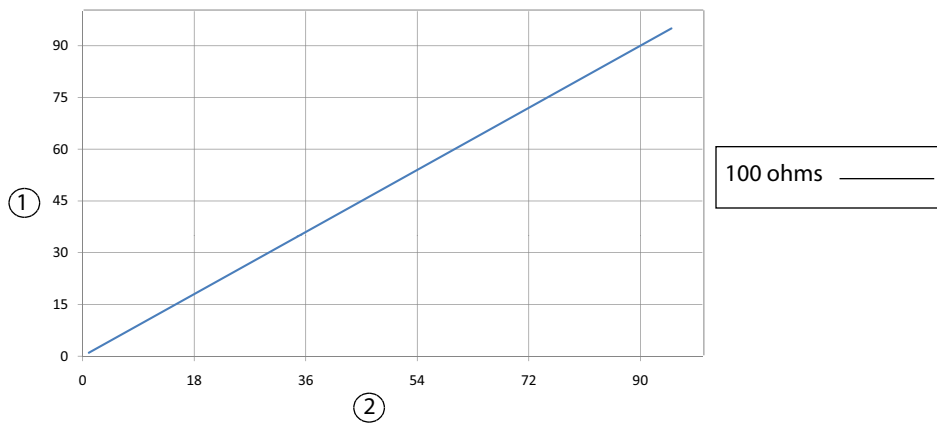Output power versus power setting for Bipolar Standard power



- ① Output power (watts)
- ② Power setting

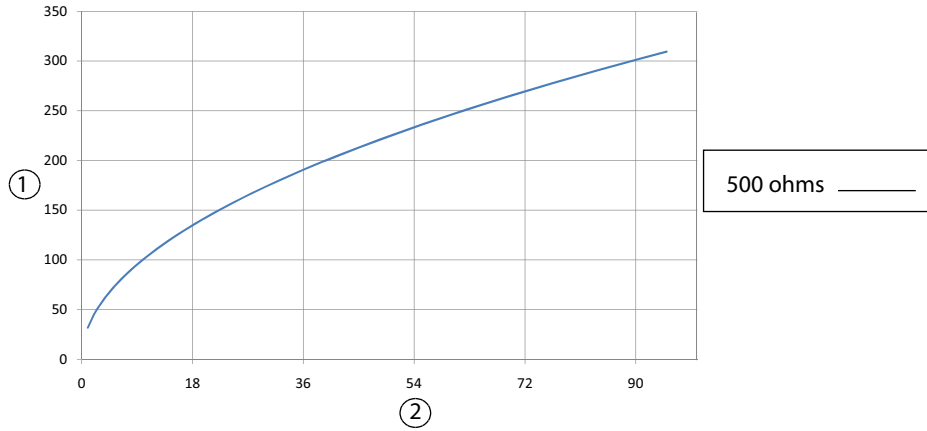Peak voltage versus power setting for Bipolar Standard power




- ① Peak voltage
- ② Power setting

Note: Maximum peak voltage in the Bipolar Standard mode occurs at 500 Ω , not open circuit.


Bipolar Macro

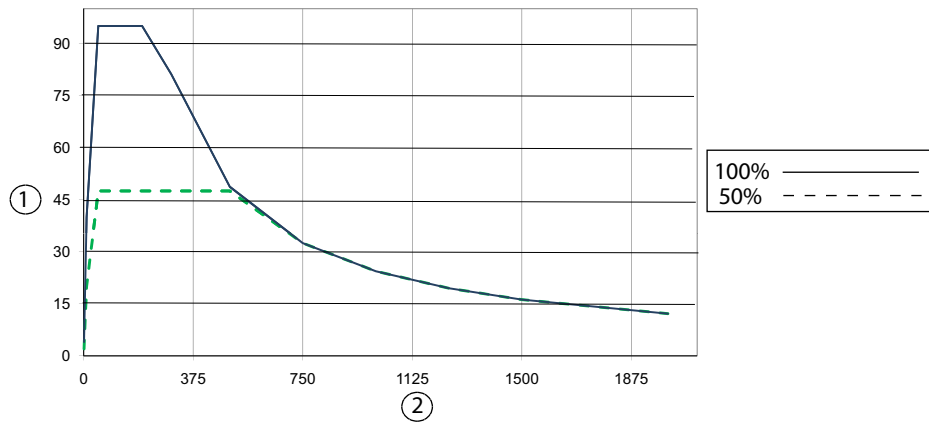Output power versus impedance for Bipolar Macro power



- ① Output power (watts)
- ② Load impedance (ohms)

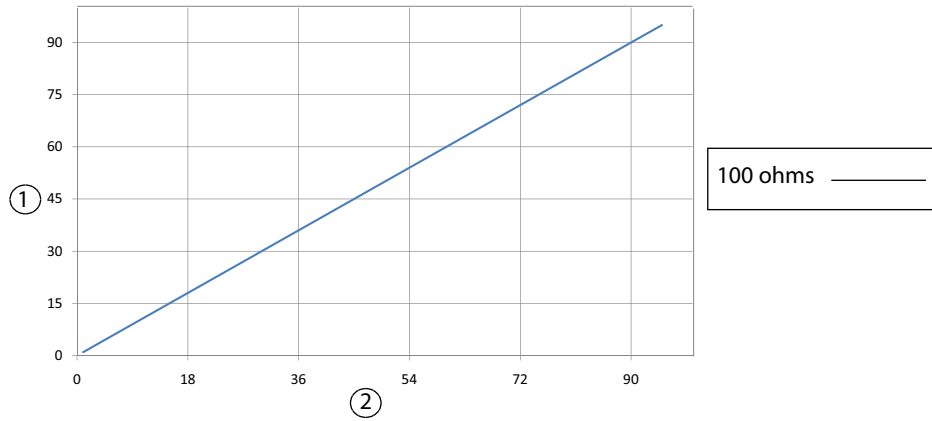
Output power versus power setting for Bipolar Macro power



① Output power (watts)

② Power setting

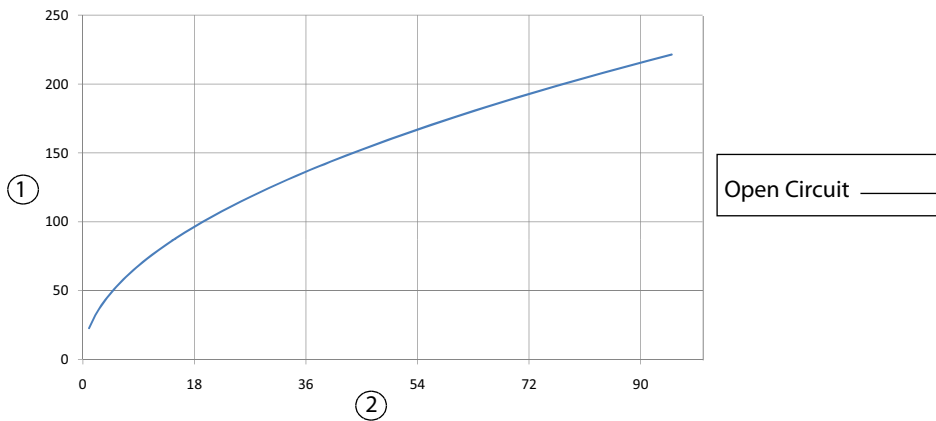Peak voltage versus power setting for Bipolar Macro power




① Peak voltage

② Power setting


LigaSure

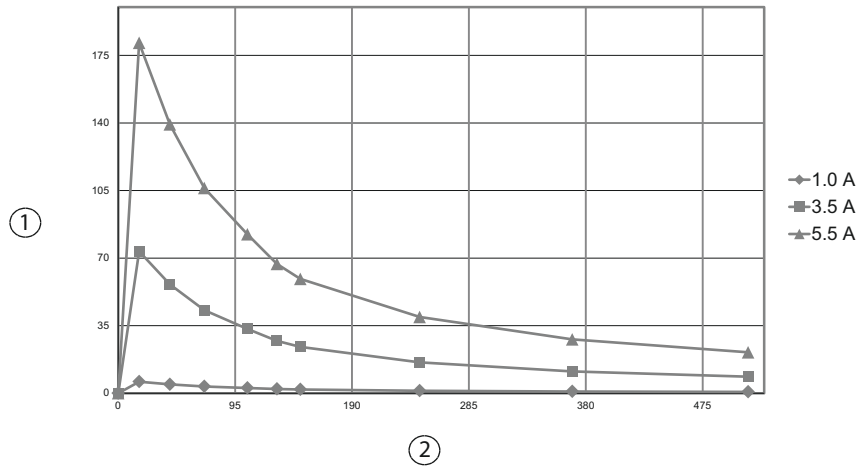Output power versus impedance for LigaSure power



- ① Output power (watts)
- ② Load impedance (ohms)

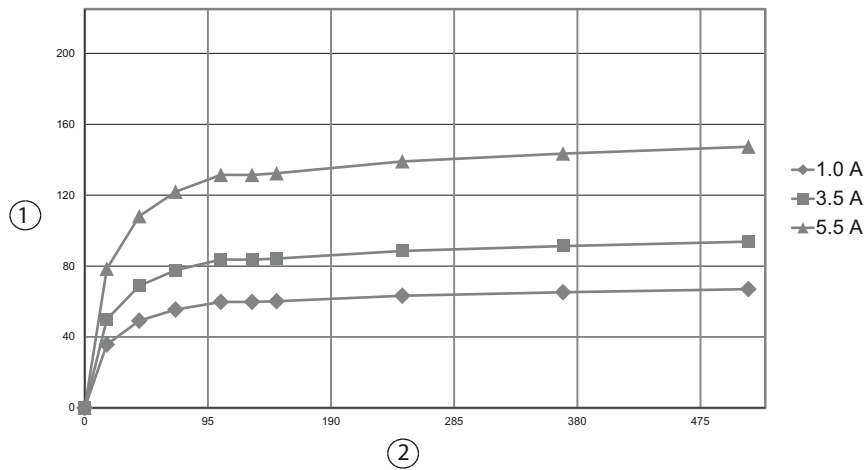Peak voltage versus impedance for LigaSure power



- ① Peak voltage
- ② Load impedance (ohms)



Rx
ONLY



COVIDIEN

Part No. 1063547

COVIDIEN, COVIDIEN with logo and Covidien logo are U.S. and internationally registered trademarks of Covidien AG.

Other brands are trademarks of a Covidien company.

Prolystica is a registered trademark of STERIS Corporation.


Enzol is a registered trademark of Johnson and Johnson Medical Inc.

U.S. Patents 5,370,645; 5,372,596; 5,772,659; 6,068,627; 6,203,541; 7,137,980; 7,222,601; 7,416,437; 7,513,896; 7,648,499; 7,651,492; 7,666,693; 7,901,400; D567,761; D574,323.

©2013 Covidien.

Made in USA. Printed in USA.

 Covidien llc, 15 Hampshire Street, Mansfield, MA 02048 USA.

 Covidien Ireland Limited, IDA Business & Technology Park, Tullamore.

www.covidien.com

1-303-530-2300 [T]

1-800-255-8522 [T]

REV 06/2013