


Osaki

TheraMedic

INSTRUCTION MANUAL

MASSAGE BED

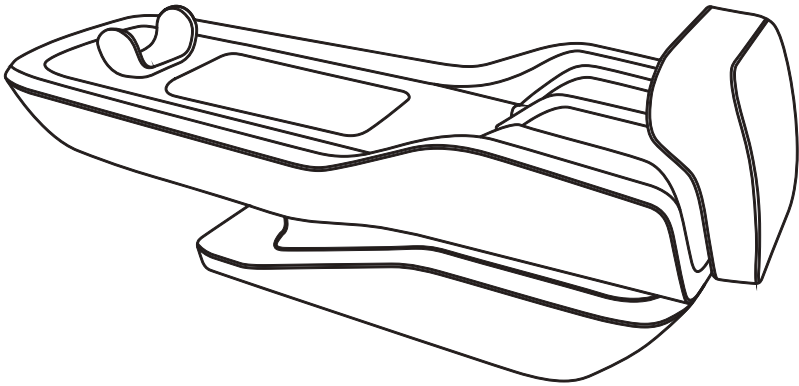

PLEASE BE SURE TO READ THIS MANUAL CAREFULLY TO LEARN THE OPERATIONS OF THE THIS PRODUCT. DO NOT DISPOSE OF THE MANUAL AND KEEP FOR FUTURE REFERENCE.

IMPORTANT SAFETY INSTRUCTIONS

When using an electrical appliance, basic precautions should always be followed, including the following: Read all instructions before using (this appliance).

DANGER--To reduce the risk of electric shock:

- Always unplug this appliance from the electrical outlet immediately after using and before cleaning.
- Never use pins or other metallic fasteners with this appliance.
- Carefully examine the covering before each use. Discard the appliance if the covering shows any sign of deterioration, such as checking, blistering, or cracking.

- **WARNING**--To reduce the risks of burns, fire, electric shock or injury to persons:
 - An appliance should never be left unattended when plugged in. Unplug from outlet when not in use, and before putting on or taking off parts.
 - Close supervision is necessary when this appliance is used by, on, or near children, invalids, or disabled persons.
 - Use this appliance only for its intended use as described in this manual. Do not use attachments not recommended by the manufacturer.
 - Never operate this appliance if it has a damaged cord or plug, if it is not working properly, if it has been dropped or damaged, or dropped into water. Return the appliance to a service center for examination and repair.
 - Do not carry this appliance by supply cord or use cord as a handle.
 - Keep the cord away from heated surfaces.
 - Never operate the appliance with the air openings blocked. Keep the air openings free of lint, hair, and the like.
 - Never drop or insert any object into any opening.
 - Do not use outdoors.
 - Do not operate where aerosol (spray) products are being used or where oxygen is being administered.
 - To disconnect, turn all controls to the off position, then remove plug from outlet.
 - Keep children away from extended foot support.
 - Connect this appliance to a properly grounded outlet only. See Grounding Instructions.
 - Use heated surfaces carefully. May cause serious burns. Do not use over insensitive skin areas or in the presence of poor circulation. The unattended use of heat by children or incapacitated persons may be dangerous.
 - Do not operate under blanket or pillow. Excessive heating can occur and cause fire, electric shock, or injury to persons.

SAVE THESE INSTRUCTIONS


GROUNDING INSTRUCTIONS

This product must be grounded. If it should malfunction or break down, grounding provides a path of least resistance for the electric current to reduce the risk of electric shock.

This product is equipped with a cord with an equipment-grounding conductor and a grounding plug.

The plug must be plugged into an appropriate outlet that is properly installed and grounded in accordance with all local codes and ordinances. Danger--Improper connection of the equipment-grounding conductor can result in a risk of electric shock.


Check with a qualified electrician or serviceman if you are in doubt as to whether the product is properly grounded. Do not modify the plug provided with the product - if it will not fit the outlet, have a proper outlet installed by a qualified electrician.

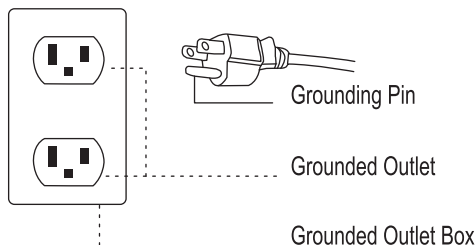

IMPORTANT SAFETY INSTRUCTIONS

OPERATING SURROUNDINGS

- OPERATING ENVIRONMENT
- Do not use or store the massage chair under damp conditions:
 - i.e. bathroom, sauna, spa, to avoid any mold development, electric shock, fire or mechanical failure.
- Don't use the massage chair outdoors
- Don't use the massage chair under dusty and temperature varying conditions.
- Don't place the massage chair near heater(s), stove or direct sunlight to avoid inner components overheating which may lead to electrical fires and shock.
- Make sure to place the massage chair on flat level ground surface
- The chair should be placed on carpet, wood or hard surface flooring
- Make sure the power plug is properly grounded to avoid shock and injury.
- Don't use the machine while smoking

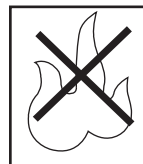

FOLLOWING CONDITIONS REQUIRING CONSULT FROM A MEDICAL DOCTORS BEFORE USE

- Suffered from malignant tumor
- Serious heart disease
- Serious skin disease
- Bone disease
- While pregnant or menstruating
- On rest by the doctor
- Spinal condition including injuries associated to the back and spine
- Hypertension and cramping
- Arthritic bone disease
- Any implanted medical device such as pacemaker
- Suffering from internal conditions (i.e. gastritis, enteritis, diarrhea, Zepatitis)
- Serious blood circulation trouble
- Consciousness issues such as vertigo and other equilibrium matters

WARNING

Children under the age of 18 years old or with physical or mental condition are not permitted to use this product unattended by an able adult and with the permission of a medical doctor. These precautions are taken to prevent injury and damage to the product

IMPORTANT SAFETY INSTRUCTIONS

For safety reasons and in order to clearly indicate potential dangers should the product be used improperly, instructions have been marked in the following manners and should be strictly noted

⊘ CAUTIONS

Non-authorized technicians are not permitted to dismantle or repair this product to avoid fire, electrical shock or serious injury



Do not clean or wipe down this product with anything wet or moist when this symbol is shown.



When storing or during maintenance, be sure to disconnect the power plug from socket. Do not detach the plug with wet hands



Make sure to turn the power switch to the "off" position and detach the power plug when cleaning or storing the product



**During cleaning or when moving the chair, make sure to detach power cord from the wall and chair to avoid fire or electrical shock
Do not place the machine on a unstable, uneven and unlevelled surface to avoid damage, insulation failure or electrical shock**



When storing the machine, make sure to disconnect power plug from the outlet and chair to avoid fire caused by electrical surges.



Clean power plug periodically to avoid buildup leading possible fire



Do not use it at any unsafe area or for any other purposes than personal use. Not for commercial use.



**Do not store the chair in areas that are not climate controlled.
Must be in room temperature with very low humidity.**



Keep the power cord unobstructed at all time, i.e. not twisted, knotted, over bent or pulled and do not place anything heavy on it



IMPORTANT SAFETY INSTRUCTIONS

⊘ CAUTIONS

Do not make any modifications or change components by an unauthorized person



Do not dismantle the machine by unauthorized person to avoid electrical shock, short-circuit or fire



Do not store it near a fire, heater or flammable substances It is one of the major reasons of causing product damage.



Should any abnormality occur(i.e. burning smell), stop using immediately and detach power plug from outlet.



please contact your local seller or authorized service agent for maintenance to avoid electrical fire and injury.



Only proper rated voltage must be used . Stop using immediately if power cord or power plug is damaged or not tightly fitted into the outlet or the base of the chair




Be sure to keep children and pets away while operating to avoid any injury Do not place the chair in direct sunlight at all times




If the supply cord is damaged,it must be replaced by the manufacturer or an authorized technician



NAME OF EACH PART

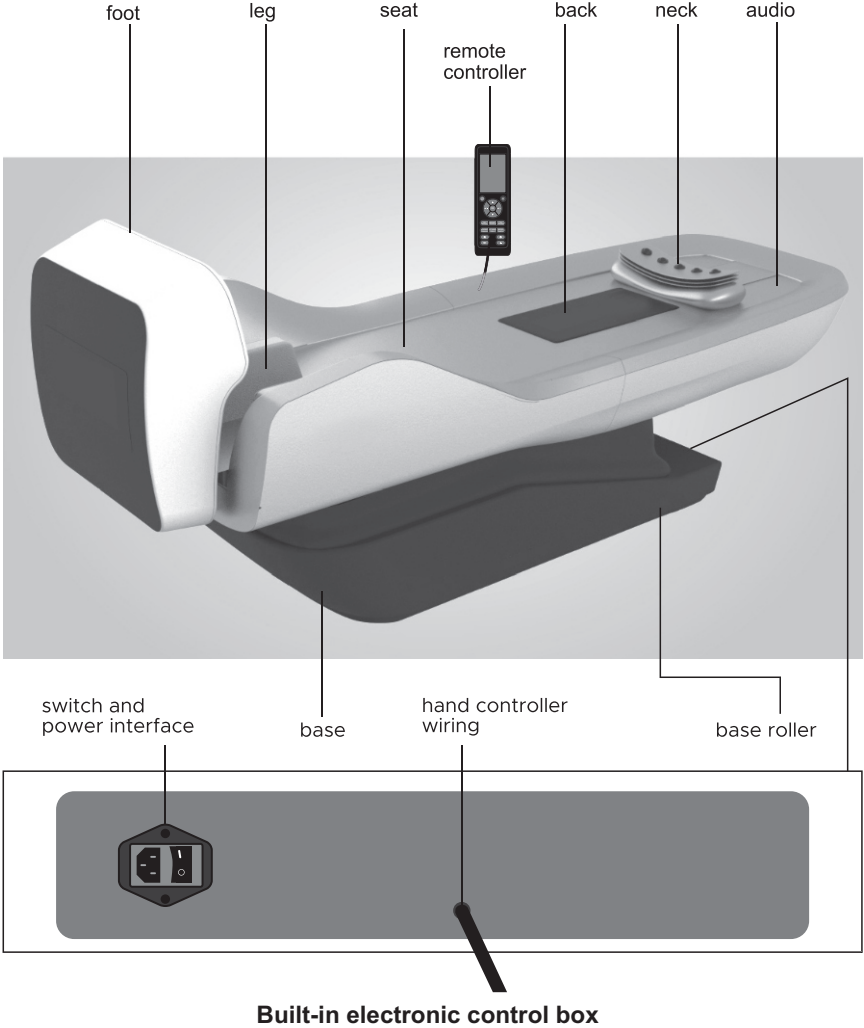

NAME OF EACH PART

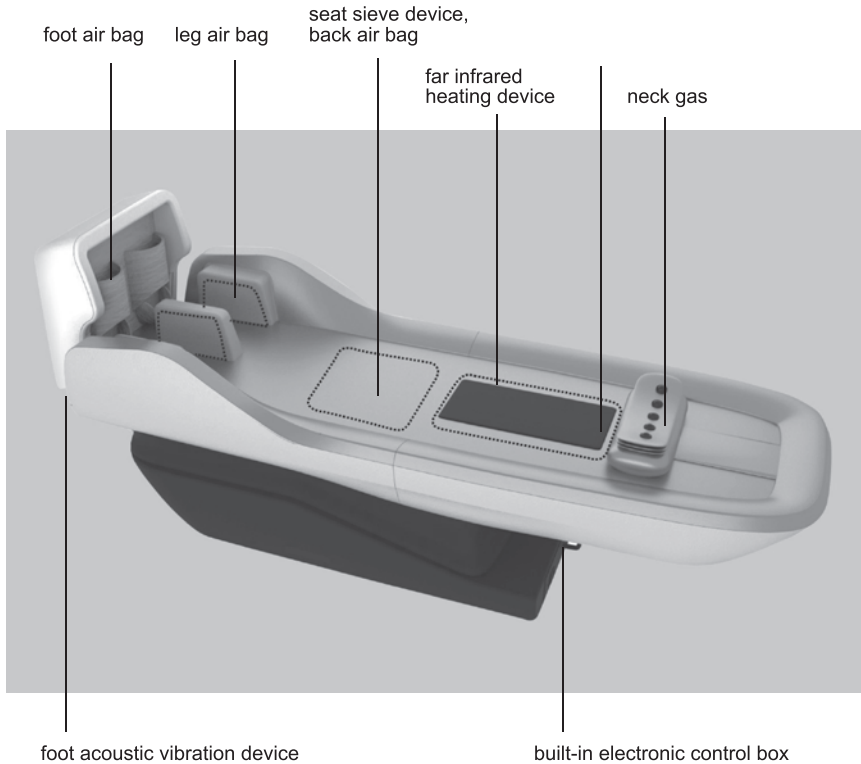

Pressure attention

Force is strictly prohibited in the ABS shell of the foot of the machine. It is not allowed to carry, break, impact, impact, etc.

It is forbidden to sit around the neck pillow.

When working on this machine, it is strictly forbidden for objects, animals or people to sit close to or under the equipment.


It is strictly forbidden to disassemble the machine, especially the electronic control part.

For leather, cloth and ABS, use of hard objects is strictly prohibited.







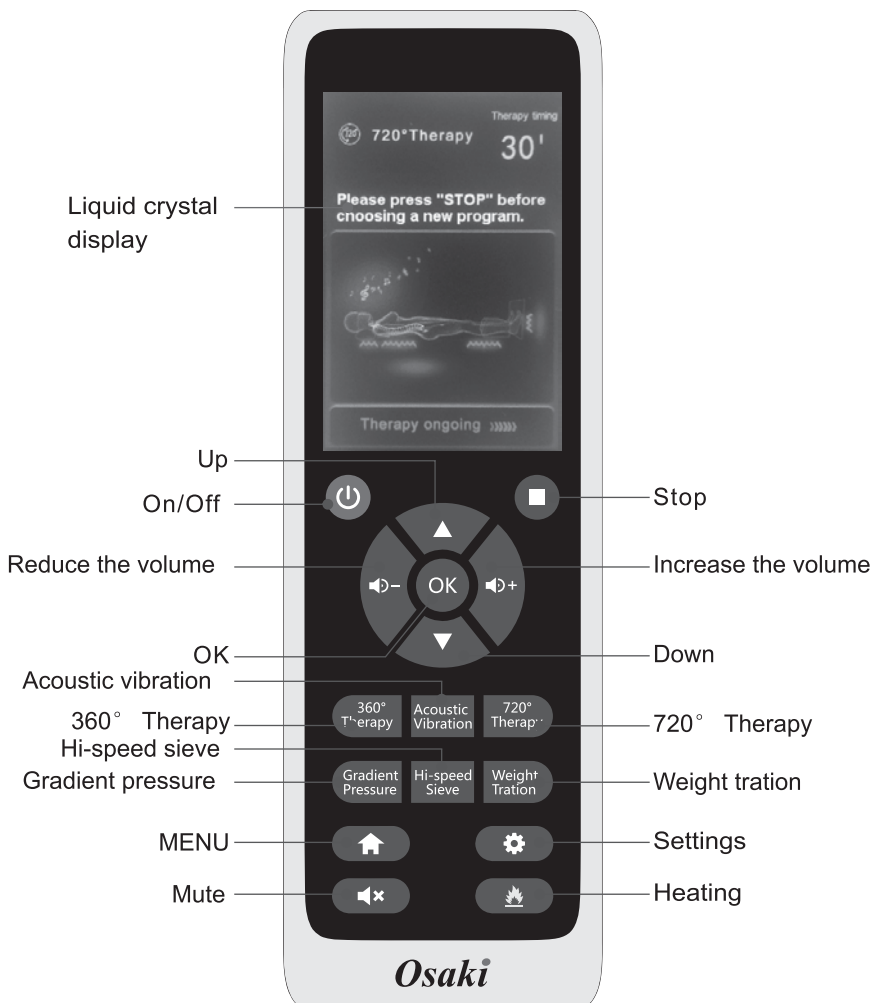








PRODUCT INTRODUCTION

720° Therapy conditioning (approximately 30 min)					
Start running, automatically start the five elements music and far infrared heating. can manually adjust the volume and turn off far infrared heating.					
First step	Step 2	Step 3	Step 4	Step 5	Step 6
Gradient pressure (about 6 min)	Acoustic vibration (about 3 min)	Weight traction (about 5 min)	Acoustic vibration (about 3 min)	Hi-speed sieve (about 4 min)	Gradient pressure (about 9 min)
360° Therapy conditioning (approximately 20 min)					
Start running, automatically start five elements music and far infrared heating, can manually adjust the volume and turn off far infrared heating.					
First step	Step 2	Step 3	Step 4	Step 5	Step 6
Gradient pressure (Approximately 5 min)	Acoustic vibration (about 2 min)	Weight traction (Approximately 3 min)	Acoustic vibration (about 2 min)	Hi-speed sieve (About 2 min)	Gradient pressure (approximately 5 min)
Gradient pressure conditioning (approximately 25 min)					
Start running, automatically start five elements music and far infrared heating, can manually adjust the volume and turn off far infrared heating.					
From neck, spine, leg to foot circulation sequence of pneumatic physiotherapy, each position four times, circulation total time of about 25 min; min; Suggest far infrared open.					
Acoustic vibration conditioning (about 10min)					
Start running, automatically start five elements music and far infrared heating, can manually adjust the volume and turn off far infrared heating.					
Horizontal position, lateral vibration, approximately 10 min.					
Gravity traction conditioning (about 10min)					
Start running, automatically start five elements music and far infrared heating, can manually adjust the volume and turn off far infrared heating.					
The whole back of the machine is inclined to 10° (foot rising, head falling), 5 air bags in the spine circulation in turn pneumatic physiotherapy; conditioning time 10min; recommended to open far infrared hyperthermia.					
Hi-speed sieve conditioning (approximately 10 min)					
Start running, automatically start five elements music and far infrared heating, can manually adjust the volume and turn off far infrared heating.					
Horizontal position, longitudinal vibration. About 10 min of conditioning.					

INSTRUCTIONS ON CONTROL PANEL



OPERATION GUIDE

	On/Off button: press this button to start massage, press again to shut down.
	Stop button: press this button to stop massage.
	720° Therapy conditioning button: press this button to enter 720° therapy conditioning mode. (For details of the model, please refer to P9. Instruction)
	360° Therapy conditioning button: press this button to enter 360° therapy conditioning mode.(For details of the model refer to the P10. of this specification)
	Gradient Pressure conditioning button: press this button to enter Gradient pressure conditioning mode. (For details of the model, please refer to P10. instruction)
	Acoustic vibration conditioning button: press this button to enter acoustic vibration conditioning mode. (For details of the model refer to the P10. of this specification)
	Weight traction conditioning button: press this button to enter self weight traction conditioning mode. (For details of the model refer to the p10. of this specification)
	Hi-speed sieve conditioning button: press this button to enter Hi-speed sieve conditioning mode. (for details of the model refer to the p10. of this specification)
	MENU button: press this button to return to the home page.
	Direction button: press these two buttons to flip the menu TAB.
	Settings: press the button to setting menu.
	Heating button: when machine starts, the heating automatically opens by default, press this button one time to turn off, and press again to turn on.
	Increase the volume button: press this button to increase the volume of music playing.
	Reduce volume button: press this button to reduce the volume of music playing.
	Mute button: press this button to mute.

TROUBLE SHOOTING

Fault phenomenon	Confirmation	Processing methods
The machine doesn't work when the power is on	The power switch of this machine is not turned on	Turn on the switch
	No native function selected	Select native function
	Poor contact with plug and socket	Make the plug in good contact with the socket
	Power cord or plug damaged	Please entrust the manufacturer, designated repairman or qualified
	Fuse damage	Replacement of fuses of the same maintenance personnel to check and replace specifications
	Failure of internal control circuit	Entrusting the designated maintenance department to check and repair
This machine is working "dā cǎo" or the buzzing	Sound produced when air pumps, motors and mechanical components work	Normal phenomena without any treatment
The sound suddenly got bigger at work	Is it too long	Cut off the power and let the machine rest for more than half an hour before restarting it
	The wear of mechanism caused by longtime use is a normal phenomenon	Please entrust the designated maintenance department to check and repair
The machine suddenly stopped	Is there a sudden power outage	Re-energize power
	May press the stop button	Re-select conditioning mode
Abnormal heating of power cord or plug	Possible overload	Cut off the power and let the machine rest for more than half an hour before restarting it
	There may be other anomalies	Stop use immediately and entrust the designated maintenance department to inspect or repair

SPECIFICATION

Product name	Massage bed		
Product model	TheraMedic		
Rated voltage	120V~ 60Hz		
Rated power	80W	Net/Gross weight	116 kg/145kg
Color	Pearl Grey & Brown	Product size	213x74x83cm
Safety design	Category 1		