User Manual

3-IN-1 BODY SLIMMING DEVICE

3-IN-1 BODY SLIMMING DEVICE

Index

- Specification.....P1
- Safety Precautions..... P2
- Important Notice.....P3
- Components.....P4
- Quick User Guide.....P5
- Instructions for using the Ultrasonic Wave.....P6
- Body Treatment.....P7
- Instructions for using the EMS.....P9
- Instructions for using the Infrared.....P12
- After useP12

- ◆ Thank you very much for choosing this product.
- ◆ Please read through carefully and understand these instructions before use.
- ◆ Please keep it in a safe place for future reference.
- ◆ This device is intended for home use and cosmetic purpose only.
- ◆ Using this device for treating any medical condition or illness is NOT supported.

Specification

Adaptor input power: AC100-240V
 Adaptor output power: 15V, 800mA
 Dimensions: 190*80*60mm
 Weight: 140g
 Ultrasonic frequency : 1MHz
 Appendix: two groups of contact pad (4pcs)

Functions

Infrared: The red light with a wavelength of 625nm can penetrate into the deep skin and generate heat to accelerate the circulation of blood.

Ultrasonic : use high frequency vibration to reduce fat and tighten the skin.

EMS: stimulate muscle with micro-current massage, against stubborn grease and relief the pain of muscle.

①


Safety precautions

Warning: Please read the following information carefully before using the machine.

1. Keep the device away from wet areas and do not immerse it in water to avoid damage to the electronic components or an electrical fault.
2. Do not store the device near heat sources, as this could shorten the life of the device causing deformation of its structure or melt the plastic parts.
3. Do not drop the product or subject it to vibration as it can cause malfunction
4. Never attempt to disassemble the product yourself. All the product components can only be maintained and repaired by professional technicians.
5. Do not use electrical adapters that are not included in the package. The use of incompatible adapters could cause damage and will void the validity of the warranty.
6. Do not tug the cable when you unplug the power. In order to avoid any bad contact of the equipment, the best way is to hold the plug and unplug it from the socket.
7. Put all the wires together tidily. Do not wind the wires, break the wires surface or cause any other damages.
8. Any people with the following conditions must avoid using this product: patients with malignant tumor, women during her pregnancy or period, patients with dermatosis focus on the hemorrhage, urgent patient with unidentified disease, and children and elders.
9. Keep the product out of the reach of children.

②

Components

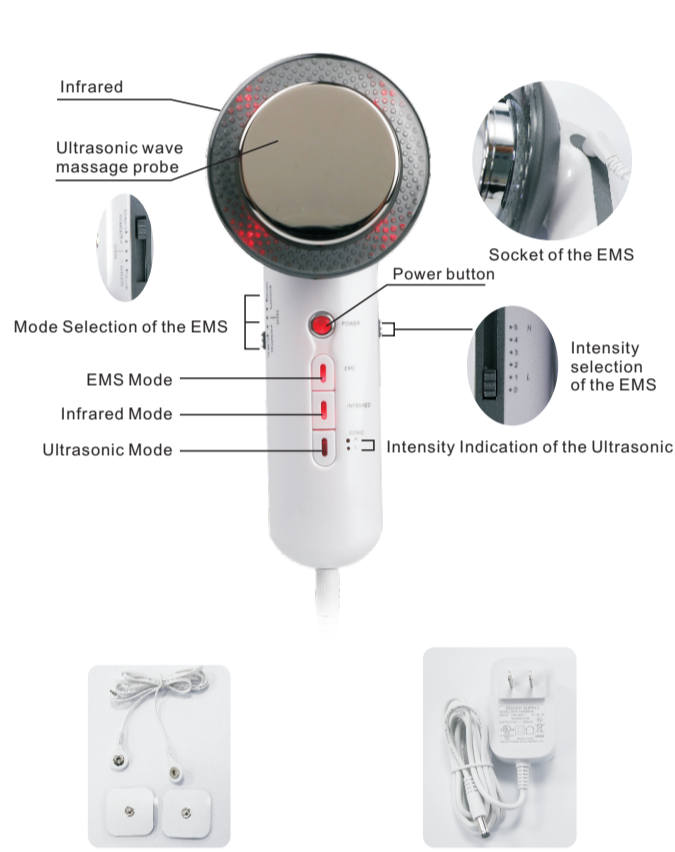③

④

Important Notice

1. You can change functions by directly pressing the corresponding button. To turn off the device, simply press the POWER button.
2. To ensure that the Ultrasonic function is working properly, press the Sonic and after 60 seconds drop a few drops of water onto the surface of the metal head and little bubbles should form. If that is the case, it means the device is functioning properly (EMS and Infra red light functions excluded).
3. When in use, place the metal head in full contact with the skin. Insufficient contact will cause a big difference in the resistance and energy dispersed on the surface of your skin.
4. Use the device under normal conditions and do not use it without applying body cream/gel.
5. The ideal length of time to use it is 5 to 15 minutes per session. Do not use it for more than 30 minutes. DO NOT apply to the same area for too long. Slide it over your skin to avoid the focus effect of the energy wave.
6. This device has a built-in safety device. When the temperature becomes too high, whether the device is being used or not, the temperature will automatically drop as a safety measure.
7. This device has a built-in time control device that will automatically shut it off when it is in idle mode, whether it is being used or not. Restart the machine if you want to use it again.
8. Remember to turn it off and unplug the cable when massage is over to avoid overheating problems caused by air vibration. This will also prevent possible burns and shortening the life of the device. Clean the centre-left side of the metal head. A soft cloth is recommended to avoid scratching the surface.
9. Never use water or chemicals, facial cleansers or liquid detergents to clean the metal head, as these could cause an electrical fault and shorten the life of the product.
10. In the event of an electrical fault, do not try to manipulate or fix the product yourself and contact your retailer.

⑤

Instructions for using the Ultrasonic

When using the Ultrasonic function, it is very important that you apply body cream/gel. Not only does it act as a conductor, but it also acts as an important component to your protection. The use of a body cream/gel is therefore mandatory when using this mode.

- Step 1:**
Clean the skin before using the device.
- Step 2:**
Apply body cream or slimming gel evenly over the skin. Remember to keep your skin hydrated when using the device.

(Note: Gel-type body creams is recommended for better conductivity. Apply the cream/gel to one area at a time to avoid it from air-drying and thus not being absorbed.)
- Step 3:**
Plug in the device, turn on the power button and the red light will start flashing.
- Step 4:**
Press the SONIC button-once for the low intensity (LOW), twice for the high intensity (HIGH). Use the low intensity on the face. Use high intensity on other parts of the body like arms, legs.
- Step 5:**
To begin the massage, refer to the details described and follow the recommended instructions you will find in the following pages of this Manual.

※**Note**
Do not use the device on joints: knees, elbows, wrists and ankles.

⑥

Body Treatment:

Arms massage:

Relax your arm muscles. Start from the top of the wrist (without going over it) to the armpit. Massage from one point to the other on both the outside and inside of the arm.

Follow the same procedure to massage the skin with repeated circular movements on the inside and outside of the arm.

Waist massage /lower abdomen

Make clockwise movements just below the stomach, by sliding the device over the skin from the waist and stomach area to the navel.

Then apply the device from the lower back to the side by moving it backwards and forwards.

Repeat this entire procedure several times.

⑦

Foot/legs massage

Starting at the lower leg above the ankle, slide the appliance up over the skin to the top of the thigh.


Use both hands to massage the lower legs upward in a circular motion.

Apply the device by sliding it over the skin on either side of the thighs from the inside to the centre.

⑧

Instructions for using EMS

- 1) Wipe off any dirt, grease or cream/gel from the parts to be massaged with a damp towel.
- 2) Connect the special cable to the electrode pads and place the adhesive electrodes on the part of the body that you want to stimulate (need to use 2 pads together).
- 3) Move the intensity selection key to position 0.
- 4) Connect the cable of the device to the mains outlet.
- 5) Turn on the device (the red light will on).
- 6) Select the mode you want to use:

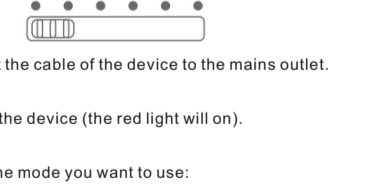- 7) Press the EMS button.
- 8) Set the desired intensity (up to 5 positions). The unit will stop working automatically after 10 minutes. If you want to repeat the operation or try other EMS mode, turn it on and start again.

At the beginning, we suggest treatment time is not more than 15 minutes, and then increase it gradually to 30 minutes each time.

After finishing the treatment, adjust the intensity selection key to the "0".

⑨

Recommend Chinese acupuncture points for EMS

- 1) When use TAPPING function, the recommended positions are as follows:
 - BELLY TIGHTENING**
 - 1. SHUIDAO (20)
 - 2. TAIYU (7)
 - 3. GUANYUAN (18)
 - 4. ZHISHI (53)
 - HIP TIGHTENING**
 - 1. HUANTIAO (54)
 - 2. CHENGFU (62)
- 2) When use Massage function, the recommended positions are as follows:
 - WAIST BEAUTY**
 - 1. TIANSHU (10)
 - 2. SHUFEN (6)
 - 3. DAIMAI (17)
 - 4. SHENYU (47)
- 3) When use Knead function, the recommended positions are as follows:
 - BELLY BEAUTY**
 - 1. GUANYUAN (18)
 - 2. DACHANGYU (48)
 - 3. ZHONGYUAN (8)
- 4) When use Scrapping function, the recommended positions are as follow:
 - SKIN BEAUTY**
 - 1. FEIYU (44)
 - 2. GANGYU (45)
 - 3. ROUMEN (5)

⑩

Instructions for using the Infrared function

When using the Infrared function, it is very important that you apply body cream/gel. Not only does it act as a conductor, but it also acts as an important component to your protection. The use of a body cream/gel is therefore mandatory when using this mode.

- Step 1**
Plug in the power
- Step 2**
Turn on the power (POWER) and the red light will blink.
- Step 3**
Press the INFRARED button.
- Step 4**
When the red light appears on the massage head, it means that the infrared wave is now working.
- Step 5**
Aim the infrared at the desired areas and the energy wave will be released to the skin.

After Use:

1. Make sure that all the buttons are off.
2. Turn off the power.
3. Put the electrode pads away.
4. Clean the head of the device and the adhesive pads with a damp cloth.
5. Never use solvents, petrol or other abrasive liquids to clean the device.
6. Put the device away in its box.

⑪

⑫

Weight Reduction

- 5) When use Slimming function, the recommended positions are as follows:

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