

A5纸张大小 要求： **105g**铜版纸 黑白印刷 装订成册

此页无需打印

coco

Heating Pad



- FOR PERSONAL USE ONLY -

IMPORTANT NOTE

Read the entire manual carefully before using this device.

IMPORTANT SAFETY INSTRUCTIONS DO NOT DESTROY

When using this or any other electrical appliance, always follow these basic safety precautions:



! DANGER-- To reduce the risk of burns, electric shock and fire, this product must be used in accordance with the following instructions.

IMPORTANT SAFETY INSTRUCTIONS

DO NOT DESTROY

DANGER:

To reduce risks of burns, electric shock, fire, and personal injury, this product must be used in accordance with the following

INSTRUCTIONS:

- 1.Read all instructions.
- 2.Do not use while sleeping.
- 3.Do not use on an infant or on an animal.
- 4.This pad is not to be used on or by an invalid, aperson, a person with

diabetes, or a person with poor blood circulation.

5.Do not use on areas of sensitive skin.

6.Burns can occur regardless of controller setting.Check skin under pad frequently.

7.Do not use in an oxygen-rich atmosphere.

8.Do not use pins or other metallic means to fasten this pad in place.

9.Do not bunch, sit on, or against, or crush pad - avoid sharp folds.

10.Never pull this pad by the supply cord.

11.Do not use the cord as a handle.

12.This unit should not be used by children without adult supervision.

13.Place the pad on top of, and not under, the part of the body needing heat.

14.Do not use with generators, power converters, or inverters.Use this pad only on a 110-120 volt AC circuit,unplug when not in use.

15.Do not tamper with this pad in any way.There are no user-serviceable parts; if for any reason this pad does not function satisfactorily, see warranty for consumer service information.

16.Do not use this pad with liniment, salve, or ointment preparations, especially ones that contain heat-producing ingredients. Skin burns could result.

17.Do not bend or pinch cord and do not use if cord is damaged.

18.Never leave the appliance unattended, especially if children are present.

19.Loosely loop cord when storing. Tight wrapping may damage cord and internal parts.

20.With not in use, store the heating pad in the original packaging, in a dark, cool, dry place.

21.Save these instructions.


This heating pad has a polarized plug (one blade is wider than the other) as a safety feature. This plug will fit into a polarized outlet only one way. If the plug does not fit fully into the outlet, reverse the plug. If it still does not fit, contact a qualified electrician.



DO NOT ATTEMPT TO MODIFY THIS PLUG OR DEFEAT THIS SAFETY FEATURE IN ANY WAY.

DO NOT ATTEMPT TO DEFEAT THIS SAFETY FEATURE.

Parts

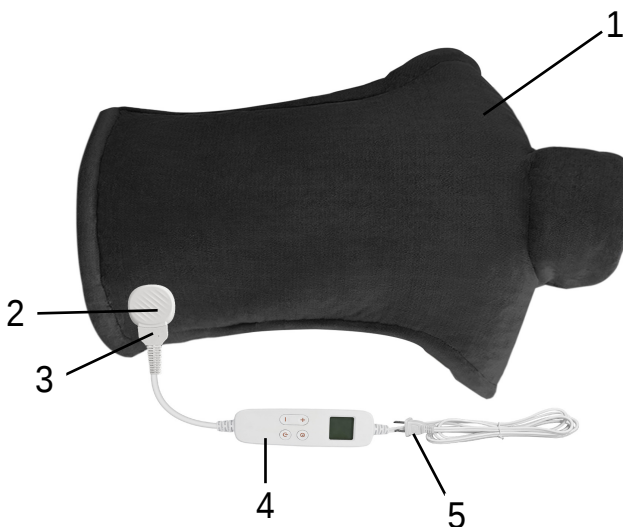


- 1.Heating Pad
- 2.Pad Connector
- 3.Controller Connector
- 4.Controller
- 5.Power Cord/Plug

How To Use Your Heating Pad

Before each use, always inspect the product to ensure there is no damage to the power cord, the Controller, or the Pad.

The Pad Connector and Controller Connector are designed to connect only one way. Align the inner prongs and then push them together until they securely snap into place.

- 1. secure (see Fig. 1).
- 2. Plug the product into a 120-volt AC electrical outlet.

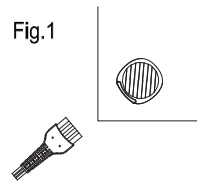

NOTE: If there is no power, make sure the AC electrical outlet is functional and the connection in step 1 is secure.

WARNING: Step 1 and Step 2 must not be reversed; otherwise, it will lead to F2 failure (see Trouble Shooting on page 9 for details).

- 3. Use the Controller to set the heat at the desired level (see Fig.2).




Display of heat level, display as heat level 6, display as heat level 3.

Corresponding temperature for heat level below:


- Heat Level 1---100°F;
 - Heat Level 2---110°F;
 - Heat Level 3---120°F;
 - Heat Level 4---130°F;
 - Heat Level 5---140°F;
 - Heat Level 6---150°F;
- Selectable timer settings: 30/60/ 90 minutes.

- 4. Press the POWER Button on the Controller to turn on the pad. Default to heat level 4 and 60 min timer upon startup. The display will illuminate, the all icon will display about 1.5 seconds, and then "130°F" "60m" will display, "60m" indicating the controller is setting 60-minute countdown to automatically shut off. (see Fig. 3).

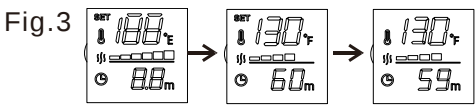

5. To change the heat level, press the “+” Button for a higher heat level or the “-” Button for a lower heat level. The desired heat level and the heat timer icon **SET** will appear for about 3 seconds, and then display the selected heat level, the real temperature and the remaining time on timer.


Fig.4

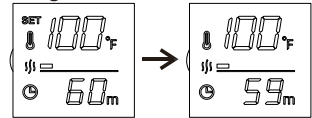

6. After turning on the Pad and selecting the desired heat setting (steps 4 and 5 above), press the TIMER Button “⌚” until the desired amount of time is displayed.

The display will illuminate the desired amount of time (for example “60m”), Time numerical display flickers for 3 seconds, and then the desired time setting will display steadily (see Fig. 5) indicates successful timer setup.

Fig.5

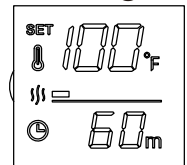

NOTE: *If the Control shuts off, make sure the Pad and Control connectors are securely connected and re-start the product by pressing “*

7. Place the Pad on top of the body. If the Pad does not evenly heat, make sure there are no folds or bunched up areas.
8. Turn the product off by pressing the POWER Button “⏻” and unplug it when it is not in use.

Auto Off Function:

This product will turn off after Timer finished.

WARNING:

Burns can occur regardless of the control setting; check skin under Pad frequently. Unplug the product when it is not in use.

Care and Maintenance

1. Disconnect the power cord from the electrical outlet.
 2. Make sure the Pad Connector is disconnected from the Control Connector by pressing the tabs on the sides of the Control Connector while pulling the connectors apart.
 3. Machine or hand wash the Pad using a mild detergent, cool water , and a gentle/delicate cycle.
 4. Line dry or lay the Pad flat to dry. Do not tumble dry.
 5. Make sure the product is completely dry before the next use.
 6. To store the product, loosely roll the Heating Pad toward the Storage Straps, loop the Straps over the rolled ends, and then store it in a cool, dry, dark location after each use.Avoid sharply folding the Heating Pad.
-

Electrical Specifications

Voltage	Frequency	Wattage	Amperage
120 VAC	60 Hz	135 W	1.13A

Frequently Asked Questions

1. When should a heating pad be used?

When ever a hot form of temporary relief is desired for personal comfort

2. What temperature does this product reach?

Maximum temperature does not exceed the Underwriters Laboratories temperature limit of 176 °F (80 °C) when set to “6”. Temperature varies by heat setting; you may need to use a lower heat setting depending on the specific application and your personal preference.

Trouble Shooting

Series NO.	Fault symptoms	Cause of the fault	Solutions
1	Show F1	Component failure	It must be replaced or repaired by professional technicians
2	Show F2	The controller cannot detect the heating pad	1.the connector is plugged in properly
			2.It needs to be replaced by professional technicians
			3.During use, please avoid pulling or exerting force
			4.Do not place the heating pad near heat sources or pets
3	Show F3	Frequency fault	Please check if the voltage and frequency of the mains power meet the requirements of the product
4	Show F4	Chip failure	It must be replaced or repaired by professional technicians
5	Show F5	Folding during use resulted in a short circuit in the heating wire	1.Please do not fold it
			2.It needs to be replaced or repaired by professional technicians
6	Show F6	Controller overheating fault	1.The controller cannot be used while wrapped
			2.The controller cannot be used near heat sources
			3.The controller can be used again after cooling down

FCC INFORMATION

WARNING: Changes or modifications not expressly approved by the party responsible for compliance could void the user's authority to operate the equipment.

NOTE: This equipment has been tested and found to comply with the limits for a Class B digital device, pursuant to Part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation. This unit generates, uses and can radiate radio frequency energy and if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation. If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try and correct the interference by one or more of the following measures:

- Reorient or relocate the receiving antenna. I
- Increase the separation between the equipment and receiver.
- Connect the equipment to an outlet on a circuit different from that to which the receiver is connected.
- Consult the dealer or an experienced radio/TV technician for help.

This device complies with Part 15 of the FCC Rules. Operation is subject to the following two conditions:

1. This device may not cause harmful interference.
2. This device must accept any interference received including interference that may cause undesired operation.

WARRANTY

If the product(s) is damaged or defective within a reasonable time, Coococo offers the warranty as below.

Period	Return for refund	Return for replacement	Who paid shipping?
Within 30 Days	Yes	Yes if undamaged	Coococo
31-90 Days	Yes after being justified	Yes after being justified	Coococo
91Days-1 Year	No	Yes after being justified	Buyer
1 Year-3 Years	No	Only to buyers who have received extended warranty	Buyer

Coococo warrants the product to be free from defects in workmanship and materials, under normal use, for a period of 12 months from the original invoice date, except as noted below. Coococo products can be repaired or replaced within warranty period or returned for refund within 90 days if the products break down in the period of proper use.

The warranty extends only to consumers, NOT to any retailers or resellers, This warranty is effective only if the product is purchased and operated in the country in which the product is purchased. Warranty claims for product defects, unless extended warranty is registered, expire 12 months from the purchase day. Warranty claims CANNOT be processed for items that have exceeded the warranty period. The warranty does NOT cover the replaced, refurbished product. The warranty does have other limitations as below "Warranty Limitations". To obtain warranty service on Coococo products, follow "Warranty Claim Procedures" below.

Warranty Claim Procedures

1. Any product the customer believes is defective and seeks for warranty, the customer must provide to Customer Support:

- a) Valid proof of purchase (see below the definition of valid proof of purchase);
- b) Pictures and/or description are required to claim defective merchandise;
- c) Model number and Series Number label attached to each product.

2. If the warranty claim is justified by Coococo Customer Support, customer(s) will obtain from Coococo a Warranty/Return Authorization Number (RMA Number).

3. Coococo reserves the right to specify that items be returned to the designated warehouse for inspection or be inspected by our representative in the field or be

retained by customer

4. Any claim for defective merchandise returns must be packed in original packaging.

5. The item(s) will be repaired or replaced or a refund credit will be issued once the item(s) are delivered. When returning items with RMA# and a prepaid shipping label provided by Coococo, Coococo takes responsibility for any damage or loss incurred in transit.

Valid proof of purchase under this Warranty

Coococo's direct-authorized resellers, currently only at Amazon/eBay/Walmart. For product(s) obtained from other Channels, please contact the customer service of the Channel directly for appropriate solution.

Warranty Limitations (this warranty does not include)

Normal wear and tear; Any condition resulting from other than ordinary residential wear or any use for which the product was not intended, such as use in rental or contract trade or commercial use; Any condition resulting from incorrect or inadequate maintenance or care; Damage resulting from misuse, abuse, negligence, accidents or shipping damage, exposed to environmental conditions temperatures, water, improper operation, improper installation, improper use of electrical supply; transportation damage; theft, vandalism; Damage resulting from the attachment or modification of any unauthorized accessories; Fail to show the valid proof of purchase; Lost or stolen products; Items that have expired their warranty period; No quality-related issues (after 30 days of purchase); Free product/Gift offer by Coococo; unauthorized repairs or modification; Replaced or refurbished product. Dissatisfaction due to buyer's remorse (after 30 days of purchase).

Refund Procedures:

(Refer to <https://coococous.com/pages/warranty-return-policy>)

Warranty Disclaimer

THE WARRANTY PROVIDED HEREIN SHALL BE THE SOLE AND EXCLUSIVE WARRANTY. COOCOCO SHALL HAVE NO LIABILITY FOR, AND EXPRESSLY DISCLAIMS ANY WARRANTY OR AFFIRMATION OF FACT, EXPRESS OR IMPLIED, OTHER THAN AS SET FORTH IN THIS WARRANTY STATEMENT, INCLUDING, WITHOUT LIMITATION:

- (1) THE IMPLIED WARRANTIES OF MERCHANTABILITY AND FITNESS FOR A PARTICULAR PURPOSE;
- (2) ANY WARRANTY OR AFFIRMATION OF FACT RELATED TO MISUSE, IMPROPER SELECTION, RECOMMENDATION, OR MISAPPLICATION OF ANY PRODUCT: AND
- (3) ANY WARRANTY OR AFFIRMATION OF FACT THAT THE CATALOGS, LITERATURE AND WEBSITES IT PROVIDES ACCURATELY ILLUSTRATE AND DESCRIBE PRODUCTS.

Coococo reserves all the right for final interpretation of these Terms of Warranty and Return Policy above and the right to change, modify, add or remove portions of these terms at any time without prior notification.

Extend Warranty for Free

Email us at service@coococous.com with "Warranty" as subject line to extend your warranty from 1 year to 3 years.



<https://www.facebook.com/uscoococo/>

OR 2) Email us at service @coococous.com with "Warranty" as subject line to extend your warranty from 1 year to 3 years.



Address: 573 BELLEVUE RD NEWARK, DE 19713 USA

Tel: (248) 819-2623 (Monday-Friday 9:00AM-4:30PM)

Facebook: www.facebook.com/uscoococo/

Email: service@coococous.com

Web: www.coococous.com

*Please provide your order number before contacting customer support