
STYLSMILE™

Invented With Love

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Toothbrush - 2 minutes twice per day

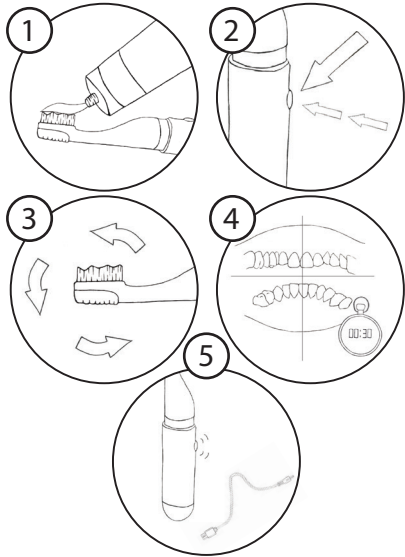(EN) Toothbrush - 2 minutes twice per day

1. Apply the **STYLSMILE** toothpaste to the full length of the toothbrush bristles. 2. Turn on the toothbrush, press once for a high speed or twice for a slower speed. 3. Brush teeth in circular motions for 2 minutes, please note the toothbrush will pulse every 30 seconds. 4. Use the 30 second pulses as a guide to brush the four sections of your mouth, move the toothbrush and brush the next section after each pulse. The toothbrush will automatically turn off after 2 minutes. 5. A flashing blue light around the power button indicates the **STYLSMILE** toothbrush needs to be recharged. Please charge the toothbrush using the USB power lead provided. **Rinse toothbrush and mouth after use. Do not swallow toothpaste. Repeat the brushing process twice per day. If the toothbrush bristles are too soft, please remove the brush head and use the FIRM brush head provided.**

Accelerator tray - 15 minutes twice per week

Before first use remove the plastic battery tab from the battery compartment of the LED Accelerator, and rinse the mouth tray in warm water.

1. Connect the LED Accelerator to the mouth tray by gently easing it on and ensuring the silicone lip is over the plastic area on all sides. 2. Apply Toothpaste to the upper and lower inside of the mouth tray, (for ease it is suggested you hold the tray vertically). 3. Place mouth tray in mouth, ensuring each tooth is coated in toothpaste. Wipe away any excess toothpaste touching the gums or lips. Turn on the LED light and leave for 15 minutes. After 15 minutes turn off the LED light and remove the tray from mouth.

Repeat twice per week and wash the mouth guard clean after each use. Do not swallow toothpaste and rinse your mouth.

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