STYLSMILE™

Invented With Love

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Instructies / Anleitung

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
Toothbrush - 2 minutes, twice per day

TOOTHBRUSH 5 v.

Accelerator Tray - 15 minutes, every day for 6 days

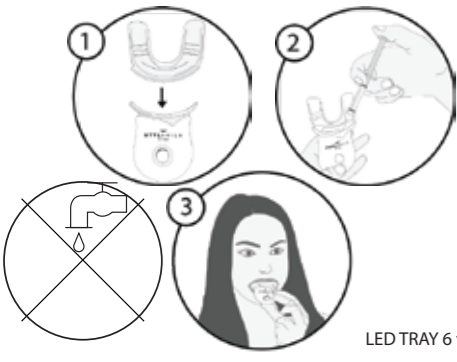LED TRAY 6 v.

EN **FR** **NL** **DE**
1 **NO** **SE** **IT** **IS**

TB01B

TB02B.06_CN

(EN)

Toothbrush - 2 minutes, twice per day

1. Apply the **STYLSMILE** toothpaste to the full length of the toothbrush bristles. **2.** Turn on the toothbrush, press once for a fast speed or twice for a slower speed. **3.** Brush teeth in circular motions for 2 minutes, please note the toothbrush will pulse every 30 seconds. **4.** Use the 30-second pulses as a guide to brush the four sections of your mouth, move the toothbrush and brush the next section after each pulse. The toothbrush will automatically turn off after 2 minutes. **5.** A flashing blue light around the power button indicates the **STYLSMILE** toothbrush needs to be recharged. Please charge the toothbrush using the USB power lead provided. **Rinse toothbrush and mouth after use. Do not swallow toothpaste. Repeat the brushing process twice per day. If the toothbrush bristles are too soft, please remove the brush head and use the FIRM brush head provided. Pull firmly to remove the head. Clean the toothbrush periodically. Rinse the handle and the brush head. Wipe surface with soft cloth. Do not use sharp objects to clean. Do not clean in dishwasher. Refill toothpaste available at www.stylsmile.com**

Accelerator tray - 15 minutes, every day for 6 days

Before first use remove the plastic battery tab from the battery compartment of the LED Accelerator, and rinse the mouth tray in warm water. **1.** Connect the LED Accelerator to the mouth tray by gently easing it on and ensuring the silicone lip is over the plastic area on all sides. **2.** Apply 1/4 Gel from syringe to the upper and 1/4 to lower inside of the mouth tray (for ease it is suggested you hold the tray vertically). **3.** Place mouth tray in mouth, ensuring each tooth is coated in gel. Wipe away any excess gel touching the gums or lips. Turn on the LED light and leave for 15 minutes. After 15 minutes turn off the LED light and remove the tray from mouth. **Repeat ever day for 6 days and wash the mouth guard clean after each use. Do not swallow gel. Rinse mouth after use. Keep LED Accelerator away from water. Wipe clean with dry cloth. Do not clean in dishwasher. Refil gel available at www.stylsmile.com**

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(EN) **TOOTHBRUSH & ACCELERATOR TRAY SAFETY PRECAUTIONS:**

THIS IS NOT A TOY. DO NOT SHINE IN EYES. THIS IS NOT INTENDED FOR USE BY PERSONS (INCLUDING CHILDREN) WITH REDUCED PHYSICAL, SENSORY OR MENTAL CAPABILITIES, OR LACK OF EXPERIENCE AND KNOWLEDGE, UNLESS THEY HAVE BEEN GIVEN SUPERVISION OR INSTRUCTION CONCERNING USE OF BY A PERSON RESPONSIBLE FOR THEIR SAFETY. CHILDREN SHALL NOT PLAY WITH THE APPLIANCE. CLEANING AND USER MAINTENANCE SHALL NOT BE MADE BY CHILDREN WITHOUT SUPERVISION. PLEASE STORE WITH CARE AND OUT OF REACH OF CHILDREN. DO NOT TAMPER WITH ANY INTERNAL COMPONENTS. THERE ARE NO USER SERVICEABLE PARTS INSIDE THE UNIT. DO NOT ATTEMPT TO REPAIR THE PRODUCT YOURSELF. DO NOT PLACE ANY PART OF THE PRODUCT NEAR HEAT SOURCES SUCH AS RADIATORS, COOKERS, ETC. DO NOT USE THIS PRODUCT FOR ANYTHING OTHER THAN ITS INTENDED USE. DO NOT OPEN THE PRODUCT OTHER THAN TO REPLACE THE ACCELERATOR BATTERIES. KEEP THE CHARGING CORD AWAY FROM WATER. DO NOT IMMERSE THE CHARGING CORD IN WATER OR ANY OTHER LIQUID. NEVER USE THE CHARGING CORD IF IT IS DAMAGED IN ANY WAY. IF THE CHARGING CORD IS DAMAGED, ALWAYS HAVE IT REPLACED WITH ONE OF THE ORIGINAL TYPE IN ORDER TO AVOID A HAZARD. IF THE APPLIANCE IS DAMAGED IN ANY WAY (BRUSH HEAD, TOOTHBRUSH HANDLE), STOP USING IT. IF THE APPLIANCE IS DAMAGED, CONTACT TEAMTOM@STYLIDEAS.COM. DO NOT USE THE CHARGING CORD OUTDOORS OR NEAR HEATED SURFACES. DO NOT CLEAN THE BRUSH HEAD, HANDLE OR CHARGING CORD IN THE DISHWASHER. DO NOT USE THE TOOTHBRUSH IN THE BATH OR SHOWER. IF YOU HAVE HAD ORAL OR GUM SURGERY IN THE PREVIOUS 2 MONTHS, CONSULT YOUR DENTIST BEFORE YOU USE THE TOOTHBRUSH. CONSULT YOUR DENTIST IF EXCESSIVE BLEEDING OCCURS AFTER USING THIS TOOTHBRUSH OR IF BLEEDING CONTINUES TO OCCUR AFTER 1 WEEK OF USE. THE **STYLSMILE** TOOTHBRUSH COMPLIES WITH THE SAFETY STANDARDS FOR ELECTROMAGNETIC DEVICES. IF YOU HAVE A PACEMAKER OR OTHER IMPLANTED DEVICE, CONTACT YOUR PHYSICIAN OR THE DEVICE MANUFACTURER OF THE IMPLANTED DEVICE BEFORE YOU USE THE **STYLSMILE**. THE **STYLSMILE** TOOTHBRUSH HAS ONLY BEEN DESIGNED FOR CLEANING TEETH, GUMS AND TONGUE. THE **STYLSMILE** ACCELERATOR HAS ONLY BEEN DESIGNED TO WHITEN TEETH. DO NOT USE EITHER FOR ANY OTHER PURPOSE. STOP USING THE APPLIANCE AND CONTACT YOUR DOCTOR IF YOU EXPERIENCE ANY DISCOMFORT OR PAIN. IF YOU HAVE MEDICAL CONCERNS, CONSULT YOUR DOCTOR BEFORE YOU USE THE **STYLSMILE**. THE **STYLSMILE** TOOTHBRUSH IS A PERSONAL CARE DEVICE AND IS NOT INTENDED FOR USE ON MULTIPLE PATIENTS IN A DENTAL PRACTICE OR INSTITUTION. STOP USING A BRUSH HEAD WITH CRUSHED OR BENT BRISTLES. REPLACE THE BRUSH HEAD EVERY 3 MONTHS OR SOONER IF SIGNS OF WEAR APPEAR. DO NOT USE OTHER BRUSH HEADS THAN THE ONES RECOMMENDED BY **STYLSMILE**. THOROUGHLY CLEAN THE BRUSH HEAD WITH WATER AFTER EACH USE. **GEL WARNINGS:** DO NOT SWALLOW GEL. AVOID CONTACT WITH THE EYES. NOT SUITABLE FOR CHILDREN UNDER THE AGE OF 7 YEARS. WHITENS NATURAL TEETH ONLY. DISCONTINUE TREATMENT IF YOU EXPERIENCE PAIN OR EXCESSIVE SENSITIVITY. THIS IS NOT A TOY. DO NOT USE THIS PRODUCT IF YOU ARE ALLERGIC TO ANY OF THE INGREDIENTS. STORE IN COOL DRY PLACE BETWEEN 5°C AND 30°C AWAY FROM HEAT AND DIRECT SUNLIGHT.

GEL INGREDIENTS: GLYCERIN, PEG-8, AQUA (WATER), PHTHALIMIDO-PEROXY-CAPROIC ACID, SODIUM POLYACRYLOYLDIMETHYL TAURATE/VP CROSSPOLYMER, MENTHA ARVENSIS LEAF OIL, SODIUM SACCHARIN, MICA, PVM/MA COPOLYMER, POLYSORBATE 20, ALOE BARBADENSIS LEAF JUICE, LECHTIN, MALTODEXTRIN, SODIUM HYDROXIDE, CI 77891. **TOOTHASTE SAFETY PRECAUTIONS:** DO NOT SWALLOW TOOTHASTE. AVOID CONTACT WITH THE EYES. CONTAINS SODIUM FLUORIDE 0.32% W/W (1450PPM F). NOT SUITABLE FOR CHILDREN UNDER THE AGE OF 7 YEARS. WHITENS NATURAL TEETH ONLY. DISCONTINUE TREATMENT IS YOU EXPERIENCE PAIN OR EXCESSIVE SENSITIVITY. DO NOT USE THIS PRODUCT IF YOU ARE ALLERGIC TO ANY OF THE INGREDIENTS. STORE IN COOL DRY PLACE BETWEEN 5°C AND 30°C AWAY FROM HEAT AND DIRECT SUNLIGHT.

INSTRUCTIONS FOR DISPOSAL (TOOTHBRUSH & TRAY): DISPOSING OF THIS PRODUCT CORRECTLY WILL HELP SAVE VALUABLE RESOURCES AND PREVENT ANY POTENTIAL NEGATIVE EFFECTS ON HUMAN HEALTH AND THE ENVIRONMENT WHICH MIGHT OTHERWISE ARISE FROM INAPPROPRIATE WASTE HANDLING. PLEASE CONTACT YOUR LOCAL AUTHORITY FOR FURTHER DETAILS OF YOUR NEAREST DESIGNATED COLLECTION POINT. PENALTIES MAY BE APPLICABLE FOR INCORRECT

APPLICABLE FOR INCORRECT DISPOSAL OF THIS WASTE, IN ACCORDANCE WITH YOUR NATIONAL LEGISLATION. USED ELECTRICAL AND ELECTRONIC EQUIPMENT (WEEE) SHOULD NOT BE MIXED WITH GENERAL HOUSEHOLD WASTE. FOR PROPER TREATMENT, RECOVERY AND RECYCLING, PLEASE TAKE THIS PRODUCT TO DESIGNATED COLLECTION POINT. PLEASE RECYCLE USED ACCELERATOR TRAY BATTERIES, USED ACCELERATOR TRAY AND **STYLSMILE** TOOTHBRUSH AT YOUR NEAREST RECYCLING CENTRE. **INSTRUCTIONS FOR DISPOSAL (TOOTHASTE / GEL):** CORRECT DISPOSAL OF THIS PRODUCT: RECYCLE PACKAGING ACCORDING TO NATIONAL LEGISLATION.

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