DR-HO'S

# TripleAction

## Back Belt

Imagine, a **good life with less pain.** Enjoy your work, play and sleep!

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# Thank you

Thank you and congratulations on the purchase of your new DR-HO'S Triple Action Back Belt! You have made a smart and important decision with regard to self-treatment and care. Many people have been pleasantly surprised with the degree of improvement they have experienced with this state-of-the-art device, and you too will discover how it helps to improve your quality of life and reduce tension.

DR-HO'S Triple Action Back Belt is an elaborate but simple-to-use device. From chronic to acute muscle tension and aches, this unique device provides quick and effective relief from many of your conditions.

Our aim is to develop unique and effective products that help people improve their quality of life. I would like to take this opportunity to express our sincere thanks to you, our customer, for your support.

DR-HO'S Triple Action Back Belt is a medical device class IIA so please read all guidelines carefully and keep it for future reference.

If you have any questions or require assistance, please contact us via [info@drhonow.eu](mailto:info@drhonow.eu) or  .

Yours in health & happiness,

Dr. Michael Ho  
Pain & Headache Specialist

## TABLE OF CONTENTS

Safety precautions.....	3	Cleaning & maintenance .....	10
Information on intended use .....	7	Problem solving.....	11
Specifications .....	8	Disposal.....	12
Treatment programs.....	8	Symbols.....	13
Batteries.....	9	Electromagnetic compatibility.....	14
Electrode pad user guidelines.....	10	Warranty & Contact.....	14

## SAFETY PRECAUTIONS

### Warnings

The long-term effects of electrical stimulation are unknown.

**PLEASE CONSULT YOUR DOCTOR BEFORE USING THIS DEVICE. THE DEVICE MAY BE DANGEROUS IF USED BY CERTAIN SUSCEPTIBLE INDIVIDUALS.**

### DO NOT USE THIS SYSTEM

- If you have a cardiac pacemaker, implanted defibrillator, or other implanted metallic or electronic device. Using the device could then cause electric shock, burns, electrical interference, or death.
- Near magnetic resonance imaging (MRI) equipment. The device is MR Unsafe.
- If you have abnormal high blood pressure; an undiagnosed pain syndrome; a tendency for internal bleeding or if hemorrhage is likely.
- Following recent surgical procedures if stimulation may disrupt the healing process.

- Near your eyes, to your head or face.
- Together with medical electronic devices such as an artificial heart or lung or respirator.
- In the presence of electronic monitoring equipment (e.g., cardiac monitors, ECG alarms), because it may cause them to malfunction.
- On open wounds, rashes, swollen, red, infected or inflamed areas or skin eruptions (e.g., phlebitis, thrombophlebitis, varicose veins) or on top of or close to cancerous lesions.
- Over insensitive areas of skin.

### DO NOT APPLY THE BODY PADS

- On the opposite sides of your head, because the effects of brain stimulation are unknown.
- To your head, face, the front and sides of your neck, your throat, or on both sides of your thorax.
- To stimulate the carotid sinus nerves, particularly in patients with a known sensitivity to the carotid sinus reflex.

## SAFETY PRECAUTIONS

- To stimulate your front neck or mouth, because this could cause severe muscle spasms resulting in closure of your airway, difficulty in breathing, or adverse effects on heart rhythm or blood pressure.
- Across your chest because this may cause cardiac arrhythmias.
- On open wounds, rashes, swollen, red, infected or inflamed areas or skin eruptions or on top of or close to cancerous lesions.
- Over insensitive areas of skin.

### DO NOT USE THE PAIN THERAPY OPERATING CONSOLE WITH THE BACK BELT OR PADS ON

- Children or infants, because the device has not been evaluated for pediatric use.
- People who have lost their mental capacity.
- Pregnant women, because it is not yet known whether it is safe to do so.

Seek advice from your doctor if you have suspected or diagnosed heart disease or epilepsy.

### IN GENERAL: IF YOU FEEL YOU WOULD BENEFIT FROM USING DR-HO'S TRIPLE ACTION BACK BELT AND WOULD LIKE TO GIVE IT A TRY, PLEASE CONSULT YOUR DOCTOR.

### DO NOT USE THE TRIPLE ACTION BACK BELT OR BODY PADS

- when in bath or shower;
- while sleeping;
- while driving, operating machinery, or any activity which might cause involuntary muscle contractions that may put the user at undue risk of injury.

### PRECAUTIONS - PAIN THERAPY OPERATING CONSOLE & BELT

- If you feel discomfort or abnormal symptoms during treatment, stop using the device and consult a doctor before resuming treatment.
- If you have any problems setting up, maintaining or using this device, contact customer service and report anything unusual.
- Keep the device out of the reach of infants, children or pets. Any inhalation or swallowing of small parts is dangerous and can be potentially fatal.
- Only use accessories that have been specified or authorized by the manufacturer. The use of unauthorized accessories may cause damage to the device or danger to the user.
- Do not use the device in areas of high humidity, such as the bathroom or in the shower, as this may damage it.
- Ensure that metallic items with conductive capabilities, such as

## SAFETY PRECAUTIONS

- belts, are kept away from the electrodes during treatment.
- If using rechargeable 1.5 V batteries, make sure you use a CE approved battery charger. Never connect the device directly to a battery charger or to other powered equipment.
- Do not use the device in the presence of a flammable anesthetic gas mixture, compressed oxygen or nitrous oxide.
- If the device is not functioning properly, you should stop using it.
- Only use the device for its intended purpose.
- Do not use the device while wearing other electronic devices, as this may cause damage.
- Do not use the Pain Therapy Device near a mobile phone, as this may cause it to malfunction.
- Do not throw the batteries into a fire because they may explode.
- Dispose of the device, batteries and components in accordance with applicable legal regulations. Unlawful disposal may cause environmental pollution.

### PRECAUTIONS - BODY PADS

- Apply pads to normal, healthy, dry, clean skin only.
- If you wish to move the pads to another location of your body during treatment, turn the device off first, then remove the pads

and resume treatment.

- Do not bend or fold the pads, as this may cause the pads to malfunction.
- Ensure that metallic items with conductive capabilities, such as necklaces, rings or belts, are kept away from the electrode pads during treatment.
- Always store the pads on the transparent protective film when not in use.
- Do not apply any solvent or any other liquids (except conductive gel or water) to the pads or to your skin because these will cause the pads not to stick and to malfunction. The pads are pre-gelled and will adhere to your skin when moist.
- Always apply a few sprays of the conductive gel onto the sticky surface of the electrodes whenever these feel dry, to ensure the pads will provide a comfortable, powerful and effective treatment.
- The simultaneous connection of a patient to high frequency surgical equipment may result in burns where the electrode pads have been attached and possibly damage the operating console.
- Do not leave pads attached to the skin after treatment.
- Placement of the electrode pads and stimulation settings should be based on the instructions of the prescribing practitioner.

## SAFETY PRECAUTIONS

- Do not share electrode pads with anyone else because this may cause skin irritation or infection. People should have their own set of pads for personal use.

### IMPORTANT:

Both of the pads in a set must be placed on the body at the same time in order to work. Make sure the pads are not touching each other because otherwise no impulses will be felt.

Stop using the pads if the individual packaging has been open for more than three months, or if the pads lose their adhesiveness.

### PRECAUTIONS - WIRES

- Do not bend or pull the end of the wires.
- Replace wires when broken or damaged.
- Long wires have the potential to strangulate. Be careful when using these in the neck and shoulder areas.
- Make sure the wires are properly plugged into the Pain Therapy Device and, when using the pads, that the pads are stuck tightly onto the skin area you wish to treat to ensure effective therapy.

### NOTE:

Do not treat one area for extended periods of time (more than 40

minutes per session, up to 3 times per day) because your muscles in that region may become exhausted and sore.

Only use the original electrode pads, gel and wires supplied by the manufacturer because other electrode pads and wires may not be electrically compatible with your device.

### ADVERSE REACTIONS

You may experience:

- skin irritation or hypersensitivity due to the electrical stimulation or Conductive Gel;
- skin irritation and burns underneath the electrode pads when applied to the skin.

### YOU SHOULD STOP USING THE DEVICE AND CONSULT YOUR DOCTOR IF YOU EXPERIENCE ANY ADVERSE REACTIONS.

### SERVICE

Do not attempt to repair or service this device, but contact Customer Service in case you have problems or questions.

### CUSTOMER CARE

If any of the instructions or warnings are not understood please contact Customer Service for assistance.

### SAVE THESE INSTRUCTIONS!

## INFORMATION ON INTENDED USE

DR-HO'S Triple Action Back Belt is an innovative back pain belt with integrated electrodes which you operate using the DR-HO'S operating console (model P.T.S.-IIB). The device has been developed to bring relief to tense or stiff muscles and muscle and nerve pain. It works by electrically stimulating the nerves on the surface of the skin and in the painful area. The Triple Action Back Belt can be used successfully in combination with other types of treatment and medication. Paramedics like physiotherapists and chiropractors have already been using neuromuscular electrical stimulation (NMES) for decades.

### How does it work?

Scientific theories state that therapy involving electrostimulation can be effective in several different ways:

- The subtle electrical impulses pass through the skin to nearby nerves, thereby preventing the pain signal from reaching the brain from the source of the pain.
- The electrical impulses stimulate the production of the body's natural painkillers, such as endorphins and serotonin.

3. Electrical stimulation is also believed to improve blood circulation. The muscles contract and relax to the rhythm of the electrostimulation and doing this repeatedly enhances blood flow and therefore stimulates your circulation.

DR-HO'S Triple Action Back Belt includes Auto-Modulating Pulse (AMP) technology, adding variety to the treatment to prevent your body becoming accustomed to it and any reduction in its effectiveness.

This product is certified as medical device in the European Union under the Medical Device Directive 93/42/EEC by SGS CE1639, exclusively for the indication(s) of relief of pain associated with sore or aching muscles. Other non-medical uses escribed to this device are not within the scope of CE certification.

Users should be aware product performance and/or safety has not been evaluated by SGS for non-medical purposes.

## SPECIFICATIONS

Power source	DC 3.0V, 2 AAA batteries
Method of line current isolation	Type BF applied part
Number of output channels	1
Number of output programs	3
Output intensity level	5 steps
Timer range	10, 20, 40 min
Weight	Pain Therapy Device: 74g, Triple Action Back Belt: approx. 80g, Small body pads: 10g, FlexTone pads: 37g
Dimensions	Pain Therapy Device: 12 x 7 x 2 cm, Triple Action Back Belt: 110 x 80 x 23 cm, Small body pads: 46 x 46 mm, FlexTone pads: 80 x 120 mm
Operating environment	Temperature 5~40°C, Humidity ≤80%RH, Atmospheric pressure 86~106kPa
Storage environment	Temperature: Pain Therapy Device -20 ~ 55°C, Electrode 5~40°C, Humidity ≤95% RH, Atmospheric pressure 50~106 kPa
Shelf-life	Electrode pads: 1 year, Pain Therapy Device: 3-10 years, Triple Action Back Belt: 5 years (Note: depending on how well these are maintained)

## TREATMENT PROGRAMS

Therapies:	A	B	C
	To temporarily increase local blood circulation and relaxation.	Symptomatic relief and management of chronic, intractable pain and of pain associated with arthritis. Stimulation of healthy muscles to improve and facilitate muscle performance.	
How long does the therapy last?	This depends on the timer setting chosen. The therapy ends automatically.		
How do I personalize the therapy?	You select program and intensity yourself.		
How long should I use the device and how often?	Each session should last 30 to 40 minutes. No more than 3 sessions per day.		
When should I stop using it?	When the pain has subsided or stopped, if treatment is uncomfortable, in the event of an adverse reaction, and if the condition worsens.		

Note: many people experience immediate relief while others require regular use to feel the benefits. The results vary and will depend on your underlying conditions and how often

you use the device. If your pain does not subside, you can try increasing the intensity or you can change the program.

## BATTERIES

### Insert 2 AAA batteries before use:

1. Remove the battery cover on the back.
2. Insert batteries correctly so that the markings correspond.
3. Refit the battery cover.

### Important Notes

Low battery power will cause the device to lose its effectiveness. Please change to new batteries whenever stimulation weakens.

If the indicator lights fade, or if the stimulation feels weak on your skin, this may be an indication of low battery power.

### WARNING!

Do not dispose of the batteries with normal household waste, but hand

them in at an official collection point in accordance with local regulations governing the disposal of such materials.

**Note:** To avoid battery corrosion, remove the batteries from the device if not used for 4 weeks or more.

If the device is used once a day for about 30 minutes per treatment, one pair of AAA Alkaline batteries should last 3 to 6 months.

We recommend the use of Alkaline batteries because they provide a steady power level until it runs out. Non-Alkaline and rechargeable batteries have a tendency to decrease in power level over time, and you might mistakenly think the device is not working right because the treatment is not as strong.

**Please refer to the Quick Guide for pad placement and user instructions.**

## ELECTRODE PAD USER GUIDELINES

- Wash skin thoroughly and then dry it before applying the electrode body pads.
- Apply the whole surface of the electrode pads firmly to the skin.
- Do not use electrode pads that do not stick properly to the skin, or only partially stick to the skin.
- If the skin looks red after treatment, you should not treat that area until the redness has disappeared.
- The electrode pads can be re-used 30 to 50 times depending on your skin type and other conditions. Whenever you notice that the pads are not sticking properly to your skin, or if the pads feel dry, you

should moisten the pads with a small amount of DR-HO'S Electro Therapy Conductive Gel using the spray bottle (do not wash them with water directly). This will help restore the adhesiveness and conductivity of the electrode pads. If moistening the pads fails to make them sticky again, they may need to be replaced.

To order new sets of electrode body pads and DR-HO'S Electro Therapy Conductive Gel, please visit the website of your supplier.

## CLEANING & MAINTENANCE

### DR-HO'S Triple Action Back Belt and Pain Therapy Device operating console:

- Use a wet cloth or neutral cleanser to wipe the surface.
- Do not use volatile liquids such as alcohol, thinner or gasoline.
- Keep well away from moisture, high temperatures, direct sunlight or water droplets.
- Please do not disassemble, repair or rebuild the device yourself because this will invalidate the warranty.

### Electrode body pads:

- Keep the electrode pads clean and never expose them to high temperature and direct sunshine.
- Keep the electrode pads at temperatures of 5~40°C.
- Gently clean the electrode pads with a moist cloth or use a brush to remove dust.
- Do not clean the pads with hot water or chemicals.
- Do not touch the gel of the pads with your fingers, as this could reduce their adhesiveness.

## PROBLEM SOLVING

Problem	Possible causes	Try this solution
The device does not appear to be working properly.	The two pads are not attached properly.	Attach both pads to the body and make sure the pads are not touching each other.
	The wires from the pads are not properly connected to the Pain Therapy Device.	Make sure the wire from the electrodes is connected to the Pain Therapy Device.
	The wire plug is not firmly connected to the output port at the top of the Pain Therapy Device.	Push the wire plug firmly into the output port of the Pain Therapy Device.
	Weakening electrical pulses.	Turn the ON/OFF intensity dial in a clockwise direction.
	Power is running out.	Replace batteries, even if you bought them recently.
	The pads are not properly in contact with the skin.	Press the electrodes firmly onto the skin.
	The gel surface of the pads is dry or contaminated.	Moisten the pads with a little water or replace with new pads.
The device stops working before treatment has ended.	Your skin is covered with oil or lotion that is preventing the pads from sticking.	Wash your skin with water and try putting the pads on again.
	The batteries have run out.	Change the batteries.
	The wire plug has not been fully inserted into the output port.	Reconnect the plug to the output port.
Three yellow lights on the Pain Therapy Device are flashing and you cannot feel any stimulation.	The wires between the electrode pads and the main wire are not connected.	Reconnect the wires from the electrodes or foot pads to the main wire leading to the Pain Therapy Device.
	One of the pads has become disconnected from your skin or a wire is unplugged.	Reconnect the body pads to your skin or reconnect the wire. Turn the Pain Therapy Device off and then turn it back on again to restart the device.
Three green lights are flashing and you cannot feel any stimulation.	The auto shutdown timer has activated and the device has powered of.	Turn the power dial until it clicks to the OFF position first. Then turn the power dial back to the ON position.

## PROBLEM SOLVING

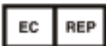
Problem	Possible causes	Try this solution
Batteries have been installed but the light does not turn on.	The battery compartment is dirty and the connectors are blocked.	Examine and clean.
	The batteries are not at full power and have not been installed correctly.	Replace the batteries or install them correctly in the compartment.
	There is a problem with the connectors in the battery compartment.	Try pulling the connectors outward slightly using a small pair of pliers.
	The power dial has been turned OFF and then ON.	Turn the power dial OFF and then ON again.
The skin has turned red and/or itches or prickles.	You are allergic to the electrode pads.	Stop using the machine until you ascertain the cause of the skin irritation.
	You have a small cut or acne.	Do not put electrodes on recently shaved skin, areas that have suffered insect bites, scrapes, acne or any open wounds.
	The adhesive surface of the pads is dirty or too dry.	Wash the adhesive surface of the pads gently with your fingertips for about 3 seconds under slow running water.
	The therapy has gone on too long.	Use the device for a maximum of 40 minutes per treatment and no more than 3 times a day.
	The 2 pads are not attached firmly and fully to the skin.	Attach the 2 pads firmly to the skin to ensure full contact.
	The surface of the pads is worn out.	Replace the set of pads.

## DISPOSAL

The Pain Therapy Device, Triple Action Back Belt and accessories should not be thrown away with household waste, but should be recycled by you or the supplier.

Dispose of packing materials appropriately in accordance with the rules and regulations in your area to prevent harm to the environment.

## SYMBOLS

Symbol	Meaning
	SERIAL NUMBER. This symbol will be accompanied by the manufacturer's serial number.
	CONSULT INSTRUCTIONS FOR USE.
	Unique Device Identifier (UDI); the Application Identifier (01) represents the device identifier (DI) ("0697074222XXXX"), the Application Identifier (11) the production date/date of manufacture ("YYMMDD", with "YY" the last two digits of the production year, "MM" the production month and "DD" the production day), the Application Identifier (21) the serial number of the device ("WWWWW-XXXXXX").
	Medical device.
	Indication of conformity with the essential health and safety requirements set out in European Directives 0120 and 1639 – Notified body identification (SGS).
	DATE OF MANUFACTURE. This symbol will be accompanied by the date of manufacture.
	MANUFACTURER. This symbol will be accompanied by the name and the address of the manufacturer.
	AUTHORIZED REPRESENTATIVE IN THE EUROPEAN COMMUNITY. This symbol will be accompanied by the name and the address of the authorized representative in the European Community, adjacent to the symbol.
	Dispose of separately from other household.
	Type BF applied part (patient isolation from electrical shock).
	IP classification.
	ATTENTION! Caution! Follow operating instructions! Failure to do so could place the patient or operator at risk.
	Batch code.
	MR UNSAFE. This device cannot be taken into the controlled access area of an MRI system.

## ELECTROMAGNETIC COMPATIBILITY

DR-HO'S Pain Therapy Device has been tested and found to comply with the electromagnetic compatibility (EMC) limits for medical devices. These limits are designed to provide reasonable protection against harmful interference in a typical medical installation.

CAUTION: do not use this device simultaneously with devices that have high EMI levels.

## WARRANTY & CONTACT

We hereby guarantee that if this product malfunctions or fails in any way within one year from the date of purchase, the purchaser will not be charged for any of the parts needed and labor costs incurred to repair the defective device, if it is used under

normal conditions and in accordance with this instruction manual. This warranty is only valid if the purchaser can provide proof of purchase.

Please contact your local supplier in case you have questions.

Dr. Ho Europe B.V.  
1261 XR - 5, the Netherlands  
info@drhonow.eu | www.drhonow.eu

For customer service:  
E: info@drhonow.eu  
Your local supplier

Are you happy with your DR-HO'S Triple Action Back Belt®? Share your posts, tips & tricks! Like & follow us @drhonow & tag your posts

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---

# DR-HO'S motionciser

Continuous passive motion exercise to recover legs & feet + relax body & mind

# DR-HO'S neck pain pro

Unique pain relief system for shoulders, neck & head incl. all body pain relief electrodes

### Fast relief and support for back pain.

**WARNING:** DR-HO'S equipment is not intended to diagnose, treat or cure any disease or condition. If you experience pain or health issues, always consult your doctor for the right diagnosis and treatment.

**PRECAUTIONARY MEASURES:** Do NOT use DR-HO'S equipment or other electrostimulators if you wear a pacemaker or are pregnant. If you have any questions, please consult your doctor before using the device.

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Guang Zhou Xinbo Electronic Co Ltd.
MedNet EC-REP GmbH