DR-HO'S

TripleAction

Back Belt

Imagine, a **good life with less pain.** Enjoy your work, play and sleep!

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DR-HO'S TRIPLE ACTION BACK BELT

1 Pain Therapy Device (model P.T.S-IIB)

This operating console is central to DR-HO'S Triple Action Back Belt system. You can use it in combination with the Triple Action Back Belt or with the individual body pads.

2 Program button

You can choose from 3 programs with different intensity intervals.

3 Timer button

You can choose from 3 treatment times, namely 10, 20 and 40 minutes. The treatment will end automatically at the end of the chosen program.

4 Electrode outlet

Connects to either the Triple Action Back Belt or body pads.

5 Intensity switch

Increase or decrease the intensity of the intended treatment by turning the wheel using your index finger.

6 Triple Action Back Belt

Ergonomic design for effective treatment or back support, whenever you need.

7 Dual-core support band

The belt can be tightened or loosened as desired to provide extra support for your core.

8 Unique built-in electrodes

Clinically proven to help relax tight muscles, improve circulation and relieve pain.

Apply DR-HO'S Electro Therapy Conductive Gel over the entire surface of the two electrodes before use and respray every 30 minutes during treatment or if you experience any slight stinging sensations.

9 Inflatable lumbar support

For extra support for your lumbar region or lower back which you can inflate as desired using the little pump.

10 Belt extender

You can use this to extend the black back pain belt if it is too short for you.

DR-HO'S Electro Therapy Conductive Gel

This patented gel is developed to moisten the electrodes and the skin, to reduce the amount of contact between the electrode surface and the skin, and to ensure more comfortable and effective treatment.

Mini Spray Bottle

To easily apply the conductive gel over the surface of the 2 electrodes.

Reusable body pads

The electrode pads provide effective treatment for specific areas of the body. These zones can be found in the overview on page 6.

IMPORTANT!

DR-HO'S Triple Action Back Belt is a medical device class IIA. Please read all guidelines carefully and keep it for future reference.

This complete system (the Triple Action Back Belt and the electrode body pads) is intended to relieve pain associated with sore or aching muscles in your (lower) back and other parts of your body, due to strain from exercise or activities of daily living, regular household and work activities.

Usage advice

Frequency of the treatment should be at least one to three times per day. However, you may treat yourself for as often as you like. Chronic problems will need more frequency and longer treatments.

Note: Do not use this system while sleeping, in the bathroom, when taking a bath or having a shower, driving or operating machinery. Treatment may cause involuntary arm movements and may make you drowsy.

During treatment you can sit straight up, lean back, or lie down.

You can use the Triple Action Back Belt with and without stimulation. Before using the belt with stimulation you should connect it to the DR-HO'S operating console. If used without stimulation the back pain belt provides excellent and ergonomic compression support for your back.

Care & maintenance

Do not immerse the Triple Action Back Belt in water or wash it in a washing machine. The electronic components may be damaged if they come into contact with water. If you wish, you can wipe the back pain belt clean with a damp cloth and then simply hang it up somewhere to dry.

To avoid battery corrosion, remove the batteries from the device if not used for 4 weeks or more.

Electrode body pads

- Do NOT use the large and small body pads under your feet.
- If the pads feel dry or irritating, spray them with a small amount of DR-HO'S Electro Therapy Conductive Gel, wait 30 seconds and re-attach the pads to your body.
- For optimal results replace the electrode pads after 30 to 50 sessions.
- Electrode pads are intended for personal use and should not be used by others.

Refer to the pad placement chart on how to position the large or small format pads, depending on the area to be treated.

If you wish to order replacement electrode body pads or conductive gel please visit the website of your supplier. If you have any questions, please contact your supplier.

1 Insert 2 AAA batteries in the back of the device. It is best to use alkaline batteries.

2 Triple Action Back Belt

A Plug the wire firmly into the top of the device.

B Both the wires must be clicked into the back belt.

C Before use you should spray DR-HO'S Electro Therapy Conductive Gel over the entire surface of the electrodes and repeat this during treatment, or as soon as you experience any light irritation or tingling.

D Put on the black back belt and make sure the electrodes are pressed against your bare skin in your lower back region.

E Use the pump to inflate the lumbar support cushion until the electrodes on the inside of the back belt are firmly in contact with your bare skin. It should feel comfortable.

Extending the belt
If the back belt is too short, you can extend it using the belt extender as shown in the picture.

2 Electrode body pads

A Insert the plug from the selected set of electrode pads into the Pain Therapy Device.

B To use the pads, remove the transparent protective film from both electrode pads to expose the self-adhesive and sticky gel layer.

C Spray a little bit of gel evenly over the surface of the pads, using the handy spray bottle supplied.


D Refer to the pad placement chart on page 6.
Note: The system ONLY works if you attach both pads from the same set to your body.

E Attach both pads to the area of the body that you wish to treat. Make sure the pads do not touch each other.

F After using the pads, stick the transparent film onto the self-adhesive gel layer again.

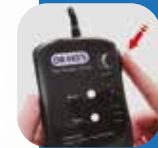3 Turn the switch to 'ON' by turning the wheel to the left with your index finger. Then increase the intensity slowly until your muscles start moving of their own accord.

4 Press to select modes A, B or C. A flashing light will indicate the selected mode.

5 Press to switch timer to OFF after your 10, 20 or 40 minutes treatment. To continue the treatment first turn the device off and then back on again.

6 Click the operating console into place on the front of the belt.

Relax and enjoy!

Extra support

You can tighten the dual support belt to provide extra support for your core, lower back or if you need to lift something.

For extra support for your lumbar region you can pump up the cushion until it feels comfortable. Keep pressing the button to deflate the inflatable pillow after use.

PAD PLACEMENT CHART

Neck and shoulders

Neck and shoulders

Elbows and forearms

Elbows and forearms

Shoulders

Wrists

Hands

Hips

Lower back

Thighs

Lower back and gluteal muscles

Gluteal muscles and thighs

Knees

Shins

Calves

Ankles

www.drhonow.eu

Fast relief and support for (lower) back pain.

Disclaimer: DR-HO'S equipment is not intended to diagnose, treat or cure any disease or condition. If you experience pain or health issues, always consult your doctor for the right diagnosis and treatment.

Precautions: Do NOT use DR-HO'S or any other electronic stimulators if you wear a pacemaker or if you are pregnant, due to unknown safety issues. If you have any questions or concerns, please consult your doctor before use.

“Imagine, a good life with less pain. Enjoy your work, play, sleep and life!”

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Guang Zhou Xinbo Electronic Co Ltd.
MedNet EC-REP GmbH