

# USER GUIDE: Malem Universal Wireless (MO15)



**MONITOR MULTIPLE PEOPLE**      **SINGLE RECEIVER**      **WIRELESS AUTOTUNING**


## Overview

- One receiver unit (Rx) can monitor up to 7 transmitters (Tx)
- When the Tx is triggered, it will sound and the light will flash
- At the same time, the Rx will sound and the corresponding LED light will flash, so that you know which Tx has been activated




## Single-tone settings:

*The MO15 is ideal for any situation when multiple patients are being monitored, e.g. school dormitories, care home, hospitals*

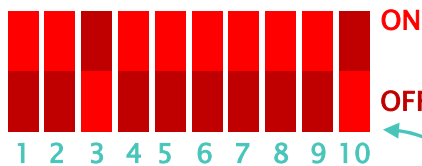

*Single tone example set for sound 3 only*

## 8-tone settings:

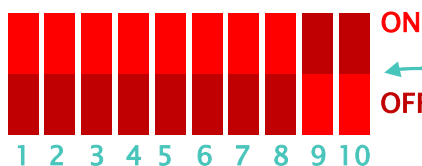

## Sound Selection

- The sound selection switch is located inside the Rx battery cover
- Choose a single sound by sliding any of the switches between 1–8 to “ON”
- To select random 8-tone setting, turn switch 9 “ON”
- Switch 10 “ON” means the Rx sounds, switch 10 “OFF” to silence it if preferred

## Receiver LED Number Selection

To select which number LED lights up on the Rx when a particular Tx is activated:

- Open the battery cover on the Tx
- Select the required number on the right hand side
- Re-tune the Tx and Rx



*Select the number that you want to light up on the Rx, e.g. 3*

*When the Tx is activated, the corresponding LED light will flash on the Rx*

