




# Ariel

## Users Manual



# Biohackn's Light therapy devices use evidence-backed Red, and Near-infrared light backed up by the science of Photobiomodulation.

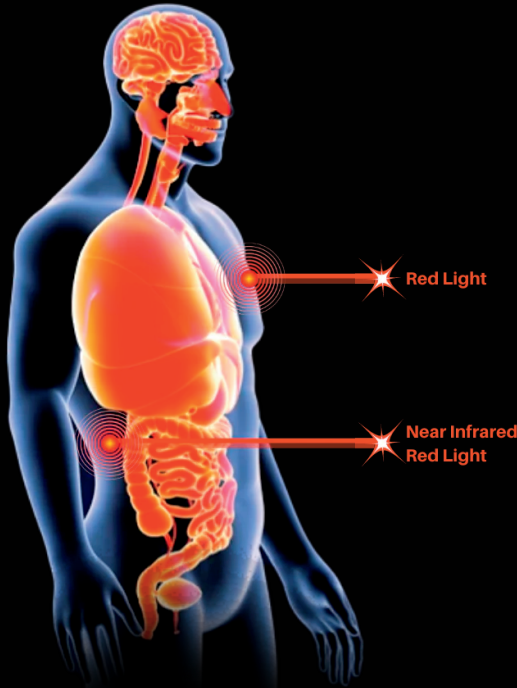

## **RED LIGHT AT 660NM**

Most of red light absorbed by your skin, leading to increase of collagen production and optimal skin health.

## **660NM: 850NM = 1:1**

Provides ultimate versatility for overall health by delivering an equal percentage of both red light at 660NM and near infrared light at 850NM.

## **NEAR INFRARED LIGHT AT 850NM**

Near infrared light penetrates deeper within your body-enhancing muscle recovery and reducing joint pain.



## How to use:

1. Clear the treated area by removing all traces of make-up, dirt, or sweat. Then, simply wash the area being treated to ensure a clean surface.
2. For best results, place the device to target bare skin not more than 6-12 inches away from the treated area.
3. When using the panel to target your face, please do not look directly into the panel. We have provided protective goggles for you to use.



## **Timing Function:**

As you turn on your device, the screen will show OFF. Switching the device on will turn on the red light only. You can then toggle between different modes.

## **Mode Setting:**

N-: 850nm.

R:660nm

NR: Red + Near Infrared



## **Precautions, warnings, and suggestive instructions:**

- Read the manual before use.
- Users with light and heat allergies should use it under the guidance of a doctor.
- If you have discomfort after using the instrument, please see a doctor immediately.
- If your device is not used for long, attention should be paid to avoiding dust accumulation.
- Please select high-quality power socket, and the power supply voltage should be consistent with the nominal voltage of the device.
- Please use socket with ground wire hole to ensure correct and good grounding.
- When not using the device, do not plug the power plug into the power socket.

- When connecting or disconnecting the power supply, be sure to hold the plug and keep your hands dry. Do not pull the wires to prevent accidental electric shock.
- It is strictly forbidden to insert fingers or other foreign matters into the protective cover of the device when it is working, preventing electric shock.
- Please do not use it in bathrooms and damp places.
- In case of power failure, please turn off the power switch and unplug the power cord.
- This product can not be used around flammable gas.
- Attention should be paid to eye protection during infrared treatment.
- The device shouldn't be used with infants and pregnant women;

## Maintenance:

1. When cleaning and maintaining the device, you must unplug the power cord plug and cut off the power supply.
2. The machine surface cannot be wiped with organic diluent or gasoline otherwise, and the machine surface will be damaged.
3. Do not disassemble, repair or change the function of the device.
4. When cleaning the device shell, please use a wet cloth dipped in neutral lotion to wipe the outer surface (when brushing, there should not be too much moisture on the fabric).

Thank you for choosing  BIOHACKN<sup>®</sup>

[www.Biohackn.com](http://www.Biohackn.com)