




The Biohackn Light Therapy Belt

Users Manual

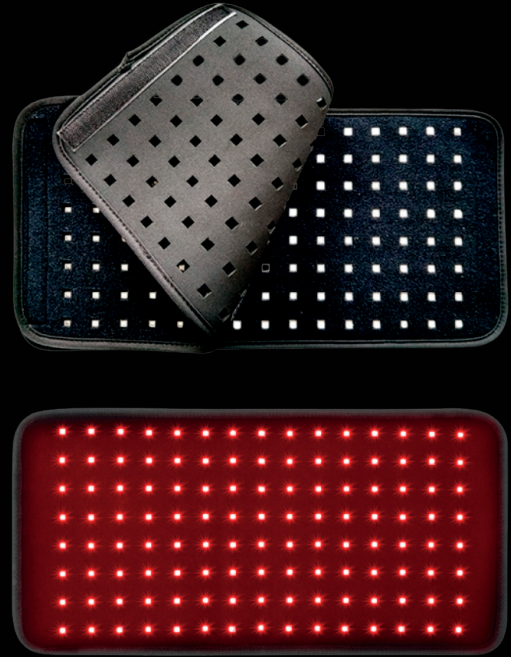

Congratulations on purchasing JUNO from Biohackn.

We are sure you are excited to use JUNO. However, we advise that you read this manual first.

| | |
|--------------|-----------------------|
| Rating Power | 180W |
| Actual Power | 48W |
| Product Size | 800*300*5MM |
| LED QTY | 360*0.5W |
| Net Weight | 1.6 KG |
| Wavelength | 1*660nm ; 2*850nm |
| EMF | 0 UT |
| Irradiance | 0 Inch: > 100 mW/c m² |
| Timer | 20 Minutes |
| Material | Fabric |
| Lifespan | Over 50,000 Hours |

Instructions:

User Instructions:

1. Plug the cable into the AC socket, and all lights should turn on. Alternatively, you can use a heavy-duty power bank (over 5V 2A) to power your device.
2. Distance: We recommend you use the belt as close to the targeted body part as possible.
3. Selection of treatment :

Press the power button to indicate green:

1. Press the button once, and All lights turn on.
2. Press twice, and only red lights (650 nm) switch on, and the near-infrared light stays off.
3. Press thrice, and only the near-infrared lights (850 nm) turn on, and the red light stays off.
4. Press 4 Times to switch both the lights off.

Pulse Recovery Mode:

Nothing in nature stays constant, even sunlight. This is why we equipped your device with a **special recovery mode**. To access this mode, press the button for more than three seconds. The button will now have a blue color.

1. Press the button once, and All lights turn on in recovery mode
2. Press twice, and only red lights (650 nm) switch to recovery mode, and the near-infrared light stays off.
3. Press thrice, and only the near-infrared lights (850 nm) turn on in recovery mode, and the red light stays off.
4. Press 4 Times to switch both the lights off.



Usage Frequency:

We suggest using your device for 20 minutes on the targetted area. After that, you can use the device up to 2 times daily. The scientific evidence and the team at Biohackn recommend that you use your device for at least 4 weeks.

After using the device, unplug it and store it in a cool, dry area.


Warning:

This device is not waterproof! Your device is not intended for use by children and pregnant women. Do not exceed 30 Min. on one same treatment.

Due to the design complexity, We highly recommend not folding the device. If the device is damaged due to folding, we will not cover it in warranty.




Folding will damage the product without warranty

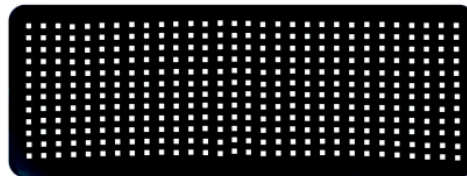

please unfold and flat when not in use

How to use:

Here is some idea on how to use your light therapy device. Get creative and share with us how you use your device.

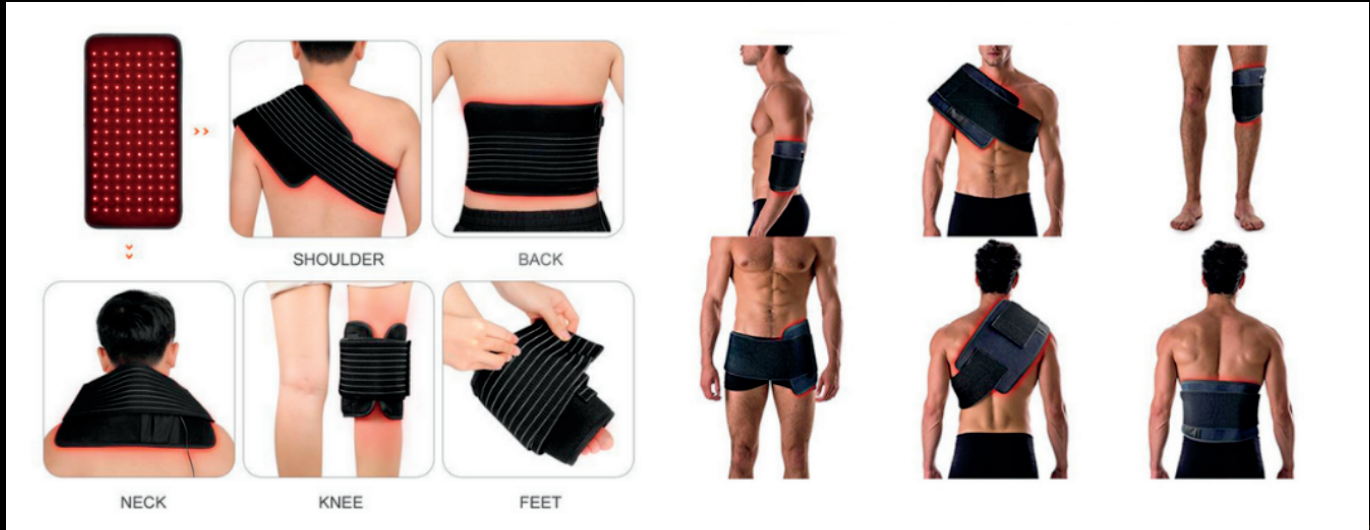

Thank you for choosing  BIOHACKN[®]

www.Biohackn.com