

USER MANUAL

Off-Check Personal Body Weighing Scale

Model : PS 141



To know about discounts & new launches connect with us on social media

@healthsenseindia



Healthsense is a proud Indian brand creating quality products idealised and conceptualised for India and manufactured by exclusive overseas partners.

www.healthsense.in










Table of Contents

Safety Information	
Safety and Usage Information	2
Your Scale and Its Environment	3
Efficient Use of Your Scale	3
Overview	
General Instructions	4
Device Components	4
Function Indication	4
List	5
LCD Display	5
Initial Start-Up	
Insert the Batteries	6
Start Measuring	
Select Measurement Unit	7
Initialising Your Scale.....	8
Start Measuring	8
Troubleshooting	
Error Prompt	9
Specifications	10
Warranty & Support	11

Safety and Usage Information

The warning signs and symbols are essential to ensure your correct and safe use of this product and protect you and others from injury. Please kindly find the meanings of the warning signs and symbols, which you may encounter in the label and user manual, as follows:

	Symbol for "THE OPERATION GUIDE MUST BE READ"		Symbol for "MANUFACTURER"
	Symbol for "DIRECT CURRENT"	 	Symbol for "ENVIRONMENT PROTECTION – Waste electrical products should not be disposed of with household waste. Please recycle where facilities exist. Check with your local authority or retailer for recycling advice"
	Symbol for "COMPLIES WITH EU REQUIREMENTS"		
	Symbol for "SERIAL NUMBER"		

Your Scale and Its Environment

To ensure your safety as well as the service life of your scale, please avoid using the scale under the following circumstances:

- Slippery floor such as tile floor
- Jumping onto the platform immediately after bath or with wet feet
- Near a cell phone or microwave oven

Avoid storage in the following locations:

- Where there is water
- Where the device may be exposed to extreme temperatures, humidity, moisture, direct sunlight, dust, or salt air
- Where there is risk of shock or drop
- Where you store chemicals or full of corrosive gases
- Where in reach of the infants or children

Efficient Use of Your Scale



To ensure the accuracy of measurement, please follow below instructions when you start measurement.

- Place the scale on a flat, hard surface. Soft surface such as carpet will affect the performance of the scale.
- Start measurement at least two hours after Getting up or Dining.
- Avoid measurement immediately after strenuous exercise, sauna or bath, drinking, and dining.
- Always start measurement in the same time slot and on the same scale located on the same flat, hard surface.
- In order to ensure the accuracy of weighing, when you moved the scale or haven't used the scale for a long time, please turn on the scale first, wait until it turns off, then you can stand on the platform again and measure your weight.
- When you don't use the scale, please place it stable. Do not put anything on it.
- Any information provided by this scale is in no way meant to treat, cure or prevent any disease or illness from happening. If in doubt, contact your physician.
- It is intended for use in the domestic setting only.

General Instructions

Healthsense Body Scale PS 141 applies Off-Check technology. It's not necessary to turn on or turn off the scale, just stand on the scale for measurement. It is delicate and trendy, convenience for using.

Device Components

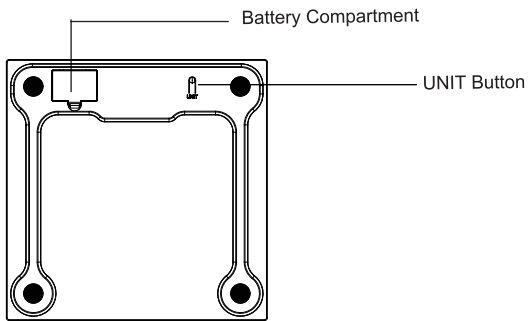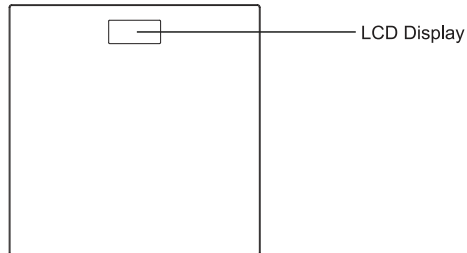


Function Indication

- ◇ UNIT Button: Press "UNIT Button" to select desired unit in weighing mode.
or
Press "UNIT Button" to turn on the scale.

In the Box

1. Body Scale PS 141
2. 3 x AAA batteries (4.5V)
3. User Manual

LCD Display

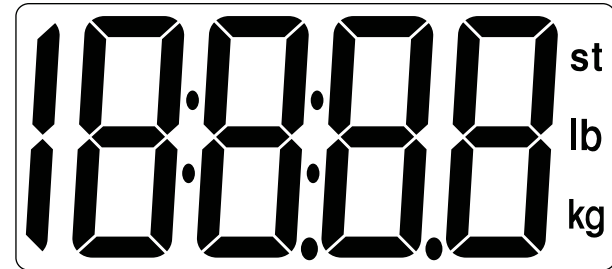


Symbol	Description
st lb	Stone
lb	Pound
kg	Kilogram

✦ Insert the Batteries

- Open the battery door in the back of the scale.
- Remove insulating strip from beneath the battery (if fitted) or insert battery (3 x AAA) observing the battery polarity signs (+ and -) inside the battery compartment.

* The digits "8888" will be shown on the LCD.

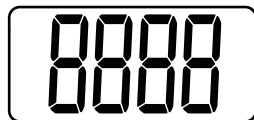

- Close the battery door and wait until the digits "0.0kg" are shown on the LCD.




⚠ CAUTION

- When the symbol "Lo" appears, the device will power off in four seconds. Then you shall replace with a new set of batteries.
- Worn batteries are hazardous waste. Do NOT dispose of them together with the household garbage. Please refer to the local ordinances and recycling instructions regarding disposal of the worn batteries.
- In order to achieve the best performance and protect your scale, please use the authorized/specified batteries.
- Risk of explosion if the battery is replaced by an incorrect type
- replacement of a battery with an incorrect type that can defeat a safeguard (for example, in the case of some lithium battery types);
- disposal of a battery into fire or a hot oven, or mechanically crushing or cutting of a battery, that can result in an explosion;
- leaving a battery in an extremely high temperature surrounding environment that can result in an explosion or the leakage of flammable liquid or gas;
- a battery subjected to extremely low air pressure that may result in an explosion or the leakage of flammable liquid or gas.

✦ Select Measurement Unit

With batteries correctly installed, press "UNIT" button in the back of the scale to select measurement unit. The default measurement unit is "Kilogram". You may press "UNIT" button to choose among kilogram, stone and pound.

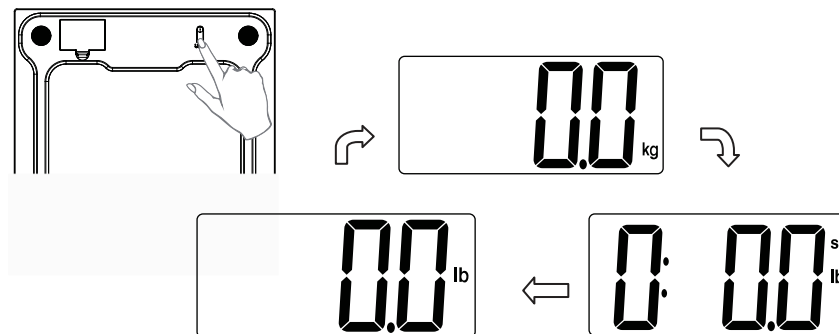

⚠ NOTE

Please put the scale on the flat hard ground when finishing the unit selection. Wait for a few seconds, then you may start measuring according to the measuring procedure below. If there is not any operation in weighing mode, the scale will turn off automatically after about 10s.

ABOUT THE ACCURACY OF THIS PRODUCT

- This product passes strict inspection before delivery and therefore its accuracy is guaranteed by the manufacturer.
- This product is specially designed for body scale as well as weight measurement. It should NOT be used by anyone during the process of transaction for verification of goods' weight.

1. When it shows 0, step on the scale.



2. Screen displaying the dynamic moving lines “----” means scale is weighing your reading.



3. When it shows “step off”, step off the scale, and the final weight will be shown.



Why is it more accurate?

Most scales fetch the zero point first before weighing and it would be used for final weight calculation.

It works fine in most conditions except that the scale fetches the zero point while it is hold in hand, and the user will get an inaccurate weight after he/she put down the scale and step on it. In this case, the shown weight will be the user's actual weight, plus the scale's weight itself.

The Off Check function reverses the weighing procedure, which is, step on and weigh first, then it fetches the zero point when the user steps off the scale, which guarantees a correct zero point on weight calculation.

Error Prompt

Error	Description	Solution
	Overload. The device will power off.	Stop using this scale for measurement.
	Low Battery. The device will power off.	Replace the battery. Please purchase the authorized batteries for replacement.

Maintenance

When carrying out usual maintenance, please ensure practice of the following Do's and Don'ts:

- DO use a dry soft cloth to wipe the dust.
- DO use a wet soft cloth, dipped into water and wrung out, to wipe the dirt. Then use a dry soft cloth to dry up the device.
- DON'T wash the device with water or immerse it in water.
- DON'T use propellant, abrasive or other chemicals to wipe the dirt in avoidance of discolor or malfunction.
- DON'T disassemble this device. If you have any problems, please contact Healthsense. *(Please refer to [Warranty](#) for contact information)*

Specifications

Product Name	Off-Check Personal Body Weighing Scale
Dimension	Scale: Approx. 300x300x22.2mm
Net Weight	Approximately 1.26kg (Excluding Batteries)
Display	Digital LCD(with backlight)
Measurement Unit	Kilogram / Stone/ Pound
Measurement Range	180kg / 397lb / 28st: 5lb
Division	0.1kg / 0.2 lb
Accuracy	50kg: ±0.3kg; 100kg: ±0.4kg; 150kg: ±0.5kg; 180kg: ±0.7kg;
Working Environment	Temperature: 5°C to 40°C Relative Humidity: ≤85% RH Atmospheric pressure: 86kPa to 106kPa
Storage Environment	Temperature: -20°C to 60°C Relative Humidity: ≤90% RH Atmospheric pressure: 50kPa to 106kPa
Power Source	4.5V (3 x AAA Battery)
Auto-ON	SENSE ON technology
Auto-OFF	About 10 seconds while showing 0.0; About 10 seconds after the weight data is locked; About 10 seconds when the weight is unstable.
Accessories	1. Three AAA Battery 2. User Manual
Mode of Operation	Continuous Operation

WARRANTY & SUPPORT

Entitle yourself for the 1 year HealthSense warranty by registering your purchase and easily reach us to submit your product related issues & queries, when you follow these hassle free ways.

SCAN - For Warranty Registration



(Must register within 30 days from date of purchase)

SCAN - To File Complaint



(Fill & submit form if you are facing any product related issues)

*** For detailed warranty terms please visit the link :**

www.healthsense.in/warranty-terms

*** For additional support :**

Mail us at : support@healthsense.in or Call us on : +91 9856 9856 36