



www.healthsense.in




Healthsense is a proud Indian brand creating quality products idealised and conceptualised for India and manufactured by exclusive overseas partners.

HEALTHsense®

USER MANUAL

GLASS-TOP PERSONAL SCALE

Model : PS 117

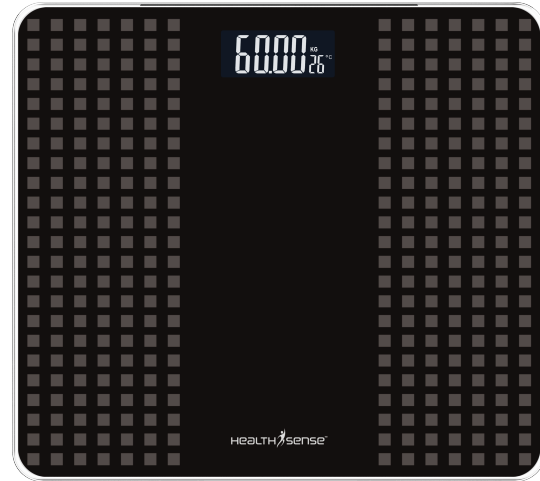

Sensing Health Everywhere

INTRODUCTION

Thank you for purchasing HealthSense Glass-Top Personal Scale. This scale has been specifically designed for your use, with an easy to use interface and superior HealthSense accuracy and design.

Your HealthSense Personal Scale comes with the following:

- 1 x Glass Personal Scale PS 117
- 1 x Instruction Sheet
- 2 x AAA batteries
- 1 x Warranty Card

Kindly read the instruction manual carefully before using your scale.

Please fill out your purchase information for future reference:

Date Purchased: _____

Place of Purchase: _____

Order ID (if any): _____

OPERATION

Battery Installation and First Time Setup

1. Locate the battery cover on the backside of the unit. Press the battery cover tabs inward and lift to remove.
2. This scale operates on 2xAAA batteries. Insert batteries or remove the insulating strip from the battery contacts before using.
3. Insert batteries as indicated by the polarity symbols (+ and -) as marked. Close the covers.
4. The scale needs to be initialized and auto-calibrated before first use or after any time you move your scale. Place the scale on a hard, flat surface. If the scale turns on while being moved, allow it time to turn off once on the ground.
5. Get on the scale and weigh yourself (please keep in mind the first reading may not be accurate and can be disregarded).
6. Get off the scale and allow it to turn off.
7. You can now step on and weigh yourself normally. If you move the scale again, you must repeat this calibration process. You will see your final weight when the display locks in on a weight and blinks two times.
8. Replace batteries when scale indicates "LO" on the display.

NOTES

1. New batteries should last for approximately one year under moderate use (less than three times daily). Batteries included with your scale are trial batteries and may last a shorter length of time.
2. Dispose of batteries within the legal restrictions of local government.

Weighing Operation

Your HealthSense Personal Scale PS 117 operates as a basic body weighing scale. No special programming is required.

1. **Select your desired measurement unit.** Your scale is set to default to kilogram at the factory, but you can change the measurement units to pound or stone by clicking the unit button inside battery chamber or at backside of your scale.
2. **Place the scale on a hard, flat surface.** Carpeted or uneven floors can negatively affect the accuracy of your readings. Optimal flooring surfaces include hardwood, ceramic tile or concrete.
3. **Initialize the scale.** If you are placing the scale down for the first time or moving your unit you must auto-calibrate the scale (see steps 4 - 7 under the "Operation" tab of this manual).
4. **Step on.** Place both feet on the center of the platform and stand still.
5. **Your weight will appear on the LCD.** Your weight will blink two times when your final weight has been calculated.
6. **Step off your scale.** It will reset to "0.0" and turn off automatically.
7. **Move your scale as infrequently as possible.** To ensure maximum accuracy and consistency digital scales generally should be moved as little as possible.

IMPORTANT: Your HealthSense scale operates with "Step-On" technology but still needs to be auto-calibrated after battery installation and when the scale is moved.

Place the scale on a hard, flat surface and weigh yourself. Ignore the first reading and get off the scale. Then without moving the scale weigh yourself again. This is your true weight.

SPECIFICATIONS

Model Name: Glass-Top (PS 117)

Capacity: 5-150kg / 11-330 lb / 0.78-22.6 st

Measurement Units: Kilograms / Pound / Stone

Graduation: 50gm / 0.2 lb / 0.01 st

Display: LCD size of 62 x 25.5 mm with Reverse LCD screen

Low Battery Indicator "LO"

Overload Indicator "Err"

Unit Dimensions: Approximately 280 x 243 x 22 mm

Unit Weight: Approximately NW: 980 g, GW: 1.2kg

Power: 2 x AAA Batteries (included)

THE DISPLAY OF AMBIENT TEMPERATURE

This product can detect and display ambient temperature automatically if starting up. And it would constantly refresh the latest ambient temperature until power OFF.

If you change service environment of the product, it is recommended to check the temperature display in 30 minutes later to ensure its accuracy.

HOW TO RECALIBRATE YOUR SCALE

Step on the scale with one foot until the digits appear on the display, then, step off. The scale will show "0.00", indicating that the calibration process is complete.

TROUBLESHOOTING

PROBLEM	SOLUTION
Scale is inaccurate or weight varies	Place feet in center of platform and stand still. Allow scale to shut off completely after weighing. Do not move scale unnecessarily.
First weight is higher than normal	Your scale must be allowed to auto-calibrate. Follow instructions under "Operation" for proper use.
Scale displays "ERR"	Place scale on hard, flat surface. Allow scale to shut off and restart.
Low battery indicated "LO" No power	Check batteries for proper installation (Polarities are correct). Replace all batteries.
Overload indicated "Err"	Too much weight has been placed on the scale and it is overcapacity. Remove weight before proceeding.
Display "C" when you step off the scale	It means this weight reading is with error. Please RE-TEST
Other conditions	Email HealthSense team for support at support@healthsense.in

CARE & MAINTENANCE

To maintain your HealthSense Personal Scale in the best possible condition follow these directions as closely as possible.

- 1. Store and operate the unit in safe dry location.** Drastic changes in temperature over time can affect the quality of weight measurement. Do not expose to direct sunlight, extreme hot/cold or humidity.
- 2. Clean with a soft, damp cloth.** Avoid abrasive cleaners and never sink or submerge in water.
- 3. Remove batteries if not in use for three months or longer.**
- 4. Use unit consistent with instructions provided in this manual.**
5. Do not use it for commercial purpose.

WARRANTY & SUPPORT

Entitle yourself for the 1 year HealthSense warranty by registering your purchase and easily reach us to submit your product related issues & queries, when you follow these hassle free ways.

At HealthSense our top priority has always been our valued customers. We very well know your HealthSense product is a part and parcel of your everyday life, which is why we are more easily accessible than earlier.

SCAN - For Warranty Registration



(Must register within 30 days from date of purchase)

SCAN - To File Complaint



(Fill & submit form if you are facing any product related issues)

* For detailed warranty terms please visit the link : www.healthsense.in/warranty-terms

* For additional support :

Mail us at : support@healthsense.in or Call us on : 080-41262836