When finished

all the mode lights

charging.

will be on

On really

low battery...

fade on and off

and the motor will

During charging...

the mode lights will

turn on and light up

from the bottom and

go up in sequence. The

base of the toothbrush


handle may get warm while charging.

the battery indicator light will

shake twice.





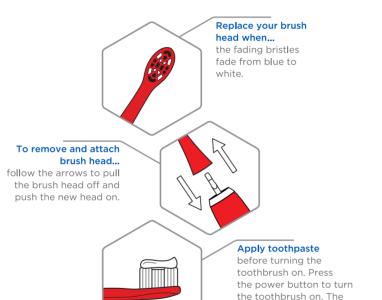

brush head vibrations will

intensify as it reaches the

full mode speed.

SPM-I Speaking of replacing brush heads... order new brush heads, toothpaste, and other oral care accessories at getevo.net RECHARGEABLE SONIC TOOTHBRUSH

=

take care of your oral health needs with our subscription service that delivers new brush head refills to your door quarterly. Changing your toothbrush head every 3 months and/or after being sick evo. helps fight nasty bacterial build up.

Access your specialty SPM-1 enamel protector brush heads at getevo.net/subscription

In case of any issues that you may have...

Let us help you...

*See Evo-Instruction Manual for warranty details.

our customer support team would love to help. Reach out by emailing support@getevo.net.

Life happens. Your toothbrush is covered with a 1 year warranty*. Register for a subscription plan and receive a lifetime warranty*.





RECHARGEABLE SONIC TOOTHBRUSH USER GUIDE - STANDARD EDITION

CVO_® stands for evolution.

Evo was born from the dream of creating futuristic products today. Conceptual products, such as concept cars, are inspiring, amazing to look at and usually something we are unable to purchase for ourselves. So we simply asked "Why?"

> This was the start of our journey. Working with some of the most creative and intelligent designers and engineers around the world, we examine every aspect of our products to make them a reality. We are dreamers, creators, and designers infusing innovative technology with great design to produce "dream products" that one may have never imagined possible.

Introducing EVO-SENSE Technology

Change the way you control your toothbrush by using EVO-SENSE. This feature allows you to control down with your finger.

Why Marvel?

technology that could be possible in our world. Looking to the future, we partnered with Marvel to create SPM-I as the next part of our collaboration.

INTRODUCING... —

SPM-I

SPM-I Enamel Protector Brush Head





Power/Mode Button ——

EVO-SENSE ----

Mode Lights ----

Standard

Sensitive

Whitening

Battery Indicator Light — -

Deep Clean

 \forall

>


V

Intensity Level Indicator ----

the intensity of the vibration by simply swiping up or

Marvel has always inspired us with the thought of futuristic

EVO-SENSE TECHNOLOGY ——

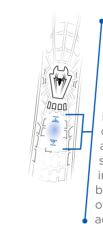

SPM-I Standard

Charging Base

with USB

USB Adapter Plug

STANDARD

EDITION

power cord

• Let us introduce to you our latest feature: EVO-SENSE

EVO-SENSE allows you to control the intensity level of the vibration by swiping up and/or down with your finger.

How does EVO-SENSE work?

When selecting a mode, you have 8 seconds to select the preferred intensity level for the current mode. Swiping up increases intensity. Swiping down decreases intensity. We created a locking system on EVO-SENSE, that automatically activates 8 seconds after selecting a mode. 8 seconds after selecting a mode, EVO-SENSE will lock, Reset if desired intensity was not chosen. To reset: turn the brush off then back on. You will have 8 seconds to select the desired level of intensity. We created this locking mechanism to avoid accidentally changing the intensity while brushing.

Change the level of Intensity

by swiping up or down on the **EVO-SENSE** area with your finger.




Automatic lock intensity while brushing.

EVO-SENSE will lock in the intensity after 8 seconds, to prevent accidentally changing the

GETTING STARTED ————

brush head while brushing.

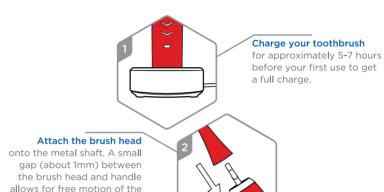

Apply toothpaste

before turning the toothbrush on, Press the power button to turn the toothbrush on. The brush head vibrations will intensify as it reaches the full mode speed.

Select your mode



by swiping up to increase intensity and swiping down to decrease intensity. You have 8 seconds to select the intensity level before it locks

in place.

Use EVO-SENSE

Switch quadrants

every 30 seconds when the interval timer pauses. The toothbrush will automatically turn off after 2 minutes. Rinse and wipe dry.

bristles to stay on the surface



vour teeth

by moving

gently along

and gum line,

allowing the

vour teeth

the brush head