

ABOUT EVO®

EVO. stands for evolution.

Evo® was born from the dream of creating futuristic products today. Concept products and concept cars are amazing to look at but almost never something we can purchase for ourselves. So we asked why?


This was the start of our challenge, thinking up creative ways to turn a dream into a reality. Working with some of the best designers and smartest engineers around the world, we examine every aspect of our products to make them a reality. We are dreamers, creators, and designers currently working on connected technology projects that you never imagined would be possible.

Why Marvel?




Marvel has always inspired us with the thought of futuristic technology that could be possible in our world. Looking to the future, we partnered with Marvel to create IRM-I as the first of our collaboration.

INTRODUCING...



IRM-I

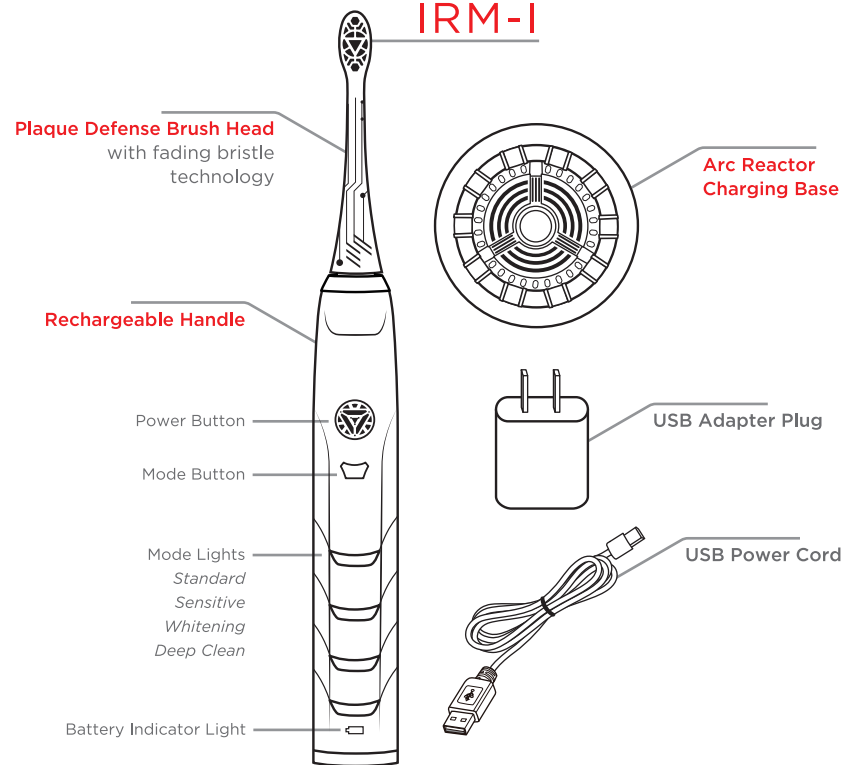

GETTING STARTED

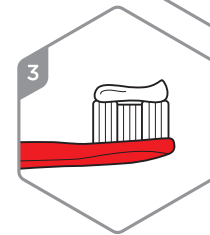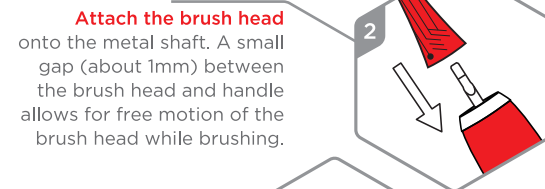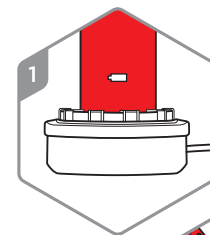

HOW TO USE

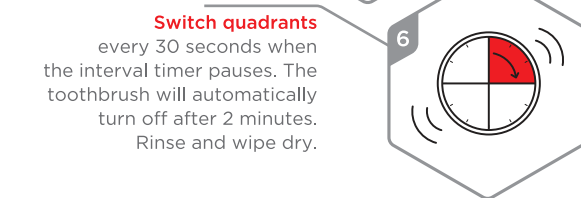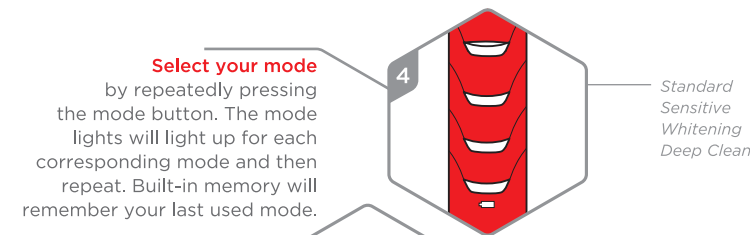

BATTERY & CHARGING

On low battery...
the battery indicator light will turn on.



On really low battery...
the battery indicator light will blink for 8 seconds and the toothbrush will not start.



During charging...
the battery indicator light will turn on and the charging base light will have a pulsating effect (light fades on and off).

The toothbrush may get warm during charging.


When finished charging...
the battery indicator light and charging base light will turn off.

REPLACING BRUSH HEADS

When it's time to replace your brush head...
the power button light will turn red.



Replace your brush head when...
the fading bristles fade from blue to white.

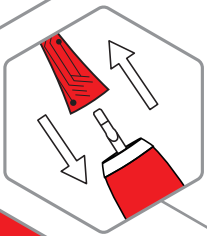


To remove and attach a brush head...
follow the arrows to pull the brush head off and push the new head on.

Once your new brush head is attached...
press the power button for 5 seconds and the light will reset to white.




STAY CONNECTED

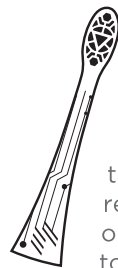

Speaking of replacing brush heads...
pick up brush heads, toothpaste, and other oral care accessories at getevo.net

Let us help you...

take care of your oral health with brush head refills, delivered to your door quarterly with our subscription service. Changing your toothbrush head every 3 months and after being sick helps fight bacterial build up.

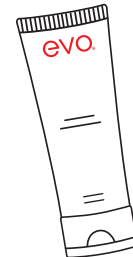

Access your specialty plaque defense brush heads at getevo.net/subscription

In case of issues...

our customer support team would love to help. Reach out by emailing support@getevo.net.

Life happens. Your toothbrush is covered with a 1 year warranty*. Register for a subscription plan and receive lifetime warranty.

*See Evo® Instruction Manual for warranty details.

evo®



RECHARGEABLE SONIC TOOTHBRUSH
USER GUIDE