




Welcome to izzo.™

It's time to get your smile on.

Congratulations!

You're about to change your smile forever with this unique 4-in-1 oral care system. We know that sounds technical, and that's because it is. We are a 100+ year-old dental manufacturing company that knows how to create great smiles and great oral hygiene.

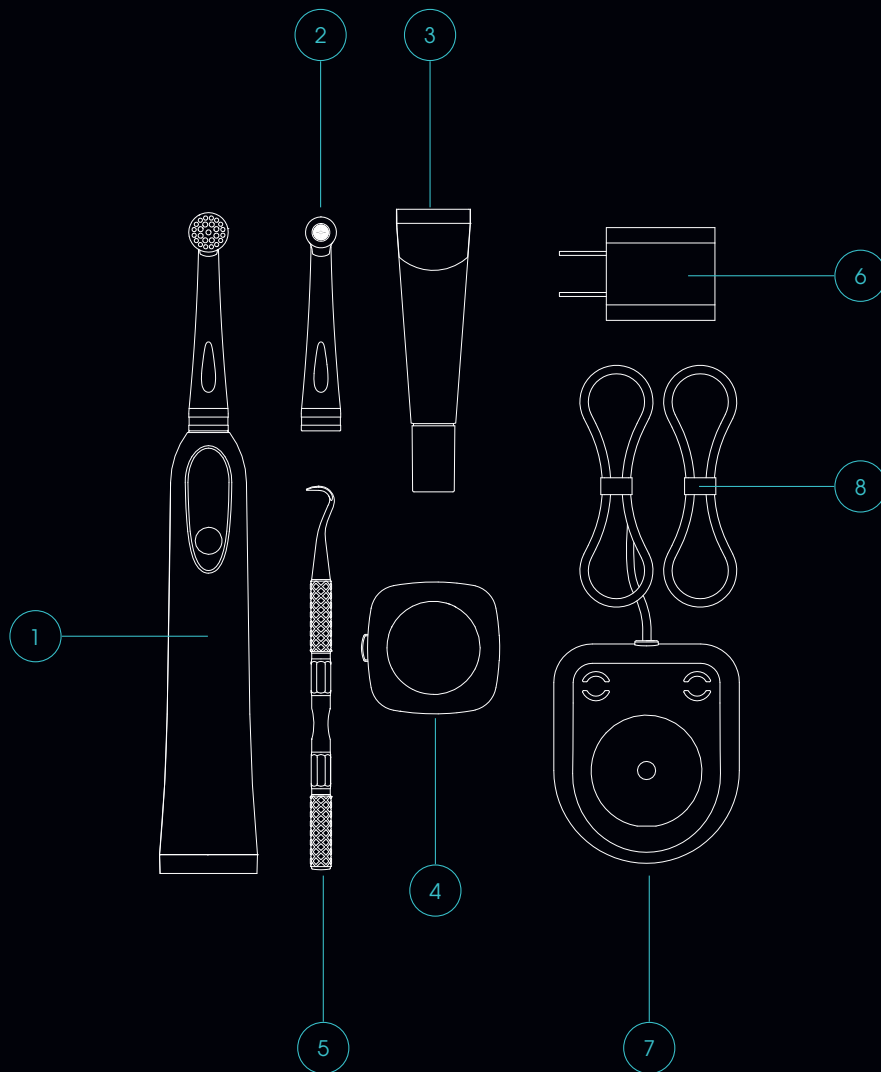

What izzo Is. What izzo Does.

izzo is professionally inspired by top dentists and hygienists, and everything about it is made for you. To get the most out of your next-level izzo features, read through this simple guide and you'll be on your way to a complete clean and a happy smile.

That's daily dentistry delivered!

1 Power Handle with Brush Head

- High-efficiency oscillating Brush Head
- Rechargeable battery for 2 weeks of use
- 3 speed settings
- "Soft Start" feature reduces splatter by steadily increasing to desired speed over 2 seconds.
- Interchangeable Brush and Polishing Heads
- 2-minute timer

2 Polishing Cup Head

- Fits Power Handle
- Provides safe and effective dental polishing for at-home use
- Soft articulating cup adapts to the shape of your teeth

3 Enamel Polishing Paste

- Professional Enamel Polishing Paste optimized for home use
- Use exclusively with the Polishing Cup Head to help remove surface stains for whiter, brighter teeth beyond brushing alone

4 UVC Sanitizer

- Sanitizes your Brush Head and Polishing Cup Head with ultraviolet light in 3 minutes
- Built-in rechargeable battery holds charge for 2 weeks

5 Scaler

- Professionally inspired dental instrument
- Helps clean debris and plaque from hard-to-reach areas in your mouth
- High performance medical-grade polymer won't scratch your teeth

6 USB Wall Charger

- Provides 2 USB ports
- Use to charge the Power Handle and UVC Sanitizer

7 Charging Base

- Inductive Charging Base for the Power Handle
- Two posts provide a place to store either the Brush Head, Polishing Cup Head or Scaler


8 Cable

- Cable to connect the UVC Sanitizer to the USB Wall Charger



Get Connected

It's hard to take charge of your oral care if you don't charge your izzo. We recommend a full 12-hour charge before your first use.

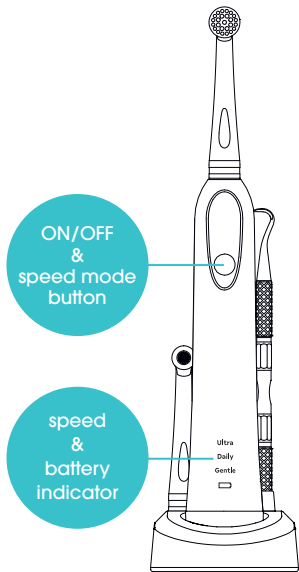

Power Handle Charging Instructions

- Plug the USB Wall Charger into a GFCI-protected 120-volt wall outlet. Then, connect the USB cable from the Charging Base into the USB Wall Charger.
- Place Power Handle into the Charging Base. The battery charge indicator will slowly pulse while charging.
- Charging will stop and the battery indicator light will turn off automatically once the battery is fully charged.
- When the Power Handle is not in the Charging Base and the battery is low, the battery indicator will begin to blink.

Notes


It's normal for the Power Handle to become warm while charging.

It's okay to keep your izzo in the charging stand when not in use—you won't harm the battery by overcharging.

A 12-hour charge will be needed if you haven't used your izzo for 3 months.

Is there a sequence I should follow when using my izzo system?

It depends on your personal preference and individual needs. We recommend: Floss, Scale, Brush, Polish, Sanitize to keep you smiling bright. The important thing is to use all components as directed. "Rinse and repeat as needed!"

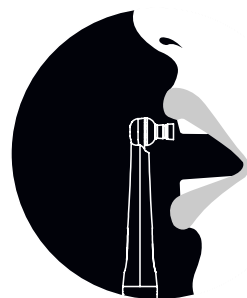

Power Handle Speed Modes

Ultra (Highest):


For polishing or when you need more brushing power

Daily (Medium):


For daily brushing

Gentle (Lowest):


For gentler cleaning or sensitive gums

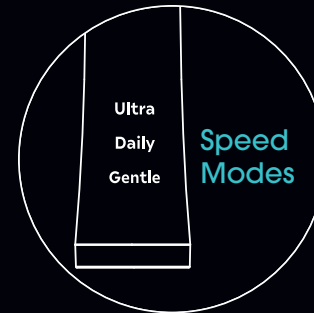


- To cycle between modes, repeatedly press the ON/OFF/Speed Mode button within the first 5 seconds of turning it on to select your desired speed. After 5 seconds, pressing that same button will turn the Power Handle off.
- To help stick to your new routine, the Power Handle restarts in the same mode from its prior use.

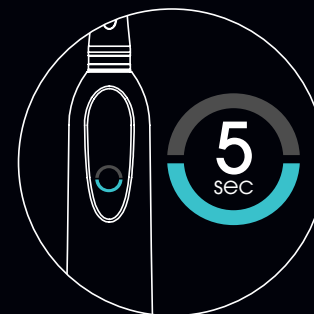

Get Brushing

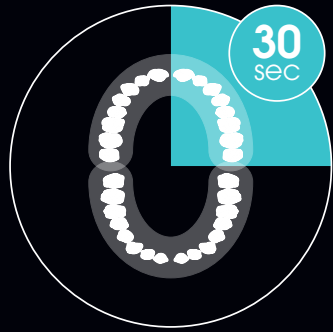


- First, rinse your mouth, wet the Brush, and add your favorite toothpaste.
- We recommend putting the Brush Head in your mouth before starting it up. Why? You don't want a mess!
- Move the Brush Head in a circular motion around your mouth. Be sure to cover the outsides, the insides, and the chewing surfaces. The bristles should just lightly touch your teeth and gums.
- Brush for 2 minutes. It's easy with izzo's built-in timer, which pulses every 30 seconds to remind you to switch from one quadrant to the next. Plus, the Brush will automatically turn off after 2 minutes.
- To clean your Brush Head, turn the Power Handle on, point the Brush Head downward 45°, and rinse under a faucet. To clean the inside of the Brush Head, run water over the hole on the back for a few seconds. After you remove your izzo from under the water, keep your Brush running for a few more seconds to expel any excess water. Wipe down as needed and sanitize in your UVC Case before your next use (see page 10 for details).

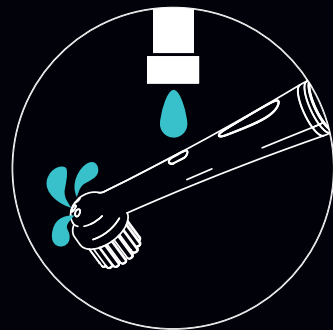


Changing the Brush Head / Polishing Cup Head

It's all about the **TWIST!**

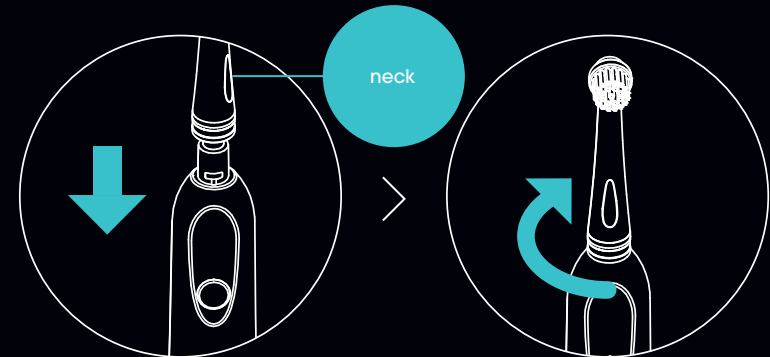

To attach your Brush or Polishing Head:

- Step 1** Place the neck of the head on the Power Handle tip
- Step 2** Twist clockwise to lock

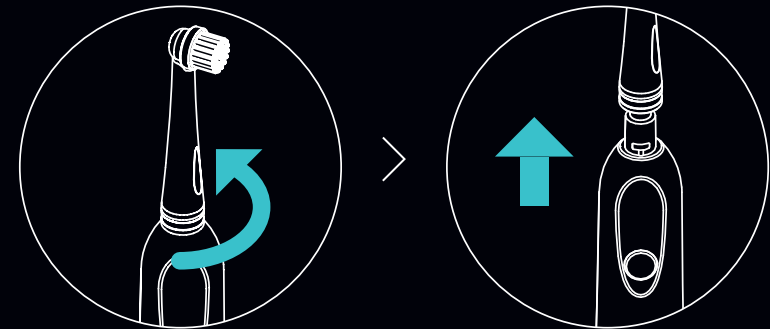

To remove your Brush or Polishing Head:


- Step 1** Hold the neck of the head (not the bristles)
- Step 2** Twist counterclockwise to unlock
- Step 3** Pull the head off

Note

Replacing your Brush Head and Polishing Cup Head every 3 months is clinically proven to remove more plaque and will give you a top-notch cleaning.

Get Polishing

It's time to take your oral care to the next level. Hello whiter brighter teeth. Goodbye surface stains!

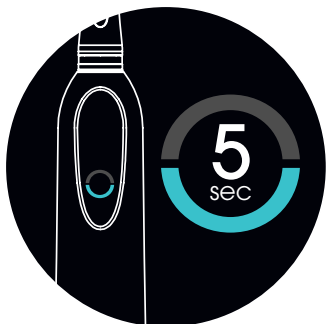


Step 1:

First, set your Power Handle mode on Ultra for the best results. Run it for 5 seconds to lock in the speed before turning off.

Fill the Polishing Cup Head with a pea-sized amount of Enamel Polishing Paste. Be careful not to overfill; a little dab will do it.

Step 2:

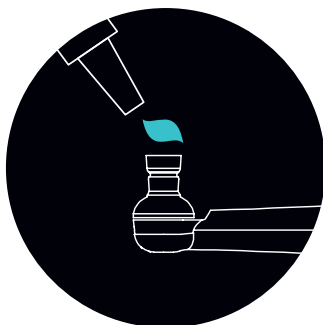

Step 3:

Position the cup at a 90 degree angle to the tooth surface and turn on. Focusing on one tooth at a time, move the cup over the entire tooth surface. For the best results, spend only 3-5 seconds per tooth.


Put another dab of Enamel Polishing Paste in the Cup Head when it runs out. You may need to do this every 3-5 teeth.

Notes

We recommend polishing once per week.

Focus on the front surfaces of teeth that are visible when you smile. To get technical, that means the vertical front (facial), back (lingual), and horizontal biting surfaces of your teeth. When you become a "pro," you can try polishing other areas as well.


For more information about polishing go to izzosmile.com

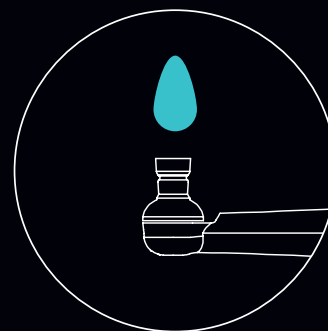

Step 4:

Now that you have a clean mouth, let's make sure you clean and care for your Polishing Cup Head just as well.

Follow Brush Head cleaning instructions as outlined on page 6. Don't forget to use your UVC Sanitizer!

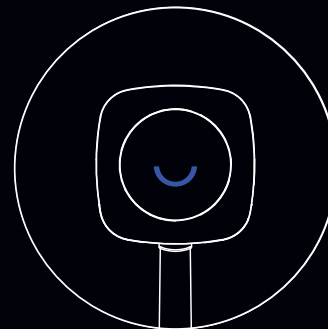

Get Germ-Free Clean

It's time to take cleaning to the next level with our UVC sanitizing case that kills 99.9% of germs. In just 3 minutes, your Brush Head and Polishing Cup Head will be cleaner than clean.



Sanitizing 101


Like most tech, the Sanitizer doesn't love water. Be sure to remove excess water from the Brush Head or Polishing Cup Head before it goes in the Sanitizer.



indicator light

Note

It's ok to leave the UVC sanitizer case on the brush or polishing head while in the stand in between uses.

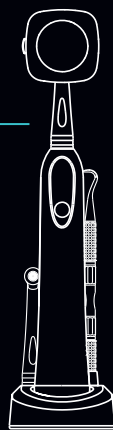

Press the latch button to open the case and place the Brush Head or Polishing Head facing the source of light. The 3-minute sanitization will begin automatically when you close the case—just listen for a click! The blue "smile" lets you know its working.

Note


Opening the case during its cycle will stop the sanitization process.

You can even store your Brush Head or Polishing Cup Head in the UVC Case in between uses. It makes a great travel case too!

When you see the indicator light flashing red, it's time to recharge via the USB-C cable. The case's indicator light will turn green when fully charged.

Note

The battery charge lasts about 2 weeks and takes about 2 hours to charge.

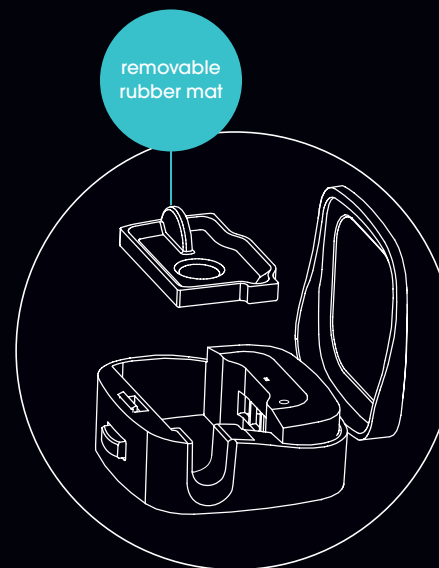


removable rubber mat

Yes, even Sanitizers need cleaning. To do so, remove and wash the silicone mat with water and then towel dry. For the interior and exterior case, all you'll need is a gentle wipe with a damp towel.

Get Scaling

Ready to learn how to scale your teeth? YES! Good, because izzo helps provide a deeper and more satisfying clean.

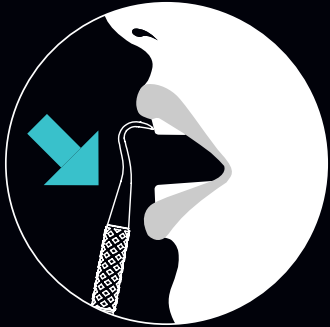


Scaling 101

Before you begin, remember the scaler should only be used to remove visible food debris and/or plaque from between the teeth. The flat side of the scaler can be used on the front and backside of teeth and along the visible border of the gums.



Note

Apply gentle pressure and be careful not to insert the pointed tip below the gumline or around loose fillings or crowns.

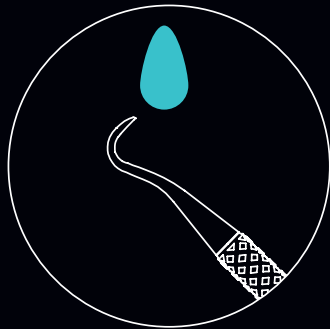

When it comes to cleaning, you've got options. Here are 3 easy ways to keep your Scaler on point:

- Rinse with warm water and wipe dry.
- Wipe with alcohol.
- Wash in a dishwasher.

A Clean izzo Is the Best izzo

We've already addressed routine cleaning, but here are a few other things to keep your izzo in tip-top shape!

Weekly: Detach the Brush Head or Polishing Head from the Power Handle at least once a week to clean the Power Handle thoroughly. Point the Power Handle downward and rinse the Head connection area under running water while the unit is on. Remove from water, and keep it running (in a downward position) for a few seconds to allow the "pumping action" to push out any fluid.

Note

Do not submerge the Power Handle in water. Allow it to dry completely before charging.

As needed: Clean the Charging Base with a damp soft cloth.

Note

Unplug the Charging Base before cleaning and allow it to dry completely before reconnecting to the USB Wall Charger.

The Planet Deserves to Be Clean Too!

Should you at some point want to say goodbye to izzo, please recycle responsibly.



- The Power Handle contains a nickel-metal hydride battery and the UVC Sanitizer contains a lithium-ion battery. Neither belong in your household waste.
- izzo must be recycled or disposed of properly, as per local regulations. Contact www.call2recycle.org for local recycling locations for rechargeable batteries.

Warranty

2-Year Warranty Included OR Upgrade to a Lifetime Warranty with Subscribe & Save.

To make sure you don't miss a beat (or a brush!), izzo includes a full 2-year "bumper-to-bumper" warranty (excluding refill components). Or upgrade to a Lifetime Warranty with an active Subscribe & Save plan (once you've completed 3 consecutive refill cycles), and enjoy savings plus the convenience of having your refills sent right to your door. It's easy and you can cancel anytime. It's the biggest no-brainer in the history of dentistry!

Have questions? Or is there a specific issue with your izzo? Check our FAQs at izzosmile.com or contact us through our online help form, chat, at customercare@izzosmile.com, or call Customer Service @ **1-888-SMILEON**.

Important Safeguards

READ ALL INSTRUCTIONS BEFORE USING

When using electrical products, especially when children are present, basic safety precautions should always be followed, including the following:

DANGER – To reduce the risk of electrocution:

1. Always unplug this product immediately after using.
2. Do not use while bathing.
3. Do not place or store product where it can fall or be pulled into a tub or sink.
4. Do not place in or drop into water or other liquid.
5. Do not reach for a product that has fallen into water. Unplug immediately.
6. Do not incinerate this appliance even if it is severely damaged. The batteries can explode in a fire.

WARNING – To reduce the risk of burns, electrocution, fire, or injury to persons:

1. A product should never be left unattended when plugged in.
2. This appliance is not intended for use by persons (including children) with reduced physical, sensory or mental capabilities, or lack of experience and knowledge, unless they have been given supervision or instruction concerning use of the appliance by a person responsible for their safety. Children should be supervised to ensure that they do not play with the appliance. Close supervision is necessary when this product is used by, on, or near children or invalids.
3. Use this product only for its intended use as described in this manual. Do not use attachments not recommended by the manufacturer.
4. Never operate this product if it has a damaged cord or plug, if it is not working properly, if it has been dropped or damaged, or dropped into water. Return the product to a service center for examination and repair.
5. Keep the cord away from heated surfaces..
6. Never block the air openings of the product or place it on a soft surface, such as a bed or couch, where the air openings may be blocked. Keep the air openings free of lint, hair, and the like.
7. Never use while sleeping or drowsy.
8. Never drop or insert any object into any opening or hose.
9. Do not use outdoors or operate where aerosol (spray) products are being used or where oxygen is being administered.
10. Connect this product to a properly grounded outlet only. See Grounding Instructions.
11. Do not charge or store the appliance or battery outside of the temperature range 32-104°F (0-40°C). Charging improperly or at temperatures outside of the specified range may damage the battery and increase the risk of fire. The recommended ambient temperature range for charging system during charging is 32-104°F (0-40°C).
12. The USB can only connect to Listed/Certified Class 2 Power Unit for charging, don't use unapproved sources.

Power Handle and Heads

- Do not use Brush Heads or Polishing Cup Heads that are not designed for the izzo system. They can damage the izzo Power Handle.
- Do not use a Brush Head or Polishing Cup Head that is damaged or worn.
- The Power Handle contains a rechargeable nickel-metal hydride battery. Dispose of the product according to your local recycling regulations.

Charging Base and USB Wall Charger

- Do not allow water, cleaners or other liquids to enter the plug of the Charging Base or the socket of the USB Wall Charger.

UVC Sanitizer

- Store away from combustible materials.
- The UVC Sanitizer contains lithium-ion rechargeable battery. Dispose of the product according to your local recycling regulations.

Enamel Polishing Paste

- To avoid dental abrasion, use the Enamel Polishing Paste only as directed with the Polishing Cup Head. Do not use the Enamel Polishing Paste with the izzo Brush Head or any other type of bristle brush.
- Enamel Polishing Paste is not intended for daily brushing or polishing; it's only intended for once-per-week use with the Polishing Cup Head. Overuse or use other than as directed may result in dentin and/or enamel abrasion.
- Use the Polishing Head on each tooth for no more than 3-5 seconds to avoid gum irritation or tooth injury.

Scaler

- Be careful not to insert the pointed tip below the gumline or around loose fillings or crowns.

SAVE THESE INSTRUCTIONS

FCC Statement: This equipment has been tested and found to comply with the limits for a Class B digital device, pursuant to part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation. This equipment generates, uses and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation. If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures:

- Reorient or relocate the receiving antenna.
- Increase the separation between the equipment and receiver.
- Connect the equipment into an outlet on a circuit different from that to which the receiver is connected.
- Consult the dealer or an experienced radio/TV technician for help.

This device complies with Part 18 of the FCC Rules.

Scaler

- To avoid irritation or trauma to your gumline, do not use the Scaler to clean under the gumline or in areas of the mouth that you cannot see.
- Do not use if any part of the Scaler is damaged.

Medical Safeguards

- Please consult with your dentist and dental hygienist to ensure there are no dental conditions that may prevent you from using izzo. Examples may include but are not limited to: recent oral surgery, severe periodontal disease, severe gingival sensitivity, local irritation or ulceration, active intraoral infection, or swelling.
- When first using a powered toothbrush, you may experience mild gum bleeding and irritation. If this persists for two weeks or more, discontinue use and contact your dental professional.

Symbols



Electromagnetic interference from the device is under limits approved by the Federal Communications Commission.



Do not rinse under running water, not waterproof.



Protection Class IPX7 – IPX7 Classification of ingress of water for Power Handle - can be submerged up to 1 meter in water for 30 minutes.



Dispose of in accordance with the Waste Electrical and Electronic Equipment Directive 2002/96/EC for product and accessory disposal.



Direct Current (DC) supply, 5 Volts, 1 Ampere



Manufacturing Date



Battery



Meets electrical safety requirements

Delivering Daily Dentistry



3 No-Brainer Reasons to **Subscribe & Save** to izzo Refills:

1. Achieve the Best Results: Because dentists recommend replacing your Brush Head, Polishing Cup Head, and Enamel Polishing Paste every 3 months, that's what we deliver.

2. Score the Most Savings & Convenience: Save big on your refills when you sign up for **Subscribe & Save**. Set it and forget it or easily modify or cancel anytime.

3. Upgrade to Lifetime Warranty: For a lifetime of clean, upgrade to a Lifetime Warranty*—yours for as long as you maintain an active **Subscribe & Save** refill plan.

For more information, check out the FAQs at [izzosmile.com](https://www.izzosmile.com).

*Requires active subscription for the past three consecutive refill cycles.

Smile On™



izzo™

1710 Romano Drive, Plymouth Meeting, PA 19462 U.S.A.

izzosmile.com | Patent Pending

Z-1002 0122045 REV3