




SWEETSPOT


MINI MESSAGE GUN™

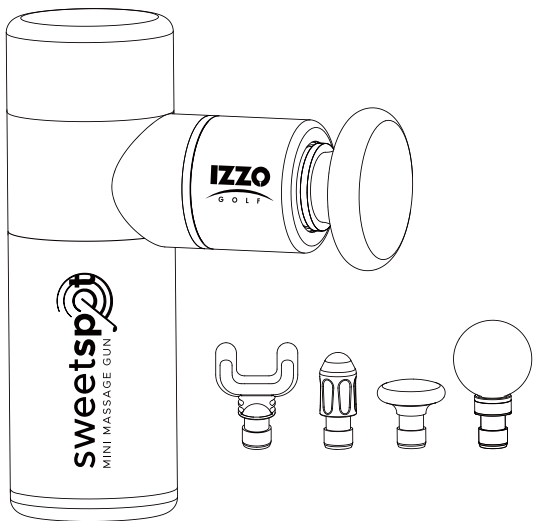


USER MANUAL

Overview

Main Body

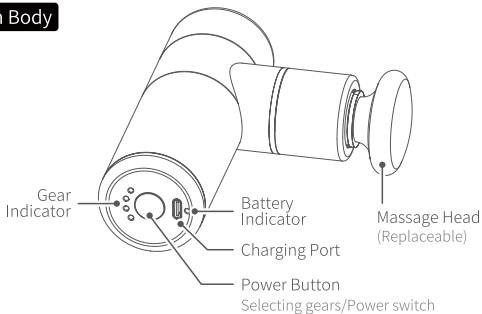

Massage Head

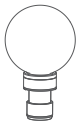

Spherical Head

For overall use, relax soft tissues & sensitive muscles




Flat Head

Massage flat major muscle groups



Cone Head

Pinpoint muscle treatment.
For Joints, Deep Tissue and Trigger Points

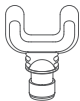


U-shaped Head

Massage Neck, Spine and Achilles Tendon

Start to Use

Install/ Replace Massage Head

Install a massage head into the hole in front of the device according to your need.

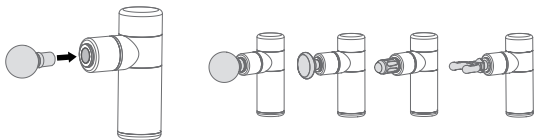

Power ON/OFF


PRESS & HOLD the Power Button to power ON/OFF.

- ※ Use it after powering on for 5s.
- ! Do not use it while charging.

Selecting Gears

Single tap Power Button to select gears after powering on.

Cycle through first gear/ second gear/ third gear and fourth gear.

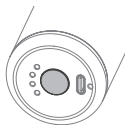

Indicator

Battery Indicator


Light	Status	Description
Red	■ ■ ■ ■ ■ ■ ■ ■ ■ ■	Charging while power OFF
Green	■ ■ ■ ■ ■ ■ ■ ■ ■ ■	Charge complete

Gear Indicator

Light	Status	Description
Blue	● ○ ○ ○	First gear
Blue	● ● ○ ○	Second gear
Blue	● ● ● ○	Third gear
Blue	● ● ● ●	Fourth gear
Blue	● ○ ○ ○	Low power and auto power OFF

Charging


! Do not use it while charging

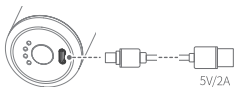


Specifications

Motor Rated Power	30W
Theoretical Battery Life	8 Hours
Weight (Reference value)	495g (Main Body)
Size	92x137x46mm

List

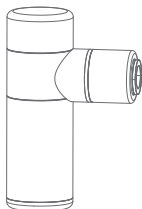

x1

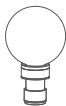

x1



x1




x1



x1



x1

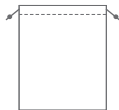

x1

Disclaimer

Thanks for using the Sweet Spot Mini Massage Gun. The information in this document affects your safety and your legal rights and responsibilities. Read the entire document carefully to ensure proper configuration before use. Failure to read & follow instructions & warnings in this document may result in serious injury to yourself or others, or damage to your products or to other objects in the vicinity.

By using this product, you hereby signify that you have read this disclaimer & warning carefully and that you understand & agree to abide by the terms and conditions herein. You agree that you are solely responsible for your own conduct while using this product, and for any consequences thereof. You agree to use this product only for purposes that are proper and in accordance with all applicable laws, rules, regulations, and all terms, precautions, practices, policies & guidelines our company has made & make available. Our company accepts no liability for damage, injury or any legal responsibility incurred directly or indirectly from the use of this product.

Safety Instructions

1. This device is not recommended to directly massage the skin. To avoid harm to the skin, customer should massage body after wearing dry clothes or placing a clean and dry towel on massage part, gently press and move the device.
2. High speed and high-pressure massaging can only be used on soft tissues without feeling pain or maladjustments. Do not massage head or any skeletal sites with high speed & high pressure.
3. Massaging the same part of your body with high speed & high pressure for a long time will cause bruises. Please stop using immediately after feeling pain or maladjustments.
4. Please keep your fingers, hair or any other body hair away from the extendable part of the device to prevent pinching.
5. Do not discard or misuse this device.
6. Do not modify this device by yourself without permission.
7. Do not operate or charge the device when nobody is around.
8. This device is for adults only.
9. People with contraindicated symptoms should use this device with doctor's permission.

Notice

1. Make sure massage head is not blocked by external force when device is powered on. If it is blocked, it will pause for a while and try to restart.
2. DO NOT contact the product with water or other liquid for the Percussion Massage Device is not waterproof.
3. DO NOT disassemble the product except where marked detachable.
4. Prolonged continuous operation may cause the motor surface temperature to rise, please operate carefully.
5. DO NOT drop or strike the product.

Storage & Maintenance

1. Keep the product out of the reach of children and pets.
2. DO NOT leave the product near heat sources such as a furnace or heater. DO NOT leave the product inside of a vehicle on hot days
3. Please store the product in dry environment.
4. DO NOT overcharge or overuse the battery, otherwise it will cause damage to the battery core.
5. Never use the product when the temperature is too high or too low.



1635 Commons Parkway
Macedon, NY 14502

CONTACT INFORMATION

PHONE: 800-284-1220

EMAIL: SALES@IZZO.COM

WEB: WWW.IZZO.COM