**Westinghouse**

INNOVATION YOU CAN BE SURE OF

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Thank you for purchasing this product. Please read this instruction manual carefully before operating the hand massager.

Please pay attention to the important safety warnings and keep this manual for future use.

Note: We reserve the right to revise the design and description of this product without notice.

## Contents

1. Safety warnings
2. Environment
3. Restrictions
4. Maintenance
5. Troubleshooting
6. Components
7. Operation
8. Product specifications
9. Warranty

## Français

### 1. Safety warnings

- This product is not recommended for children.
- Disconnect the massage unit before cleaning.
- Please operate this product according to the instructions in this manual.
- This product is designed for indoor use only.
- Do not exert excessive pressure on the massager. This may significantly slow down or stop the motor, potentially cause overheating and subsequent damage to the massager.
- Press the “stop” button to terminate the massager if discomfort is experienced.
- The total recommended usage time of the massager is 20 minutes.
- Please avoid use if the massage unit is damaged and call your Westinghouse™ technician or dealer.
- Do not use wet hands to pull the plug.
- Do not spill liquids onto this product to avoid electric shock or cause damage to this product.
- In the event of a power failure disconnect the massager.

### 2. Environment

The Westinghouse™ wireless hand-held massager should be used only in well ventilated, normally dry and dust-free environments.

### 3. Restrictions

People with osteoporosis, heart disease and heart pacemakers, pregnancies, injuries or any other health-related issues should consult their doctor before using this massage unit.

Stop immediately the use of the massager if you feel unwell and consult your doctor.

### 4. Maintenance

- This product should only be maintained by Westinghouse™ certified technicians; otherwise, the warranty will be void.
- Make sure the wall outlet is in good working order before use.
- If the massager will not be used for an extended period, store in a dry and dustless environment.
- Avoid product exposure to high temperatures and direct sunlight.
- Clean this product with a dry cloth. Do not use thinner, benzene or alcohol.


- We recommend using the massager intermittently and not for a long time.

## 5. Troubleshooting

- It is normal for the motor to make noise while in use.
- If the massager doesn't normally work, please check that the plug and socket are properly connected, and the switch is turned on.
- After several uses, the massager might power off because of overheating. It can be used again after half an hour of cooling time.

## 6. Components

### Structure description
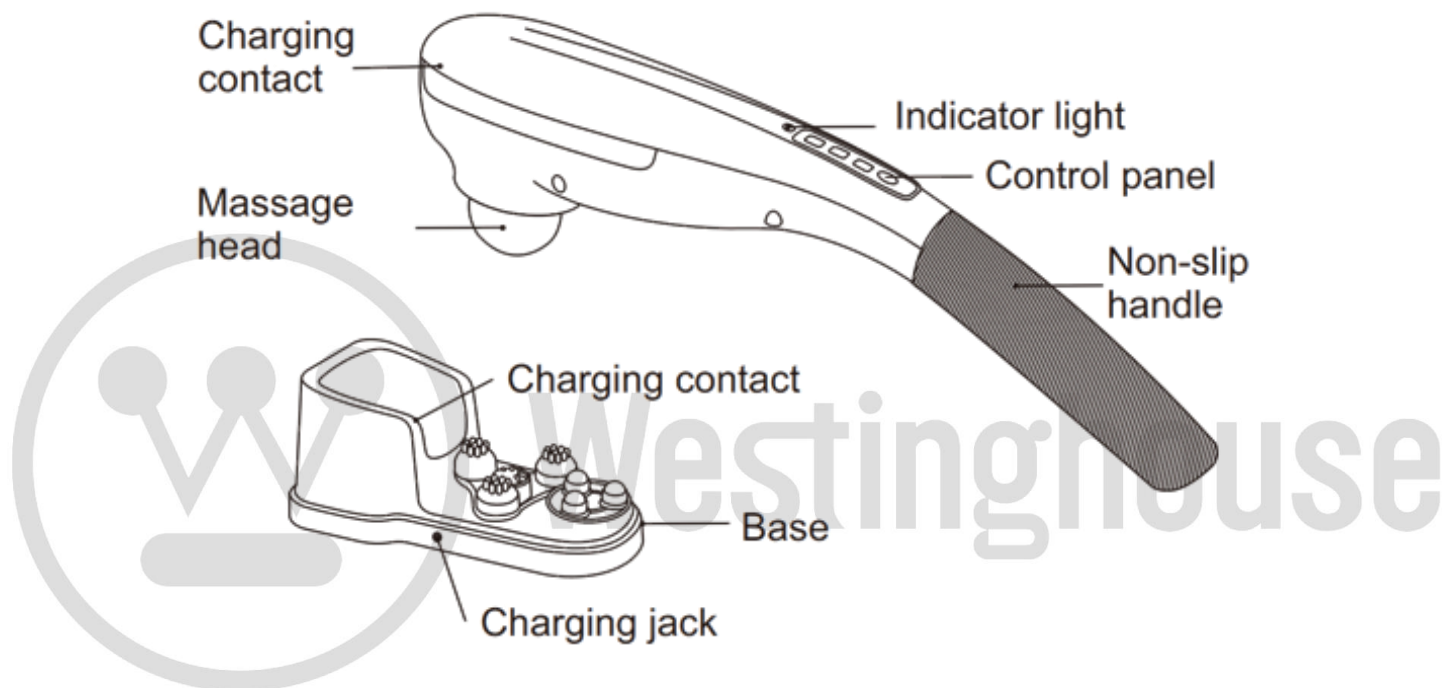
## 7. Operation





- Connect the massager to an AC wall outlet.
- Push the power button on the unit's control panel to turn the unit on or off.
- Press the MODE button to select the type of massage. There are 4 types, continuous, slow bursts, medium burst and long burst.
- Press on the up or down arrow to select the burst speed. There are three speeds.

The WES42-0916 hand massager comes with 2 pairs of massage head accessories. These vary the massage intensity and personality.

To add these massage head accessories, simply screw them onto the massage head. Unscrew to remove.

You can also enjoy the WES42-0916 and its relaxing effects on the neck, shoulder, the legs, abdomen and other body parts.
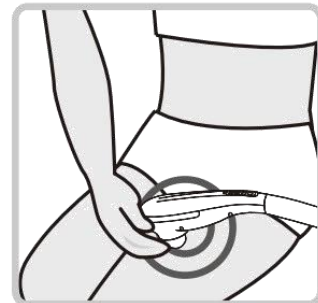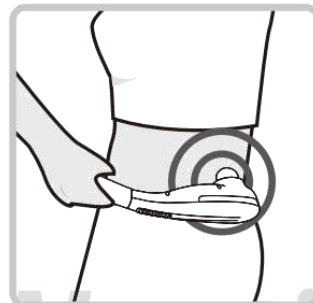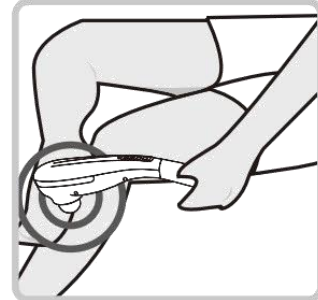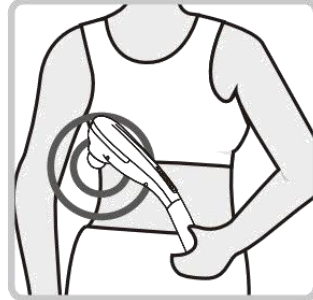
## 8. Product specifications

Input Voltage: AC110 V - 60 Hz

## 9. Warranty

### LIMITED ONE (1) YEAR WARRANTY

This product is guaranteed against all manufactured faults and defects for a period of one (1) year from the date of purchase.

For repairs or product defects, please call our customer service department at:

1-866-967-7333. Please have your proof of purchase on hand. If we are not able to replace your defective part, we will send you a service technician for repairs.

Important notes: Any damages caused by an accident or an abusive usage are NOT covered under warranty. We also decline all responsibility concerning indirect or consecutive damages. In certain Canadian provinces, it is not permitted to limit indirect or consecutive damages. Therefore, it is possible these limitations are not applicable to you. This warranty gives you specific legal rights. You may also have other legal rights which vary from one province to another.

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