

Oral stretching  Waving gear

S A L I O R A
b e a u t é

Massage gums with near infrared power ®

SALIORA beauté

Instruction Manual

Please read the instruction manual carefully
before using this product, and use it correctly.
In addition, keep the instruction manual in a safe place.

TRUSTLEX
with your life

Safety Precautions

Danger

Matters that are likely to cause death or serious injury

- Never use the product near electrical medical equipment or electric therapy apparatus for household use as described below.
- It may cause the malfunction of implantable medical devices such as pacemakers, medical equipment for life support such as heart-lung machines, wearable medical electrical equipment such as electrocardiographs and droppers, and electric therapy apparatus for household use.
- This product uses dry batteries. Do not throw it into fire, heat it, or use it or leave it at high temperatures, as it may cause overheating, ignition, or rupture.

Warning

Matters that may cause death or serious injury

- Do not disassemble, modify or repair this product, except to dispose of it. Doing so may cause fire, electric shock or injury. Consult the distributor for repairs.
- Do not replace the batteries with wet hands. Doing so may cause electric shock or injury.
- Do not place this product within reach of infants.
- Stop use immediately in the event of abnormality or malfunction. Continued use may cause fire, electric shock or injury.
- If you feel ill after using, discontinue use. It may cause an accident or poor health.
- This product is designed for gum massage. Do not use it for any purpose other than that described in this manual. If you cannot use this product properly, stop using it immediately. If you experience pain or discomfort, consult a physician immediately.
- This product should not be used alone by people who cannot express their intentions (including infants) or people with physical disabilities. It should not be used by people who have weak oral sensation. It may cause accident or injury.

Caution

Matters that may cause personal injury and cause only property damage

- Use by children should be under the supervision of a guardian.
- Do not apply strongly to teeth or gums, or apply to one area for too long. It may damage teeth or gums. Consult a dentist if you have severe periodontal disease.
- If you have bothersome symptoms in the oral cavity, such as of teeth during treatment, do not use in that area. Doing so may cause injury.
- Do not wash or dry this product in a dishwasher.
- Do not use thinner, benzene, or alcohol. Doing so may cause malfunction or damage.

Precautions for Use

- This product is for personal care, and is not for use with multiple patients in a dental clinic or medical institution.
- If you have had oral treatment or gum surgery in the past two months, please consult with a dentist before using this product, the same as when you use a general oral hygiene product for the first time.
- If severe bleeding occurs as a result of using this product, or if bleeding is observed with one week after using this product, discontinue use and consult a dentist.
- If you have any health concerns, consult your doctor before use.
- Press the button in the center of the button with the pad of your finger. Pressing with your fingernail may cause a malfunction.

Method of Operation and Usage

● Power On/Off

When the power button is turned on, the ultrasonic vibration and infrared output will start, and will automatically turn off after 3 minutes. To turn them off during operation, press the power button again.

① Exchangeable head cap: Pink




② Exchangeable head cap: Pink




③ Exchangeable head cap: Blue

Massage the gums

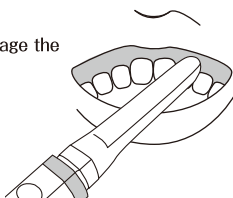


④ Exchangeable head cap: Blue

Massage the lips

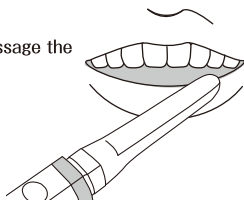

● Head caps

Remove the cap of the main body and attach an exchangeable head cap.

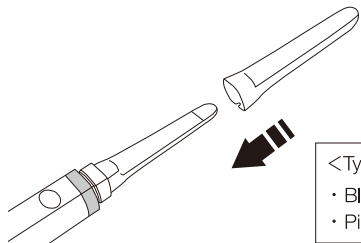


<Type of exchangeable head cap>

- Blue: For gums
- Pink: For beauty

● How to change the batteries

Turn the lower part of the main body counterclockwise and remove the lid. Insert two batteries in a row, with the negative poles toward the head, and turn the lid clockwise to close.

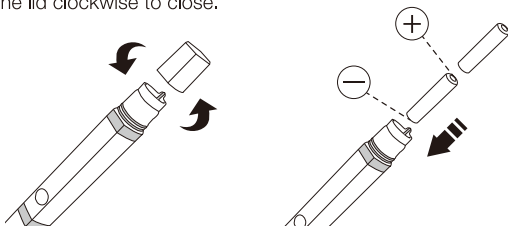


* Never replace the batteries with wet hands.


- * Number of uses (with new batteries): 150 to 200 times with 3 minutes of use.

After use, rinse only the head cap with water,
and wash it to keep it clean.

- * Since the main body is not waterproof,
do not expose it to water.

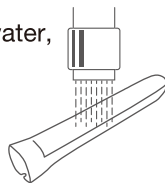

Contents of the Set and Name of Each Part

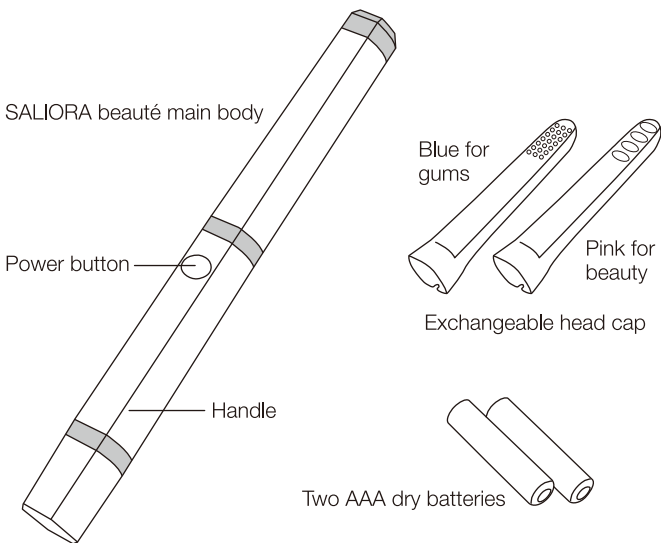

Product Specifications

Product name: SALIORA beauté
Body size: Length 220 mm
Handle diameter 23.5 mm
Main body weight: 53 g
Power supply: 2 AAA dry batteries
Operating temperature: 0°C to 45°C

Total usage time: 10 hours or more
Power consumption: about 0.5 W
Infrared output: 640 mW
Ultrasonic vibration: 240 kHz
Warranty: 1 year from purchase

Warranty

- Be sure to fill out the warranty card, including the date of purchase and name of the distributor, and keep it in a safe place after reading the contents carefully.
- The warranty period is one year from the date of purchase. Contact the distributor if a problem occurs during use.

Product No.	
Date of purchase	
Retailer	
Sale destination	

This product will be repaired free of charge up to one year from the date of purchase, if handling has been normal.

- When applying for repairs, submit the product and this warranty card.
- If any of the following falls within the warranty period, it will be beyond the scope of free repair, and materials will be charged separately.
- Failure or damage due to incorrect operation
- Failure or damage due to improper repair or modification
- Failure or damage due to transport or dropping after purchase
- Failure or damage due to fire, earthquake, flood, lightning or other weather
- Material and service fees will be charged separately for repairs after the warranty period.
- This warranty will not be reissued. Please keep it in a safe place.

Manufacturer **TRUSTLEX INC.**

2F Dai-2 Kimi Bldg., 2-11 Toyotsu-cho, Suita, Osaka 564-0051 Japan

Tel: +81-6-6378-1611 <http://www.trustlex.co.jp>