

User Guide

AXXZIA MATE FOR EYES



beauty-nj.com

Product Specifications:

Dimensions: Approximately 35mm (W) × 35mm (D) × 175mm (H)

Weight: Approximately 120g

Power Source: Rechargeable

Input Voltage: AC100-240V 50/60Hz

Power Consumption: 1.6W ±10%

Power Supply: Maximum output mode approximately 1 hour on full charge

Output Levels: 3 levels (*Start with the low-frequency EMS for testing)

Accessories: Main unit/Charging stand/Charging cord (75cm)

Dedicated AC adapter/Instruction manual

Battery Charging Time: Approximately 5 hours (Red light indicates charging, green light indicates fully charged).

General Information

AXXZIA MATE FOR EYES is the renewed eye-focused beauty device that targets stimulation of the orbicularis oculi muscle that plays a crucial role in the appearance of the eyes and surrounding areas. Designed for effective anti-aging care and effortless home eye lift care without downtime.

Renewed AXXZIA MATE FOR EYES has evolved from its initial release in January 2019 to meet the diverse needs of modern skincare. With enhanced EMS power, advanced RF technology, and a focus on user convenience, this device represents a significant leap in at-home eye care.

- **Key Features:**

- Strengthens orbicularis oculi muscle for a natural eye lift effect
- Three evolving modes for customized care: Purple (EMS Low), Orange (EMS Medium + LED), Red (EMS High + LED)
- New RF technology for deep skin warming and enhanced skincare absorption.

- **Benefits:**

- Promotes collagen production for firmer, youthful-looking skin
- Reduces signs of aging like sagging and dryness around the eyes
- Compatible with your favorite skincare products for enhanced results

Product Features:

- **Smooth flat head:** Easy to use over eye sheets or sheet masks with a smooth flat head. Directly apply serums or gels to the head.
- **Simple operation:** A single button for power ON/OFF and switching EMS modes.
- **Cordless charging:** Place it on the charging stand for approximately 4 hours for easy charging.

Precautions

! Not recommended to use for the following people:

- People with heart disease.
- People with facial nerve disorders.
- People with diabetes.
- People with bleeding disorders, epilepsy, tuberculosis, or immunosuppressive diseases.
- People with blood abnormalities such as hemophilia.
- People with high blood pressure.
- People with retinal detachment or a history of it.
- People who have undergone surgery to implant metal, silicone, plastic, or similar materials in their face or body due to fractures or orthopedic procedures.
- Infants, children in their growth period, and those unable to express their will.
- People with head or ophthalmic sequelae from accidents, etc.
- People who have undergone eye surgery.
- People with particularly sensitive skin.
- People with skin conditions (melanosis, atopic dermatitis, or other skin diseases) on their face.
- People experiencing capillary dilation due to long-term use of steroid hormones or liver function impairment.
- People with allergies.
- People during menstruation, pregnancy, or breastfeeding.
- People with abnormalities such as trauma or burns on their face.
- People with anemia or high blood pressure (180mmHg or higher).

! Do not use in the following cases or areas:

- **When feeling unwell or severely fatigued.**
- **When having a fever.**
- **While consuming alcohol or medication.**
- **After surgical procedures, laser treatments, or deep chemical peels in the application area within the past 3 months.**
- **When having skin abnormalities such as dermatitis or excessive sunburn.**
- **On wounds, around the ears, on areas where hair removal was done on the same day, or on areas where plastic surgery was performed.**
- **On the upper eyelids or the throat.**

Handling Precautions

- ! - Do not unplug or plug the adapter with wet hands.
- Do not immerse the main unit or AC adapter in water or wash it with water.
- Do not use or store in bathrooms, humid environments, or places prone to water exposure.
- Do not use if the AC adapter is damaged or the plug is loose in the socket.
- Do not use damaged cord.
- Do not use beyond the rated capacity of the socket or wiring devices.
- Use only the dedicated AC adapter.
- Do not use the dedicated AC adapter with other products.
- Always use within the range of AC 100-240V.
- Insert the plug of the adapter or device securely to the base.
- Always unplug the adapter when cleaning the device.
- Regularly clean the power plug and device plug to prevent dust accumulation.

! For safety reasons, do not use the following in combination with or near this product:

- Implantable medical electronic devices such as pacemakers.
- Life-support medical electronic devices such as artificial hearts.
- Wearable medical electronic devices such as electrocardiographs or infusion pumps. (These may cause malfunction or accidents due to medical electronic devices or home treatment devices.)
- **Do not use while bathing. (May cause strong electric shock, electric shock, or burns.)**
- Do not leave the product powered on and unattended during use. This can cause accidents.
- **Remove accessories, watches, etc. from your body before use.**

Handling Precautions

- ! -Do not throw the device into fire or heat.
- Do not charge, use, or leave near heat-producing devices such as stoves or hot air machines.
- Do not use or store within reach of infants or toddlers, or allow them to use it.
- Do not modify or repair it yourself.
- In case of abnormalities or malfunctions, immediately stop using it.

<Abnormalities or Malfunctions Examples>

The main unit, adapter, or cord:

- Feels unusually hot
- Shows signs of deformation
- Has a burnt smell
- Makes abnormal noises

! Attention:

- Do not use continuously for long periods on the same area.
- Do not fix the head part to the same area during use.
- Do not stare directly at the indicator light.
- When unplugging the adapter or device plug, always hold the adapter or device plug without touching the cord.
- Unplug the adapter from the socket when not charging.

EMS (Electrical Muscle Stimulation)

EMS Technology: The device utilizes Electrical Muscle Stimulation (EMS) technology to deliver gentle electrical impulses to the orbicularis oculi muscle. These impulses mimic the natural contraction and relaxation of muscles, effectively toning and strengthening the muscle fibers.

Modes:

1. Purple (EMS Weak Mode):

- Ideal for first-time users. Start with the weak mode indicated by the purple LED.
- Gentle EMS strength targets fatigued muscles around the eyes.

2. Orange (EMS Medium Mode):

- The orange LED targets the skin, providing firmness and a brighter complexion.

3. Red (EMS Strong Mode):

- Specifically targets areas where firmness and dryness are concerns.
- Tighten the skin with high-output EMS while receiving treatment from the red LED.

RF (Radio Frequency)

AXXZIA MATE FOR EYES incorporates Radiofrequency (RF) technology, which generates heat energy to penetrate deeper layers of the skin. This heat stimulation not only enhances muscle activation but also promotes collagen production, leading to firmer and more resilient skin around the eyes.

- Comfortably warms the skin up to the stratum corneum to enhance the penetration of beauty ingredients.
- Transforms dull skin around the eyes into a brighter appearance.
- Provides soothing warm care for eyes exposed to PCs or smartphones for extended periods.

LED

LED Therapy: The device features LED lights with specific wavelengths (e.g., orange and red) that work synergistically with EMS and RF functions. These LED lights contribute to skin rejuvenation, targeting various skin concerns such as fine lines, wrinkles, and uneven skin tone.

- LED wavelengths are available in **Medium Mode (orange) and Strong Mode (red)**.
- Use this mode in combination with eye care cosmetics (eye cream, serum, gel or eye sheet masks) to enhance the penetration of beauty ingredients.

How to use

Recommended treatment method for eye and facial care:

1. Remove makeup and cleanse your face to ensure your skin is clean.
2. When caring for the eye area, always remove contact lenses beforehand.
3. Take an appropriate amount of serum or gel and apply it to the areas you want to treat.
4. **Press and hold the power button** at the bottom of the device to turn it on, then **select the output level**. (Upon startup, it will automatically set to low with a purple light.

Pressing the power button again switches to medium with an orange light, and pressing it once more switches to high with a red light. For initial use, low level is recommended.)

- **Recommended usage frequency:** Approximately 5-10 minutes once a day. (Adjust based on your skin's response).


- Start from the side of the nose and move towards the outer corner of the eye, drawing a circle.

Usage Tips:


- Always use this product in combination with a serum, gel, or eye care seed.
- If device doesn't slide well on your skin apply more cream/serum/gel.
- Do not use on the throat, eyeballs, or mucous membranes.

Recommended Products for Eye Care:

1. AXXZIA Beauty Eyes Routine Care Essence Premium
2. AXXZIA Beauty Eyes Essence Sheet Premium Plus

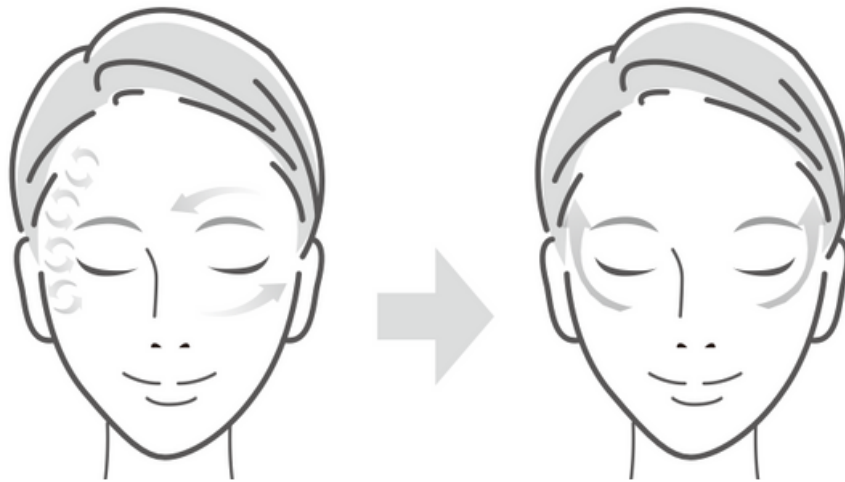

Draw circles around the corners of your eyes and temples and gently lift the skin from bottom to top.

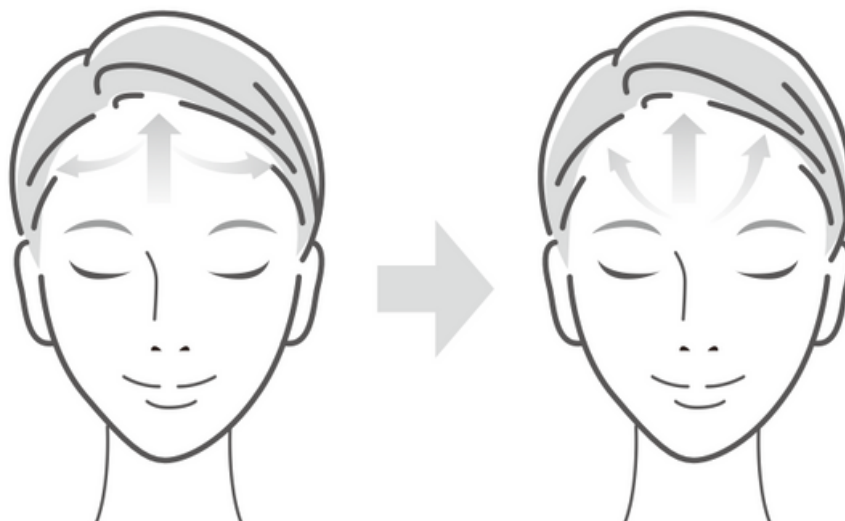

Move the device up from eyebrows and then down to both sides of the forehead. Move from eyebrows to the hairline.

This product will **automatically turn off 10 minutes after starting operation**. If you need to turn off the power during use, please press and hold the power button.

Please Note:

- When starting to use, you may not immediately feel noticeable warmth, but after a few minutes, you should gradually feel warmth from within the skin. If you experience abnormal heat on the skin during use, stop using it immediately and contact the sales store.

Recommended Usage:

- Use once a day for about 5-10 minutes per session. Adjust based on your skin's response.
- Move the device from bottom to top to lift the skin.

Daily Care and Storage:

After use, wipe the head part with a tissue to remove dirt. If there is noticeable dirt, use a cloth or cotton soaked in water or lukewarm water to wipe it clean. **Never wash the device under the water as it can cause malfunction.**

- Avoid using paint thinner, gasoline, acetone, chlorine-based, alkaline, volatile substances (such as nail polish remover, alcohol, etc.), or acidic detergents as they can cause damage, part breakage, or discoloration.
- Store in a clean, well-ventilated place away from stoves, heaters, direct sunlight, high temperatures, or humidity.

Safety Tips:

- If the device has not been used for a long time, ensure it functions correctly before use.
- Do not use if you notice a burnt smell or hear abnormal sounds.
- Clean any dirt or dust from the plug socket or head part before use.
- **This product is not fully waterproof. Avoid water exposure or immersion.**