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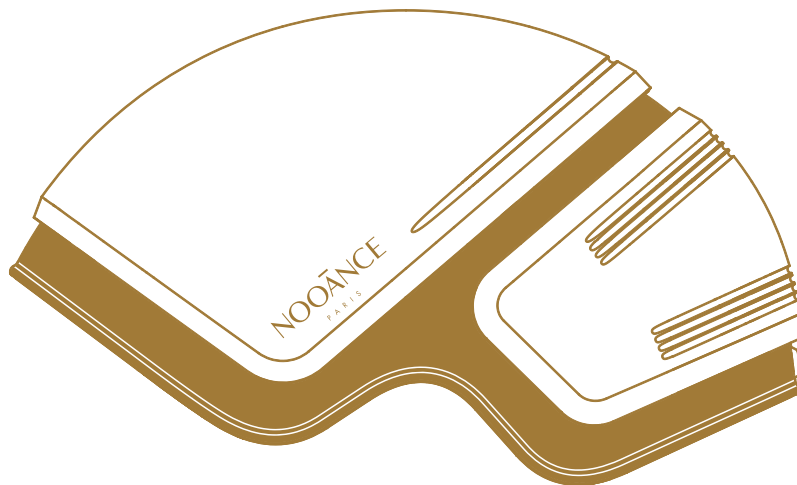
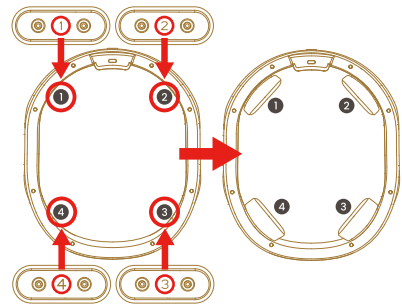
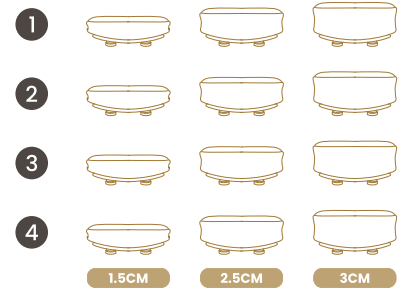
How To Use The Device (Two ways to use)

Preparation: Make sure your hair is dry before using the NOOANCE helmet.

Please select the appropriate sponge pad size to fit inside the NOOANCE helmet.

Connect the controller to use

Step 1

Connect the NOOANCE helmet to the controller using a USB cable.

Step 2

Position the NOOANCE helmet in its upright position on the top of your head.

Step 3

Sit down and position the dome comfortably. Long press the ON/OFF Button to turn on the NOOANCE helmet.

(Note: If the start/pause button is not pressed within five minute, the NOOANCE helmet will turn off automatically.)

Step 4

Tap the ON/OFF Button to select the preferred mode. (Mode 1 is continuity; Mode 2 is the pulse.)

Step 5

Press the start/pause button to begin the 25-minute treatment. A "beep" indicates the treatment has started. Once the treatment is complete, the NOOANCE helmet will turn off automatically with two "beep" sounds.

Step 6

To turn off the NOOANCE helmet, hold the button for 2 seconds.

Connect the adapter to use

Step 1

Connect the NOOANCE helmet to the adapter.

Step 2

Long press the On/Off button to turn on the NOOANCE helmet default start mode 1, and the indicator light on the NOOANCE helmet will always be white.

Step 3

Press the On/Off button again to switch to mode 2, and the indicator light on the NOOANCE helmet will flash white.

Step 4

Place the NOOANCE helmet vertically on your head after selecting your preferred mode. Once the treatment is complete, the NOOANCE helmet will turn off automatically with two "beep" sounds.

Step 5

To turn off the NOOANCE helmet, hold the button for 2 seconds.

Light Sensitivity Test

We highly recommend users perform a light sensitivity test before their first use of the NOOANCE helmet.

Please refer to the "Light Sensitivity Test" section of the manual for specific test steps.

Installation of bearings

There are 3 sizes to choose from. At the bottom of each sponge is the corresponding number. Please put the correct number in the specified position.