
Product manual of **PELCAS**
Nanometer Electric Micro Needles Machine


USE MANUAL

USE MANUAL

1. Use methods of stoste or essence liquid:

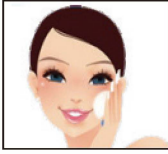
Clean skin

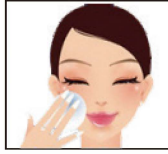Coate with stoste or essence liquid


Use nanometer needles tipsto act vertically on the skin about 10 to 15 minutes

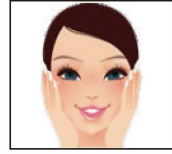Extra facial essence can be gently pressed to make it naturally absorbed

2. Use methods of mask paste or eye paste mask:

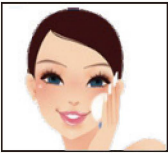Clean skin

Apply Mask Paste or Eye paste mask

Application of nursing at the same time use nanometer needles tips to act vertically on the skin about 10 to 15 minutes

Extra facial essence can be gently pressed to make it naturally absorbed

COMMON PROBLEMS BEFORE USE

1. How to clean?

Face and hands need fully clean. Please use hand washing or toilet soap to clean hands, use facial cleanser or liquid face wash to clean face.

2. Do the nanometer needles tips need to be disinfected?

Nedles nanometer tips do not need to be disinfected because we use aseptic packaging to pack them.

3. How many times do you need import essence to Nanometer Electric Micro Needles Machine in a week?

Usually twice a week, thin horny once a week.

4. What products can be Imported to Nanometer Electric Micro Needles Machine?

We advise import independent packaging and components of a single products to the machine, such as, lyophilized powder, stoste (hyaluronic acid), water needle and hydrating mask.

COMMON PROBLEMS IN USE

1. How to operate the import methods?

All nanometer needles tips need point lift and slide. If your skin tolerance is weak, you had better point lift. If your skin tolerance is strong and thick horny, you should slide. However, eyes, forehead, mouth, cheekbones, neck need point lift because their less fat. We advise point lift when you use the machine for the first time. You could gently slide on your face with step by step methods and intensity the second time if you do have no discomfort. You need find yourown appropriate methods and intensity previous period. You could operate according to the your experience later period.

2. How to master degree scale and intensity when you import stoste or essence liquid to

Nanometer Electric Micro Needles Machine?

Please adjust degree scale to under 0. 25mm for the first time. You could adjust degree scale according to your own situation later period. First please contact with the skin lightly and slowly, according to the skin feeling. enlarge the intensity when you have no discomfort.

3. How to import various products to machine?

Essence products: after applying the essence to the face, the essence can be vertically imported into the face with nanometer micro needles and can be applied to the face at one time. You could also apply the essence at the same time import into the face.

Mask products (wash-free mask):after using mask to the face, the essence liquid can be vertically imported into the face with nanometer micro needles about 5 to 10 minutes.

4. How many products could you import to machine one time?

We advise importing no more than two products one time. If your skin tolerance is weak, you had better import one product one time.

COMMON PROBLEMS AFTER USE

1. Do you need use other skin care products after importing products to machine?

You had better use moisturizing mask or repair skin mask after importing products to machine. You could use moisturizing cream if your skin is dry

2. Could your face touch the water after importing products to the machine?

We advise importing products to the machine at night, your face does not touch water at that night.

The nex day, you could use normal skin care products if your face has no discomfort (red and swollen or itches). Otherwise you need use moisturizing skin care products.

3. Could you bask in the sunshine, eat spicy food, seafood and alcoholic drink?

We advise you not exposing your skin to the sun, eating spicy foods. seafood and alcoholic drink and other stimulating foods within 3 days.

4. Could the Needles nanometer tips reuse?

Needles nanometer tips is disposable, they can oxidize when you open the packing. Do not reuse needles nanometer tips because the used tips may breed bacteria.

5. How to deal with if your face becomes red and swollen or itches after importing products to the machine?

You could use moisturizing mask or repair skin mask within three consecutive days. It usually could improve the next day.

HOW TO USE AND IMPORT LYOPHILIZED POWDER?

Lyophilized powder does not with alcoholic products or machines using. A pair of lyophilized powder include a bottle of solvent liquid and a bottle of dry powder. Please mix solvent liquid with dry powder before use. Then they could dissolve to freeze dried essence. In addition to the use of importing lyophilized powder, the daily care need smear. A well dissolved bottle of freeze-dried essence need run out within three days. importing a bottle one time. The activity of lyophilized powder is relatively high. If your skin is easy to be allergic, you smear some part of your face. You could import lyophilized powder without allergy after 2days. If there is blain on your face, do not import into full face, only using to the part with pockmark and pox pit

HOW TO IMPORT WATER NEEDLES TO MACHINE?

Please import a piece of water needle once a week. After importing water needle, you need paste moisturizing mask within seven consecutive days.

If you did not import water needles, you had better use smear type water needle (most of ingredient is hyaluronic acid) twice a week (once 1 to two pieces). After importing smear type water needle, please paste moisturizing mask.

HOW TO USE VARIOUS PRODUCTS FOR DIFFERENT PROBLEMATIC SKINS?

Dry or sensitive skin:

Please import replenishment water products (hyaluronic acid or moisturizing mask)first three times before using efficacy products. You could import efficacy products when your skin is not hydroperic. Spotted skin : we advise freckle cream with replenishment water products(for example freckle whitening essence) using.

Blain blain skin :

You could use removal blain products to import. Blain blain skin has inflammation, so please only import part of blain blain skin so as to avoid infecting other parts.

Notes of importing products to blain blain and pockmark:

You could use removal blain products with lyophilized powder imported. Firstly, please use of lyophilized powder to import into part of pockmark at the same time steering clear of blain blain, avoiding infection. Finally, removal blain products could be vertically imported into part of blain blain.

Concrete operation methods to pox pits:

After cleaning face, please use machine point lift on part of pox pits. Then you need smear lyophilized powder on pox pits. Then please use machine point lift until pox pits become red or slight pain. Please import removal blain products twice a week, lyophilized powder smeared at daily morning and night, it has obvious effect that it is continuous use more than 3 months.

OPERATION METHODS

Point lift: nanometer needles tips fixed to nanometer electric micro needles machine, then use nanometer needles tips act vertically on the skin. When you could fell slight sunken feeling. please point lift with gentle, speed, and concentrated methods. Then you need adjust intensity according to touching skin.

Slide : nanometer needles tips fixed to nanometer electric micro needles machine, then use nanometer needles tips act vertically on the skin. When you could fell slight sunken feeling, you need slide with straight line, at the same time. please stay 3 to 5 seconds on the same skin and not beyond twice.

NANOMETER ELECTRIC MICRO NEEDLES MACHINE RANGE OF APPLICATION

1. Skin care area:

Apply to daily nurse, you could coordinate with skin care products, in order to promote efficiently absorbing nutrient substance.

2. Skin repair area:

Apply to wheelk. spot, scar and slow death, you could coordinate with, as to promote efficiently absorbing products, in order that your skin could get better repair effect.

USE METHODS

1. Clean skin: take off formal dress and ornaments and clean face, you could use toner to secondary clean.
2. Get out nanometer needles tips, then fix to the top of Nanometer Electric Micro Needles Machine.
3. Drip moderate stoste or essence liquid on the top of nanometer needles tips or smear on the face.
4. Please long press power switch key with 3 seconds to turn on Nanometer Electric Micro Needles Machine, then machine acts vertically on the surface of skin, every area of skin does not be stayed on 3 to 5 seconds.
5. When you finish operating machine, please turn off switch and take down nanometer needles tips, scar daily skin care products(essence liquid)
6. Do not reuse nanometer needles tips, because they are disposable. Please drop them after you use.

KEEP IN STORAGE

1. Please put machine on well-ventilated surroundings, keep clean and dry place and away from direct sunlight, damp, and weight press.
2. Please do not put the place that children could touch.
3. Please take down nanometer needles tips when you store machine.

PAY ATTENTION TO

Notes of Nanometer Electric Micro Needles Machine:


1. Do not unload or modify charger and machine, do not use charger under power supply damage, otherwise it will lead to touch power shock, fire and switch damage.
2. Do not touch charger with wet hands. Please cut off the electricity supply once charger touch fire or water liquid.
3. Do not use short-circuit charger. Do not mechanically shake or rock charger, do not be direct sunlight. Do not use charger in the place that the high temperature (such as bathroom). Do not get wet in the rain or become damp.
4. Please do not use the charger near the TV set or radio.
5. Pay attention to charger with inserting the plug, the power line damages is you pull the power cord. Nanometer Electric Micro Needles Machine and charger must be put on a temperature of minus 10 degrees to plus 55 degrees well-ventilated place. It is dangerous that you use charge without identification charger, at the same time, it goes against
6. Please do not put the charger to the place that children could touch.

Notes of skin care products:

1. Please choose the skin care products that they are safe, valid and no harm for you.

Notes of nanometer needles tips:

1. Do not use if the disposable packing is damage or broken. Because it may cause skin contamination.
2. Please clean your hands and face before operating.
3. Please make a skin test combining with other skin care products in your hands for the first time. Do not use if your skin appears itches etc allergic phenomenon.
4. The nanometer needles tip is the role of moderate intensity, not to weight press, there is no need to raise wrist to nanometer needles tip lead to touch the skin inhomogeneously.
5. It is normal phenomenon that some people's skin becomes red after using nanometer needles tips within 10 to 20 minutes. It will disappear naturally within 24 hours.
7. The same nanometer needles tip is used to only one person at a time. Do not share the same nanometer needles tip with others, do not reuse many times with the same nanometer needles tip, otherwise it must lead to skin sensibility, inflammation, redness and spots and so on bad result.
8. Saving nanometer needles tips, please drop used nanometer needles tips and put unopened nanometer needles tips and machine on well-ventilated surroundings and away from direct sunlight, damp, weight press and children.
9. Nanometer needles tips need to continuously moved, the same area of skin does not be stayed on 3 to 5 seconds and not beyond twice.

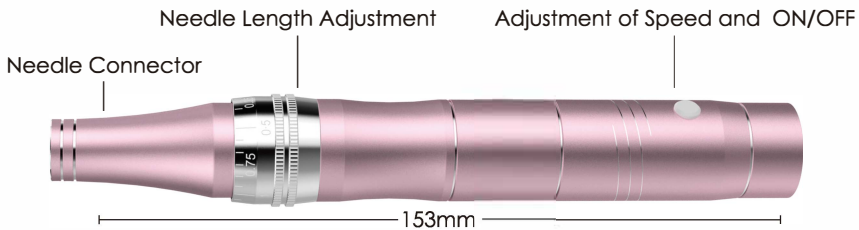level—Blue level (low speed)

level—Green level (intermediate speed)

level—Red level (highspeed)

Friendly remind:

Micro Needles Machine is a physical way to open the skin channels to promote product absorption. Solving the skin problems also need to coordinate with the corresponding products.

Nanometer needles round tips:

It mainly be suitable for daily care skin, importing nutrients, promoting absorbing various essence care skin products

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