

# PSORIASIS ON THE BACK

A patient from Switzerland, who presented chronic psoriasis symptoms, but from time to time there were severe flare-ups, when the whole back area was affected.

During the last cycle, the patient treated herself with Bioptron Pro1. Bioptron was used as a monotherapy (no other treatment was applied), once a day at a 10 cm distance for 60 minutes, in 10-minute treatments on 6 areas - in order to cover the whole surface of the back.

Before the treatment with Bioptron Hyperlight Therapy, the whole surface of her back was covered with red patches. Moreover, there were pustules on the skin and it was very rough.

After 9 days of Bioptron Hyperlight Therapy treatment, the red patches turned into milder pink, which was a marker of inflammation decrease. The number of pustules was highly reduced and the skin was smoother.

After additional 6 days of Therapy, the redness further decreased, while the pustules almost disappeared.

After 17 days in total, the patient's back has recovered drastically. The redness almost completely vanished, there are no pustules, and the repaired skin is smooth and without inflammation.



Start of treatment



3 days of treatment



6 days of treatment




9 days of treatment



12 days of treatment



15 days of treatment

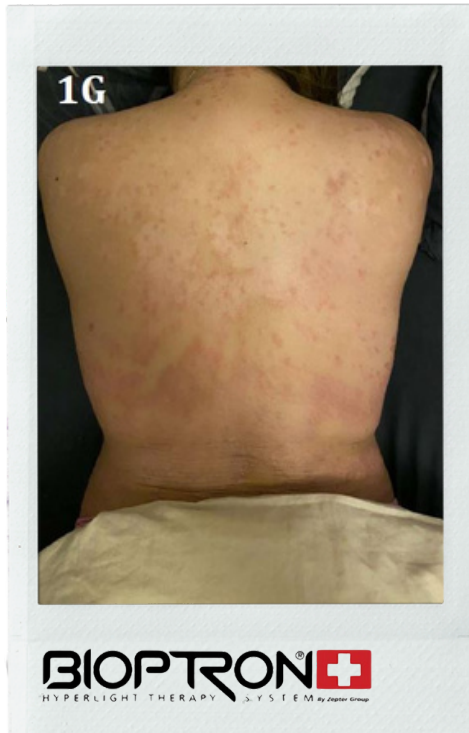

17 days of treatment