



SKIN
FOR ALL
£
HAIR
FOR MEN

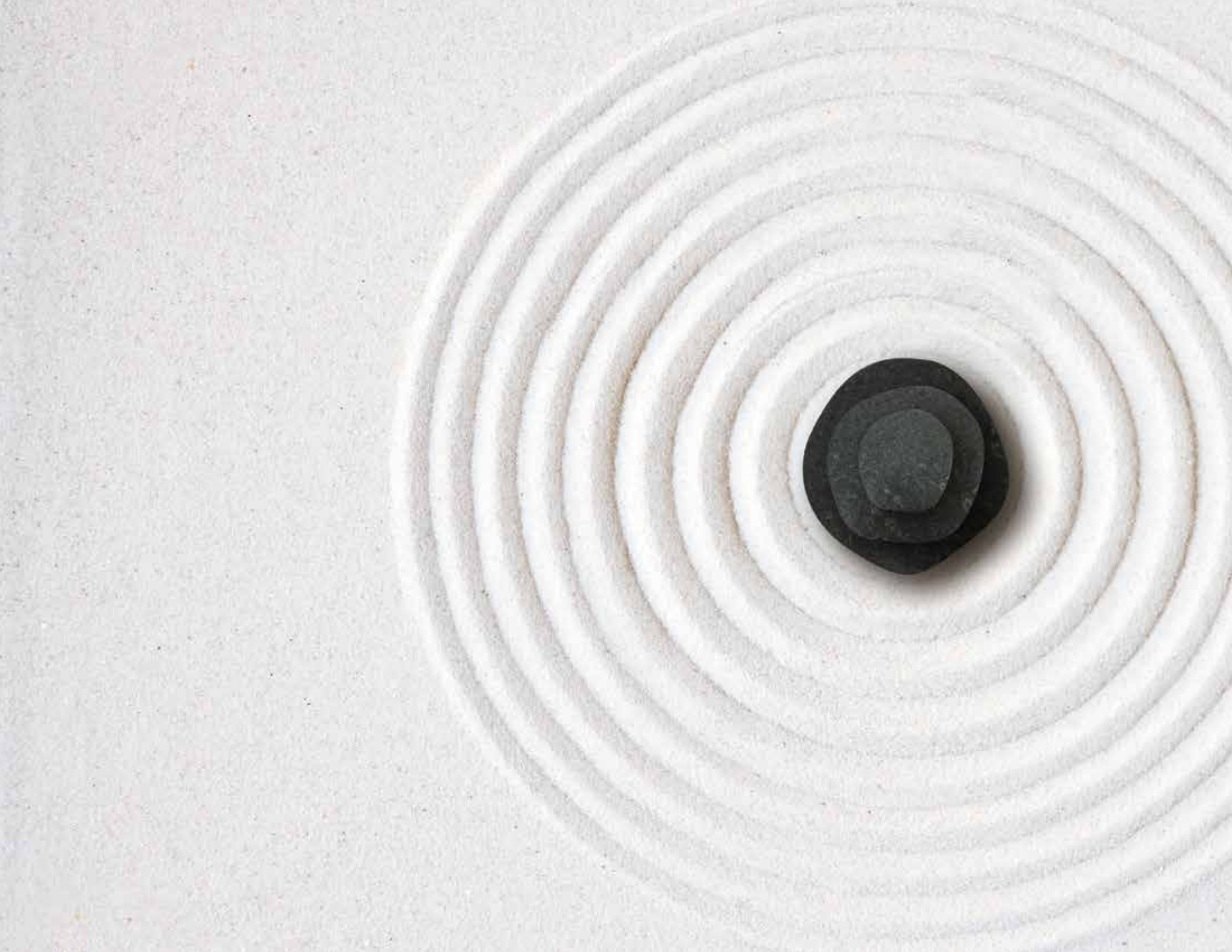




TABLE OF CONTENTS

A. Introduction	03
B. Luxe-Skin Revitalizing Treatment by Myufull	05
C. Features	07
D. Facial Muscles & Contours	09
E. Hands	10
F. At-home Instructions	11
G. Safety Guidelines	13
H. Product Specifications	14



A. Introduction

0o Skin

Facial muscles are the muscles underneath your skin. Typically around 30% of your facial muscles are used on a daily basis. One of the main causes of wrinkles and weakening or sagging of your skin is due to the limited use of your facial muscles.



Firm, tight, and healthy-looking skin is a result of the balance between the firmness of your facial muscles, collagen, and moderate fat, all which support each other. When any one of these factors are out of balance, the skin will begin to lose its firmness and develop fine lines at a rapid rate. The typical areas that tend to experience this are the cheeks, mouth, lines around the chin and eyes.

The multiple different forms of currents will assist in lifting, sculpting, and firming the skin - which are the immediate results - however, studies have shown that these currents also can increase elastin fibers and collagen thickness.

This is a composite photo of a 55-year-old mother and her 30-year-old daughter. Notice the change in the positioning and shape of the chin, jaw line, and around the lips and eyes.



B. Luxe-Skin Revitalizing Treatment by Myufull

Oo Skin

The Luxe-Skin Revitalizing Treatment by Myufull is an anti-ageing device that utilizes advanced EMS, waveform, and AEON technology to reduce wrinkles and lift facial contours. The look and feel of the face is supported by the facial muscles that stretch throughout the entire face. Just as the muscles in your arms and legs will weaken over time without proper exercise. The Luxe-Skin Revitalizing Treatment by Myufull is designed to rebuild your facial muscles, stimulate collagen production, and assist in restoring elasticity. The waveform from the device also increases ATP production within your skin cells. This increase in ATP production is necessary for your skin cells to produce collagen and elastin. Both of which decrease as you age. The Luxe-Skin Revitalizing Treatment by Myufull is one of the only devices on the market with EMS technology, AEON technology as well as waveform. By utilizing these 3 advanced technologies, you can help reduce fine lines and wrinkles, tighten and lift your skin, as well as reduce breakouts and scarring. All at the comfort of your home or at the Myufull Facial Treatment centre within Oo Spa.

EMS (Electronical Muscle Stimulation)

An exercise method that promotes muscle activity by passing a pulse current. Muscles repeatedly contract and relax when a pulse current is transmitted. By taking advantage of this characteristic, you can effectively train your facial muscles without moving your body.





5 Simple Steps

- 1.** Prior to use on the face, turn on the Device. Then place it on the back of your hands and rotate it in a circular motion for up to 1 minute per hand.
- 2.** Now place the Device on the back and side of your neck and rotate it in a circular motion then move it up along your neck. Repeat the same movements on your shoulders and move up to your neck. Do this for up to 1 minute on each shoulder.



- 3.** Now begin by placing the Device on your left cheek beside your mouth and move up your cheek towards your temple.



- 4.** Now repeat the same motion on the right cheek.



- 5.** Now repeat the same motion on both sides evenly for up to 10 minutes on each side of your face.

We recommend following these 5 simple steps every day to maintain the results



C. Features

0o Skin





Technologically more advanced than any other machine on the market. Featuring the ability to not only target your face and hands but multiple other body parts simultaneously

Featuring EMS technology, AEON technology, as well as waveform technology. These features provide the ability to not only target your face, hands, and neck, but also multiple other parts of your body simultaneously

The low to high frequency intervals also mimic the effects of soft massaging.


In addition, the battery lasts longer than standard devices by reducing power consumption. A set of batteries can be used up to 4 weeks with 30-minute daily use.



D. Facial Muscles and Contours

0o Skin

There are five key areas that the Luxe-Skin Revitalizing Treatment by Myufull provides benefits to: the corner of the eyes, the length and depth of the nasolabial fold, the corners of the mouth, and the cheeks.

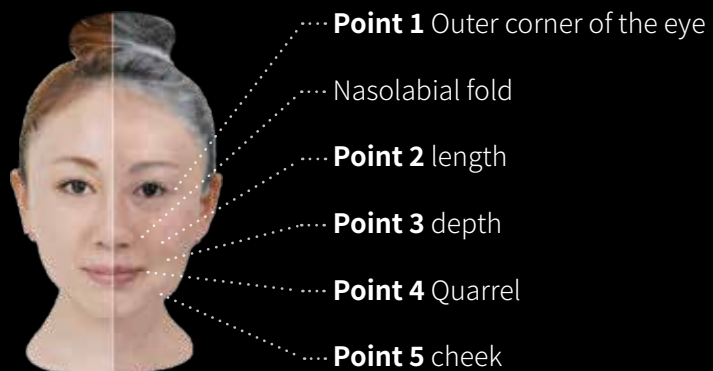




E. Hands

Oo Skin

Connection between the hands and organs of the body

Similar to your feet, your hands have various points which directly connect to certain organs. By using the Luxe-Skin Revitalizing Treatment by Myufull device on hands, you stimulate the brain and organs.

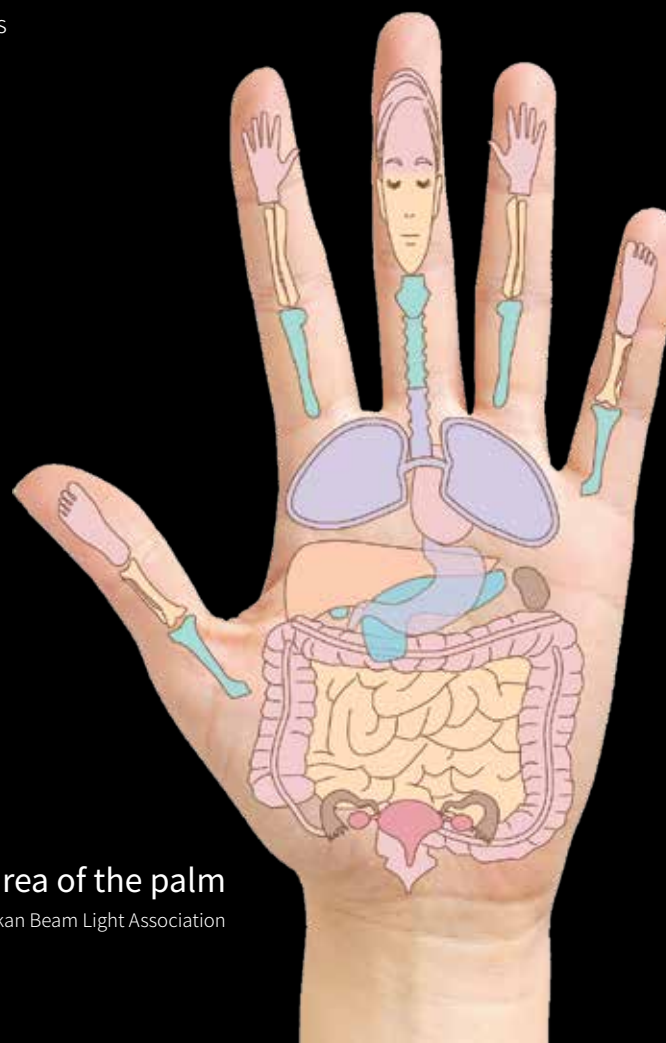

Reflective area of the palm

Source: Seishi Kaikan Beam Light Association

F. At home instructions

Oo Skin



Face, Neck, & Shoulders

Begin by placing the Luxe-Skin Revitalizing Treatment device on the shoulder, one side at a time, then rotate in a circular motion. As you rotate the Luxe-Skin Revitalizing Treatment device, move it up the back of your neck to behind your ears. Follow through by placing the device near your mouth and move it up the side of your cheeks to your temple. Repeat these steps for 10 minutes on each side of your face.



EMS & Hair Treatments


By placing the Luxe-Skin Revitalizing Treatment device on your scalp and rotating clockwise as well as moving up the sides of your head, you are providing your hair with negative ions which stimulate hair growth and also aids with thinning hair.





Acne & Skin Complexion

Skin complexion, years of acne and red spots can be naturally reversed and the skin can be enhanced with the Luxe-Skin Revitalizing Treatment device. By increasing blood flow and stimulating the production of ATP and collagen, you can have youthful skin.

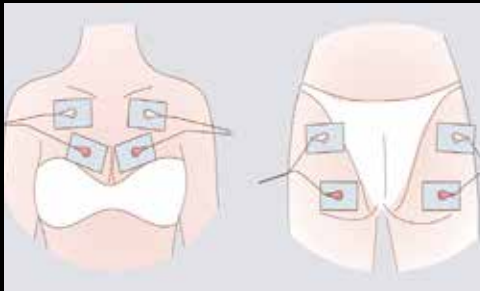

Upper Chest & Hips


Place the EMS patches that are connected to the device below your collar bone on each side of your chest. Then place the patches on your hips (front and back).



Benefits of whole body treatments

Negative Ions

Whole body treatments with the Luxe-Skin Revitalizing Treatment device provides the entire body with Negative Ions. The benefits of negative ions for your body range from lower inflammation, energizing your cells, tightening and firming your skin and muscles and much more!

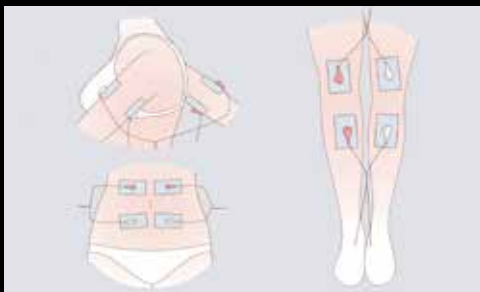

Whole Body Use

Apply the patches on multiple parts of your body. Arms, stomach, inner thighs, and calves are just some of the spots you can target with the patches.

G. Safety Guidelines

Oo Skin

- ⚠ Danger** Do not use the Myuffull EMS device in combination with the following devices as they may cause malfunction.
- 1) Implantable medical electronic devices such as pacemakers.
 - 2) Life-supporting medical electronic devices such as heart-lung machines.
 - 3) Wearable medical electronic devices such as electrocardiographs.

- ⚠ Caveat** Please consultate with your Medical Doctor if any of the following conditions apply, prior to use.
- | | |
|-------------------------|----------------------------|
| 1) Acute health illness | 5) Heart conditions |
| 2) Malignant tumours | 6) Sick or high fever |
| 3) Infectious diseases | 7) Abnormal blood pressure |
| 4) Pregnant | 8) Skin abnormalities |

Do not use near the heart, inside the mouth, near genital area, or areas where the skin is damaged, etc.
Do not use this device with any other electrical devices at the same time.
Do not use for any purpose other than low frequency treatment.
Do not repair or attempt to modify the Device. * Abnormal operation may result in injury.

- ⚠ Note** If you feel any abnormalities on your skin or body, stop using the device immediately.
* There is a risk of accidents and poor physical condition.



H. Product Specifications

Oo Skin

Power supply AC adapter (switching method) [AC-100V 50 / 60Hz output voltage: DC3V / 0.3A]

Battery 1.5V (AA) x 2

Output terminal Stereo mono jack ... 1 piece

Weight About 170g

External dimen- 213 (W) x 55 (H) x 49 (D) mm

Accessories 4-pole pad cable x 1
Single pole cable x 1
Portable pouch x 1
AA alkaline batteries (1.5V) x 2
Conductive gel pad x 4 small / 4 medium / 4 large
Pin holder x 1
AC adapter x 1
AC adapter connection cable x 1



B.I.D. International
Arch of Europe Gold Award
(2010 in Frankfurt, Germany)



OTHERWAYS
Platinum Technology Award
(2010 in Geneva, Switzerland)





Oo Japanese Spa
105-1428 Granville St.
Vancouver, BC V6Z 1N2
604 848 8470
info@oojapanesespa.com



@oojapanesespa