




**AAOCARE.COM**

22732 Granite Way  
Suite B  
Laguna hills CA 92653  
support@aaocare.com

Register your device today at [aaocare.com](http://aaocare.com)  
Receive our one year extended warranty.

Warranty will cover manufacturing defects and  
we will repair or replace your device at no cost.

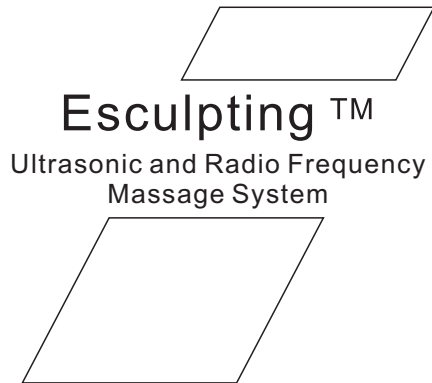

## 1.PRODUCT STRUCTURE




## 2.FEATURES

- 1.40000 times/s ultrasonic low frequency vibration. Instantly improve cell metabolism, promote fat block decomposition and dredges lymph nodes.
- 2.Powerful radio frequency shaping function with a heat effect on the dermal layer of the skin which enhances fat-burning.
- 3.Wavelength 640nm energy red light with high purity, powerful photosynthetic energy and uniform energy density characteristics that activates skin cells and promotes cell metabolism.
- 4.Mechanical vibration to induce collagen fiber contraction and accelerate the rate of fat decomposition which is effective in reducing stubborn fat.
- 5.Salon quality, zinc alloy massage head that feels extraordinary and is easy to operate.

## 3.SPECIFICATION


Ultrasonic Frequency	40KHz
Radio Frequency	10MHz
Adapter Input	110V or 220V
Operating Current	0.15-2A
Power	16W
Continuous use time	60minutes

## 4. HOW TO OPERATE

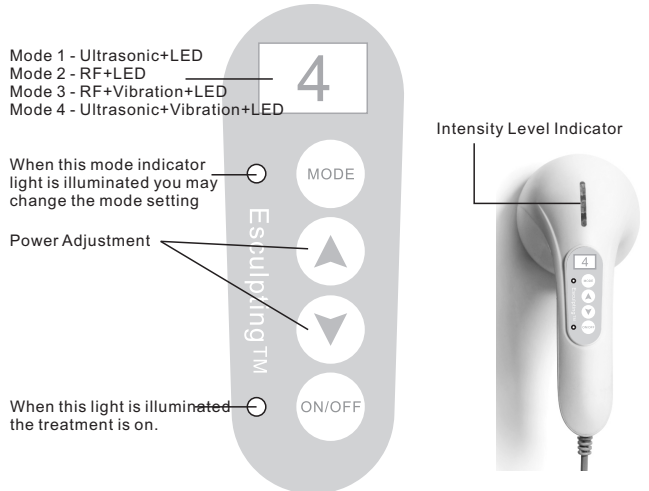

1. Plug cord into power supply. The default mode is mode 1.
2. Press mode to change to the desired treatment type. The mode indicator light is on the handle.
3. Press ON/OFF to start desired program.
4. After pressing the ON/OFF button initially, you can then change the level of intensity, by using the arrows.

## 5. DEVICE HEAD SURFACE

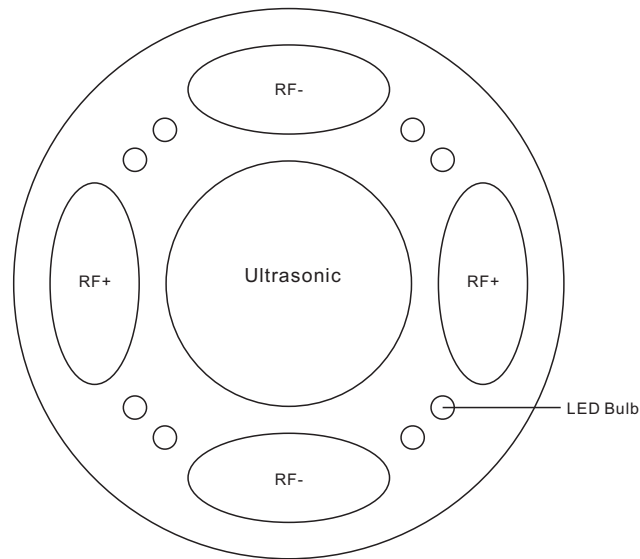


## 6.HOW TO USE ESCULPTING

1. Apply ultrasound gel or massage oil to the treatment area and spread evenly.
2. Select the appropriate treatment mode, Press On/Off button, Place head on desired area.
3. Hold the handle firmly on the treatment area and follow the specific instructions under treatment directions for the specific area.
4. Use the device for 5 to 10 minutes per treatment area.
5. Only use the device on the body. Do not use on the face, chin or chest.
6. Directly after each use, clean the treatment head with rubbing alcohol and a soft cloth.
7. For best results, use 2 to 3 times per week and expect results within 6 to 8 weeks.


## 7.TREATMENT DIRECTION

### 7-1.Abdomen central

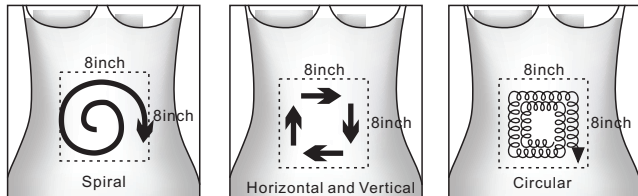


**Warning, do not to apply Probe to the navel, breast or the thorax area!**


### 7-2.Body Sides

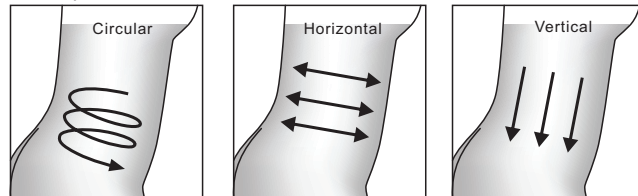


### 7-3.Hips Area

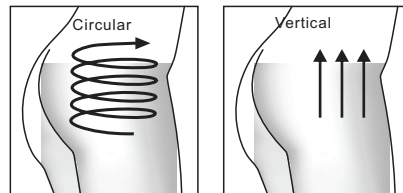


7-4. Buttocks

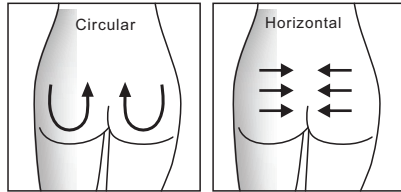


7-5. Inner thigh

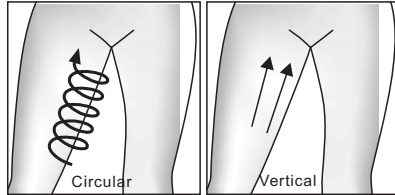

7-6. Back of the knee

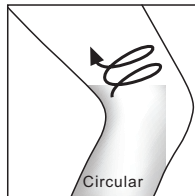

7-7. Arms

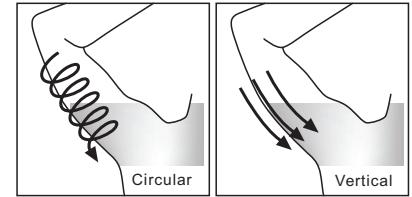

## 8.MODE INTRODUCTION

Do not exceed a 60 minute treatment time for whole treatment.

### 8.1.Mode 1:Fat Burn

Red light+Ultrasonic, 1-3 power level.

### 8.2.Mode 2:Skin tightening

Red light + Radio Frequency, 1-3 power level.

### 8.3.Mode 3:Fat burst

Red light+Radio Frequency+Mechanical vibration, 1-3 power level.

### 8.4.Mode 4:Shaping

Mechanical vibration +Ultrasonic, 1-3 power level.

## **9.TREATMENT TIPS**

- 1.If you have a lot of hair on treatment area, shave first, then use the Esculpting Device.
- 2.Always use ultrasound gel or massage oil with the device.
- 3.Do not exceed 4 treatments on the same area in one week.
- 4.Hydrate before treatment for best results.

## **10.WARNINGS AND PRECAUTIONS**

- 1.Do not use this product if you are pregnant or have a skin disease.
- 2.Avoid using the device within an hour after eating.
- 3.Do not overuse this device, overheating can cause machine failure.
4. Keep the device away from water at all times.
- 5.If skin appears abnormal after treatment, immediately stop using device and consult your doctor.
- 6.Store device in a dry, ventilated location.
- 7.Keep out of reach of children.
- 8.Do not use on: head, neck, collar bone, elbow, wrist, knee, ankle, foot or spine.