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GENIUS LIGHT

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GENIUS LIGHT

Multi-programme Hi-tech Care
8 care programmes
on demand.

INSTRUCTIONS FOR USE

Join the TALIKA community

www.talika.com

THE BIRTH OF A UNIQUE BRAND

TALIKA was created in 1948. During the Second World War, Danielle Roches, a nurse in Paris hospitals at the Department of Ophthalmology, formulated a plant-based cream to treat wounded eyes of burns victims. To her great astonishment, the balm's side effect was to boost lash growth of her patients giving rise to **LIPOCILS**.

EXPLORE - DISCOVER SHARE

Since this foundational milestone, **TALIKA** has continued to create unique skincare products, drawing on revelations from the living world, the most visionary scientific discoveries and the secrets of the Ancients. These skincare products allow women to realise and reveal all of their beauty potential.

HOW IT ALL BEGAN

"Who has never dreamed, like Aladdin, of having a magic lamp with their very own genie? » Today, **TALIKA** has just created such a lamp, employing its technological discoveries and know-how to enable women to realise their dreams with the new **GENIUS LIGHT**; the 1st multifunctional LED/electrostimulation* face mask inspired by aerospace research. A genuine hi-tech device, **GENIUS LIGHT** is the answer to most women's skin concerns, granting, at their request - as all good genies should - wishes for younger-looking, firmer, purified, refreshed and glowing skin, calmed from the harmful effects of sun and pollution. Step into the light!"

Alexis de Brosse
President & Founder

*by Talika

PRODUCT DESCRIPTION

For the first time ever,
TALIKA has brought together
8 Light Therapy® programmes combined
with electrostimulation/micro-current
functions, in a single hi-tech mask.

TECHNOLOGIES

8 care programmes on demand, combining
2 frontline technologies:

1. Light Therapy®

Using work done by NASA, we discovered that each wavelength has a specific cosmetic benefit which has been proven and identified.

- **Pulsing orange light (590 nm): anti-wrinkles**
Reduces fine lines and deep wrinkles (crow's feet, frown lines), firms the fine skin in the eye area.
- **Red light (630 nm): anti-redness**
Soothes and reduces redness, helps regenerate the skin.
- **Blue light (440nm): anti-blemishes**
Reduces blemishes and regulates the production of sebum.
- **Green light (525nm): anti-dark spots**
Reduces the appearance of dark spots and evens out skin tone.

2. Electrostimulation

Gentle transmission of electrical impulses ("lifting" or soothing effects) to stimulate the sub-dermal tissues.

Iontotherapy

- **Ionophoresis micro-current:**
an imperceptible micro-current which increases skin permeability to improve the penetration and efficacy of the active ingredients in cosmetic treatments.
- **Inverse ionophoresis micro-current:**
helps purifying the epidermis.

DESCRIPTIONS OF PROGRAMMES

1. Light (only)

| DAILY CURES once a day | | | |
|-----------------------------------|----------------------------------|--|--|
| | PRO- GRAMME NAMES | 1. GLOBAL ANTI- AGING | 2. FIRST WRINKLES |
| LIGHT THERAPY® | LED 630nm | ● | ● |
| | LED 590nm | ● | ● |
| | LED 525nm | ● | |
| | LED 440nm | | |
| | Ionto- phoresis | ● | ● |
| | Total duration | 3 min 45 | 1 min 45 |
| | Benefits | Firms skin Smooths wrinkles Boosts the effectiveness of cosmetic products | Firms skin Reduces wrinkles Boosts the effective- ness of cosmetic products |

1. Light (only)

| DAILY CURES once a day | | | |
|-----------------------------------|---|---|-------------------------------|
| PRO- GRAMME NAMES | | 3. ANTI- DARK SPOTS | 4. ANTI- BLEMISHES |
| LIGHT THERAPY® | LED 630nm | ● | ● |
| | LED 590nm | | |
| | LED 525nm | ● | |
| | LED 440nm | | ● |
| | Ionto- phoresis | ● | ● |
| | Total duration | 3 min | 6 min |
| Benefits | Reduces the appearance of dark spots Boosts the effec- tiveness of cosmetic products | Regulates sebum levels Improves skin texture Boosts the effectiveness of cosmetic products | |

| FLASH EFFECT once a day | | | |
|----------------------------|----------------|---|---|
| PRO-GRAMME NAMES | 5. DETOX | 6. ANTI-STRESS | |
| LIGHT THERAPY® | LED 630nm | ● | |
| | LED 590nm | | |
| | LED 525nm | ● | |
| | LED 440nm | | |
| | Iontophoresis | ● | ● |
| | Total duration | 2 min | 3 min |
| | Benefits | Helps to brighten the complexion Boosts the effectiveness of cosmetic products | Anti-redness Boosts the effectiveness of cosmetic products |

**FLASH EFFECT
once a day**

| | | | | FLASH EFFECT once a day | | |
|-----------------------|-----------|--|---|---|---------------------------|------------------|
| | | | | PRO- GRAMME NAMES | 7. RA- DIANCE BOOST | 8. AFTER- SUN |
| LIGHT THERAPY® | LED 630nm | | | ● | ● | |
| | LED 590nm | | | ● | | |
| | LED 525nm | | | | ● | |
| | LED 440nm | | | | | |
| Iontophoresis | | | ● | ● | | |
| Total duration | | | 3 min | 4 min | | |
| Benefits | | | <p>Firms skin, restores ra- diance and regenerates the skin Boosts the effective- ness of cosmetic products</p> | <p>Anti-redness Boosts the effectiveness of cosmetic products</p> | | |

TECHNOLOGIE

2. Light + Electrostimulation

| DAILY CURES once a day | | | |
|----------------------------------|-------------------------|---|--|
| PRO-GRAMME NAMES | | 1. GLOBAL ANTI-AGING | 2. FIRST WRINKLES |
| LIGHT THERAPY® | LED 630nm | ● | ● |
| | LED 590nm | ● | ● |
| | LED 525nm | ● | |
| | LED 440nm | | |
| ELECTROSTIMULATION | “Lifting” effect | ● | ● |
| | Soothing effect | | |
| | Total duration | 3 min 45 | 1 min 45 |
| | Benefits | Smooths deep wrinkles Firms skin Gives the complexion a pinky glow Reduces pore size Anti-dark spots Anti-redness boosts the effectiveness of cosmetic products | Reduces wrinkles and fine lines Firms the skin Reduces pore size and refines skin texture Boosts the effectiveness of cosmetic products |

2. Light + Electrostimulation

| DAILY CURES once a day | | | |
|---------------------------|------------------|--|--|
| PRO-GRAMME NAMES | | 3. ANTI-DARK SPOTS | 4. ANTI-BLEMISH |
| LIGHT THERAPY® | LED 630nm | ● | ● |
| | LED 590nm | | |
| | LED 525nm | ● | |
| | LED 440nm | | ● |
| ELECTROSTIMULATION | "Lifting" effect | | |
| | Soothing effect | ● | ● |
| | Total duration | 3 min | 6 min |
| | Benefits | Improves radiance and evenness of skin tone Reduces the appearance of dark spots Boosts the effectiveness of cosmetic products | Regulates sebum levels Improves skin texture Detox and anti-pollution effect |

FLASH EFFECT
once a day

| | | PRO-GRAMME NAMES | 5. DETOX | 6. ANTI-STRESS |
|---------------------------|------------------|------------------|--|---|
| LIGHT THERAPY® | LED 630nm | | | ● |
| | LED 590nm | | | |
| | LED 525nm | | ● | |
| | LED 440nm | | | |
| ELECTROSTIMULATION | "Lifting" effect | | | |
| | Soothing effect | | ● | ● |
| | Total duration | | 2 min | 3 min |
| | Benefits | | Detox and anti-pollution effects Helps brighten the complexion Improves skin texture | Soothing, anti-redness Improves well-being |

**FLASH EFFECT
once a day**

**PRO-
GRAMME
NAMES**

**7. RADIANCE
BOOST**

**8. AFTER-
SUN**

LIGHT THERAPY®

**LED
630nm**

**LED
590nm**

**LED
525nm**

**LED
440nm**

ELECTROSTIMULATION

**“Lifting”
effect**

**Soothing
effect**

**Total
duration**

3 min



4 min

Benefits

Gives skin a healthy pink glow
“Anti-fatigue” effect
Reduces the deepest wrinkles
Overall skin firming
Reduces pore size.

Soothing, anti-redness
Improves well-being
Helps prevent the appearance of brown spots

DEVICE OVERVIEW

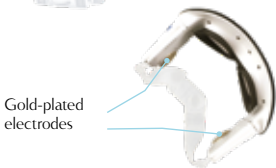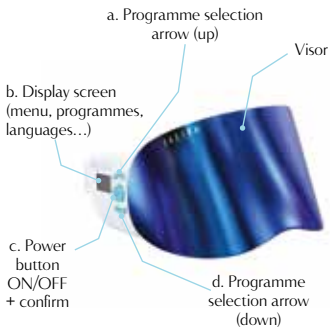DIRECTIONS FOR USE

Use the device on clean free of make-up skin and ensure that hair is tied back.

1. Pre-fit the **GENIUS LIGHT** mask to your face: put it on as you would glasses and check that the gold pads located on the sides of the mask are in contact with your skin. The pads must be positioned as close as possible to your cheekbones. If they feel too tight, remove the detachable supports (gold electrodes and white silicone) or alternatively if the mask is too big, add the clips to tighten it. Adjust the strap to your head. Remove the mask.

2. Turn **GENIUS LIGHT** on by pressing button (c) and select your chosen language using arrows (a) and (d), then confirm using button (c).

3. Select the **LIGHT** mode (using button (a)) if you don't want simultaneous electrostimulation, or **LIGHT + ELECTROSTIMULATION** if you want a programme with simultaneous electrostimulation, then confirm using button (c).

4. Choose the care programme of your dreams from the list below (using arrows (a) and (d) then confirming with button (c)), and the programme is ready:

1. GLOBAL ANTI-AGING

2. FIRST WRINKLES

3. ANTI-DARK SPOTS

4. ANTI-BLEMISHES

5. DETOX

6. ANTI-STRESS

7. RADIANCE BOOST

8. AFTER-SUN

5. Wet the area from the temples to the cheekbones (the area where the gold pads are in contact with the skin) with water, a gel or your usual lotion, to help conduction. The vibrations will feel strongest with water, as this is the most conductive element.

6. Fit **GENIUS LIGHT** to your face and fasten the strap.

7. The programme will start.

8. If you chose **LIGHT + ELECTRO-STIMULATION**, select the intensity that you would like using arrows (a) and (d) (intensity levels ranging from 1 to 20, a beep will sound each time you change level). If you chose **LIGHT**, activate the programme by confirming with button (c).

9. Carry on with a normal activity (computer, reading...); the programme will stop automatically. A beep will sound when the programme is finished.

During the electrostimulation programmes:

- In the event that you feel tingling sensations, reduce the intensity or stop the programme immediately or use a **LIGHT** mode only programme.
- In the event that you feel sensations in your mouth or teeth (those with fillings...), reduce the intensity or stop the programme immediately or use a **LIGHT** mode only programme.

NOTES

- Place the electrodes as close as possible to your cheekbones, by sliding the pads into the forward position. The electrodes are gold-plated.
- To go back, just select the back arrow shown on the screen. Select this arrow using buttons (a) and (c) to confirm.
- A beep will sound at the end of each programme to let you know that it has finished.
- Use the device on skin which is free of make-up to increase the effects of your **GENIUS LIGHT**.
- After using the device, you may continue to feel vibrations on your face for a few minutes and see slight redness in the areas where the pads have been in contact with the skin.
- We recommend that you do one programme per day. We advise a maximum of two “lifting effect” electrostimulation programmes, with an interval of eight hours.
- We advise you to alternate the **LIGHT / LIGHT + ELECTROSTIMULATION** programmes.
- When you have reached the maximum intensity level, a “triple beep” will sound.
- Remove any jewellery from the treatment area before use (piercings).
- If the device hits a hard surface, it automatically switches off.

CARE AND STORAGE

- Keep the device in its sleeve.
- Clean the bezel using a microfibre cloth. Clean the electrodes with a standard disinfectant.

RECHARGING

- To recharge the device: take the cable supplied with the mask. Connect the mini USB part of the cable to the mask and the USB part to a USB charger, which you plug into the mains or directly into the USB port of your computer.
- Once charged, the battery lasts for between 4 and 6 hours.

PACK CONTENTS

- 1 mask (built-in rechargeable lithium battery)
- 1 USB recharging cable
- 1 pouch
- 2 clips/electrodes
- 1 instructions leaflet
- 1 warranty card

WARNINGS

DO NOT USE THE DEVICE UNDER ANY OF THE FOLLOWING CIRCUMSTANCES:

- Not suitable for people having a metal implant at the level of the head
- If you are pregnant, breast feeding or if you have had a baby in the last 6 months.
- Following an episode of sunstroke (2 weeks)
- On skin which is "unwell" (irritation, cuts, allergy, swollen skin, infected skin...)
- On recent scars (less than 12 months)
- Sufferers of epilepsy
- Individuals who wear a pacemaker
- Individuals who wear braces (for electrostimulation programmes only)

In case of doubt, ask your doctor for advice.

WARNINGS

- Only use the device for the purposes described in the instructions leaflet.
- Keep the device out of reach of children.
- Avoid using the device on particularly sensitive or reactive skin. In case of doubt or skin conditions, seek medical advice before using and consult the Talika customer service team.
- In the event of a skin reaction during or just after use of the device, stop use immediately and seek medical advice.
- Protect the device from any violent knocks or from being shaken or dropped.
- Avoid prolonged exposure of the device to a very cold environment, direct sunlight or UV rays.
- Ensure the device does not get wet. Do not use the device in a damp environment (near a
 - bath, shower, swimming pool...).
- Do not use the device in an environment in which there is radiation being emitted from electromagnetic devices. This may cause the device to malfunction.
- Do not use the device if it is damaged.
- If the device is damaged, do not attempt to repair it. Contact the after-sales service, its approved repair service or any other person designated by the manufacturer.
- Do not use the device if you are allergic to metals.

WARRANTY

Please refer to the warranty card.

No replacement or refund will be given if the product is damaged by the user:

- Damage caused by improper use of the device by the user (impact, water exposure, over-use...).
- Damage caused by repairs not covered by the after-sales guarantee (see relevant leaflet).
- Damage caused intentionally by the user.

This device is designed for private, not commercial or professional use.

ENVIRONMENTAL PROTECTION

Please dispose of or recycle the device or battery in accordance with the local regulations which apply.

CUSTOMER SERVICE

customerservice@talika.com
www.talika.com

ALSO INCLUDED

- Day cream
- Serum
- Eye contour treatment

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